tria.®

ACNE CLEARING BLUE LIGHT

Quick Start Guide

tria.


TRIA ACNE CLEARING BLUE LIGHT ENJOY HEALTHY, CLEAR SKIN WITH FEWER BREAKOUTS

The clinically-proven Tria Acne Clearing Blue Light gives you the same high-energy, non-UV blue light as professional in-office treatments. Unlike traditional acne creams and medications, Tria's gentle blue light penetrates deep within the skin to naturally eliminate breakout-causing bacteria without irritation to clear existing blemishes and help prevent future breakouts before they begin.

Simply glide the Tria Blue Light over your face, chest, back or neck for 6 minutes a day, and start enjoying consistently clear skin.

ELIMINATE EXISTING BLEMISHES IN AS FEW AS 2 WEEKS
AND HELP PREVENT FUTURE BREAKOUTS

GET TO KNOW YOUR ACNE CLEARING BLUE LIGHT

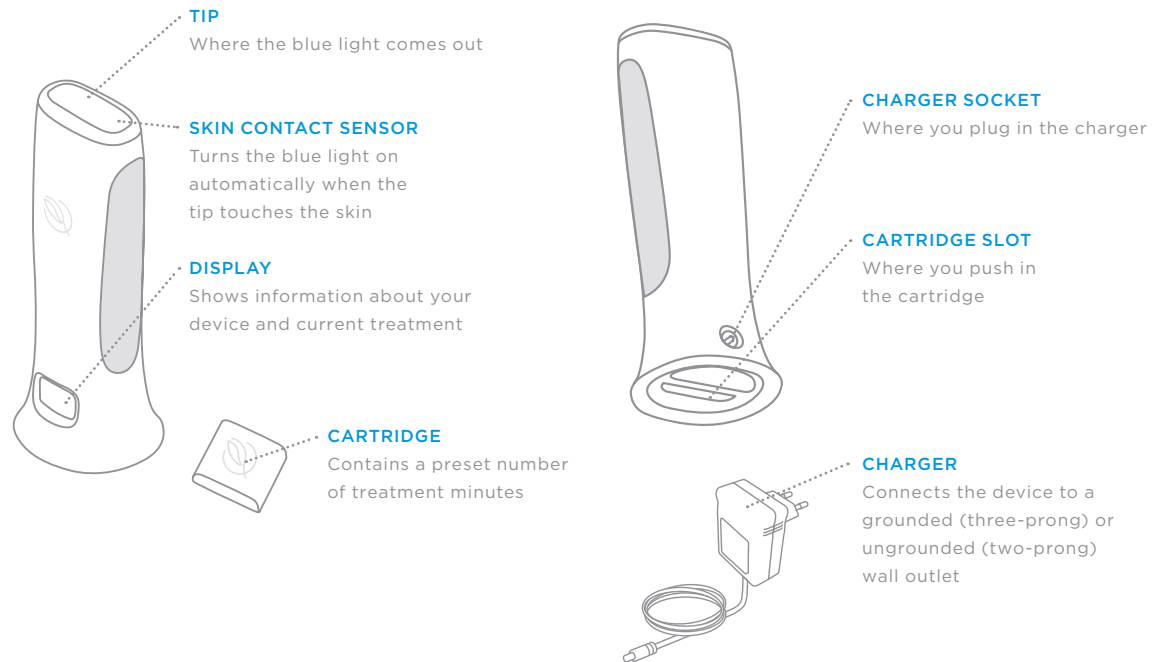DISPLAY

BATTERY INDICATOR
Shows how charged the battery is: blinks when it is completely out of power and you need to recharge it.

ELAPSED TIME
The large numbers show how long the current treatment has lasted (in minutes and seconds).

MINUTES REMAINING
The small numbers show how many minutes of treatment time are left on the cartridge.

When there is no more time on the cartridge, or if there is no cartridge in the device, the treatment time disappears, the picture of the cartridge blinks, and the display says "REPLACE CARTRIDGE."

GET STARTED

CHARGE AND LOAD YOUR DEVICE

To ready the device for first use, simply plug the charger into the wall and into the device.

After the device is fully charged, unplug the device from the charger. Then, insert the cartridge into the bottom of the device. Treatment time will illuminate. Your device is ready to use.

1

CLEANSE SKIN

To ensure maximum blue light effectiveness, cleanse your face thoroughly, removing all traces of makeup and sunscreen.

2

TREAT

To begin, press the tip of the device to your face. This will automatically turn on the device. Glide the tip over the treatment areas for 6 minutes once daily or 3 minutes twice daily.


For your convenience you can treat your skin in two ways:

1. Touch the Tria Blue Light to your skin and glide it over your face, chest, back or neck as if you are "painting" your skin with the light (see fig. 1).
2. Or, concentrate on those areas of most concern, such as the T-zone, chin or jaw.

3

NOURISH

Follow your Tria Blue Light treatment with a protective, antioxidant-rich facial serum.

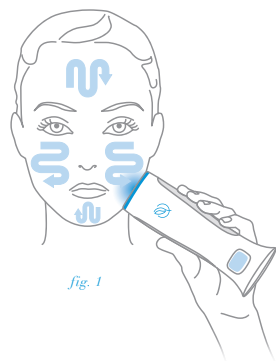Visit triabeauty.co.uk for a foam cleanser, specially formulated to exfoliate and prepare skin for the Tria Blue Light, and a serum, specially formulated to calm inflammation and reduce discolouration.

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TRIA BEAUTY CUSTOMER SERVICE / TRIA BEAUTY KUNDENDIENST

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customerservice@tribeauty.co.uk
tribeauty.co.uk
tribeauty.co.uk/instructions

Germany/Deutschland

customerservice@tribeauty.de
tribeauty.de
tribeauty.de/instructions

Austria/Österreich

customerservice@tribeauty.at
tribeauty.at
tribeauty.at/instructions

tria.

Please refer to the Instructions for Use booklet for complete information about the Tria Acne Clearing Blue Light, including warnings & indications for use.
©2014 Tria Beauty, Inc. 4160 Dublin Blvd., Suite 200, Dublin, CA 94568. All rights reserved. Tria, Tria Beauty, and the Tria logo are trademarks of Tria Beauty, Inc.
Protected by patent(s). See tribeauty.com/legal

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