tria[®]

AGE-DEFYING LASER

Instructions for Use

EN

WELCOME

WELCOME TO LASER SKIN REJUVENATION WITH TRIA BEAUTY. ENJOY PROFESSIONAL ANTI-AGEING RESULTS AT HOME.

With the Tria Age-Defying Laser you can treat multiple signs of ageing, restoring your skin's natural radiance. Your new device is a true revolution in personal skincare because it gives you professional anti-ageing laser treatment results through short, simple, daily treatments in the comfort of your own home. The laser light stimulates your body's own natural renewal process, regenerating cells and producing new collagen.

The result is more youthful-looking skin with fewer fine lines and wrinkles, a more even colour, and a smoother texture. Proven over a decade of clinical use and favoured by dermatologists, laser skin rejuvenation treats multiple signs of ageing safely and with better results than traditional anti-ageing creams and lotions. You'll have visible results in as few as two weeks.

Before you use your Tria Age-Defying Laser, please read these *Instructions for Use* and watch the *How to Use* video available on our website, for important safety information and to ensure best results. You can also find the *Quick Start Guide* on our website (see back cover).

TABLE OF CONTENTS


| | |
|--|----|
| GETTING TO KNOW YOUR TRIA AGE-DEFYING LASER | 2 |
| BEFORE USING YOUR TRIA AGE-DEFYING LASER FOR THE FIRST TIME | 4 |
| HOW TO USE YOUR TRIA AGE-DEFYING LASER | 4 |
| WHAT TO EXPECT | 7 |
| HELPFUL TIPS FOR BEST RESULTS | 8 |
| HOW OFTEN TO USE YOUR TRIA AGE-DEFYING LASER | 8 |
| TAKING CARE OF YOUR TRIA AGE-DEFYING LASER | 9 |
| TECHNICAL SPECIFICATIONS | 10 |
| TROUBLESHOOTING | 11 |
| IMPORTANT SAFETY INFORMATION | 14 |

GETTING TO KNOW YOUR TRIA AGE-DEFYING LASER


WHAT'S INCLUDED:

- TRIA AGE-DEFYING LASER
- CHARGING CRADLE
- POWER SUPPLY AND CORD
- INSTRUCTIONS FOR USE

IMPORTANT FEATURES:

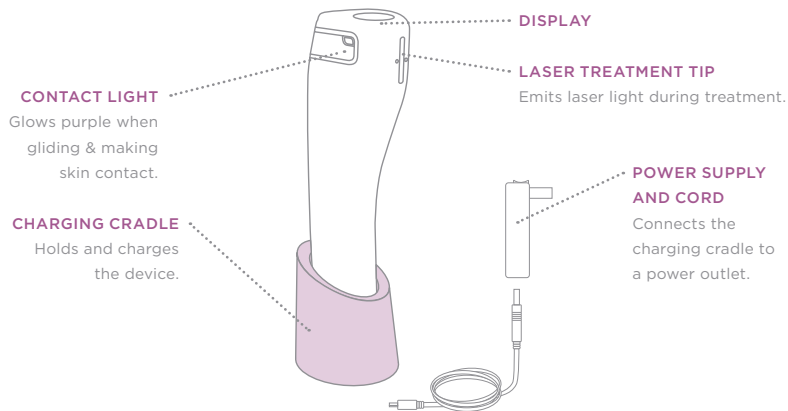
TRIA AGE-DEFYING LASER

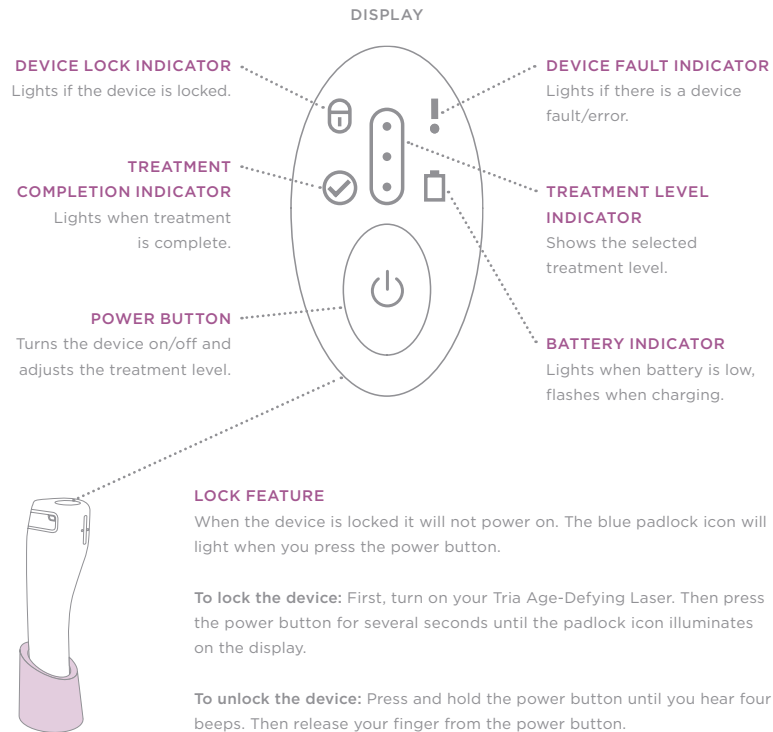BEFORE USING YOUR TRIA AGE-DEFYING LASER FOR THE FIRST TIME

CHARGE YOUR TRIA AGE-DEFYING LASER

- Insert the cord into the power supply and the charging cradle. Then plug the power supply into an outlet.
- Place your Tria Age-Defying Laser into the charging cradle.
- When your Tria Age-Defying Laser is charging, the battery indicator will flash. When fully charged (about 2½ hours), the battery indicator will turn off.
- For your convenience, your Tria Age-Defying Laser may be stored in its charging cradle between treatments. This ensures that it is always fully charged and ready for treatment.

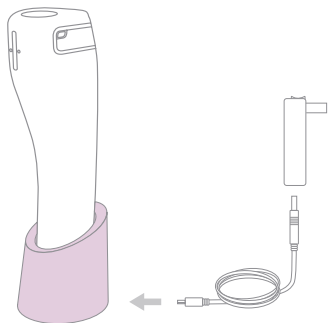WATCH THE HOW TO USE VIDEO

Visit our website for helpful information and techniques for using your Tria Age-Defying Laser. (see back cover).

HOW TO USE YOUR TRIA AGE-DEFYING LASER

THE TRIA AGE-DEFYING LASER IS INDICATED FOR THE TREATMENT OF WRINKLES, SKIN DISCOLOURATION AND ROUGH SKIN.



THE INSTRUCTIONS BELOW DESCRIBE HOW TO USE YOUR TRIA AGE-DEFYING LASER CORRECTLY.

1. PREPARE YOUR SKIN

- Wash and dry your face.
- For best results, use the Tria Priming Cleanser, which is specially formulated to remove makeup, sunscreen, and other products that can block the laser light and make it less effective.

2. SELECT YOUR TREATMENT LEVEL

- Remove your Tria Age-Defying Laser from its charging cradle and press the power button to turn it on.
- Use the power button to select your desired treatment level (Low, Medium or High). For convenience, your Tria Age-Defying Laser remembers the treatment level you last used.

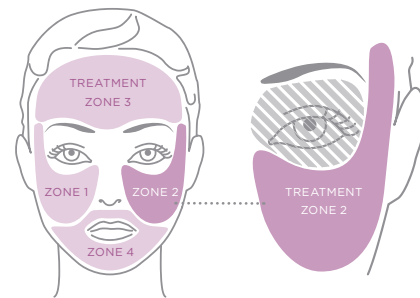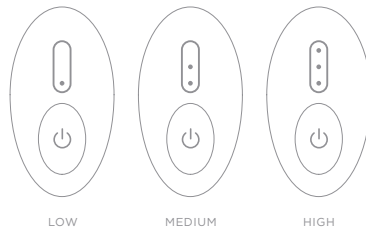
3. TREAT YOUR SKIN

- When treating, divide your face into 4 zones as shown in the picture below to ensure even treatment over the entire face. Treat zone 1 until the device beeps, then move on to zones 2, 3, and 4.


⚠ WARNING: Do not use to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary decrease in vision or permanent blindness. Treat around the eyes only where you can feel bone beneath the skin.

⚠ WARNING: Do not treat each zone more than once per day (once every 24 hours). Doing this could result in increased risk of skin irritation.


- To treat each zone, place the entire laser tip flat against your skin, and glide it slowly from side-to-side to cover the entire zone as shown in the pictures below. Keep the device moving and in good skin contact to keep the purple contact light glowing.

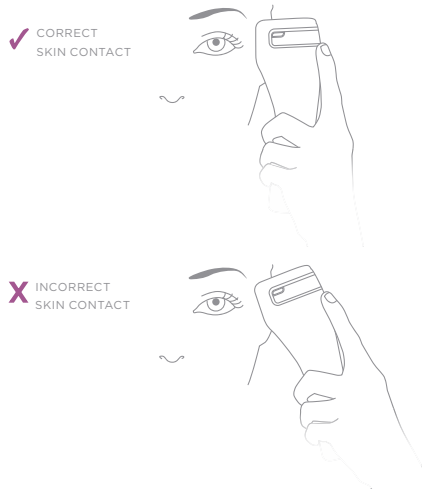- When treating your forehead you may find it helpful to hold your device horizontally and move it up and down. Remember to keep the device moving and in good skin contact to keep the purple contact light glowing.

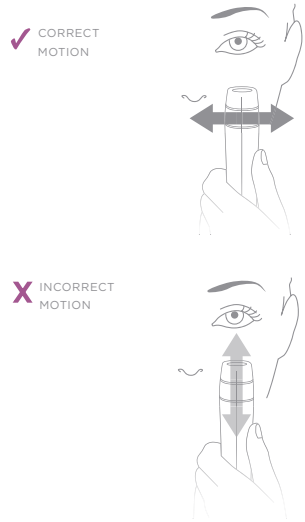- Treat each zone until your Tria Age-Defying Laser beeps. This will take about 30 seconds for each zone on Low, 1 minute for each zone on Medium and 2.5 minutes for each zone on High.

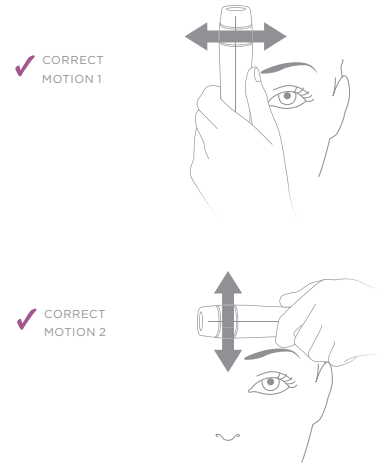- Your Tria Age-Defying Laser will beep once after the first zone is treated, twice after the second zone is treated and three times after the third zone is treated.
- After you treat all four zones, the treatment completion indicator (checkmark) will light,

a special tone will sound, and the Tria Age-Defying Laser will automatically turn off.

NOTE: Moving to the next zone before the beep may cause inadequate treatment of that zone. For best results, do not move to the next zone until you hear the beep.

4. NOURISH YOUR SKIN

- For best results, use the Tria Finishing Serum, which is specially formulated to improve your laser results by:
 - o Soothing and hydrating treated skin
 - o Reducing the appearance of fine lines and wrinkles
 - o Improving firmness and increasing elasticity

WHAT TO EXPECT

- We recommend making the Tria Age-Defying Laser treatment part of your evening skincare routine. This allows the skin renewal to begin and side effects to diminish overnight.
- Users typically notice an improvement in their skin in as few as 2 weeks, with the best results occurring a few weeks after the 8 week course of treatment.
- During or immediately after treatment, you may feel mild discomfort such as stinging, itching, or prickling. There can also be a warming or heat sensation similar to that caused by sunburn.

These sensations are normal and generally subside within minutes or hours and lessen with continued use.

- During or immediately after treatment, you may also develop temporary redness on your skin. This is normal and typically resolves overnight or within 3 days and lessens with continued use.
- Other side effects that may occur in normal use include:
 - o Skin swelling. This generally diminishes overnight.
 - o Flaking or dryness.
 - o Temporary small brown dots on the skin.
- Apply a soothing post-treatment serum or moisturiser to help alleviate skin discomfort after treatment. For best results try Tria's Finishing Serum.

NOTE: If you have more significant or persistent side effects, discontinue treatment and contact Tria Beauty Customer Service or a health care professional.

HELPFUL TIPS FOR BEST RESULTS

WATCH THE HOW TO USE VIDEO

Visit our website to watch our How to Use video and see how best to hold and move your Tria Age-Defying Laser across your skin (see back cover).

PRACTICE IN A MIRROR

Watch yourself in a mirror the first few times you treat. You'll see the glow of the purple contact light when you keep the device flat against your skin and continually moving.

TREAT AT NIGHT

Treat in the evening before bedtime to allow any redness or swelling to diminish overnight.

INCREASE YOUR TREATMENT LEVEL

It may take a week or so for your skin to adjust to the laser treatment. During this period you may want to use the Low treatment level. Once your skin has adjusted, use a higher treatment level for better results.

PROTECT YOUR SKIN FROM SUNLIGHT

Apply a premium quality sunscreen every morning. With daily sunscreen, you'll protect the benefits from your Tria Age-Defying Laser treatments, maintain results longer, and help prevent new sun damage.

HOW OFTEN TO USE YOUR TRIA AGE-DEFYING LASER

- For the first 8 weeks, use your Tria Age-Defying Laser once per day.

- Stop using your Tria Age-Defying Laser for the next 4 weeks to allow maximum cellular and collagen regeneration.
- After this 4 week period, you may begin another 8 week treatment cycle.

TAKING CARE OF YOUR TRIA AGE-DEFYING LASER

CLEANING

For best results, keep the laser treatment tip clean. Oils, lotions, makeup, or other debris can block the laser light or keep the device from functioning properly.

If the laser window or device tip does become visibly dirty, clean it with a lint-free cloth.

If other parts of the device or charging cradle become dirty, clean as needed with a soft cloth moistened with rubbing alcohol.

Make sure the charging cradle is not plugged in and do not get water or other liquid inside any part of the device or the charging cradle.

MAINTENANCE

The device does not need regular or special maintenance. If you think your device or charging cradle is not working properly or is damaged, do

not use the device. Contact Tria Beauty Customer Service for assistance.

STORAGE

Keep the device and charging cradle away from water, dampness, or debris that may get inside the device or the charging cradle.

DISPOSAL

This device contains lithium batteries. Please follow the regulations where you live for recycling or disposing of electrical devices.

TRAVEL

If you travel with your Tria Age-Defying Laser, pack the device and charging cradle in its original packaging to prevent damage.

The charging cradle will operate properly at all common household voltages (any supply between 100 and 240 VAC at 50/60 Hz). You may need a plug adapter so that the charging cradle can plug into the wall socket.

LIFE OF THE DEVICE

The device is designed to work properly for at least 3 years of normal use. There are no parts that need regular replacement. The device is working properly if it functions as described in this booklet.

TECHNICAL SPECIFICATIONS

MODEL: SRL

LASER TYPE: Diode laser

LASER CLASS: 1, in accordance with IEC 60825-1:2014

NOTE: This Class 1 laser device emits invisible laser radiation. Do not treat the eye, eyelids, or point directly at the eye. This device is safe to eyes and skin if used in accordance with these *Instructions for Use*.

OUTPUT WAVELENGTH: 1440 nm

OUTPUT LASER ENERGY: 5-12 mJ/pulse

BATTERY CHARGER INPUT:
100-240 V - 50/60 Hz, 0.15 A

BATTERY CHARGER OUTPUT:
5 V DC, 1.0 A, alternatively USB connection

PHYSICAL

OPERATING TEMPERATURE:

5°C to 30°C

OPERATING RELATIVE HUMIDITY:

Non-condensing 15-95%

OPERATING ATMOSPHERIC PRESSURE:

700 hPa to 1060 hPa

TRANSPORT AND STORAGE

TEMPERATURES:

-20°C to 50°C

TRANSPORT RELATIVE HUMIDITY:

Non-condensing 0-95%

TROUBLESHOOTING

Use this chart to help solve any problems that come up. If you still have problems, please contact Tria Beauty Customer Service.

| YOUR QUESTION | IT'S POSSIBLE THAT | HOW TO TROUBLESHOOT |
|--|--|--|
| The device fault indicator (!) is lit. | Your Tria Age-Defying Laser isn't working properly. | Turn your Tria Age-Defying Laser off, then back on. If the indicator remains lit, contact Tria Beauty Customer Service. |
| My Tria Age-Defying Laser doesn't charge when I place it in the charging cradle. | Your Tria Age-Defying Laser isn't connected properly. | Check that the device is sitting correctly in the charging cradle. Then check that the cable is properly connected at both ends and fully plugged into the wall outlet. (see page 4) |
| | The outlet in your home isn't working. | Try a different outlet. |
| | Your Tria Age-Defying Laser isn't sitting in the charging cradle properly. | Re-position your Tria Age-Defying Laser in the charging cradle, making sure it rests against the bottom of the cradle. Check that the battery indicator is flashing. |

| YOUR QUESTION | IT'S POSSIBLE THAT | HOW TO TROUBLESHOOT |
|---|---|---|
| When I am using my Tria Age-Defying Laser the contact light doesn't glow purple or stops completely during treatment. | Your Tria Age-Defying Laser isn't charged. | Charge your Tria Age-Defying Laser by placing it in the charging cradle. |
| | The laser treatment tip isn't making good contact with your skin. | Make a slow, steady gliding motion with good skin contact. (see page 6-7) |
| | Your Tria Age-Defying Laser isn't working properly. | Charge your Tria Age-Defying Laser and restart your treatment session. If the problem persists, contact Tria Beauty Customer Service. |
| When I am using my Tria Age-Defying Laser, the contact light illuminates intermittently. | The laser treatment tip isn't making good contact with your skin. | Make a slow, steady gliding motion with good skin contact. (see page 6-7) |
| | Dirt or debris is blocking the laser treatment tip. | Clean the laser treatment tip with a dry lint-free cloth. (see page 9) |

| YOUR QUESTION | IT'S POSSIBLE THAT | HOW TO TROUBLESHOOT |
|--|---|---|
| My Tria Age-Defying Laser won't turn off. | You accidentally released the power button too quickly. | Hold down the power button for 1-2 seconds, then release. |
| My Tria Age-Defying Laser won't turn on. | Your Tria Age-Defying Laser is locked. | Unlock your Tria Age-Defying Laser. (see page 3) |
| My Tria Age-Defying Laser won't lock. | Your Tria Age-Defying Laser isn't turned on. | First, your Tria Age-Defying Laser must be on. Then press the power button for at least four seconds until the padlock icon illuminates on the display. |
| My Tria Age-Defying Laser is locked/the device lock indicator is lit. | You accidentally locked your Tria Age-Defying Laser. | Unlock your Tria Age-Defying Laser. (see page 3) |

IMPORTANT SAFETY INFORMATION

Read the following important safety information before using the device. If you have any questions contact Tria Beauty Customer Service or your doctor.

WARNINGS

- Do not use to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as permanent or temporary decrease in vision or permanent blindness. Treat around the eyes only where you can feel bone beneath the skin.
- Do not treat skin that is infected, burned, or cut. Doing this could cause discomfort or irritation. Wait for these conditions to heal before using the device.
- Do not treat tattooed skin or areas with permanent makeup. Doing this could result in a damaged tattoo or increased risk of skin injury.
- Do not treat a zone on your face more than once per day (once every 24 hours). Doing this could result in increased risk of skin irritation.
- Do not treat if your skin is significantly red or swollen from a previous treatment. Wait until the redness and swelling have disappeared before treating again. If you have skin redness or swelling that lasts for more than a week, stop using the device and contact your doctor.

- Do not use the device if you experience skin colour darkening. Failure to follow these instructions may cause further skin colour darkening that may take more than 3 months to diminish.
- Do not use the device to treat skin which is painful to touch, raised (higher than surrounding skin), crusting or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- Do not use if you are pregnant. The device has not been tested for safety on pregnant women and the risks are unknown.
- Do not use if you are under the age of 18. The use of the device on or by anyone under the age of 18 has not been tested and could result in serious skin or eye injury.
- Do not use if the device or charging cradle is damaged.
- Do not use the device if it is wet, or if it or the charger could fall into water while it is plugged in. If the device or charger has fallen into water, do not touch until you first unplug it. Touching the device while it is wet could result in possibly fatal electric shock.
- Invisible laser light. Do not view directly with optical instruments. Class 1 Laser Product.

Failure to follow these warnings may result in potential skin injury.

PRECAUTIONS

To avoid potential device damage or injury:

- Do not treat body areas other than the face. Using this device on other areas of the body has not been tested and the risks are unknown.
- Applying treatment applications in the exact same spot without moving the device will increase redness and swelling and could potentially lead to blisters or scabs.
- If the area of skin for which you are using the device to treat wrinkles appears to be worsening, you should cease use and check with your doctor to make sure you do not have a more serious condition.
- Only store the device in a dry room at a temperature between -20°C and +50°C. Use the device only in a normal room temperature environment.
- Do not drop the device on hard surfaces, as this could damage it. Using a damaged device could result in serious skin or eye injury.
- Only use the provided charging cradle and power supply to recharge the device. Using different equipment could damage the device.

- Do not attempt to open or repair the device, charging cradle, or power supply.
- Do not use if you are undergoing professional cosmetic treatments, such as a skin peel, laser resurfacing or having toxins (e.g., Botox) or fillers injected.
- Do not use if you are taking a medicine, such as isotretinoin or other oral retinoid (e.g. Accutane), that requires you to avoid sunlight or that contains a steroid, or you have a medical condition that increases sensitivity to light.
- Do not use if you have a significant medical skin condition, your skin is unusually sensitive to light, or your skin easily develops a rash or an allergic reaction.

Failure to follow these precautions could result in skin injury, device damage or injury, including vision loss or potentially fatal electrical shock.

SYMBOL GLOSSARY

Warning Symbol

Laser Warning Symbol

Type BF Symbol

Serial Number Symbol

Manufacturer Symbol

Power On/Off Symbol

Accompanying Documentation
Must be Consulted Symbol

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Protected by patent(s). See tribeauty.com/legal. All rights reserved. Tria, Tria Beauty, and the Tria logo are trademarks of Tria Beauty. / Protegido por patente(s). Ver tribeauty.com/legal. Todos los derechos reservados. Tria, Tria Beauty y el logotipo de Tria son marcas registradas de Tria Beauty. / Protégé par un ou plusieurs brevets. Voir tribeauty.com/legal. Tous droits réservés. Tria, Tria Beauty et le logo Tria sont des marques de commerce de Tria Beauty. / Protetto da uno o più brevetti. Vedi tribeauty.com/legal. Tutti i diritti riservati. Tria, Tria Beauty, e il logo Tria sono marchi di fabbrica di Tria Beauty. / Geschützt durch Patent(e). Siehe tribeauty.com/legal. Alle Rechte vorbehalten. Tria, Tria Beauty und das Tria-Logo sind Marken von Tria Beauty. / Chronione patentem/-ami. Patrz: tribeauty.com/legal. Wszystkie prawa zastrzeżone. Tria, Tria Beauty i logo Tria są znakami towarowymi Tria Beauty. / Beschermd door octrooi(en). Zie tribeauty.com/legal. Alle rechten voorbehouden. Tria, Tria Beauty, en het Tria-logo zijn handelsmerken van Tria Beauty.
MDSS GmbH, Schiffgraben 41
■ Tria Beauty, a company of Channel Investments, LLC