tria[®]


AGE-DEFYING LASER

Quick Start Guide


EN

GET TO KNOW YOUR AGE-DEFYING LASER

Age-Defying Laser

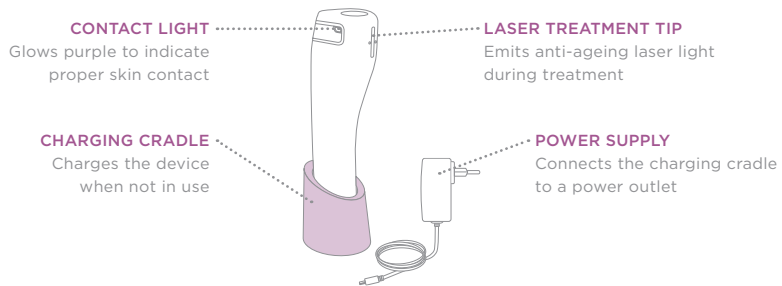Display

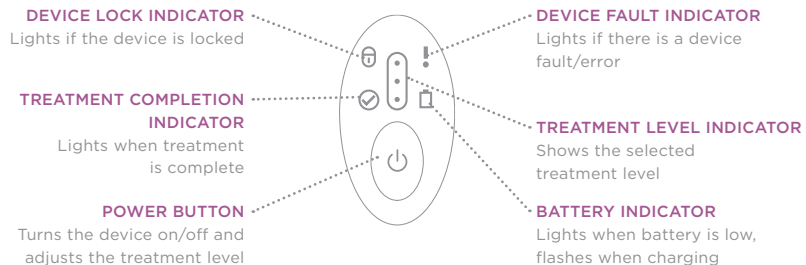TAKE OUR ADVICE

PRACTICE IN A MIRROR

Perfect your technique by using your device in front of a mirror. The purple light will glow when you are treating correctly.

INCREASE YOUR TREATMENT LEVEL FOR BEST RESULTS

Some women feel a slight stinging or warming sensation during treatment. This is normal and will diminish. To help your skin adjust, start on Low. Increase to Medium and then High for best results.

MAKE IT A NIGHTLY BEAUTY RITUAL

For best results, treat skin in the evening to allow any redness to diminish.

WEAR SUNSCREEN

Make it a routine to apply a broad-spectrum sunscreen in the morning and reapply during the day.

BEFORE FIRST TREATMENT

For detailed directions and safety information, read the complete *Instructions for Use*.

CHARGE: Charge the laser until the battery light turns off, indicating the device is fully charged (~2.5hrs).

HOW TO TREAT

Treat nightly by first cleansing and drying your skin. Treat with your laser, and then apply a post-treatment serum. For best results, use the Tria Priming Cleanser and Finishing Serum.


1 TURN ON DEVICE

Remove the laser from its charging cradle and press the power button.


2 SET THE TREATMENT LEVEL

Press the power button to adjust levels (Low, Medium or High). Your Age-Defying Laser defaults to the level you last used.

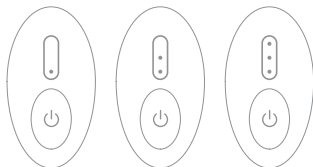3 START YOUR TREATMENT

- Divide the face into 4 zones. Start in Zone 1 by placing the laser tip flat against your skin and gliding it slowly from side-to-side.
- Cover the entire zone and maintain proper skin contact. The purple light will glow if you are treating correctly.
- Treat until you hear a beep and then progress to the next zone. The device will automatically shut off after you treat the fourth zone.

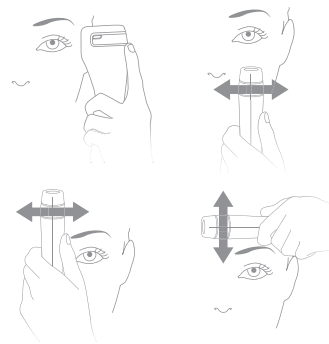
Increase levels for best results

As your skin adjusts, increase the treatment level for best results.

⚠ DO NOT use your Age-Defying Laser to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary decrease in vision or permanent blindness. Treat around the eyes only where you can feel bone underneath the skin.

TREAT THE FULL FACE BY ZONE


Treatment Zones

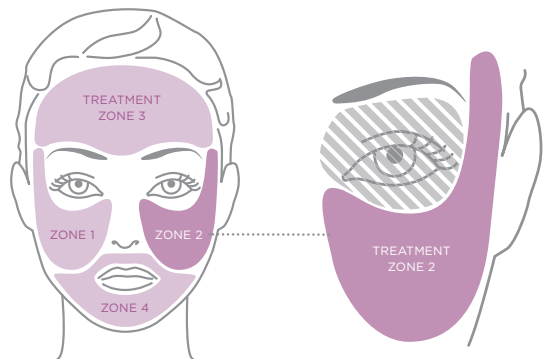Approximate Time Per Setting

Low	30 seconds per zone
Medium	1 minute per zone
High	2½ minutes per zone

Treatment Levels

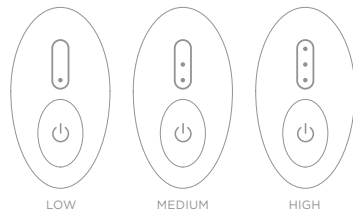GIVE YOUR ANTI-AGEING RESULTS A LIFT

*Enhance your Age-Defying Laser Experience
with our Cleanser and Serum.*

PRIMING CLEANSER

Prepare skin for laser treatment with this luxurious cleanser. Formulated with Willow Bark Extract to promote exfoliation and Birch Tree Sugars to help preserve moisture, our cleanser will leave your skin feeling refreshed and hydrated.

FINISHING SERUM

Our nourishing serum is formulated to be the perfect finish to your laser treatment. Sea Fennel helps to soothe and hydrate laser-treated skin, while Retinol and Biopeptides promote elasticity and firmness and help to improve the appearance of fine lines and wrinkles.

For more information on Tria's anti-ageing skincare products visit our website.

TRIA BEAUTY
CUSTOMER SERVICE

United Kingdom
customerservice@triaebeauty.co.uk
triaebeauty.co.uk
triaebeauty.co.uk/instructions

tria.

Refer to the *Instructions for Use* booklet for complete information about the Tria Age-Defying Laser, including warnings and indications for use.

©2014 Tria Beauty, Inc.

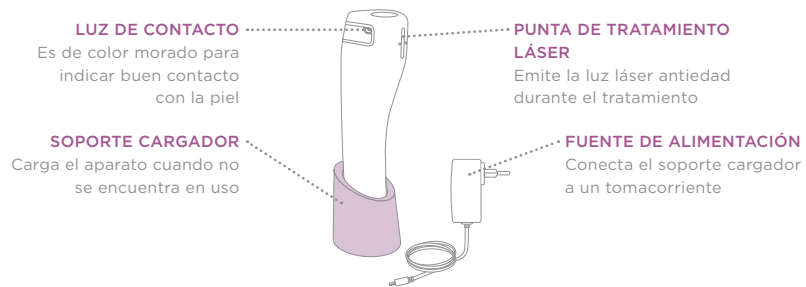
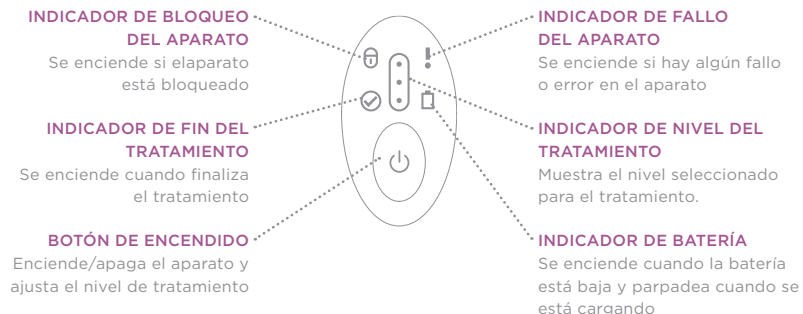
All rights reserved. Tria, Tria Beauty and the Tria logo are registered trademarks of Tria Beauty, Inc. Protected by patent(s). See triaebeauty.com/legal.

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