tria®

**AGE-DEFYING  
EYE WRINKLE CORRECTING LASER**

Instructions for Use

# TABLE OF CONTENTS

<b>GETTING TO KNOW YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER</b> .....	3
<b>WHO SHOULD USE THE AGE-DEFYING EYE WRINKLE CORRECTING LASER (INDICATIONS FOR USE)</b> ..	5
<b>WHEN NOT TO USE THE AGE-DEFYING EYE WRINKLE CORRECTING LASER (CONTRAINDICATIONS)</b> ...	6
<b>IMPORTANT SAFETY INFORMATION</b> .....	6
Warnings .....	7
Precautions .....	9
Benefits and Risks .....	10
<b>BEFORE USING YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER FOR THE FIRST TIME</b> .....	13
<b>HOW TO USE YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER</b> .....	14


<b>HOW OFTEN TO USE YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER</b> .....	17
<b>WHAT TO EXPECT</b> .....	17
<b>HELPFUL TIPS FOR BEST RESULTS</b> .....	18
<b>CARE INSTRUCTIONS</b> .....	19
<b>TROUBLESHOOTING</b> .....	21
<b>LIMITED WARRANTY</b> .....	25
<b>TECHNICAL SPECIFICATIONS</b> .....	31

## GETTING TO KNOW YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER

### WHAT'S INCLUDED:

- AGE-DEFYING EYE WRINKLE CORRECTING LASER
- INSTRUCTIONS FOR USE
- POWER SUPPLY AND CORD
- QUICK START GUIDE

### IMPORTANT FEATURES:

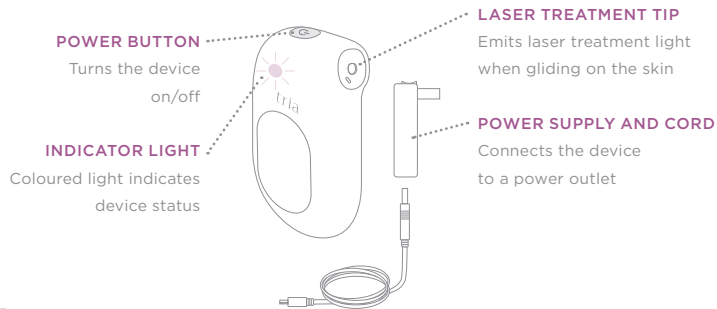See table below for important device indicators:

WHEN THE INDICATOR LIGHT IS . . .	IT MEANS . . .
Green	Steady green: Device is ready to use.
	Slowly flashing green: Device needs charging.
	Rapidly flashing green: Device is charging.
Purple (with vibration)	Device is gliding on skin and emitting laser treatment light.
Blue	Device is locked and must be unlocked for use.
No Light	Treatment is complete. Device has turned off.
Red	Device has an internal error. See troubleshooting guide.

### LOCK FEATURE

A lock feature is available on your device to prevent unintentional product use. When the device is locked, it will not power on. The indicator light will flash blue when you press the power button.

### To lock or unlock the device:

Press and hold the power button and release after four blue flashes. The power button must be released before the fifth blue flash to lock or unlock the device.

## WHO SHOULD USE THE AGE-DEFYING EYE WRINKLE CORRECTING LASER (INDICATIONS FOR USE)

The Age-Defying Eye Wrinkle Correcting Laser is indicated for the treatment of periorbital wrinkles (“crow’s feet”)

Use this chart to see if the Age-Defying Eye Wrinkle Correcting Laser is right for you. If you are not sure whether the device is right for you, do not use it until you have contacted Tria Beauty Customer Care or a health care professional for more information.

SKIN CONDITION	USE EYE WRINKLE CORRECTING LASER	DESCRIPTION	EXAMPLES
Periorbital wrinkles	✓ OK to use	Wrinkles around the eye (also called “crow’s feet”).	

## WHEN NOT TO USE THE AGE-DEFYING EYE WRINKLE CORRECTING LASER (CONTRAINDICATIONS)

Do not use Age-Defying Eye Wrinkle Correcting Laser on skin that is infected, burned or cut, as this could cause discomfort or irritation. Wait until the infection, burn or cut has healed before using the laser.

## IMPORTANT SAFETY INFORMATION

Read the following important safety information before using the device. If you have any questions, contact Tria Beauty Customer Care or a healthcare professional.

## WARNINGS

⚠ Do not use to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary decrease in vision or permanent blindness. Treat round the eyes only where you can feel bone beneath the skin.

⚠ Do not use on tattooed skin. Doing this could result in a damaged tattoo or increased risk of skin injury.

⚠ Do not treat each eye area more than once per day (once every 24 hours). Doing this could result in increased risk of skin irritation.

⚠ After treatment is complete and the device automatically turns off, wait at least 24 hours to begin a new treatment session. Failure to do so may increase skin irritation.

⚠ Do not use the device if you are under the age of 18. The use of the device on or by anyone under the age of 18 has not been tested and could result in serious skin or eye injury.

⚠ Do not use if you are pregnant. The device has not been tested for safety on pregnant women and the risks are unknown.

⚠ Do not use the device if it is wet or if it could fall into water while it is plugged in. If the device has fallen into water, do not touch until you first unplug it. Touching the device while it is wet could result in a possibly fatal electric shock.

⚠ Do not use the device if you experience skin colour darkening. Failure to follow these instructions may cause further skin colour darkening that may take more than 3 months to diminish.

⚠ Do not use the device to treat skin which is painful to touch, raised (higher than surrounding skin), crusting or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.

⚠ Do not treat if your skin is red or swollen from a previous treatment. Wait until the redness and swelling have disappeared before treating again. If you have skin redness or swelling that lasts for more than a week, stop using the device and contact a health care professional.

⚠ Invisible laser light. Do not view directly with optical instruments. Class 1 Laser Product.

⚠ No modification of this equipment is allowed.

## PRECAUTIONS

- Do not use the device on any area of the body other than round the eyes. Using the device on other areas of the body has not been tested and the risks are unknown.
- The device has only been clinically tested on periorbital wrinkles. Using the device on other wrinkle areas may not improve their appearance.
- Applying treatment applications in exactly the same spot without moving the device will increase redness and swelling and could potentially lead to blisters or scabs.
- Do not use the device if you have a skin condition or are recovering from professional cosmetic treatments, such as a skin peel, skin resurfacing, injectable toxins or fillers.
- If the area of skin for which you are using the device to treat wrinkles appears to be worsening, you should cease use and check with your doctor to make sure you do not have a more serious condition.
- Only store the device in a dry room at a temperature between -25°C to 70°C (-13°F to 158°F). Use the device only in a normal room temperature environment.
- Keep away from pets and children.
- Do not drop the device on hard surfaces, as this could damage it. Using a damaged device could result in serious skin or eye injury.
- Only use the provided power supply and cord to recharge the device. Using different equipment could damage it.
- Do not attempt to open or repair the device or power supply.

## BENEFITS AND RISKS

The benefits and risks of the Age-Defying Eye Wrinkle Correcting Laser were tested in a clinical study involving 45 participants, ages 35-70. Participants treated the skin round their eyes 7 days a week for 8 weeks and were evaluated throughout the study for wrinkle improvement and side effects.

### **BENEFITS**

In the study, the Age-Defying Eye Wrinkle Correcting Laser was evaluated for improving periorbital wrinkles in 45 men and women. Improvement was defined as at least a 1-point change on a 9-point scale, as judged by at least 2 of the 3 grading doctors. Participants also evaluated their own level of improvement.

Improvements in periorbital wrinkles were measured at 4 and 12 weeks after the last treatment by three independent doctors not involved in the clinical study. Blinded evaluations were performed by reviewing before and after photographs. Doctors did not know when the photographs had been taken. Four weeks after the last treatment, 34 of 45 participants (76%) had improvement in periorbital wrinkles. Twelve weeks after the last treatment, 33 of 45 participants (73%) had improvement in periorbital wrinkles. Please note that the observed improvement could be in part due to other non-device-specific effects. Based on subject satisfaction results, 35 of 41 participants (85%) noted improvement in their periorbital wrinkles 4 weeks after the last treatment and 34 of 41 participants (83%) noted improvement in their periorbital wrinkles 12 weeks after the last treatment.

Your results may vary because of how you use your Age-Defying Eye Wrinkle Correcting Laser or how your skin responds to treatment.

#### **RISKS**

The side effects that appeared most frequently after treatment in the clinical study included: redness (44 of 45 subjects), stinging or prickling sensation (42 of 45 subjects), warmth sensation (30 of 45 subjects) and burning sensation (28 of 45 subjects). Clinical results show that these side effects were usually resolved within several minutes or hours after a treatment. When using the Age-Defying Eye Wrinkle Correcting Laser daily, these side effects may diminish and resolve before your next daily treatment; however, depending on how your skin responds

to treatment, some may not fully resolve before your next daily treatment. Do not perform another treatment until all side effects have fully resolved.

Other side effects reported in the clinical study included: bronzing (temporary small brown dots on the skin) (14 of 45 subjects), dryness (5 of 45 subjects) and oedema (skin swelling) (3 of 45 subjects). Other side effects that were reported in one subject each included hyperpigmentation (skin colour darkening), acne, itching, tightness, stinging eyes, and watery eyes. In the clinical study, all these side effects were resolved within 3 months after stopping treatment. Not all of these side effects may have been caused by the Age-Defying Eye Wrinkle Correcting Laser.

Users with darker skin tones are more prone to persistent redness that can result in skin darkening. If persistent redness or skin darkening occurs, you should discontinue use of the Age-Defying Eye Wrinkle Correcting Laser. Sun exposure may worsen the condition. Apply a broad-spectrum sunscreen of at least SPF 30 daily to protect your skin from sunlight.

When you are treating your skin with the Age-Defying Eye Wrinkle Correcting Laser, you may feel discomfort that many participants describe as stinging or prickling. There can also be a warming or sunburn-like sensation that can last several minutes or hours after treatment. Many study participants reported these sensations lessened with continued use.


## BEFORE USING YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER FOR THE FIRST TIME

### CHARGE YOUR DEVICE

- Insert the cord into the power supply and the device. Then plug into a power point.
- While your device is charging, the indicator will rapidly flash green. When fully charged (up to 2 hours), the indicator light will turn steady green. Your device may arrive fully charged and not require additional charging before first use.

**NOTE:** The indicator light will slowly flash green when the battery needs to be charged. The indicator light will rapidly flash green while the Age-Defying Eye Wrinkle Correcting Laser is charging and turn steady green when fully charged.
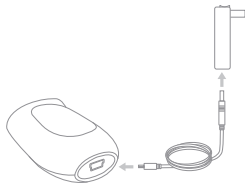
### TEST THE LASER ON AN AREA OF SKIN

If using the Age-Defying Eye Wrinkle Correcting Laser for the first time, we recommend you test the device on your inner forearm. This will help you determine what the treatment feels like. Most users report becoming used to the sensation after several days of treatment. If an adverse reaction such as blistering or scabbing occurs, do not use the device. This type of reaction is very rare.

## HOW TO USE YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER

THE INSTRUCTIONS BELOW DESCRIBE HOW TO USE YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER CORRECTLY.

### 1. PREPARE YOUR SKIN

Wash and dry your eye area.

### 2. CLEAN BEFORE USE

Your device is most effective when clean. Ensure optimal results by wiping the laser treatment tip with a dry cotton bud or cloth prior to each use.

### 3. TURN ON YOUR DEVICE

**Press and release the power button:** the indicator light will turn green to indicate that the Age-Defying Eye Wrinkle Correcting Laser is ready to

use. Do not hold the power button down while treating.

**NOTE:** You must unplug the device to perform a treatment.


### 4. TREAT YOUR FIRST EYE AREA

Identify the treatment area round your eye as shown below.

- **Slide the device on your skin:** Place the laser treatment tip flat against your skin and glide it slowly in small circular motions, performing multiple passes throughout your treatment to cover the entire area.

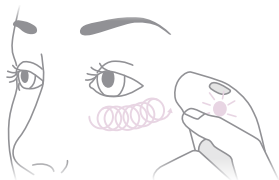- **Keep the purple light on:** The indicator light glows purple and the device vibrates when laser treatment is being delivered. Keep the device moving and in good skin contact to keep the purple light glowing.
- **Continue to treat until device turns off:** Treatment will complete after about 1 minute. The indicator light will go out, the device will vibrate on and off several times, and the device will automatically turn off.

**NOTE:** Moving to the other eye area before the device turns off may result in inadequate treatment of the first eye area.

**NOTE:** If the indicator light begins flashing green during treatment, complete your treatment before charging the device.

## 5. REPEAT ON YOUR OTHER EYE AREA

(opposite side of the face)

## 6. APPLY SUNSCREEN

Daily use of a broad-spectrum sunscreen of at least SPF 30 is recommended. Follow the sunscreen instructions.

**! WARNING:** Do not use to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary decrease in vision or permanent blindness. Treat round the eyes only where you can feel bone beneath the skin.

**! WARNING:** Do not treat each eye area more than once per day (once every 24 hours). Doing this could result in increased risk of skin irritation.

**NOTE:** It is normal to feel a slight stinging or warming sensation during treatment.

**NOTE:** If for any reason you need to turn off the device, press and hold the power button for 1-2 seconds to manually turn off the device.

**NOTE:** If the indicator light is flashing green before treatment, please charge the device, as this indicates the battery is low.

## HOW OFTEN TO USE YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER

Use your Age-Defying Eye Wrinkle Correcting Laser 7 days a week for 8 weeks. We recommend making the treatment part of your evening skin care routine. This allows the skin renewal process to begin and side effects to diminish overnight. Do not use this device for longer than an 8-week course of treatment.

## WHAT TO EXPECT

- Users typically notice an improvement in their skin in as few as 2 weeks, with the best results occurring a few weeks after the 8-week course of treatment.
- During or immediately after treatment, you may feel mild discomfort, such as stinging, itching or prickling. There can also be a warming or heat sensation similar to that caused by sunburn. These sensations are normal and generally subside within minutes or hours and lessen with continued use.
- During or immediately after treatment, you may also develop temporary redness or mild swelling on your skin. This is normal and will generally be resolved before your next treatment.

- Other side effects that may occur in normal use including flaking or dryness and temporary small brown dots on the skin.

**NOTE:** If you have more significant or persistent side effects, discontinue treatment and contact Tria Beauty Customer Care or a health care professional.

## HELPFUL TIPS FOR BEST RESULTS

### SOOTHE AFTER TREATMENT

Apply a post-treatment eye cream to help soothe skin after treatment.

### KEEP SKIN HYDRATED

During the 8-week course of treatment, apply a moisturiser as needed to keep skin well hydrated.

### PRACTISE IN A MIRROR

Perfect your technique by watching yourself in a mirror the first few times you apply the treatment. You'll see the glow of the purple contact light when you keep the device flat against your skin and continually moving.

## CARE INSTRUCTIONS

### TREAT AT NIGHT

Apply the treatment in the evening before bedtime to allow any redness or swelling to diminish overnight.

### PROTECT YOUR SKIN FROM SUNLIGHT

Apply a broad-spectrum sunscreen of at least SPF 30 every morning. With daily sunscreen, you'll protect the benefits from your Age-Defying Eye Wrinkle Correcting Laser treatments, maintain results longer, and help prevent new photodamage.

### MAINTENANCE

The device does not need regular or special maintenance. If you think your device is not working properly or is damaged in any way, do not use the device. Contact Tria Beauty Customer Care for assistance.

### CLEANING

For best results, keep the laser treatment tip clean. Oils, lotions, makeup or other debris can block the laser light or keep the device from functioning properly.

If the laser treatment tip does become visibly dirty, clean it with a dry cotton bud. Be sure to turn off the device before cleaning. If other parts of the device

become dirty, clean as needed with a soft cloth moistened with rubbing alcohol. Avoid getting water or other liquid inside any part of the device

### STORAGE

Keep the device away from water, dampness, or debris that may get inside the device.

### DISPOSAL

The device contains a lithium battery. Please follow the regulations where you live for recycling or disposing of electrical products.

### TRAVEL

If you travel with your device, pack the product in its original packaging to prevent damage. The power supply will operate properly at all common global voltages (any supply between 100 and 240 VAC at 50-60 Hz). You may need a plug adapter so the power supply can plug into the wall socket.

### LIFE OF PRODUCT

Age-Defying Eye Wrinkle Correcting Laser is designed to last for 3 years or more. There are no parts that need regular replacement. The device is working properly if it functions as described in this booklet.

## TROUBLESHOOTING

Use this chart to help solve any problems that come up. If you still have problems, please contact Tria Beauty Customer Care on 0800 0484 322 or [customerservice@triabeauty.co.uk](mailto:customerservice@triabeauty.co.uk).

YOUR QUESTION	IT'S POSSIBLE THAT	HOW TO TROUBLESHOOT
My Age-Defying Eye Wrinkle Correcting Laser won't turn on.	Your device isn't charged.	Charge your device (see page 13).
	Your device is plugged in.	Unplug the device before attempting to apply treatment.
My Age-Defying Eye Wrinkle Correcting Laser won't stay on and the indicator light flashes blue while the power button is depressed.	Your device is locked.	Unlock your device. (see page 4).
My Age-Defying Eye Wrinkle Correcting Laser won't turn off.	You released the power button before the green light turned off.	Hold down the power button until the green light goes off, then release.

YOUR QUESTION	IT'S POSSIBLE THAT	HOW TO TROUBLESHOOT
My Age-Defying Eye Wrinkle Correcting Laser won't lock.	You didn't hold the power button down long enough to see four blue flashes.	Repeat the locking sequence (see page 4).
	You didn't release the power button before the fifth blue flash.	
The indicator light is slowly flashing green.	Your device needs charging.	Charge your device (see page 13).
		If the indicator light begins flashing green during treatment, complete your treatment before charging the device.
The indicator light is red.	Your device isn't working properly.	Turn your device off, then back on. If the indicator light remains red, contact Tria Beauty Customer Care.

<b>YOUR QUESTION</b>	<b>IT'S POSSIBLE THAT</b>	<b>HOW TO TROUBLESHOOT</b>
<b>My Age-Defying Eye Wrinkle Correcting Laser doesn't charge.</b>	Your device isn't connected properly.	Check that the cord is properly connected at both ends and the power supply is plugged into a working power point. (see page 13).
<b>When I am using my Age-Defying Eye Wrinkle Correcting Laser, the contact light doesn't glow purple or stops completely during treatment.</b>	Your device isn't charged.	Charge your device (see page 13).
	The laser treatment tip isn't making good contact with your skin.	Make a slow, steady gliding motion with good skin contact. (see page 15).
	Your device isn't working properly.	Charge your device and restart your treatment session (see page 13).

<b>YOUR QUESTION</b>	<b>IT'S POSSIBLE THAT</b>	<b>HOW TO TROUBLESHOOT</b>
<b>When I am using my Age-Defying Eye Wrinkle Correcting Laser, the contact light illuminates intermittently.</b>	The laser treatment tip isn't making good contact with your skin.	Make a slow steady gliding motion with good skin contact. (see page 15).
	Dirt or debris is blocking the laser treatment tip.	Clean the laser treatment tip with a dry cotton bud. (see page 19).
<b>I am using my Age-Defying Eye Wrinkle Correcting Laser and the treatment sensation greatly diminished.</b>	Dirt or debris is blocking the laser treatment tip.	Clean the laser treatment tip with a dry cotton bud. (see page 19).

## LIMITED WARRANTY

The Age-Defying Eye Wrinkle Correcting Laser comes with a 1-year limited warranty detailed below.

### WHAT IS COVERED BY THE LIMITED WARRANTY?

This 1-year limited warranty covers defects in materials and workmanship in your Age-Defying Eye Wrinkle Correcting Laser if you purchased it from an authorised partner or direct from the manufacturer (Tria Beauty). Your 1-year limited warranty is valid for you alone, the original purchaser, and cannot be transferred.

### WHAT IS NOT COVERED BY THIS LIMITED WARRANTY?

This limited warranty is effective only if the Age-Defying Eye Wrinkle Correcting Laser is purchased and used in Europe. This limited warranty does not cover

damage, defects, or failure caused by abuse or misuse; improper care; impact or insertion of foreign objects; vandalism; improper storage; improper or incorrectly performed maintenance or repair; application of harmful chemicals; Acts of God, fire, or water; abrasives; negligence; failure to follow manufacturer's instructions; loss of the Age-Defying Eye Wrinkle Correcting Laser; normal wear and tear; alteration; over-use; professional use; or commercial use.

### HOW LONG DOES THIS LIMITED WARRANTY LAST?

The limited warranty begins on the date of purchase and expires one year later. The warranty period is not extended for any period during which the Age-Defying Eye Wrinkle Correcting Laser (or any parts thereof) is in repair or if your Age-

Defying Eye Wrinkle Correcting Laser is replaced by Tria Beauty during the warranty period.

### WHAT WILL TRIA BEAUTY DO IF YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER IS DEFECTIVE?

If you discover a covered defect in your Age-Defying Eye Wrinkle Correcting Laser before the 1-year warranty period expires, contact Customer Care for assistance.

If we confirm a defect covered by this limited warranty, it will be Tria Beauty's decision as to whether we will: (i) replace your Age-Defying Eye Wrinkle Correcting Laser with either a new or a refurbished Age-Defying Eye Wrinkle Correcting Laser that is the same as or similar to the Age-

Defying Eye Wrinkle Correcting Laser you purchased; or (ii) repair your Age-Defying Eye Wrinkle Correcting Laser using either new or refurbished parts.

If we find that the defect in your Age-Defying Eye Wrinkle Correcting Laser is not covered under this limited warranty, we will notify you about the service alternatives that are available on a fee basis.

If we find that your Age-Defying Eye Wrinkle Correcting Laser is not defective, we will notify you and return your Age-Defying Eye Wrinkle Correcting Laser to you. All Limited Warranty services must be authorised and performed by Tria Beauty or a service provider authorised by Tria Beauty.

### WARRANTY LIMITATIONS AND DISCLAIMERS

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND IS IN ADDITION TO YOUR STATUTORY RIGHTS.

YOUR REMEDIES FOR THE BREACH OF THIS LIMITED WARRANTY ARE LIMITED SOLELY TO THE REPAIR AND REPLACEMENT OF YOUR AGE-DEFYING EYE WRINKLE CORRECTING LASER AS SET FORTH IN THIS WARRANTY STATEMENT.

TO THE MAXIMUM EXTENT PERMITTED BY APPLICABLE LAW, TRIA BEAUTY DISCLAIMS ALL IMPLIED WARRANTIES. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF IMPLIED WARRANTIES ARE HEREBY EXPRESSLY LIMITED

TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH ABOVE. IN NO EVENT SHALL TRIA BEAUTY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, TORT (INCLUDING NEGLIGENCE OR STRICT PRODUCT LIABILITY), OR OTHERWISE, EVEN IF TRIA BEAUTY IS ADVISED OR AWARE THAT SUCH DAMAGES MIGHT ARISE. SUCH DAMAGES INCLUDE, BUT ARE NOT LIMITED TO, LOSS OF USE OF THE AGE-DEFYING EYE WRINKLE CORRECTING LASER, COST OF ANY SUBSTITUTE EQUIPMENT, THE CLAIMS OF THIRD PARTIES, OR DAMAGE TO PROPERTY.

TO THE MAXIMUM EXTENT PERMITTED BY

APPLICABLE LAW, IN NO EVENT SHALL TRIA BEAUTY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, STRICT PRODUCT LIABILITY, TORT OR OTHER LEGAL THEORY, EVEN IF TRIA BEAUTY IS ADVISED OR AWARE THAT SUCH DAMAGES MIGHT ARISE. SUCH DAMAGES INCLUDE, BUT ARE NOT LIMITED TO, LOSS OF USE OF THE AGE-DEFYING EYE WRINKLE CORRECTING LASER, COST OF ANY SUBSTITUTE EQUIPMENT, THE CLAIMS OF THIRD PARTIES, OR DAMAGE TO PROPERTY. IF THE APPLICABLE LAW DOES NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS, THE LIMITS IN

THIS PARAGRAPH AND THE DISCLAIMER OR EXCEPTIONS FROM WARRANTY MAY NOT APPLY. THIS LIMITATION DOES NOT APPLY TO CLAIMS FOR DEATH OR PERSONAL INJURY CAUSED BY A DEFECT IN THE AGE-DEFYING EYE WRINKLE CORRECTING LASER OR ANY OTHER LIABILITY THAT CANNOT BE EXCLUDED OR LIMITED BY LAW.

THESE WARRANTY TERMS, CONDITIONS, AND LIMITATIONS ALLOCATE THE RISKS OF PRODUCT FAILURE BETWEEN TRIA BEAUTY AND YOU, THE PURCHASER OF THE AGE-DEFYING EYE WRINKLE CORRECTING LASER. THIS ALLOCATION IS RECOGNISED BY BOTH PARTIES AND IS REFLECTED IN THE PRICE OF THE AGE-DEFYING EYE WRINKLE CORRECTING LASER. BY USING THE

AGE-DEFYING EYE WRINKLE CORRECTING LASER  
YOU ACKNOWLEDGE THAT YOU HAVE READ  
THESE TERMS, CONDITIONS, AND LIMITATIONS,  
UNDERSTAND THEM AND ARE BOUND BY THEM.

#### **PRODUCT CHANGES**

Tria Beauty reserves the right to make changes and improvements to its products without incurring any obligation to incorporate such improvements in products previously sold.

#### **DISPUTES**

This warranty and any dispute or claim arising out of or in connection with it or its subject matter or formation (including non-contractual disputes or claims) shall be governed by and construed in

accordance with the law of the country of purchase.

You agree that the courts of the country where you purchased the Age-Defying Eye Wrinkle Correcting Laser shall have non-exclusive jurisdiction to settle any dispute or claim that arises out of or in connection with this warranty or its subject matter or formation (including non-contractual disputes or claims).

#### **TIME FOR FILING CLAIMS**

Any claim for breach of this limited warranty must be commenced within twelve (12) months following the original purchase of the Age-Defying Eye Wrinkle Correcting Laser.

#### **NO OTHER WARRANTIES**

Unless modified in a writing signed by both parties, these warranty terms, conditions, and limitations are understood to be the complete and exclusive agreement between the parties, superseding all prior agreements and representations, oral or written, and all other communications between the parties relating to the subject matter of these warranty terms. No employee of Tria Beauty or any other party (including without limitation distributors, dealers, and retailers) is authorised to make any warranty in addition to those made herein.

#### **MISCELLANEOUS TERMS**

As used in this contract, “we” and “us” refer to Tria Beauty. “You” refers to the original purchaser of the Age-Defying Eye Wrinkle Correcting Laser. A waiver

of strict compliance with any term herein shall not be deemed a continuing waiver or waive any other terms or conditions. We shall be excused from performance to the extent of any delay or inability to perform arising due to weather events, disasters, acts of terror, governmental action or regulation, labor or materials shortages or disputes, or other circumstances not under our control. If any term herein is illegal, invalid, or unenforceable, then such term shall be ignored, but the remaining terms and conditions shall remain in effect. herein is illegal, invalid, or unenforceable, then such term shall be ignored, but the remaining terms and conditions shall remain in effect.

## TECHNICAL SPECIFICATIONS

**NAME:** Age-Defying Eye Wrinkle Correcting Laser

**MODEL:** FANp

**LASER TYPE:** Diode laser

**LASER CLASS:** 1 per IEC 60825-1:2014

**NOTE:** This Class 1 laser emits invisible laser radiation. Do not treat the eye, eyelids, or point directly at the eye. This product is safe to eyes and skin if used in accordance with these Instructions for Use.

**OUTPUT WAVELENGTH:** 1440 nm

**OUTPUT LASER ENERGY:** 10 mJ/pulse

**POWER SUPPLY INPUT:** 100-240 V-, 50-60 Hz, 0.15A

**POWER SUPPLY OUTPUT:** 5 V DC, 1.0 A,  
alternatively USB connection

### PHYSICAL

**OPERATING TEMPERATURE:**  
5 °C to 30 °C

**OPERATING RELATIVE HUMIDITY:**  
Non-condensing 15-95%

**OPERATING ATMOSPHERIC  
PRESSURE:**  
700 hPa to 1060 hPa

### TRANSPORT AND STORAGE

**TEMPERATURES:** -25 °C to 70 °C

**TRANSPORT RELATIVE HUMIDITY:**  
Non-condensing 0-93%

**TRANSPORT ATMOSPHERIC PRESSURE:**  
200 hPa to 1060 hPa

There are no parts that need regular replacement. The Age-Defying Eye Wrinkle Correcting Laser is working properly if it functions as described in this booklet.

This device complies with Part 15 of the US FCC Rules. Operation is subject to the following two conditions:  
(1) this device may not cause harmful interference, and  
(2) this device must accept any interference received, including interference that may cause undesired operation.

## SYMBOL GLOSSARY

Warning Symbol

Laser Warning Symbol

Type BF Symbol

Serial Number Symbol

Manufacturer Symbol

Power On/Off Symbol

Accompanying Documentation Must be Consulted Symbol

Temperature Limitation Symbol

WEEE Symbol

CE Mark

Class II Equipment

Date of Manufacture

Part Number

Class 1 Laser Product

Keep Dry

RoHS Compliance Symbol

Recycling Symbol

**TRIA BEAUTY**  
**CUSTOMER CARE**

0800 0484 322  
customerservice@triabeauty.co.uk  
triabeauty.co.uk

Complies with 21 CFR 1040.10 and 1040.11 except deviations pursuant to Laser Notice No. 50. Dated June 24, 2007.  
The device may potentially pose electromagnetic interference with other equipment. The device has been tested to be safe in accordance with IEC 60601-1-2:2007.

Protected by patent(s). See [triabeauty.com/legal](http://triabeauty.com/legal). ©2021 Tria Beauty. All rights reserved. Tria, Tria Beauty, and the Tria logo are registered trademarks of Tria Beauty.

 Tria Beauty, a company of Channel Investments, LLC  
14796 Wicks Blvd., San Leandro, CA 94577 USA

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MDSS GmbH