
tria®

HAIR REMOVAL LASER PRECISION

Quick Start Guide

EN


GET TO KNOW YOUR TRIA HAIR REMOVAL LASER PRECISION

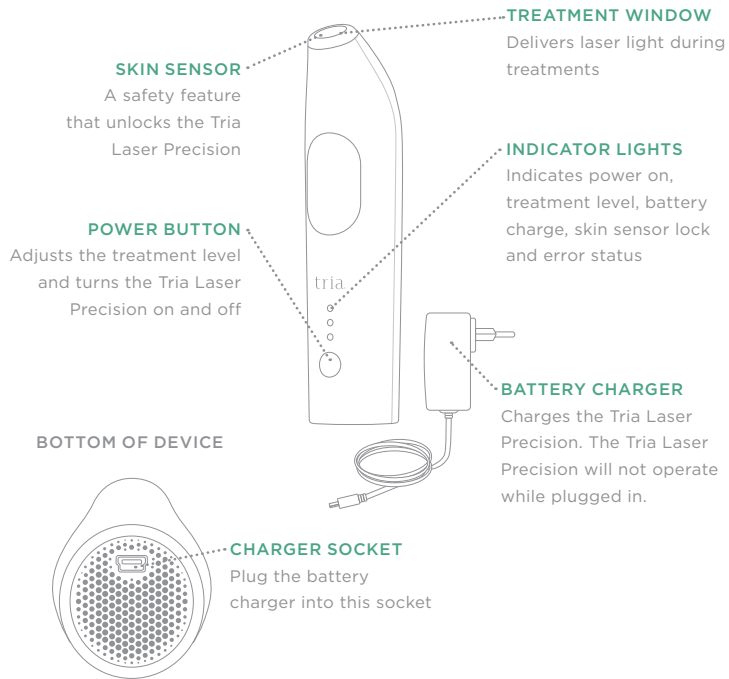
INDICATOR LIGHTS GUIDE

Treatment Levels

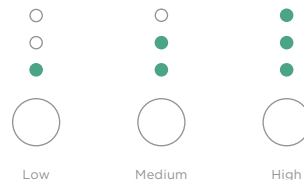Other Indicators

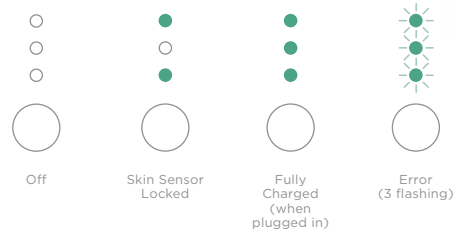BEFORE FIRST TREATMENT

For detailed directions and important safety warnings, please read the complete *Instructions for Use* guide before use.

CHARGE: You will also need to charge the Tria Laser Precision for at least 2 hours before you use it or anytime the battery is empty. Each battery charge lasts approximately 15-20 minutes.

HOW TO TREAT

Before every treatment, shave, clean and dry the area you wish to treat. Before your first full treatment on any new area, perform a patch test in that area. Repeat treatments every 2 weeks for up to 3 months or until you're satisfied with results.


1 TURN ON DEVICE

Unplug the Tria Laser Precision from the battery charger and press the power button.


2 UNLOCK

Unlock by holding the device tip with the built-in skin sensor against the skin you wish to treat. If the Tria Laser Precision is safe for your skin tone, the laser will emit a series of beeps and unlock.

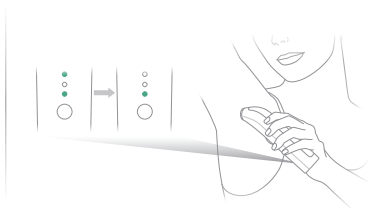
3 SET THE TREATMENT LEVEL

Press the power button to adjust levels. For best results, treat at the highest level that you can tolerate. In a clinical study, higher levels were found to be 2-3 times more effective.

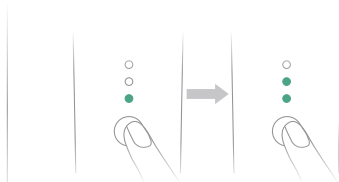4 TREAT. HOLD FOR ONE BEEP.

Touch and hold the treatment window to the skin until you hear a beep. Then, lift and place the laser back on the skin so that it overlaps the previous spot by about 1/4 inch, or 5 mm.

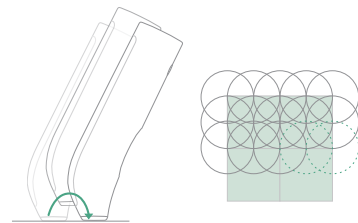
NOTE: It is normal for time between beeps to lengthen as you progress through each treatment. This should not affect your results.

TREATMENT GUIDE

Treat in an overlapping pattern until you have completely covered the area you want to treat. Use this guide as an example of how many minutes you may need to treat each area when treating on Medium. Treatment times will be shorter when treating on Low and longer when treating on High.

Always follow the treatment instructions and warnings provided in the *Instructions for Use* booklet.

- Each battery charge lasts for about 15-20 minutes. You may need to recharge if treating multiple areas.
- Facial hair is only cleared for women, cheek line downward.

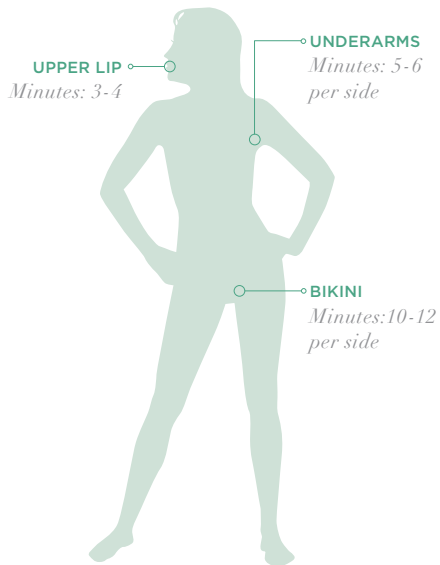TREATMENT AREA

The Tria Laser Precision is specially designed to treat small, sensitive areas, such as your underarms and bikini line. Say good-bye to constant shaving, waxing and skin irritations and say hello to silky, smooth skin.

CHECK YOUR SKIN TONE

The Tria Laser Precision is designed to work on fair to light brown skin tones and light brown to black body hair. It does not work on white, grey, blonde or red hair. If your skin tone matches any of the skin tones labeled "Not Safe" for treatment, then **DO NOT USE** the Tria Laser Precision.

SKIN TONE CHART

INTENDED SKIN TONES:

NOT SAFE:

**TRIA BEAUTY
CUSTOMER SERVICE**

United Kingdom

customerservice@triabeauty.co.uk

triabeauty.co.uk

triabeauty.co.uk/instructions

tria®

©2014 Tria Beauty, Inc. All rights reserved.

Tria, Tria Beauty, and the Tria logo are trademarks of Tria Beauty, Inc.
Protected by patent(s). See triabeauty.com/legal

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©2014 Tria Beauty, Inc. Alle rechten voorbehouden.
Tria, Tria Beauty, en het Tria-logo zijn handelsmerken van Tria Beauty, Inc.