tria[®]

**POSITIVELY CLEAR
ACNE CLEARING BLUE LIGHT**

Instructions for Use

WELCOME AND CONGRATULATIONS!

Congratulations on your new Positively Clear Acne Clearing Blue Light. This handheld device uses blue light to treat mild to moderate acne.

For important information about safety and how to use the device, read all the instructions in this booklet before you use the Positively Clear Acne Clearing Blue Light.

We want you to achieve the best results. If you have any questions or want more information, please call Customer Service at 1-877-321-8742.

TABLE OF CONTENTS

GLOSSARY	4	HOW TO USE YOUR POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT	14
ABOUT ACNE	5	Step-By-Step Instructions	
ABOUT POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT	6	How Often to Use Your Device	16
Device and Battery Indicators	7	How Long to Continue Treatments	
WHO SHOULD USE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT (INDICATIONS FOR USE)	8	TAKING CARE OF THE DEVICE	17
WARNINGS	9	TROUBLESHOOTING	18
RISKS AND BENEFITS	10	PRODUCT GUARANTEE	19
BEFORE USING THE DEVICE FOR THE FIRST TIME	12	WARRANTY	20
		TECHNICAL SPECIFICATIONS	23

GLOSSARY

BACTERIA: *Propionibacterium acnes* or *P. acnes* that causes acne breakouts. Lives on sebum secreted by the follicle.

BLACKHEAD: A small bump on the skin with a black center, not inflamed.

COMEDONAL ACNE: Blackheads and whiteheads.

CYST: A closed, bag-like structure under the skin, containing liquid. Acne cysts are present only in severe acne.

INFLAMMATORY/INFLAMED: Red and sensitive to the touch.

NODULE: A solid mass, like a knot, either raised or felt under the skin. Acne nodules are present only in severe acne.

PIMPLE: A red and sensitive (inflamed) bump on the skin.

PUSTULE: A red and sensitive (inflamed) bump on the skin filled with white or yellow pus.

WHITEHEAD: A small bump on the skin with a white center, not inflamed.

ABOUT ACNE

Acne starts with the buildup of oil under the skin, which causes a type of bacteria called *P. acnes* to multiply. Pores become clogged with this bacteria, excess oil, and dead skin cells, and the skin becomes inflamed. The result is acne blemishes, such as pimples and pustules. Blue light treats acne blemishes by killing the bacteria that cause them.¹

There are several kinds of acne. Comedonal acne consists of whiteheads and blackheads. Mild or moderate acne involves red and sensitive pimples and pustules. Nodules and cysts are signs of more severe acne. Positively Clear Acne Clearing Blue Light is intended only for mild or moderate acne (pimples and pustules).

¹ Kawada Y et al. Acne phototherapy with a high-intensity, enhanced, narrow-band blue light source: an open study and in vitro investigation. *Journal of Derm Sci* 30 (2002) 129-135.

ABOUT POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT

Your Positively Clear Acne Clearing Blue Light is a dermatologist-recommended device that uses blue light to treat acne.

TREATMENT TIP

Emits blue light during treatment

SKIN CONTACT SENSOR

Ensures the device emits blue light only when in contact with the skin

POWER BUTTON

Turns the device on and off

DEVICE INDICATORS


Indicate device status

BATTERY INDICATOR

Indicate battery status

CHARGING PORT

Connects the device to the USB Cable for charging

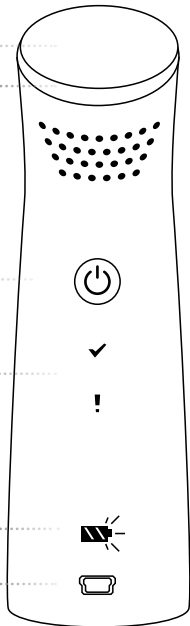TREATMENT GUIDE




Frames the treatment area on your skin

POWER SUPPLY AND USB CABLE

Connects the device to a power outlet

PACKAGE CONTENTS

- 1 Positively Clear Acne Clearing Blue Light
- 12 treatment guides
- 1 Power Supply
- 1 USB Cable
- Instructions for Use and Quick Start Guide

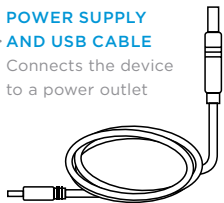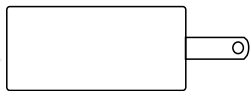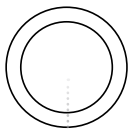IMPORTANT DEVICE INDICATORS

TREATMENT COMPLETION INDICATOR

Lights when treatment of the selected area is complete (5 mins). The device will also chime and then automatically turn off.

DEVICE ERROR INDICATOR

Lights if there is a device error. The device will also make a buzzing tone.

LOW BATTERY INDICATOR

Battery indicator flashes. The battery needs to be recharged. When the device is plugged in and charging, the battery indicator will show the battery filling up.

FULL BATTERY INDICATOR

Battery indicator shows a full battery (3 bars). The device is fully charged and ready to use.

WHO SHOULD USE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT (INDICATIONS FOR USE)

Positively Clear Acne Clearing Blue Light is indicated for the treatment of mild to moderate inflammatory acne.

Use this chart to see if Positively Clear Acne Clearing Blue Light is right for you.

ACNE TYPE	OK TO USE?	DESCRIPTION	EXAMPLE
MILD ACNE	✓ YES	Several inflamed pimples (red and sensitive bumps)	
MODERATE ACNE	✓ YES	Many inflamed pimples and pustules (red and sensitive bumps filled with white or yellow pus)	
COMEDONAL ACNE	✗ NO	Blackheads and whiteheads (small bumps on the skin with a white or black center, not inflamed)	
SEVERE NODULAR OR CYSTIC ACNE	✗ NO	Many pimples and pustules, with several nodules (solid masses of skin, like knots, either raised or felt under the skin) or cysts (closed, bag-like structures under the skin, containing liquid)	

WARNINGS

TO REDUCE THE RISK OF POTENTIALLY SERIOUS INJURY:

⚠ **DO NOT** use Positively Clear Acne Clearing Blue Light on or around your eyes or eyelids, and **do not** aim Positively Clear Acne Clearing Blue Light at your eyes or stare directly into the blue light. This could cause serious eye injury, including temporary blindness or temporary loss of vision.

⚠ **STOP USING** Positively Clear Acne Clearing Blue Light and see a health care professional if you feel pain or your skin becomes uncomfortably hot during treatment. Continuing to use Positively Clear Acne Clearing Blue Light could result in skin injury, such as temporary redness, discoloration, or blisters.

⚠ **STOP USING** Positively Clear Acne Clearing Blue Light and see a health care professional if you have skin redness that persists more than 24 hours after treatment. Failure to follow this warning may result in skin injury, such as temporary redness, discoloration or blisters.

⚠ **DO NOT** use Positively Clear Acne Clearing Blue Light to treat severe nodular or cystic acne. Delaying proper treatment

of severe nodular or cystic acne may be harmful. You should seek proper treatment without delay. See a healthcare professional.

⚠ **DO NOT** use Positively Clear Acne Clearing Blue Light to treat moles, warts, open sores, cancerous lesions, or any skin condition other than acne. Delaying proper treatment of these conditions may be harmful. You should seek proper treatment without delay. See a healthcare professional.

⚠ **DO NOT** treat the same area of skin more than once a day (once every 24 hours). Doing this could result in increased risk of skin injury, such as temporary redness, discoloration or blisters.

⚠ **DO NOT** use Positively Clear Acne Clearing Blue Light on children under the age of 13. **Do not** let anyone under 13 use Positively Clear Acne Clearing Blue Light. Positively Clear Acne Clearing Blue Light has not been tested on children under 13 and may cause serious skin or eye injury.

⚠ **DO NOT** share Positively Clear Acne Clearing Blue Light with other people. Sharing the device could transfer bacteria from one person to another and may cause harmful infection.

WARNINGS (continued)

⚠ DO NOT use Positively Clear Acne Clearing Blue Light if the device, power supply or USB cable is visibly damaged and never attempt to open or repair the device, power supply or USB cable. Doing this could result in a serious eye or skin injury or possibly fatal electric shock.

⚠ DO NOT use Positively Clear Acne Clearing Blue Light if the device, power supply or USB cable is wet, or if the device, power supply or USB

cable could fall into water while it is plugged in. Doing this could result in possibly fatal electric shock.

⚠ DO NOT use if you are pregnant. The device has not been tested for safety on pregnant women and the risks are unknown.

⚠ CAUTION - Performance of the procedures other than those specified herein may result in hazardous radiation exposure.

RISKS AND BENEFITS

The risks and benefits of Positively Clear Acne Clearing Blue Light treatments were tested in a clinical study of 40 individuals with acne.

RISKS

If you use Positively Clear Acne Clearing Blue Light as directed, the risks are very low. In the clinical study, none of the individuals treated (0 out of 40) experienced any side effects. However, the following mild side effects may occur with normal use:

SKIN WARMTH: The skin may feel warm during treatment.

TEMPORARY SKIN REDNESS:

The skin may develop a slight redness that fades within 24 hours.

TEMPORARY SKIN DISCOLORATION:

The skin may develop a slight discoloration that lasts more than 24 hours. If this occurs, stop using the device.

TEMPORARY SKIN DRYNESS:

The skin may develop slight dryness that goes away during the first 1-2 weeks of treatment.

If the device is used improperly or it malfunctions, there are other risks. The following chart lists some of these risks:

POTENTIAL ADVERSE EVENT	LIKELIHOOD OF POTENTIAL ADVERSE EVENT	POTENTIAL ADVERSE EFFECT	LIKELIHOOD OF POTENTIAL ADVERSE EFFECT
General overuse of the device, such as performing many treatments per day	Unlikely	Skin injury, such as temporary redness, discoloration or blisters	Minor skin injury is unlikely; major skin injury is very unlikely.
Use of device on severe acne	Unlikely	User delays known effective treatments	Likely
Use of device on or by children under 13	Unlikely	Same as for adults, although the severity of the effects may be different.	Unknown because use on children under 13 has not been tested
Use of device on non-acne skin conditions	Unlikely	User delays known effective treatments	Likely
Use of device on or around eyes or eyelids	Very unlikely	Blindness or temporary loss of vision in the exposed eye	Very unlikely
Use of physically damaged or wet device or charger that exposes user to electrical hazard	Very unlikely	Serious eye or skin injury, shock or electrocution	Very unlikely (same as other common electrical devices)
Device malfunction with excessive blue light output	Very unlikely	Skin injury, such as temporary redness, discoloration or blisters	Minor skin injury is unlikely; major skin injury is very unlikely.
Continued use despite side effects other than warmth, redness or dryness that quickly fades	Very unlikely	Skin injury, such as temporary redness, discoloration or blisters	Minor skin injury is unlikely; major skin injury is very unlikely.
Use of device despite skin sensitivity to light	Very unlikely	Skin injury, such as temporary redness, discoloration or blisters	Minor skin injury is unlikely; major skin injury is very unlikely.
Sharing a device with others who have common skin bacteria	Unlikely	Worse acne or a minor skin infection	Unlikely
Sharing a device with others who have serious skin infections	Very unlikely	Serious skin infection	Unlikely

BENEFITS

In the clinical study, daily use of Positively Clear Acne Clearing Blue Light was shown to eliminate about 70% of mild to moderate inflammatory acne blemishes (pimples or pustules) after 2 weeks of use. This study was conducted by a physician using a prescription version of the device. Your results may differ because of differences in the device, how you use it, or how your acne responds.


Blue light immediately starts killing the bacteria that cause acne,¹ but it takes some time for the blemishes on your skin to go away. You should start to notice a visible reduction in the number of blemishes after about 1 week of daily use, with improved results after 2 weeks of daily use. Individual results may vary.

¹ Kawada Y et al. Acne phototherapy with a high-intensity, enhanced, narrow-band blue light source: an open study and in vitro investigation. *Journal of Derm Sci* 30 (2002) 129-135.

BEFORE USING THE DEVICE FOR THE FIRST TIME

CHARGE YOUR DEVICE

Insert the USB cable into the power supply and into the device charging port. Then, plug the power supply into an outlet. When your device is charging, the battery indicator will show the battery filling up until the device is fully charged (less than 6 hours). Once the device is fully charged, the battery indicator bars will

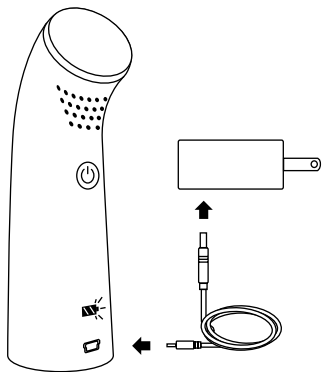stop filling up and display a full battery. Your device may arrive fully charged and not require additional charging before the first use.

TEST THE DEVICE ON AN AREA OF SKIN

Test Positively Clear Acne Clearing Blue Light on your forearm before you begin treatments. This is important because your skin may be sensitive to light due to medications you may have taken or other medical conditions. If this test shows that the device is not for you, do not use the device.

1. TREAT AN AREA OF SKIN ON YOUR FOREARM

Follow the step-by-step instructions in the "How to Use Your Positively Clear Acne Clearing Blue Light" section (pages 14-16). Apply one treatment guide and perform a full treatment inside the guide. It is normal for your skin to feel warm, but if you feel pain or your skin feels uncomfortably hot during treatment, **stop** using the device.

The device is not for you.

2. WAIT 24 HOURS AND LOOK AGAIN AT THE TREATED SKIN

No redness or difference: If you do not notice any difference between the treated and untreated skin, you may begin Positively Clear Acne Clearing Blue Light treatments.

Redness or difference: If you notice any difference between the treated and untreated skin, such as redness, discoloration, or discomfort, see a healthcare professional and do not use Positively Clear Acne Clearing Blue Light again.

NOTE: Repeat this test if you think your skin sensitivity may have changed for any reason, for example as a result of changes in your medication or other medical conditions.

HOW TO USE YOUR POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT

The instructions below describe how to use your Positively Clear Acne Clearing Blue Light correctly:

1. CLEANSE YOUR SKIN.

Cleanse and dry the skin you want to treat.

2. APPLY THE TREATMENT GUIDE.

Take one treatment guide, and place it on the area of skin that you wish to treat. Place the guide so that it "frames" an area with pimples or pustules. The area in the middle of the guide is the treatment area.


WARNING: Do not use the device on or around your eyes or eyelids, and do not aim the device at your eyes or stare directly into the blue light. This could cause serious eye injury, including temporary blindness or temporary loss of vision.

3. TURN ON THE DEVICE.

Press and release the power button. The device will beep. The battery indicator will display the current battery level. The device is now ready to use.

NOTE: The device will not turn on while plugged in. You must unplug the device before treating.

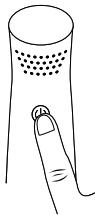
4. BEGIN YOUR TREATMENT.

Place the treatment tip against your skin inside the treatment guide. The device will now begin to emit blue light onto the treatment area. You will hear the fan whenever the device is emitting blue light.

Glide the treatment tip continuously and evenly across the entire treatment area. Maintain contact with the skin.

NOTE: It is important to keep the tip moving across the entire treatment area. It is normal to feel your skin become warm. Do not hold the device in one place long enough for your skin to become uncomfortably hot.

NOTE: If you lift the device off your skin, the device will stop emitting blue light and the fan will turn off. Place the treatment tip back onto the treatment area to continue treating.

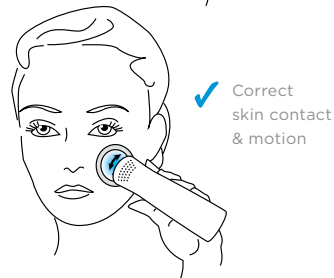NOTE: If you lift the device off your skin, the device will stop emitting blue light and the fan will turn off. Place the treatment tip back onto the treatment area to continue treating.

TAKING CARE OF THE DEVICE

5. COMPLETE YOUR TREATMENT

Continue gliding the device within the treatment area for 5 minutes. When treatment of the selected area is complete, the device will chime and the treatment completion indicator (✓) will light. The device will then automatically turn off.

WARNING: Stop using the device and see a health care professional if you feel pain or your skin becomes uncomfortably hot during treatment. Continuing to use the device could result in skin injury, such as temporary redness, discoloration, or blisters.

NOTE: If you need to turn off the device before your treatment is complete, press and hold the power button until the device beeps. Then, release. The device is now turned off.

6. TREAT OTHER AREAS AS NEEDED.

To treat another area of your skin, repeat steps 1-5

WARNING: Do not treat the same area of skin more than once a day (once every 24 hours). Doing this could result in increased risk of skin injury, such as temporary redness, discoloration or blisters.

HOW OFTEN TO USE YOUR DEVICE

Treat your skin with Positively Clear Acne Clearing Blue Light once a day (every 24 hours). For best results, do not skip treatments.

HOW LONG TO CONTINUE TREATMENTS

Continue to use for as long as you have visible pimples and pustules.

CLEANING

Clean the tip or the outside of the device whenever it appears dirty. Make sure the device is not plugged in, then moisten a soft cloth with water and wipe down the tip or the outside of the device. **DO NOT** get water or liquid inside any part of the device.

MAINTENANCE

The device does not need regular or special maintenance. If you think your device is not working properly or is damaged in any way, **do not** use the device. Call Customer Service at **1-877-321-8742** for assistance.

STORAGE

Keep the device away from water, dampness, or debris that may get inside the device. The device works best when charged at least once every six months when storing.

DISPOSAL

Please follow the regulations where you live for recycling or disposing of electrical devices.

TRAVELING WITH YOUR DEVICE

If traveling with the device, pack the device in its original packaging to prevent damage. The device will work properly at all common international voltages (any main supply between 100 and 240 VAC at 50-60 Hz). You may need an adapter so that the power supply that comes with the device can plug into a wall outlet. Use only the accompanying power supply and USB cable to charge the device.

LIFE OF THE DEVICE

The device is designed to work for 3 years or more of normal use. There are no parts that need regular replacement. The device is working properly if it functions as described in this booklet.

TROUBLESHOOTING

Use this chart to help solve any problems that come up.
If you still have problems, please call Customer Service at **1-877-321-8742**.

YOUR QUESTION	IT'S POSSIBLE THAT	HOW TO TROUBLESHOOT
My device will not turn on.	Your device is not charged.	Charge your device. (see page 12)
	Your device is plugged in.	Unplug the device before attempting to treat.
My device will not turn off.	You released the power button before the device beeped.	Hold down the power button until the device beeps, then release.
My device error indicator (!) is lit.	Your device is not working properly.	Turn your device off, then back on. If the indicator remains lit, contact Tria Beauty Customer Care.
My device does not charge.	Your device is not connected properly.	Check that the USB cable is properly connected at both ends and the power supply is plugged into a working wall outlet. (see page 12)
	Your device is already charged.	Turn on device. If the battery indicator shows a full battery (three bars lit), then your device is already charged and ready for use.
When I am using my device, the fan and blue light do not turn on or they stop completely during treatment.	Your device is not charged.	Charge your device. (see page 12)
	Your device is not turned on.	Press and release the power button. The device will beep. The battery indicator will display the current battery level.
	The treatment tip is not making good contact with your skin.	Make sure the treatment tip is in full contact with the skin when gliding the treatment tip across the treatment area. (see page 15)
		Lift treatment tip off skin and then place it back on skin in treatment area.
	If you heard the device chime and saw the treatment completion indicator (✓) light, your 5 minute treatment is complete.	Wait 24 hours to treat the same area of skin. To treat another area of skin, repeat steps 1-5. (see pages 14-16)
Your device is not working properly.	Charge your device and restart your treatment session. If the problem persists, contact Tria Beauty Customer Care.	
When I am using my device, the fan and blue light turn on and off intermittently.	The treatment tip is not making good contact with your skin.	Make sure the treatment tip is in full contact with the skin when gliding the treatment tip across the treatment area. (see page 15)

PRODUCT GUARANTEE

The Positively Clear Acne Clearing Blue Light comes with a 1-year limited warranty detailed below.

If you find that the Positively Clear Acne Clearing Blue Light is not right for you, or are unsatisfied with your device for any other reason, you can return it for a full refund within 90 days of purchase. Please call Customer Service at **1-877-321-8742**.

What is covered by the limited warranty?

This 1-year limited warranty covers defects in materials and workmanship in your Positively Clear Acne Clearing Blue Light if you purchased it from an authorized partner or directly from the manufacturer (Tria Beauty). Your 1-year limited warranty is valid for you alone, the original purchaser, and cannot be transferred.

What is not covered by this limited warranty?

This limited warranty is effective only if the Positively Clear Acne Clearing Blue Light is purchased and used in the United States or Canada. This limited warranty does not cover

damage, defects, or failure caused by abuse or misuse; improper care; impact or insertion of foreign objects; vandalism; improper storage; improper or incorrectly performed maintenance or repair; application of harmful chemicals; Acts of God, fire, or water; abrasives; negligence; failure to follow manufacturer's instructions; loss of the Positively Clear Acne Clearing Blue Light; normal wear and tear; alteration; overuse; professional use; or commercial use.

How long does this limited warranty last?

The limited warranty begins on the date of purchase and expires one year later. The warranty period is not extended for any period during which the Positively Clear Acne Clearing Blue Light (or any parts thereof) is in repair or if your Positively Clear Acne Clearing Blue Light is replaced by Tria Beauty during the warranty period.

What will Tria Beauty do if your Positively Clear Acne Clearing Blue Light?

If you discover a covered defect in your Positively Clear Acne Clearing Blue Light

before the 1-year warranty period expires, contact Customer Care for assistance.

If we confirm a defect covered by this limited warranty, it will be Tria Beauty's decision as to whether we will: (i) replace your Positively Clear Acne Clearing Blue Light with either a new or a refurbished Positively Clear Acne Clearing Blue Light that is the same as or similar to the Positively Clear Acne Clearing Blue Light you purchased; or (ii) repair your Positively Clear Acne Clearing Blue Light using either new or refurbished parts.

If we find that the defect in your Positively Clear Acne Clearing Blue Light is not covered under this limited warranty, we will notify you about the service alternatives that are available on a fee basis.

If we find that your Positively Clear Acne Clearing Blue Light is not defective, we will notify you and return your Positively Clear Acne Clearing Blue Light to you.

All Limited Warranty services must be authorized and performed by Tria Beauty or a service provider authorized by Tria Beauty.

WARRANTY LIMITATIONS AND DISCLAIMERS

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND IS IN ADDITION TO YOUR STATUTORY RIGHTS.

YOUR REMEDIES FOR THE BREACH OF THIS LIMITED WARRANTY ARE LIMITED SOLELY TO THE REPAIR AND REPLACEMENT OF YOUR POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT AS SET FORTH IN THIS WARRANTY STATEMENT.

TO THE MAXIMUM EXTENT PERMITTED BY APPLICABLE LAW, TRIA BEAUTY DISCLAIMS ALL IMPLIED WARRANTIES. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF IMPLIED WARRANTIES ARE HEREBY EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH ABOVE. IN NO EVENT SHALL TRIA BEAUTY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, TORT (INCLUDING NEGLIGENCE OR STRICT PRODUCT LIABILITY), OR OTHERWISE, EVEN IF

TRIA BEAUTY IS ADVISED OR AWARE THAT SUCH DAMAGES MIGHT ARISE. SUCH DAMAGES INCLUDE, BUT ARE NOT LIMITED TO, LOSS OF USE OF THE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT, COST OF ANY SUBSTITUTE EQUIPMENT, THE CLAIMS OF THIRD PARTIES, OR DAMAGE TO PROPERTY.

TO THE MAXIMUM EXTENT PERMITTED BY APPLICABLE LAW, IN NO EVENT SHALL TRIA BEAUTY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, STRICT PRODUCT LIABILITY, TORT OR OTHER LEGAL THEORY, EVEN IF TRIA BEAUTY IS ADVISED OR AWARE THAT SUCH DAMAGES MIGHT ARISE. SUCH DAMAGES INCLUDE, BUT ARE NOT LIMITED TO, LOSS OF USE OF THE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT, COST OF ANY SUBSTITUTE EQUIPMENT, THE CLAIMS OF THIRD PARTIES, OR DAMAGE TO PROPERTY. IF THE APPLICABLE LAW DOES NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS, THE LIMITS IN THIS PARAGRAPH AND THE DISCLAIMER OR EXCEPTIONS FROM WARRANTY MAY

NOT APPLY. THIS LIMITATION DOES NOT APPLY TO CLAIMS FOR DEATH OR PERSONAL INJURY CAUSED BY A DEFECT IN THE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT OR ANY OTHER LIABILITY THAT CANNOT BE EXCLUDED OR LIMITED BY LAW.

THESE WARRANTY TERMS, CONDITIONS, AND LIMITATIONS ALLOCATE THE RISKS OF PRODUCT FAILURE BETWEEN TRIA BEAUTY AND YOU, THE PURCHASER OF THE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT. THIS ALLOCATION IS RECOGNIZED BY BOTH PARTIES AND IS REFLECTED IN THE PRICE OF THE AGE-DEFYING EYE WRINKLE CORRECTING LASER. BY USING THE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT YOU ACKNOWLEDGE THAT YOU HAVE READ THESE TERMS, CONDITIONS, AND LIMITATIONS, UNDERSTAND THEM AND ARE BOUND BY THEM.

PRODUCT CHANGES

Tria Beauty reserves the right to make changes and improvements to its products without incurring any obligation to incorporate such improvements in products previously sold.

TECHNICAL SPECIFICATIONS

DISPUTES

This warranty and any dispute or claim arising out of or in connection with it or its subject matter or formation (including non-contractual disputes or claims) shall be governed by and construed in accordance with the law of the country of purchase.

You agree that the courts of the country where you purchased the Positively Clear Acne Clearing Blue Light shall have non-exclusive jurisdiction to settle any dispute or claim that arises out of or in connection with this warranty or its subject matter or formation (including non-contractual disputes or claims).

TIME FOR FILING CLAIMS

Any claim for breach of this limited warranty must be commenced within twelve (12) months following the original purchase of the Positively Clear Acne Clearing Blue Light.

NO OTHER WARRANTIES

Unless modified in a writing signed by both parties, these warranty terms, conditions, and limitations are understood to be the complete and

exclusive agreement between the parties, superseding all prior agreements and representations, oral or written, and all other communications between the parties relating to the subject matter of these warranty terms. No employee of Tria Beauty or any other party (including without limitation distributors, dealers, and retailers) is authorized to make any warranty in addition to those made herein.

MISCELLANEOUS TERMS

As used in this contract, "we" and "us" refer to Tria Beauty. "You" refers to the original purchaser of the Positively Clear Acne Clearing Blue Light. A waiver of strict compliance with any term herein shall not be deemed a continuing waiver or waive any other terms or conditions. We shall be excused from performance to the extent of any delay or inability to perform arising due to weather events, disasters, acts of terror, governmental action or regulation, labor or materials shortages or disputes, or other circumstances not under our control. If any term herein is illegal, invalid, or unenforceable, then such term shall be ignored, but the remaining terms and conditions shall remain in effect.

MODEL: BL2

ELECTRICAL INPUT

VOLTAGE: 100 - 240 V - at 50 - 60 Hz

CURRENT: 0.2 A max

CHARGER ELECTRICAL OUTPUT

VOLTAGE: 5 VDC

CURRENT: 1.0 A max

LIGHT OUTPUT

WAVELENGTH: 415 nm


PHYSICAL

SIZE: 4cm x 4cm x 14cm

WEIGHT: 227g

THIS DEVICE IS:

1. A Class II ME Equipment
 2. Internally powered in normal operation
 3. IP22 rated to protect particles >12.5mm entering the device and from dripping water.
 4. Classified in Risk Group 1 in accordance with IEC 60601-2-57: 2011. The charger provides isolation from supply mains.
- Caution - UV emitted from this device.

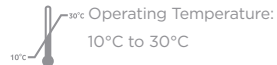
Operating Relative Humidity:

Up to 85%, Non-condensing

Operating Atmospheric Pressure:

700 hPa - 1060 hPa

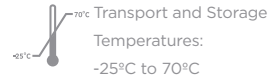Transport Relative Humidity:

Up to 90%, Non-condensing

 This device is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxygen.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

SYMBOL GLOSSARY

Warning Symbol

WEEE Symbol

Type BF Symbol

Serial Number Symbol

Date of Manufacture Symbol

Manufacturer Symbol

CE Mark Symbol

Temperature Limitation Symbol

Accompanying Documentation
Must be Consulted Symbol

FCC Declaration of Conformity

Class II equipment

RoHS Compliant

IP22 Ingress Protection Rating

Keep Dry Symbol

Recycling Symbol

Card widely recycled

TABLE 1 - GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC EMISSIONS

The Positively Clear Acne Clearing Blue Light is intended for use in the electromagnetic environment specified below. The customer or the user of the Sapphire should assure that it is used in such an environment.

EMISSIONS TEST	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
RF emissions CISPR 11	Group 1	The Positively Clear Acne Clearing Blue Light must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected when RF is energized.
RF emissions CISPR 11	Class B*	The Positively Clear Acne Clearing Blue Light is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not Applicable	The Positively Clear Acne Clearing Blue Light is suitable for use in all establishments, other than domestic and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not Applicable	

*Meets Class B limits but is not intended for use in domestic establishments.

TABLE 2 - GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY

The Positively Clear Acne Clearing Blue Light is intended for use in the electromagnetic environment specified below. The customer or the user of the Positively Clear Acne Clearing Blue Light should assure that it is used in such an environment.

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV Contact ± 8 kV Air	± 6 kV Contact ± 8 kV Air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV Mains ± 1 kV I/Os	± 2 kV Mains ± 1 kV I/Os	Mains power quality should be that of a typical home environment.
Surge IEC 61000-4-5	± 1kV Differential ± 2kV Common	± 1kV Differential ± 2 kV Common	Mains power quality should be that of a typical home environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	>95% Dip for 0.5 Cycle 60% Dip for 5 Cycles 30% Dip for 25 Cycles >95% Dip for 5 Seconds	>95% Dip for 0.5 Cycle 60% Dip for 5 Cycles 30% Dip for 25 Cycles >95% Dip for 5 Seconds	Mains power quality should be that of a typical home environment. If the user of the Positively Clear Acne Clearing Blue Light requires continued operation during power mains interruptions, it is recommended that the Positively Clear Acne Clearing Blue Light be powered from a battery.
Power frequency (60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical home environment.

TABLE 3 - GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY

The Positively Clear Acne Clearing Blue Light is intended for use in the electromagnetic environment specified below. The customer or the user of the Positively Clear Acne Clearing Blue Light should assure that it is used in such an environment.

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	(V1)=3Vrms	Portable and mobile RF communications equipment should be separated from the Positively Clear Acne Clearing Blue Light by no less than the distances calculated/listed below.
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	(E1)=3V/m	$D=(3.5/V1)(\text{Sqrt } P)$ 150kHz to 80MHz $D=(3.5/E1)(\text{Sqrt } P)$ 80 to 800 MHz $D=(7/E1)(\text{Sqrt } P)$ 800 MHz to 2.5 GHz where P is the max power in watts and D is the recommended separation distance in meters. Field strengths from fixed transmitters, as determined by an electromagnetic site survey, should be less than the compliance levels (V1 and E1). Interference may occur in the vicinity of equipment containing a transmitter.

tria®

TABLE 4 - RECOMMENDED SEPARATION DISTANCES BETWEEN PORTABLE AND MOBILE RF COMMUNICATIONS EQUIPMENT AND THE POSITIVELY CLEAR ACNE CLEARING BLUE LIGHT

The Positively Clear Acne Clearing Blue Light is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Positively Clear Acne Clearing Blue Light can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Positively Clear Acne Clearing Blue Light as recommended below, according to the maximum output power of the communications equipment.

MAXIMUM OUTPUT POWER W	SEPARATION (M) 150 KHZ TO 80 MHZ $D = (3.5/V1)(\text{SQRT } P)$	SEPARATION (M) 80 MHZ TO 800 MHZ $D = (3.5/E1)(\text{SQRT } P)$	SEPARATION (M) 800 MHZ TO 2,5 GHZ $D = (7/E1)(\text{SQRT } P)$
0,01	0.11667	0.11667	0.23333
0,1	0.36894	0.36894	0.73785
1	1.1667	1.1667	2.3333
10	3.6894	3.6894	7.3785
100	11.667	11.667	23.333

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© 2019. All rights reserved. Tria, Tria Beauty, and the Tria logo are trademarks of Tria Beauty.
Channel Investments, LLC
MDSS GmbH, Schiffgraben 41