


# Quick Start Guide

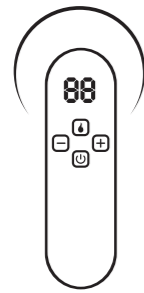

## Fat Burning Iron F666

### Iron off stubborn fat and stretch marks

- Thank you for using the fat burning iron F666. Please read Quick Start Guide before using, and follow its instruction.
- In order to use this product safely, please read the "Important Safety Tips" in the Quick Start Guide carefully.
- After reading, please keep manual in a safe place for reference at any time.
- Please keep the product out of the reach of children.
- This product belongs a non-medical device.

### Important Safety Tips

Tips	
Users Disabled	<ul style="list-style-type: none"> <li>• Patients who use pacemakers or have serious arrhythmias.</li> <li>• Patients with cognitive impairment and skin sensory paralysis.</li> <li>• Babies, pregnant women, parturients (60 days after surgery).</li> </ul>
Parts Forbidden	<ul style="list-style-type: none"> <li>• Implanted sites such as metal and plastic.</li> <li>• Skin ulcers.</li> <li>• The heart, both carotid sinuses and the front of the throat.</li> </ul>
Action Prohibited	<ul style="list-style-type: none"> <li>• Do not immerse this product in water.</li> <li>• It is forbidden to use it while operating the other machine or driving.</li> <li>• It is prohibited to disassemble, modify and repair this product by yourself.</li> </ul>
Instructions Followed	<ul style="list-style-type: none"> <li>• This product should not be near strong magnets.</li> <li>• Do not use it on the same part of the body for a long time, it may cause skin inflammation.</li> <li>• The charging cable cannot be replaced by yourself. If it is damaged, please contact the customer for confirmation before replacing it.</li> </ul>

### Instructions For Use


Please follow the steps below as first-time usage:

1. Remove the protective film on the electrode sheet and handle on the device.
2. Wipe off the remaining glue on the electrode pad and the protective film on the handle with a wet paper towel firstly, then wipe off the water marks with a dry paper towel.
3. Power on and confirm whether each function key works normally.

#### Tips:


Please do not dispose of the product randomly after the end of its service life. The built-in battery may contain environmentally harmful ingredient. Please dispose of the product at local disposal bins and make sure that the internal power of the battery is exhausted.

### Product Specification

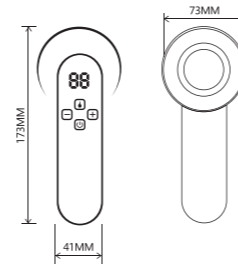



Battery Capacity	3.7V/1200mAh	Power Supply	AC220V/50Hz
Charge	5V/2A	RF Frequency	1M HZ
Bio-Wave Frequency	30-300HZ	Red Wavelength	630(±10)nm
Temperature Control	37-42℃	Material	ABS, PC, Stainless Steel
Product Weight	160g	Package Weight	230g

### Accessories Details

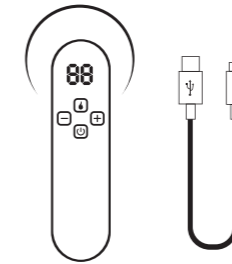

Device	USB charging cable
1	1

### Component Function

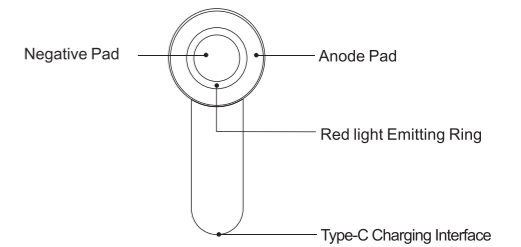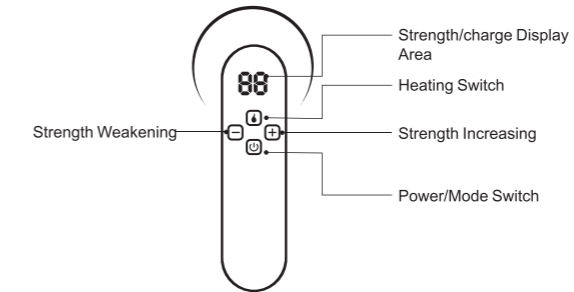

### Charge

- Please use up the remaining power of the device before the first charging, and use the device after fully charged. For later usage, you can charge as needed. But it is recommended that you discharge the battery fully per month. This can extend the battery life.
- If you haven't used it for more than 6 months, please charge it for 10-12 hours before using it.
- Insert the Type-C connector that comes with the device into the charging port of the device, then plug the USB port into an adapter that is not higher than 5V 2A, and then plug the adapter into the 220V power port to charge.
- During the charging process of the device, the charging display area will display a dynamic rectangular white light. There will display "FU" after fully charged.

### Maintenance

#### Clean the Electrode Pad

- Wipe the residual hydrogel on the electrode pad with a wet paper towel after usage, then wipe the water stains with a dry paper towel.
- Do not wipe with rough objects to avoid scratching the electrode pads.
- If the hydrogel between the red light emitting ring and the electrode sheet cannot be wiped clean, wipe it with a cotton swab and then wipe with a wet tissue. Do not immerse device in water for cleaning.
- Do not immerse device in water for cleaning.

### Clean the Fuselage

- Wipe the device with a tissue or soft cloth.
- Do not immerse it in water for cleaning.
- Do not clean with dirty wipes to avoid scratching the fuselage.

### Working Mechanism

#### 4 Types of Technology

- Bio-Wave technology (30-300 Hz)

Bio-Wave Penetration technology is like your skins weekly workout helping lift and firm while retaining skin moisture to drastically stimulate collagen and elastin production. Wave Penetrating Cell Stimulationion, Wave Penetration technology is like your skins weekly workout helping lift and firm.

- RF technology (1M Hz)

The RF waves penetrates into the dermis and converts electrical energy into heat energy, increasing the temperature inside the skin. This helps burn fat, tighten the skin, and regenerate collagen.

- Red light technology (630(±10)nm)

It can promote blood circulation and rejuvenate the skin, restore skin elasticity, and minimize fine lines, wrinkles and stretch marks.

- Constant temperature technology (37-42℃)

Constant temperature technology quickly increases the temperature of the active part and improves the effect of other technologies.

#### 3 Kinds of Working Modes

- Dispel Lines Mode

Use 200Hz Bio-Wave and 630(±10)nm red light to remove stretch marks, body fine lines and scars. The active part has a feeling of pulling outward.

- Burn Fat Mode

Use 30Hz Bio-Wave, 1M Hz RF and 630(±10)nm red light to remove stubborn fat in various parts of the body. There is a sense of hammering and heat at the action part.

- Lift Skin Mode

Use 300Hz Bio-Wave and 630(±10)nm red light for lifting and tightening loose muscles in various parts of the body. The action part has a feeling of pulling inward.

#### Intensity Levels

Each mode is designed with 10 intensity levels, suitable for different parts of the body (do not use it on the face).

#### Heating Function

The heating function can be turned on and off by the heating switch key according to needs. When the heating function is turned on, the mode is in red, when the heating function is turned off, the mode is in white.

#### Operating Hours

The working time is 15 minutes for each operation, and then the device automatically shuts down (the working mode and switching intensity will not affect the working time after the startup).

### Operation Steps

Step	Icon	Action/Notes
Hydrogel		<p>Squeeze the hydrogel from the bottle and place it in the center of the middle electrode pad, and then hold the handle to smear the hydrogel evenly.</p> <p>※ It is prohibited to use of essential oils.</p> <p>※ The hydrogel should not be smeared too much, just size of heating switch button.</p>
Boot up		Press and hold the power key until the intensity and mode display lights up, then release the key.
Select Mode		<p>Short press the power key to select the demand mode.</p> <p>※ Due to the size limitation of the display interface, the mode names displayed on the device are all abbreviations. The corresponding explanation are as follows:</p> <p>D-Lines - Dispel Lines B-Fat - Burn Fat L-Skin - Lift Skin</p>

Select Gear		<p>Short press "+" and "-" keys to select intensity needed.</p> <p>※ Because different parts have different perceptions, please select the required intensity according to your needs.</p>
Massage		<p>Hold the handle, put the electrode pads close to the massage site, and move gently and repeatedly after feeling the appropriate strength.</p> <p>※ It is not advisable to squeeze vigorously during massage, and the electrode pads should be close to the skin to cause irritation to skin.</p> <p>※ The intensity can be adjusted at any time during massage to ensure comfort.</p> <p>※ During the massage, if you feel increased friction, weakened irritation or tingling, please smear appropriate amount of hydrogel.</p>
Switch /Heating Function		<p>Press the heating switch key shortly to turn on the heating function, then the mode is displayed in red;</p> <p>Press the heating button shortly to turn off the heating function, then the mode is displayed in white.</p>
Shut Down		Press and hold the power button until all displays are off.
Clean		<p>Wipe the electrode surface where the hydrogel is smeared with a wet paper towel until the residual hydrogel is wiped clean, then wipe water stains with a dry paper towel.</p> <p>Wipe the hydrogel smeared to the body with a wet tissue or other clean item.</p>

### When in Doubt

Anomalies	Reason	Handling Method
Does not boot	Battery runs out	Charge
Can not be charged	The charging cable is broken	Replace the charging cable
	The charging port is not connected successfully	Reconnect
No effect or effect is not obvious	Unadjusted intensity	Use high gear

No effect or effect is not obvious	Excessive skin resistance	Smear hydrogel
	Device failure	Contact Customer Service
Strong irritation or tingling	Dry skin	Smear hydrogel
	Acts on skin ulcers	Replace the active part
	Use high gear	Lower gear
	Continuous use time is too long	Stop using

### Harmful Substance Declaration

Part Name	Pb	Hg	Cd	Cr6+	PBB	PBDE
Plastic Parts	○	○	○	○	○	○
Circuit Board	○	○	○	○	○	○
Battery	○	○	○	○	○	○
Charging Cable	○	○	○	○	○	○
Packaged Goods	○	○	○	○	○	○

This form is compiled in accordance with SJ/t11364.

○ : Indicates that it meets the limit requirements of GB/T 26572-2011 or less.

X : Indicates that it exceeds the limit requirement specified in GB/T 26572-2011.

## Warranty Card

Thank you for purchasing our products!

In order to use this product safely and satisfactorily, please read the instruction manual carefully before usage.

### Warranty Record

User name: \_\_\_\_\_ Product number: \_\_\_\_\_

Name/Model: \_\_\_\_\_ Purchase date: \_\_\_\_\_

Distributor: \_\_\_\_\_

Date	Repair Content	Repairman	Remarks

#### Warranty:

- ※ Within three months from the sale date of the product, if there is a performance failure, the product itself and the outer packaging must be intact and free of scratches to replace the same type of product, (subject to inspection) but does not include man-made damage.
- ※ Commodities are guaranteed for one year from the date of sale, and lifetime maintenance, accessories are not covered by the warranty.
- ※ The warranty is valid only under normal use.
- ※ All external damages, such as man-made damage, self-dismantling, unpacking, improper usage, etc., are not covered by the warranty. The warranty must provide this card. If the card is not provided or the card is altered privately, we have the right to do non-warranty treatment.