

RF RADIO FREQUENCY BEAUTY INSTRUMENT

USE INSTRUCTIONS



- Thank you very much for your purchase.
- Please read the manual carefully for normal operation before use. After reading it, please keep it for future reference.
- The use of this beauty instrument is premised on the care of healthy skin, for the purpose of beauty, and its use or sale for the purpose of treatment is prohibited.


1、Function introduction

RF radio frequency acts on fat dissolution, lymphatic detoxification, skin tightening, and enhances skin elasticity. It can reach the fat layer directly. It has the function of directional radio frequency output, which accelerates the movement of fat cells to generate heat energy. Fat and toxins are excreted from the body, so as to dissolve fat; When acting on the skin, the collagen fibers in the dermis layer can be heated to 45 degrees to 65 degrees, the collagen fibers contract, causing the loose skin wrinkles to be tightened, and at the same time, the proliferation of stimulated collagen in the dermis layer promotes the increase in the thickness and density of the dermis, So as to remove wrinkles, remove scars, restore skin elasticity and gloss to achieve the purpose of firming and wrinkle removal;

2、Instructions

1. Before use, wash away the sebum and dirt accumulated on the facial skin surface with warm water, so that the skin with smooth pores can better absorb nutrients.
2. Connect the handle to the corresponding hole, connect the power supply and start by default, press the Mode key to select the mode, then press the Timer key to adjust the working time, and then press Level + or Level- to adjust the appropriate working intensity and start working. When the countdown is over, the work will stop automatically.
3. A separate handle for the body probe.
4. The face probe and the eye probe share the same handle. Unscrew and replace the probe.

3、Part name

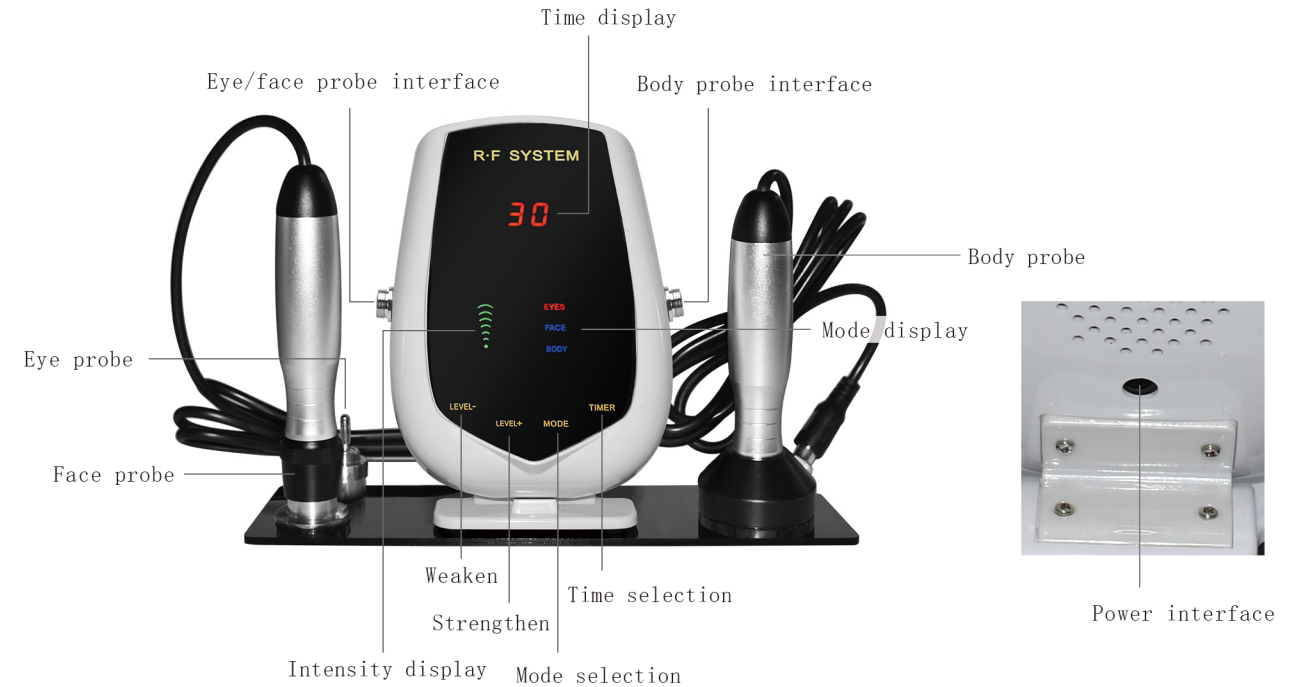

4、Instrument parameters

Power input voltage: 100V-240V	Screen: Led digital display
Power output voltage: 24V	Gears: 7 gears
Machine power: 50W	3 modes: body/face/eyes
RF output frequency: 5M	Time: 5、10、15、20、25、30 minutes
Number of handles: 2	

5、Special attention

- The product can only be used under rated voltage and frequency;
- Please do not cut off the power supply with wet hands and forcibly pulling the power cord;
- Please do not move the product by pulling the power cord;
- Use this product only in the manner described;
- If you are not using the product, please turn off the power;
- When the product is damaged, especially when the power cord and casing are damaged, please do not use or stop using this product;
- Please do not immerse the product in the liquid;
- If the power cord is damaged, it must be replaced by the manufacturer or agent or similar professional person to avoid injury;
- The surface of the vent of the product should not be blocked with newspapers, tablecloths, curtains and other items;
- This product has been used in normal operation and has not been dismantled since the date of purchase. It has been confirmed by the warranty personnel that it is a product quality problem and enjoys free warranty service within one year.