
EN
PT
ES
LV
RO
EL
ET
FR
HE
DE
RU

POWER BACK

The ultimate at-home back massage

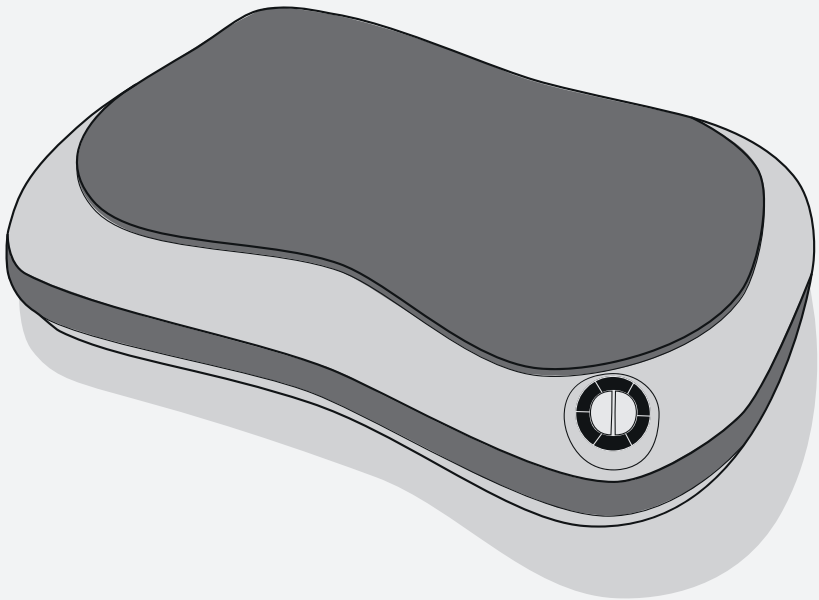POWER BACK


The ultimate at-home back massage

ENGLISH	2
PORTUGUESE	7
ESPAÑOL	12
LATVIAN	17
ROMANIAN	22
GREEK	27
ESTONIAN	32
FRENCH	37
HEBREW	42
GERMAN	47
RUSSIAN	52

POWER BACK MASSAGER

Innovative design with four heating massage heads can produce kneading massage effect. It can deeply massage fatigue muscle tissue. Accelerate circulation of blood. Ease the pain, relax mood, and make you more comfortable to fully enjoy the fun of massage.

DESCRIPTION OF POWER BACK

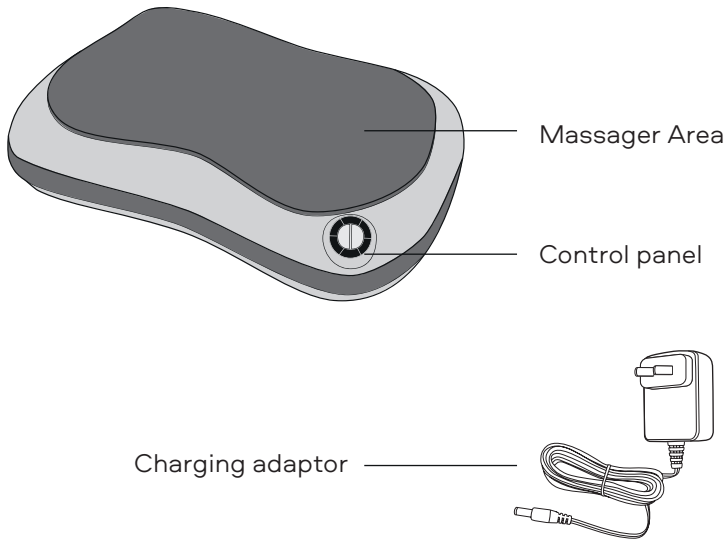SAFETY PRECAUTIONS

1. Before using, please consult a doctor if:

- you suffer from heart diseases or high blood pressure
- you suffer from infectious diseases
- you suffer from a cancer
- you have a diagnosis of schizophrenia
- you have an open wound, cut, redness, infections, burn or bruise on the skin
- you are pregnant, in the period of menstruation and incessant fever
- you have a sensitive skin, angiectasia or were hospitalized
- you had a neck surgery
- your joints and muscles can be described as slow to warm up

2. Do not use the Massager when driving a car.

3. Do not disassemble the device yourself. If necessary, contact the service center.

4. Do not use, charge or store the device in hot or humid environment or in direct sunlight in order to avoid failure and malfunction.

5. Do not allow children and disabled people to use the device unattended.

6. Do not leave the massager under water or in rooms that provoke corrosion.

7. In case of malfunction stop using the device immediately.

8. Do not drop or hit the device.

9. Do not put heavy objects on the massager.

10. The massager is not a replacement for treatment or doctor's recommendations.

11. The massager shall not be used for diagnosis, treatment or for commercial purposes.

12. Do not use the massager during charging.

13. Do not use the Massager at the same time as other gadgets with a heating function. heat function, such as an electric blanket.

14. The Massager is not a substitute for medical treatment or medical advice. In case of questions, consult your physician.

15. The massager may not be used for diagnosis, treatment or other commercial purposes. other commercial purposes.

MASSAGER OPERATING INSTRUCTIONS

- Check the availability of the device and components.
- Fully charge the massager before use.
- Do not start the device if the power and power supply voltage are out of range.

HOW TO CHARGE THE DEVICE

- Lithium battery inside the massager, full charge before first using.
- The indicator on the adaptor will be red when charging, and it will turn to be green when full charged. Unplug the adaptor to use the massager.
- In daily use, please charge it when the massager can't turn on or the massage strength turns weak.

OPERATING PROCEDURES

Rating massage time: 15mins, automatically power off after 15mins use.

1. Press the On/Off button  again if want to continue.
2. Press the On/Off button  the first time, start massage with heating function.
3. Press the On/Off button  the second time, reverse rotation.
4. Press the On/Off button  the third time, turn off the heating function but keep massage.
5. Press the On/Off button  the forth time, turn off the device, stop massage.
6. Stay at the strong strength when turn on the massager, press the button  can turn to weaker strength.

WARRANTY

- The mechanical components of the product are guaranteed for 12 months from the date of purchase. from the date of purchase.
- The product cover is not covered by the warranty.
- The warranty is valid if the product is faulty due to the manufacturer's fault, of the manufacturer.
- The warranty does not apply in case of mechanical damage and if there are signs of repair by unauthorised persons.

TROUBLESHOOTING

If there is any problem when you use the device, please follow the instructions to check by yourself first. If still can't solve the problem, please cut off the power and contact the local service center for help. Do not disassemble or fix the device by yourselves.

Symptom	Possible Cause	Solution
The product is not working	Please check if you have press the power on switch button	Turn on the power switch
	Check if the battery has been run out	Charging using the original adaptor
	Check if time to the default setting time (15 mins)	Restart the power switch
	Weighing too much on the device	Unplug the device and remove the overload things
Device over heating	Long time use	Stop using the device for 2 hours
Power cord damaged, smoke comes out from the device		Stop using and contact with the local service center

CLEAN AND CARE

AFTER USAGE:

- Turn off the massager and disconnect the charger.
- Store in a dry cool place. DO NOT store the device up in the air.

CLEAN:

- Wipe the massager with a soft dry cloth. Do not use water to clean the device.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, furniture polishes or their solutions.

TECHNICAL SPECIFICATIONS OF THE PRODUCT

Name: S8 BACK MASSAGER POWER BACK

Model No: E71-MP005-S8

Power supply: AC 100-240V/DC12/6V

Power consumption: 18 W

Lithium battery: 2000 mAh

Timer: 15 minutes

Dimension: 330 x 220 x 110 mm

Weight net: 1,5 Kg

MASSAGEADOR PARA AS COSTAS POWER BACK

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