
EN
PT
ES
LV
RO
EL
ET
FR
HE
DE
RU

POWER MOOD

*A rejuvenating head massage from
the comfort of your home*

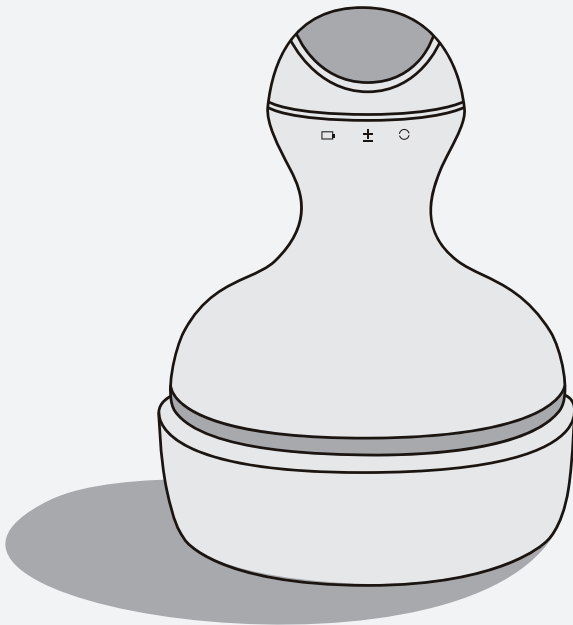
POWER MOOD

*A rejuvenating head massage from
the comfort of your home*

ENGLISH	2
PORTUGUESE	4
ESPAÑOL	6
LATVIAN	8
ROMANIAN	10
GREEK	12
ESTONIAN	14
FRENCH	16
HEBREW	18
GERMAN	20
RUSSIAN	22

Product Structure

Content: Machine, Charge base, USB Cable, User Manual

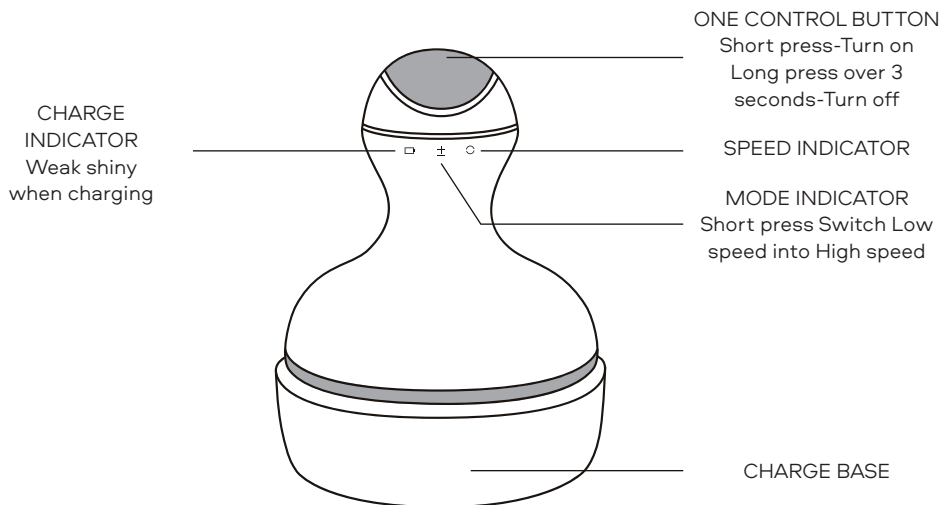Recharge detail

1. When recharging, green light indicator will be shining.
2. When full charge, the indicate light will tirn green.
3. Rechargeable time about 2 hours.
4. Machine will be stoped when charging.

Product Feature

1. Built-in Lithium battery.
2. One control button.
3. 10 minutes massager time.
4. Automatically protect function.
5. Two speed.
6. Two switches of scartching.
7. Micro-USB cable and power adapter.
8. DC 5V 1000mAh.
9. Waterproof.
10. 3D head massage.

NOTES AND MAINTENANCE

NOTES

1. It is waterproof device, you can use cleaning fluid and water to clean it (do not sink it into water for a very long time).
2. Please use dampcloth or neutral detergent to clean the product and wipe dry afterwards.
3. Do not use sharp brush, sharp tool, gasoline and diesel to clean it.

Maintenance

1. Please keep it into a box and keep in a safe, dry and cool place.
2. Do not keep it in any place with high temperature, near fire or direct sunlight, or it may faded or damaged.
3. Do not contact benzene, dilient or any other corrosive liquid.

TRANSPORTATION

This Device has simple shock proof protection, can be shipped by air, railway, road and by seaway transport.

Read the following safety instruction carefully before using this massager. If you suffer from any of the medical conditions listed below, consult with a doctor before using this relaxation massage device.

IF YOU:

- suffer from a heart condition
- have a sensory disorder or pacemaker
- suffer from osteoporosis or diabetes
- have a malignant tumor
- have neurological disorder or perceptual disorder on head
- have fractured or head surgery of your massaging area
- head skin disease and allergic dermatitis
- suffering from acute diseases or are undergoing medical treatment

We advise you to consult with your doctor before using this massager.

Name: S8 Neck Massager Power Mood
Model No: E28-SL002
Power supply: 100-240V ~ 50/60Hz/3,7V
Power consumption: 5W
Lithium battery: 1000 mAh
Adapter output: DC 5V 1.0A
Timer: 15 minutes
Dimension: 140 x 110 x 110 mm
Weight net: 0.66 Kg

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