




# LUXE K

## Hair Removal Device USER MANUAL

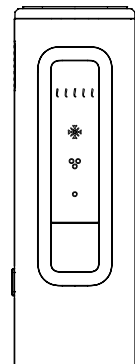


Please read this instruction manual carefully before using the product

### Product introduction:


The instrument adopts photonic beauty IPL technology, which directly irradiates skin with intense pulse light to produce photochemistry.

### Functions overview:

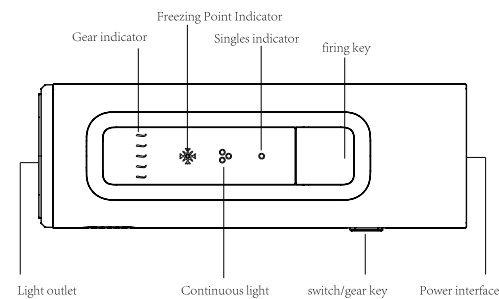


1. Long press the power button: turn the instrument on and off.
2. Short press the power button: switch gears.
3. Short press the firing key: stimulate the light.
4. Long press the firing button when not touching the skin: turn on the automatic firing mode.
5. Long press the firing button when touching the skin: equivalent to continuous firing.
6. Simultaneously press the on-off key and the firing key: turn the ice cold function on and off.

### Product Overview:

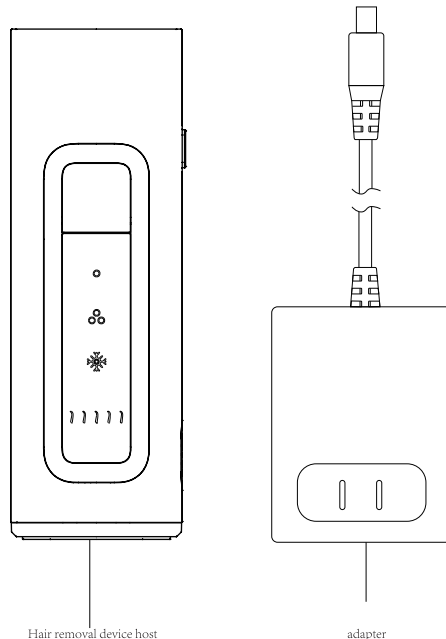


### Instructions:

1. Plug the adapter power cord into the power supply

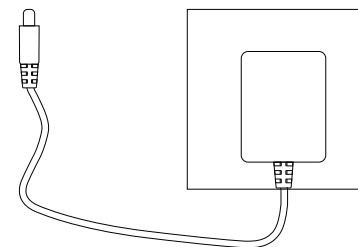


2. Insert the output end of the adapter into the bottom jack of the epilator.

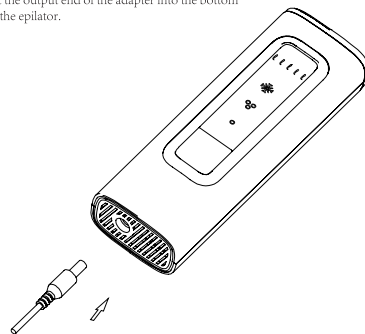

3. Standard method:


Press and hold the "On/Off button" for 3 seconds, the instrument will be turned on, and one indicator light will light up. At this time, the instrument is in the first gear state Weak, 5th gear is the strongest, one gear corresponds to one indicator light. At this time, the instrument is in the single-click mode, the single-click mode light is on, the light head and the skin are fitted at 90° and the hair removal can be started by clicking.

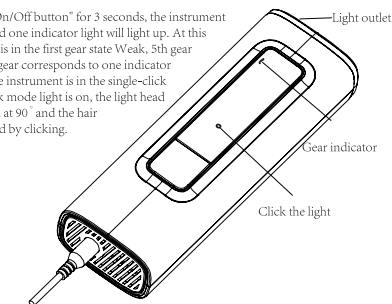

4. How to turn on the cold function: Click the power button and the firing button at the same time to turn on the cold function, and the cold indicator light is on. If the cold function is not needed, press the power button and the firing button at the same time to turn off.

5. Turn on the continuous light output mode:



When not touching the skin, long press the firing button to turn on the automatic light-emitting mode. Without touching the skin, press and hold the firing button to turn off the automatic light-emitting mode. At this time, the continuous light-emitting indicator is off, and the single-click light is on.

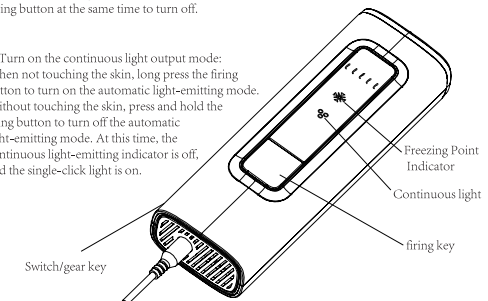

Note: If the instrument cannot be turned on, please check whether the power cord is plugged in properly.

### Precautions for use:

Do not block the air outlet during use

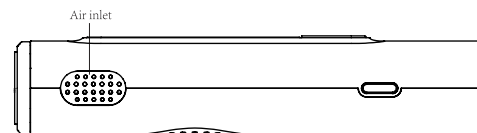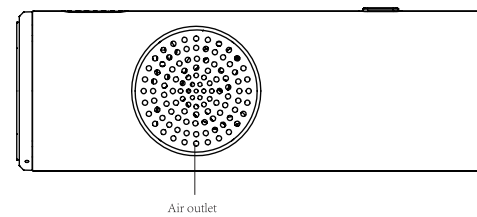

How to use ice compress function: After hair removal, you can use the ice compress function alone to cool down. Adjust the instrument to single-click mode, and then turn on the ice application function to apply ice back and forth against the skin. When the temperature of the lamp head decreases, the ice application effect will be better. (Do not press the firing key at this time).

In order to prolong the service life, please avoid using this product in a humid environment. The light needs to be close to the skin before it can work.

### Skin test:

Before using first or after tanning, perform a skin test on each treatment area. Skin testing is a necessary step to check the skin's response to the treatment and to determine the appropriate light intensity for a particular area.

1. Select an area close to the treatment area.
2. Turn on the power of the instrument, and select the 1st gear by default.
3. Place the device against the skin and press the fire button to emit a single flash.
4. After finishing the skin test, please wait 24 hours to check for any skin defects.

### Suitable for skin tone:

Scope of application: lip hair, armpit hair, limbs hair, frontal hairline and other effects of beautiful hair. Painless hair removal is suitable for black hair and dark hair, not for white, gray, golden hair, please refer to the comparison card below.

Skin Tone Hair Comparison Card

SKIN COLOR \ HAIR COLOR	WHITE	BEIGE	LIGHT BROWN	BROWN	DARK BROWN	BLACK
WHITE/GRAY	✗	✗	✗	✗	✗	✗
RED	✗	✗	✗	✗	✗	✗
LIGHT BLONDE	✗	✗	✗	✗	✗	✗
DARK BLONDE	✓	✓	✓	✓	✗	✗
BROWN	✓	✓	✓	✓	✗	✗
DARK BROWN	✓	✓	✓	✓	✗	✗
BLACK	✓	✓	✓	✓	✗	✗

### Precautions:

1. People with skin perception disorder cannot use it.
2. People with allergic dermatitis, eczema and other skin diseases cannot use it.
3. Pregnant women, women who are menstruating and breastfeeding cannot use it.
4. Those who have been exposed to sunlight in the past 2 weeks cannot use it.

### Instructions for use:

1. Shave the hair.
2. Clean your skin and make sure it is completely dry and free of any oily substances. Before using the photon epilator, it is forbidden to smear products such as coloring agent and aloe gel. Otherwise there will be no light.

### Product parameters:

Product name: IPL photorejuvenation hair removal device

Output: DC 12V/3A

Power: 36W

Number of shots: unlimited

Lighting time: 1-3 seconds

Energy gear: 5 gears

Irradiation area: 3.3cm

Product weight: about 200g

Adapter specification: input AC100-240V

Working temperature: 5 C ~ 30 C

Light source technology: pulsed light

### Steps for usage:

1. Before use, use cotton swabs or clean cloth dipped in medical alcohol to clean the light outlet and filter. Keep the light outlet dry before use, otherwise the light outlet will be affected.
2. Plug in the power plug, short press the switch button for 3 seconds to turn on the hair removal mode, the default is 1 gear and in manual mode, the cooling fan rotates.
3. The power switch key is also the gear key, double-click the "gear key" to select the appropriate gear, the display energy gear, the first gear is the smallest, the fifth gear is the largest, and the gears can be adjusted cyclically.
4. After the tap position is adjusted, the light outlet of the device is completely attached to the skin part for hair removal. At this time, press the light button to light out when the display screen flashes.
5. Press the light out button for 3 seconds to enter or exit the manual light out mode, and the screen will light up automatically.

### Open the automatic light output mode:

In step 5 of the cycle, each skin area that needs to be shed is flashed one time and moved slowly. The shed area is illuminated three times at a time, but not three times in a row in the same area. And so on to take off the whole hair removal area, and then repeat 2 times.

6. After using the hair removal device, selectively apply aloe vera gel (or not). After using the device, short press the "power switch button" for 3 seconds to shut down, and cut off the power in time. After use, the filter of the lamp head may have small black spots attached, just wipe it with alcohol.

### Treatment cycle reference table:

Time	Frequency	Effect
1-4 weeks	Use 2-3 times a week	Over 75% of the hair begins to fall out naturally
5-6 weeks	2 times a week	New hair loss and inhibition of growth
7-8 weeks	1 time a week	Hair becomes thinner and less, no longer bouncing
8+ weeks	1-2 times a week	Use selectively according to hair growth, When the hair starts to grow, (use it again until it stops growing.)

A depilation course of treatment for irradiation 3 times, the first irradiation of a depilation of the part, and then do the second and third times, not in the same place continuous light 3 times. The higher the gear, the higher the energy and the better the effect. In the early stage, it is recommended to use 1-2 gears first. After the skin adapts to the later stage, you can choose to use 3-5 gears. Depending on your personal hair growth rate, the penetration of each part will be different.

### Completion of use:

1. The skin may be slightly red, swollen or tingly or warm. These reactions are completely harmless and quick disappear.
2. After using the instrument for 2 hours, you can safely apply lotion, cream deodorant, moisturizing product or cosmetics on the treatment site: If the skin is prickly or red after the treatment, please wait for these conditions to disappear before applying any product on the skin. If the skin feels prickly after applying the product, please wash it off with clean water.

### Product maintenance:

1. Please turn off the power supply of the device immediately after use, remove the plug, and let the instrument cool down.
2. Wet a soft cleaning cloth with a little water, and then use it to clean the following parts: filter/skin contact sheet.
3. Store the instrument in a dry place with a temperature between 0 C ~ 60 C and no dust.

### Troubleshooting:

Problem Description	There may be a reason	Solution
The device/power adapter becomes warm during use	that is normal	No action required
I put the epilator on my skin and the (flash ready) light didn't come on, but the cooling fan started running.	Device needs to be reset	When the overheat protection is activated, the fan still keeps running, please do not turn off the power of the device, please wait for it to cool down for 15 minutes before continuing to use it.

Unpleasant smell when using the device	Accessories or filters are too dirty You have not completely removed the hair in the treatment area, which may burn and produce an odor.	Please clean the accessories carefully. Please remove the hair on the skin surface before disusing the device.
During the treatment, the skin felt more sensitive and the device felt unbearable pain.	The light intensity setting you are using is too high.	Please check it. Make sure you select the appropriate lighting intensity setting and lower it if necessary
	You have not completely removed the hair from the treatment area.	before using this instrument. Please remove the hair on the skin surface first.
	The UV filter at the light outlet is damaged.	Once the filter of the effective outer line is damaged, Do not continue to use this instrument.
Not satisfied with the results of the treatment	The light intensity setting you are using is too low.	Please select a higher gear setting.
	You did not shine a second flash around the previously treated area	Fire the flash as close as possible to the previous part.
	This device has no effect on your hair color.	If your hair is light blonde, gold or red. These treatments won't work. Remove all hair successfully. To discuss how you follow the recommended course of treatment, you can program shorter treatment intervals. But it should not be treated more than once a day.
	Insufficient frequency of instrument use may slow your response to IPL therapy.	Continue using the device for at least 6 months. Because during this time the hair regrows but decreases.

### Operation Notes:

1. Do not look directly at the flash or flash directly to the eyes during operation.
2. Mild redness of skin after operation is a normal phenomenon, which can be applied with ice or cooled by itself.
3. When cleaning the dirt on the glass plate inside the outlet, do not spray alcohol directly on the outlet.
4. During the operation of the instrument, it is not allowed to smear any gelatinous cosmetics such as gel on the skin, so as not to damage the instrument.
5. Try not to eat light-sensitive food (such as spinach, rape, mustard, etc.) during operation.
6. The adjustment of energy gear is from weak to strong. Under the condition of not damaging the skin, the higher the energy, the better the effect.
7. Darker skin may experience delayed skin reactions and cannot gain energy because the immediate response is not obvious.
8. Pay attention to skin moisturizing and sunscreen after flash.
9. Avoid ultraviolet radiation during use, do a good job of shading protection outside, try to avoid computer radiation, apply isolation cream, sunscreen protective layer.
10. During use, please do not block the air inlet and heat dissipation hole, so as not to affect heat dissipation.
11. During the use of the instrument, if there are hairs or stains sticking to the plastic edge of the light outlet, the plastic edge of the light outlet will be black or produce black spots, but it does not affect the performance and normal use of the instrument.
12. The "light tube" inside the light outlet of the instrument, in the later period of the specified life, the inner wall of the tube will show different degrees of blackening phenomenon, which belongs to the normal aging process of the tube during the life of the energy attenuation rate will not exceed the specified range, can be used normally.
13. The product does not have waterproof function, please do not use it in the bathroom, in the shower and in other places that may be easily exposed to water.
14. If the shell or window of the product is damaged or the instrument is soaked carelessly, please send it back to the manufacturer for maintenance to determine whether it can be repaired before use. If it is soaked with electricity, please unplug the power supply to the contactor first to ensure your safety.
15. Avoid storage in places with strong magnetic field and serious humidity and dust or places with high temperature.
16. Do not strongly impact or hit the instrument.
17. Use of this product for purposes other than designated functions is prohibited.