lightinderm

REGENERATIVE SKIN CARE.

ASSISTANCE & SÉCURITÉ

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SAFETY AND SUPPORT

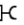
GETTING STARTED

Thank you for choosing Lightinderm to take care of your skin! Before you start using your care programme, please take a few moments to read the information below.

WARNINGS

The Lightinderm system is a care programme designed for home use. It is not, under any circumstances, designed to treat health problems.

- Choose the serum capsule adapted to your needs.
- Do not use this device if your skin is photosensitive or if you are taking medication that may cause skin hypersensitivity to light.
- Do not use this program on extremely or chronically ill, injured or irritated skin (for example, on swollen or red skin, on skin rashes caused by allergies, burns, bruises, swelling or open or healing wounds).
- If you experience an adverse reaction, please consult a doctor and the Lightinderm customer service.
- Do not use this program if you are likely to be pregnant or if you are pregnant or breastfeeding, as the risks involved are currently unknown.
- Do not use this product near bathtubs, showers, washbasins or other containers of water.
- When the device is used in a bathroom, unplug it after use, as a nearby source of water is a hazard, even if it is turned off.

- The Lightinderm device must not be used, cleaned or maintained by children.
- Children, even under supervision, should not play with the device.
- This appliance may be used by children over 8 years of age and by persons with reduced motor, sensory and mental capabilities or persons who have never used the appliance before, provided that they are supervised or follow the instructions for the safe use of the device and are aware of the potential risks.
- Assessments conducted by a certified independent service have shown that the Lightinderm device is safe for the skin and eyes. However, it is recommended not to point the device directly at the eyes.
- For additional protection, the installation of a residual-current device (RCD) rated at less than 30 mA is recommended on the bathroom electric circuit. Ask your installer for advice.
- To recharge the battery, the part number of the removable power supply unit is given below.
- Only the charger supplied by Lightinderm can be used to charge the device. Charger Part Number: W&T-AD05 W050 100 
- This device is fitted with non-replaceable batteries.

SAFETY PRECAUTIONS

- Do not allow water or liquid to enter any part of the device, as this may cause it to malfunction.
- The Lightinderm device only runs using Lightinderm serum capsules.
- Do not use the Lightinderm device if it is visibly damaged and never attempt to open or repair the device. If you have a problem with the device, please contact Lightinderm customer service.

SUPPORT

If you can't find the solution to your problem in the answers below, our customer service is available:

contact@lightiderm.com
Paris Santé Cochin - Lightiderm
27 rue du faubourg Saint-Jacques - 75014 Paris
Monday to Friday from 10am to 5pm:
+33 9 63 50 03 27

- The device does not turn on when you press the on button for more than 2 seconds:
 - Charge the device for 20 minutes and then try to turn it on again.
- The device turns on and indicates an error by beeping 3 times in a row and turning on the battery indicator 3 times:
 - Charge the device for 20 minutes and then try to turn it on again.
- The device turns on and indicates an error by beeping 3 times in a row and by turning on all the indicator lights and remaining dose indicators 3 times:
 - There is no serum capsule in the device: insert one to start your care session.
- The device turns on and indicates an error by beeping 3 times in a row and by switching on all remaining dose indicators 3 times:
 - The serum capsule is empty, it must be changed.
- You are unable to insert a serum capsule into the device:
 - Do not force, press the eject button on the back of the device and then try to reinsert the capsule. If you still cannot insert it into the device, please contact Lightiderm customer service.
- You cannot remove the serum capsule from the device:
 - Press the eject button firmly and then gently pull the capsule out.
- We recommend not having more than 2 consecutive sessions.

GUARANTEE

- The Lightiderm device is guaranteed for 2 years from the date of purchase. If an anomaly due to a material or manufacturing defect occurs during the guarantee period, Lightiderm will repair or replace the products at its own expense.
- Lightiderm will only pay for the repair or replacement of products if proof of purchase covering the guarantee period is provided.
- The guarantee is not valid if the device is defective due to misuse or improper maintenance, or if it has been handled or repaired by a person not authorised by Lightiderm.
- To use the Lightiderm device properly, it is important to strictly follow the instructions in the user guide.

STORAGE AND MAINTENANCE

- Maintenance: The Lightiderm device does not require any special or regular maintenance. No modifications to the product are authorised.
- Cleaning: Unplug the Lightiderm device before cleaning it. Clean the outside of the device using a dry cloth.
- Travel: If you travel with the Lightiderm device, slip it into the supplied carrying case.
- Storage: It is best to place the Lightiderm device on its storage base when not in use. Do not store near a heat source or on a hot surface.
- Charging: The device has an autonomy of 7 sessions after a full charge. When the battery is completely discharged, 20 minutes of charging is enough for a new session.

TECHNICAL SPECIFICATIONS

MATERIALS: Aluminium / Silicone / ABS-PC / PMMA
body safe

COLOURS: Pale pink / White

DIMENSIONS: 163 x 60 x 34 mm (without cartridge)

WEIGHT: 120 g (without cartridge)

BATTERY: 800 mAh

PORTS: USB-B USB 2.0

AUTONOMY: A minimum of 7 sessions per full charge

INTERFACE: 1 on button and 1 eject button

CHARGER PART NUMBER:

W&T-AD05 W050 100 

This symbol indicates the conformity of the device with European Union Directive 2014/35/EU.

This symbol means that, in accordance with European laws and regulations, you must dispose of your product and/or its battery separately from household waste. When this product reaches the end of its life, take it to a collection point indicated by local authorities.

Please refer to the documentation supplied with the device.