



It's different when we make it.

Deep dive research into the therapy based on the light source and electricity for years.

We are a research-oriented company that studies skin technology.

Our technologies originated from our own researches,

while our beauty devices are created by

the efforts of several Masters and Doctors.

We proudly offer you unparalleled skincare solutions.



The scientific basis of all stages for raw material and technology development.




High level of technology based on original technology and patents.

GALVANIA LED BOC

Triple Action Homecare

LIGHT THERAPY

RED(630nm) + NIR(850nm) penetrates deep into the skin!

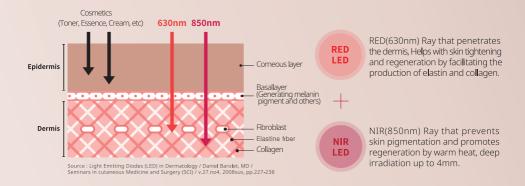


ELECTRONIC THERAPY

ALTERNATING CURRENT (FARADIC)

Essence Keratin I STIMULATION Epidermis Dermis

If you take care of your skin with a faradic current, not only will your muscles become softer, but you will also increase your cell regeneration and metabolism in addition to promoting blood circulation. It is effective in improving skin elasticity, removing wrinkles, and discharging facial wastes.

DIRECT CURRENT (GALVANIC)

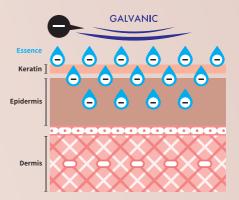

Galvanic current uses the property of repulsing the same polarity and attracting the opposite polarity when molecules are ionized and become electrically charged. It plays a role in penetrating active ingredients into the skin and discharging wastes.

Quadro-effect LED Mask U±Well



- 4 in 1, Quadruple effect LED mask (LED+IR+Microcurrent+Galvanic Together)
- Focusing on the nasolabial zone, where 80% of skin trouble such as wrinkles, pore is concentrated.
- Ergonomic design and highly condensed LED chip for best result (RED + IR 163)
- Lightweight and charging cradle, detachable battery for easy carry.
- Mobile Bluetooth control by dedicated app
- Clinically tested





DUAL EFFECT

Light and electricity two functions at once



ULTRALIGHT & COMPACT

Lightweight and slim design



DESIGNED FOR SAFETY

By smart IC circuitry design



CORDLESS DESIGN

Perfectly hands-free



ENERGY SAVING & RECHARGEABLE

Low power design allows for minimal power consumption for long hours of use and rechargeable function



FOCUS

Intensive care for small areas vulnerable to facial wrinkles



BLUETOOTH

Smart control by your mobile phone Select skin type



WATERPROOF

Can be washed lightly with water for cleanliness after use



CONTINUOUS USE EXTRA BATTERY

Replaceable 2 battery packs



Ear Loop

The mask is tightened to the face by an ear loop.
It helps staying effective and more easy to move while wearing it.



- It is recommended to use together with the Ecowell's Silver Mask Pack for best result.
- It can be used together with basic cosmetics such as essence or serum. If you use it without
- applying a mask pack or cosmetics, the effect will be reduced.





Indicator color	Microcurrent and LED strength is adjustable by each touch	Management mode (LED+Microcurrent+Galvanic)	Alarm sound	Full operation time
Green	Short Touch	Microcurrent Level 1~3	Beep sound once	Total 20min (LED 10 min /
	Long Touch	LED Level 1~3	Beep sound twice	microcurrent + galvanic 20min)


Articles	Specifications
6:	Main Device: 133mm x 143mm x 50mm
Dimension	Cradle: 140mm x 140mm x 130mm
	Main Device :
\\(\lambda \)-:- -+	38g (Before battery Install),
Weight	44g (After battery install)
	Cradle: 192g
Light source	630nm(RED) 104ea + 850nm(NIR) 38ea (Total 142ea)
Operation time	LED+Microcurrent+Galvanic/20min
Recommended operating temperature	0°C - 40°C
Recommended operating humidity	10% - 90%
Number of consecutive uses	1 time
Full charge time	Max 1 hour
Battery Specifications	3.7V / 110mAh

LED MASK Clincal test


(KC Skin Research Center)




- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask

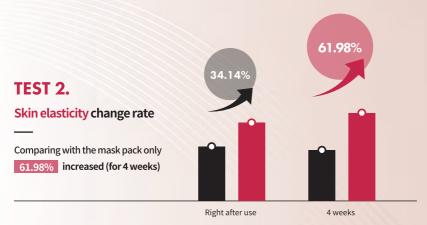

- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask

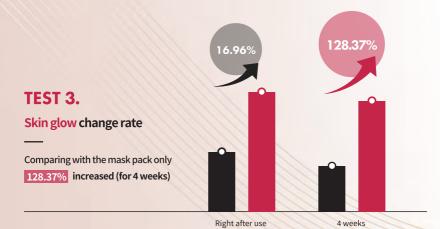

- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask

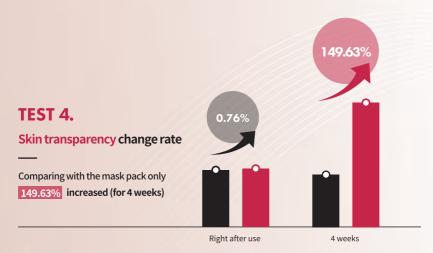

- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask

MICROCURRENT MASK PACK SKINZAP





Double effects by microcurrent treatment and excellent essence



Focusing on skin elasticity by generating collagen deeply in your face



Incredible effect by ingredients of essence (Adenosine, Niacinamide, 3 patented ingredients)



Easy to use: All you have to do is attach the device on your forehead

SILVER MASK



BLACK MASK (Graphene & CNT)



Silver mask

- Ecowell's best-selling mask pack
- Essence composition with whitening and wrinkle improvement effects / Use of vitamin C coated sheets
- True microcurrent effect with silver-containing silver pattern

Black mask

- Ecowell's new mask pack
- Essence composition focused on soothing and moisturizing
- Excellent microcurrent effect with graphenecontaining black pattern

Pick Points

Microcurrent + galvanic simultaneous action

-One device simultaneously implements microcurrent, which helps skin elasticity, and galvanic function, which helps absorb active ingredients in cosmetics, helping with efficient skin care.

Real Microcurrent Treatment

- -Have seen others without device? Then, it can't make microcurrent.
- -Increase of collagen thickness / Increase of protein synthesis / Elevation of elastin fiber

True Home Aesthetic

-Save your time from spa and skincare shop.

Differentiated Essence on Mask Sheet

- -It fills your skin with elasticity from the inside out by double function essence.
- -Niacinamide, Adenosine, etc together with 3 patented ingredients. (Silver Mask)

One more thing you will get from this maskpack

-Make up goes so well after using this.

Automatic Power ON / OFF

- -Just leave the device on your forehead, and it will last 20 minutes.
- -Skin care is not a hassle anymore.

Microcurrent Maskpack Clinical Test

(P&K Skin Research Center, Korea)

TEST 1.

Increased skin elasticity of the cheeks

Comparing with the mask pack only

107.8% increased (for 4 weeks)



TEST 2.

Increased skin elasticity of the forehead

Comparing with the mask pack only

62.3% increased (for 4 weeks)



TEST 3.

Increased skin elasticity of the eye region

Comparing with the mask pack only

101.3% increased (for 4 weeks)



TEST 4.


Increased skin density

Comparing with the mask pack only

53.2% increased (for 4 weeks)





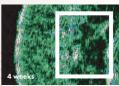

MC:WELL Mask pack + Device test group, Female of 40's with dry skin type

TEST 5.

Enhanced skin tone (Brightening)

Comparing with the mask pack only

171.7% increased (for 4 weeks)







MC:WELL Mask pack + Device test group, Female of 40's with normal skin type

03

PERFECT COMBINATION OF BEAUTY DEVICE AND CREAM Intensive care around eyes and lips

WRINKLELESS PERFECTOR & GALVANIC WRINKLELESS CREAM



4 in 1 Home Aesthetic Device



Eyebags and Facial Line Killer



Easy to use : Automatic power ON / OFF



Just grab the device and roll it around your face and neck

Pick Points

Galvanic Treatment

Boosting the penetration of active ingredients (i.e., Niacinamide, Vitamin-C, hyaluronic acid)

- RED LED (630nm) Light Therapy
 Irradiating (630-660 nm) to aid collagen production
- NIR (Near Infrared / 850nm) Light Therapy
- Micro Vibrations

Lifting effect with approximately 12,000 micro vibrations per minute

Diamond Kneading

Optimizing skin massage by 42-face

Articles Specifications

Battery LR44 / 2EA

Continuity of Applicator Usage (Actual Using)

Motor Cylindric Vibrator

Function Vibration, Red LED, Galvanic, NIR

80 min (120 min)

ON / OFF Biometric Switch



Compact Size 3 in 1 GALVANIC LED BOOSTER



Renew your skin with 3-minute massage every day



Handy and small enough to carry even when you're out of home



Easy to use: All you have to do is attach the device to your forehead

Pick Points

Galvanic Therapy

-It's so called, iontophoresis, which boosts penetration of cosmetic ingredients compared with rubbing your hands.

-It is effective when used with essence containing niacinamide, vitamin C, and hyaluronic acid, etc.

Red LED Treatment (630~600nm)

-With the close LED irradiation, it helps collagen generation inside of your skin.

• NIR (Near Infrared / 850nm) Light Therapy

-It promotes regeneration of skin by warm heat.

• Automatic Power ON / OFF

- -The device is automatically turned on/off using biometric current.
- -Skin care is not a hassle anymore.

Articles	Specifications
Battery	CR2032 / 2EA
Continuity of Applicator Usage (Actual Using)	120 min (180 min)
Motor	Coin Type Vibrator
Function	Vibration, Red LED, Galvanic
ON / OFF	Biometric Switch



Compact size 5-in-1 elasticity improver ELB (EMS:Low Frequency Therapy +microcurrent+LED Booster)



- 5-in-1 smart beauty device with EMS+microcurrent+vibration+LED (RED+IR) care
- Enhanced elasticity care by EMS and microcurrent
- Red + IR LED provides deep intensive care.
- Compact size achieved by adopting smart micro technology (62mm x 45mm /18g)
- · Easy to carry and use
- Various designs and colors (OEM)

Articles	Specifications
Battery	CR2032 / 2EA
Continuity of Applicator Usage (Actual Using)	100 min (160 min)
Motor	Coin Type Vibrator
Function	Vibration, Red & RGB LED, EMS
ON / OFF	Biometric Switch



Compact Silicone Vibrating Brush bubblux



Vibrating cleansing device, excellent in removing dirt and wastes in pores without irritation



Easy to use by lightweight, waterproof (IPX8)

- Exceptional Convenience
 - -Slide your fingers into Flexible Finger Band (No ON / OFF switch)
- Microvibration
 - -It delivers MORE cleansing power and LESS skin irritation approximately 12,000 RPM
- Ease of Using
 - -Hand-to-hand convenience in compact design
- No Hassle of Recharging
 - -Built-in battery case (replaceable CR2032 battery)

Articles	Specifications
Battery	CR2032 / 1EA
Continuity of Applicator Usage (Actual Using)	50 min (75 min)
Motor	Bar Type Vibrator
Function	Vibration
ON / OFF	Magnetic Switch
Waterproof level	IPX8




HIFU multi-care device O'WELL

- Multi-beauty device with built-in HIFU (focused ultrasound function)
- Radio Frequency (RF)
- Low frequency (EMS)
- Microcurrent
- Galvanic iontophoresis
- Hot SPA
- LED (R,G,Y)
- IR infrared
- Rotatable head
- Beautiful charging stand

Specifications Articles Battery Rechargeable **Continuity of Applicator Usage (Actual Using)** 60 min (80 min) Function MODE1 (GALVANIC+MICROCURRENT) MODE2 (RF+EMS+NIR) MODE3 (HIFU) ON / OFF **Button Type**

Skin makeup by micro-vibration Multi-function Gua Sha Massage Device

- In addition to the effects of micro-vibration and traditional massage, galvanic, microcurrent, and thermal functions are implemented at the same time.
- The effect of Guasha massage
 - Removal of stagnant toxins in the skin of the face
 - Lifting effect of facial skin
 - Relieves skin swelling and improves blood circulation
- Easily change modes with one Start/Stop button
 - -Green: galvanic + micro-vibration mode
 - -Blue: microcurrent + micro-vibration mode
 - -Red: thermotherapy + micro-vibration mode
- Operates by holding the device without a separate switch and touching the other side of the device to your face
- Auto-off function: Auto-off after 1 minute when not in use
- Built-in USB C-Type charging port (1-hour continuous use when fully charged for 2 hours)

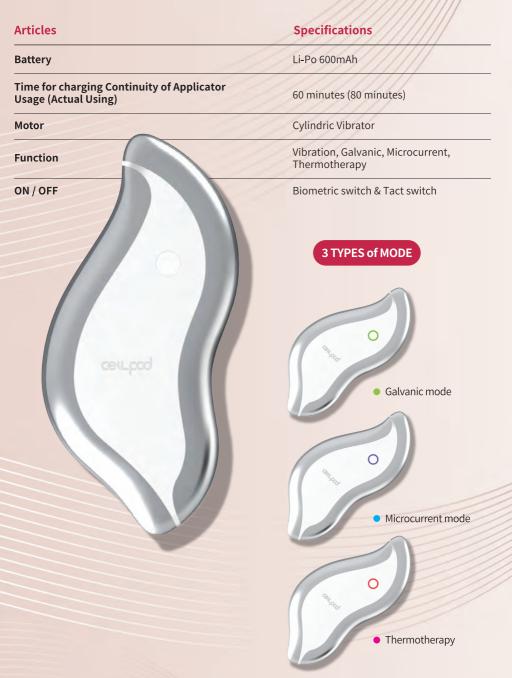





SKINZAP

与CELLPOD的相遇,注定与众不同。

为进入皮肤深处进行多年研究, CELLPOD 是研究皮肤科学的研究中心型企业。

于自主研究所中开发多种原创技术及专利, 由硕-博士组成的研究团队承诺深层皮肤研究

> 可提供现有护肤设备无法体验到的 CELLPOD 独家皮肤护理。



已认证的成分, 科学技术开发的协同组合




基于现有技术及专利的 高端技术产品

GALVANI LED BOC

Triple Action Homecare

光疗法

RED LED光(630nm)及近红外线(850nm)光线,可对皮肤深处进行照射。

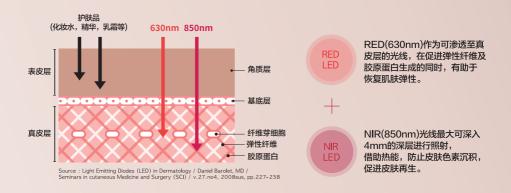


电流疗法

交流(感应电流)

機电流 角质 I 表皮层 A STIMULATION 表皮层

使用微电流对皮肤进行管理,肌肉就会变得松弛柔软, 当然不仅如此,微电流还可促进血液循环, 促进细胞再生及新陈代谢,于改善皮肤弹性,缓解皱纹 以及排出面部老废物质方面都具有明显功效。

直流(Galvanic)

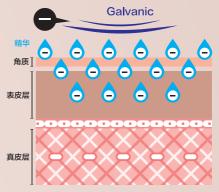

当分子离子化时,就会携带有电荷 输出相同的极性,Galvanic电流就是借用 相同极性输出这一特点,将有效成分渗透至皮肤 内部,同时排出皮肤中堆积的老废物质。

四重功效的LED面膜 LU±Well



- 四重功效的LED面膜(LED + IR + 微电流図+ Galvanic复合效果)
- 集中护理汇集皱纹,毛孔等80%皮肤问题的蝴蝶区。
- 为实现最佳效果,采用满足人体工学的设计以及高密度的LED芯片(RED + IR LED 163)
- 重量轻,使用充电式支架,分离式电池,携带方便。
- 安装专用程序,通过蓝牙进行控制。
- 完成临床试验





DUAL EFFECT

光疗及电流的 双重同步刺激效果



ULTRALIGHT & COMPACT

轻便,纤细的设计



DESIGNED FOR SAFETY

智能IC回路设计,确保安



CORDLESS DESIGN

解放双手的设计, 使用方便



ENERGY SAVING & RECHARGEABLE

耗电最低化 长时间使用及充电功能



FOCUS

集中护理皱纹脆弱部位



BLUETOOTH

根据皮肤类型,利用智能手机 进行阶段调整



WATERPROOF

使用后,可用水轻柔清洗



CONTINUOUS USE EXTRA BATTERY

提供两节电池



Ear Loop

使用耳挂将面膜紧紧贴合于面部 提升产品使用效果的同时,方便活动



- 建议与Ecowell银色面膜一起使用,以获得最佳效果。
- 可与精华等基础化妆品一起使用。
- 如果你使用它没有使用面膜或基础化妆品,效果会降低。



LED启动指示灯(绿色)




LED启动指示灯	可调整微电流及 RED LED的强度	管理模式 (LED+微电流+Galvanic)	提示	启动时间
绿色	轻触	微电流强度变更1~3	Beep声音1次	LED10分钟/电流20分钟
	长按	LED强度变更1~3	Beep声音2次	ここの10万钟/电流20万钟

项目	详细事项
	设备: 133mmx143mm x 50mm
尺寸	支架: 140mm x 140mm x 130mm
	设备:
壬巳	38g (装备电池前)
重量	44g (装备电池后)
	支架: 192g
光源	630nm(RED) 104个 + 850nm(NIR) 38个 (总计142个)
运行时间	LED+微电流+Galvanic/总计20分钟
建议运行温度	0°C − 40°C
建议运行湿度	10% -90%
连续使用次数	1次
充电完成时间	Max 1小时
电池规格	3.7 V / 110mAh

LED MASK 临床试验 (KC Skin Research Center)




- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask

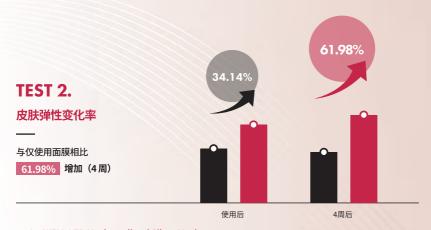

- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask



- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask

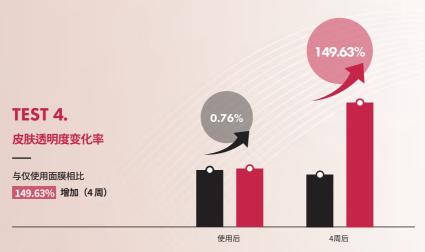

- Lu:WELL LED Mask + cellpod Silver Mask
- cellpod Silver Mask

O2 ^{微电流面膜} SkinZAP





通过微电流及优质精华 达到双重效果



深入面部皮肤, 生成胶原蛋白的同时 集中改善皮肤弹性



精华成分的卓越 效果(腺苷,烟酸 酰胺,三大专利成分)



简单的使用方法: 只需将产品 贴放在额头部位即可。

• 银色面膜



黑色面膜(石墨烯&CNT)



银色面膜

- ECOWELL的畅销面膜
- 由兼具美白及皱纹改善功能的精华组成/使用维他命C 涂层贴膜
- 添加含有金属银的图案,达到真正的微电流效果

黑色面膜

- ECOWELL全新推出的面膜
- 由兼具镇定及保湿为主的精华组成
- 添加石墨烯的黑色图案,达到卓越的微电流效果

Pick Points

• 微电流+GALVANIC同时作用

- 一台设备同时实现有助于皮肤弹性的微电流和有助于吸收化妆品活性成分的GALVANIC功能, 有助于高效护肤。

• 真正的微电流疗法

- -严选专用设备的真正微电流效果
- 增加胶原蛋白厚度/增加蛋白质合成/增加弹性纤维

• 真正的家庭式美容

- 节省了前往美容院及护肤品专柜的时间。

• 不同种类的精华(面膜)

-双重功能精华,有助于皮肤内外充满弹性。 含有烟酰胺,腺苷等三种专利成分(银)

• 具有改善肌肤纹理的效果,有助于后续上妆。

自动电源ON/OFF

- 将设备放在额头上, 持续使用20分钟。
- 护肤已不再麻烦。

微电流面膜临床试验结果 (P&K 皮肤临床研究中心)

TEST 1.

增加脸颊部位的皮肤弹性

相较于只使用面膜

效果增加 107.8% (四周期间)



TEST 2.

增加额头部位的皮肤弹性

相较于只使用面膜

效果增加 62.3% (四周期间)



TEST 3.

增加眼周部位的皮肤弹性

相较于只使用面膜


效果增加 101.3% (四周期间)



TEST 4.

皮肤紧致度增加

相较于只使用面膜效果增加 53.2% (四周期间)

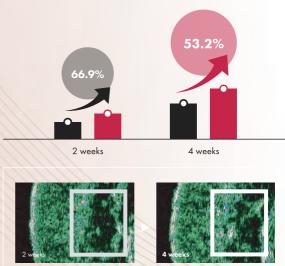


面膜+设备试验组,40岁+,女性,干性皮肤

TEST 5.

肤色改善(美白效果)

相较于只使用面膜效果增加 171.7% (四周期间)

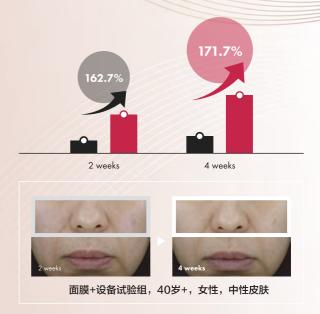

03

美肤设备与乳霜的完美组合 WRINKLELESS PERFECTOR& GALVANIC WRINKLELESS乳霜 集中护理眼部及嘴周



4 in 1 家庭式美容 设备



黑眼圈及 脸部线条护理



简单使用: 自动ON/OFF



手持产品 于脸颊及颈部进行 轻柔涂抹即可。

Pick Points

Galvanic效果

促进有效成分的渗透(i.e., 烟酰胺, 维他命-C, 玻尿酸)

• RED LED (630nm)光疗

使用光线(630 - 660nm)照射,促进胶原蛋白的生成

- NIR (近红外线 / 850nm)光疗
- 微震动

每分钟约12,000的微震动,达到提拉皮肤的功效

• 钻石切割

42面的钻石切割设计,更适合对弧线的脸部皮肤进行按摩

项目 详细事项

电池 LR44×2EA

连续使用时间(实际使用时间) 80分钟(120分钟)

电机 气缸式震动器

功能 微震动,红色LED,Galvanic,近红外线

ON / OFF 人体识别开关



小巧尺寸3 in 1 Galvanic LED助推器



尝试每日按摩三分钟,让皮肤焕然一新吧。



不管是外出还是 随身携带,小巧设计 随时方便使用。



便捷使用, 只需将产品 贴合面部即可。

- Galvanic疗法
 - -又名离子电渗疗法(iontophoresis),比起用手涂抹,护肤品有效成分的渗透力会显著增强。
 - -与添加有烟酰胺,维他命C,玻尿酸的精华等-起使用效果更佳。
- 红色LED (630nm)
 - -密集的LED照射,有助于皮肤内胶原蛋白的生成。
- NIR(近红外线/850nm)光线治疗
 - -温热触感有助于皮肤再生。
- 自动电源ON/OFF
 - -利用人体电流,设备可自动控制ON/OFF。
 - -护肤已不再麻烦。

项目	详细事项
电池	CR2032 / 2EA
连续使用时间(实际使用时间)	120分钟(180分钟)
电机	纽扣类型震动器
功能	微震动,红色LED,Galvanic
ON / OFF	人体识别开关



小巧尺寸5 in弹性改善设备 ELB EMS: 低频疗法



- 兼具EMS+微电流+震动+LED(RED+IR)护理功能的5 in 1智能美容设备
- EMS及微电流,强化弹性护理
- Red+IR LED,可进行深度集中护理。
- 采用智能面膜技术,实现小巧尺寸设计(62mm x 45mm 18g)
- 实现小巧尺寸(62mm x 45mm 18g)
- 携带及使用方便
- 多款设计及颜色(OEM)

项目	详细事项
电池	CR2032 / 2EA
连续使用时间(实际使用使用)	100分钟(160分钟)
电机	纽扣类型震动器
功能	微震动,红色&RGB LED, EMS
ON/OFF	人体识别开关



O 6 hob blux



可将毛孔中的老废物质, 于无刺激的情况下清除干净的 震动清洁设备



轻量,防水(IPX8) 使用简单

- 优秀的便捷性
 - 将手指插入柔性绷带后即可启动产品。(无另设ON/OFF开关)
- 微震动
 - 提供更强清洁力的同时,减轻对皮肤的刺激
 - 约12,000rpm
- 使用简单
 - 小巧设计,便于手持
- 无需进行麻烦的充电
 - 内置电池盒(可更换电池CR2032)

项目	详细事项
电池	CR2032 / 1EA
连续使用时间(实际使用时间)	50分钟(75分钟)
电机	条状震动器
功能	微震动
ON / OFF	电磁开关
防水等级	IPX8



7 HIFU多重护肤设备 O'WELL

- 具备HIFU功能的多重美容设备 (聚焦超声功能)図
- 高频(RF)
- 低频(EMS)
- 微电流
- Galvanic离子电渗疗法
- 热SPA
- LED(R,G,Y)
- IR红外线
- 旋转式触头
- 时尚的充电支架

项目 详细事项 电池 充电式 连续使用时间(实际使用时间) 60分钟(80分钟) 功能 MODE1 (GALVANIC+MICROCURRENT) MODE2 (RF+EMS+NIR) MODE3 (HIFU) ON/OFF 按键式开关

利用微电流为皮肤上妆 利用微电流为皮肤上妆

- 微震动,传统按摩的功效,Galvanic,微电流及加热功能兼备。
- 刮痧按摩效果
 - 清除脸部堆积的毒素
 - 面部皮肤提拉效果
 - 缓解皮肤浮肿现象, 改善血液循环
- 通过一个开始/停止按键进行轻松切换
 - 绿色: Galvanic+微震动模式 - 蓝色: 微电流+微震动模式 - 红色: 加热疗法+微震动模式
- 未另设开关,仅需手握设备,将其接触面部即可自动开启。
- 自动关闭功能: 不使用一分钟后将自动关闭
- USB C型充电端口(两小时满额充电后可连续使用一个小时)

详细事项 项目 Li-Po 600mAh 电池 连续使用时间(实际使用时间) 60分钟(80分钟) 气缸式震动器 电机 微震动,Galvanic,微电流,加热疗法 功能 ON/OFF 人体识别开关,按键式开关 三种模式 Galvanic模式 ● 微电流模式 • 加热模式



