

**cellpod<sup>®</sup>**  
**WRINKLELESS PERFECTOR &  
GALVANIC WRINKLELESS CREAM**





**It's different when cellpod makes it**

cellpod is the beauty device brand of Ecowell Co., Ltd. a research-oriented company specializing in skin technology. We have developed technologies to address customers' skin concerns through systematic research and innovation led by various experts with doctoral degrees.

With our expertise and technological prowess, we not only offer solutions under the cellpod brand but also provide customized solutions for industry partners to implement their own specialized beauty devices.

Through collaboration with industry partners, we aim to grow together and provide optimized skincare solutions to a broader range of customers.








# FOR THOSE WHO

- ✓ Want to manage increasing fine lines
- ✓ Seek to improve habitual skin wrinkles
- ✓ Desire to lift sagging facial lines with elasticity
- ✓ Aim to maintain a vibrant and radiant skin tone
- ✓ Haven't experienced the effects of various products and devices
- ✓ Wish to manage high-end skincare effects at home
- ✓ Seek noticeable skincare experiences

# Internationally recognized smart Kbeauty device

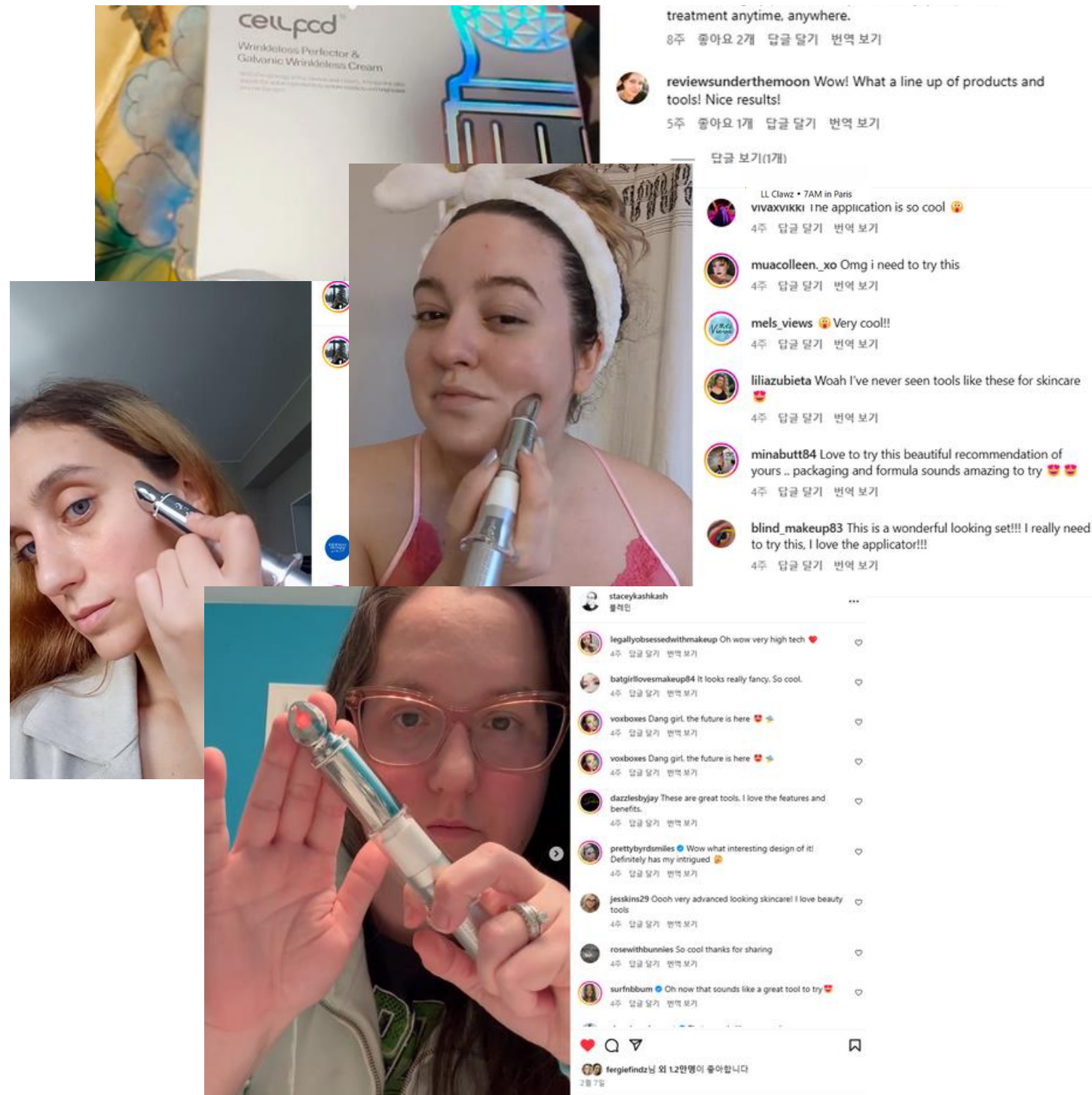

treatment anytime, anywhere.  
8주 좋아요 2개 댓글 달기 번역 보기

reviewsunderthemoon Wow! What a line up of products and tools! Nice results!  
5주 좋아요 1개 댓글 달기 번역 보기

답글 보기 (1개)

LL Clawz • 7AM in Paris  
vivaxvixki I ne application is so cool 🥰  
4주 댓글 달기 번역 보기

muacolleen\_xo Omg i need to try this  
4주 댓글 달기 번역 보기

mels\_views 🥰 Very cool!!  
4주 댓글 달기 번역 보기

liliazubieta Woah I've never seen tools like these for skincare  
4주 댓글 달기 번역 보기

minabutt84 Love to try this beautiful recommendation of yours .. packaging and formula sounds amazing to try 🥰🥰  
4주 댓글 달기 번역 보기

blind\_makeup83 This is a wonderful looking set!!! I really need to try this, I love the applicator!!!  
4주 댓글 달기 번역 보기

staceykashkash  
legallyobsessedwithmakeup Oh wow very high tech ❤️  
4주 댓글 달기 번역 보기

batgirllovesmakeup84 It looks really fancy. So cool.  
4주 댓글 달기 번역 보기

voxboses Dang girl, the future is here 🥰🥰  
4주 댓글 달기 번역 보기

voxboses Dang girl, the future is here 🥰🥰  
4주 댓글 달기 번역 보기

dazzlesbyjay These are great tools. I love the features and benefits.  
4주 댓글 달기 번역 보기

prettybydismiles Wow what interesting design of it! Definitely has my intrigued 🥰  
4주 댓글 달기 번역 보기

jeskins29 Oooh very advanced looking skincare! I love beauty tools  
4주 댓글 달기 번역 보기

rosewithbunbles So cool thanks for sharing  
4주 댓글 달기 번역 보기

surfnbbum Oh now that sounds like a great tool to try 🥰  
4주 댓글 달기 번역 보기

fergiefindz 님 외 12만명이 좋아합니다  
2월 7일

★★★★★

Kenzie Baker Verified

Verified by Judge.me

Love 🥰

Feels so calming and good on my skin! Definitely worth the buy!

★★★★★

SarahRD Verified

Works so good!

After just using this perfector and cream for just a short time my skin is looking more hydrated and my fine lines are less noticeable. I highly recommend this product.



**Wrinkleless Perfector & Galvanic Wrinkleless Cream come with galvanic, red LED, and massage vibration. Providing ease of use, they offer skin elasticity and radiance while aiding in the improvement of fine lines and wrinkles.**

★★★★★

Shimaa Verified

Must have and try!

This is the first time WRINKLELESS PERFECTOR & GALVANIC WRINKLELESS CREAM has been used with the device. It has effective red, galvanic and neon LED massage vibrations. This device helps with good absorption, regeneration, skin elasticity and radiance. It is easy to use and control. You put a tube of cream and put it in the device with a little pressure as you need, and then you start massaging the skin from the inside out. You will notice the difference after the first use, as the skin becomes plump, clear, and moisturized. Helps reduce fine lines and wrinkles. I like the experience because it is good for refreshing the skin and renewing.



**Wrinkleless Perfector & Galvanic Wrinkleless Cream provide a quadruple effect, offering easy use while brightening and rejuvenating the skin, and restoring elasticity. They leave the skin feeling moisturized and firm after massage.**

★★★★★

Pattil Verified

Relaxing

This Wrinkleless Perfector & Galvanic Wrinkleless Cream has a quadruple effect, it Galvanic Iontophoresis, Red LED, NIR (Near Infrared) and Vibration Massage. This device helps products absorb better, restores elasticity, brightens and regeneration of the skin. It is very easy to use, insert the cream container into the device and then press the cream container to dispense, to the tip, the amount you need. Then just gently massage the device on your skin from the inside out, as if smoothing out fine lines. The device is very easy to hold and control. I love how this feels on my skin, it is very relaxing. My skin looks plump, hydrated and refreshed after using this beauty device. The only con I have for this product, which isn't a deal breaker, is that it does take button cell batteries. I would prefer it to be rechargeable. Overall, I love this device and would recommend it.



★★★★★

Jadori Verified

Love it!

I love the ease of use, the cream is also non-sticky and performs wonderfully with the device! Results are amazing!






# 4-IN-1 Skincare Solution

Galvanic, Near-Infrared (NIR),  
Micro-Vibration, RED LED



## **GALVANIC**

Enhanced Penetration of Active  
Ingredients (such as Vitamin C,  
Hyaluronic Acid, etc.)

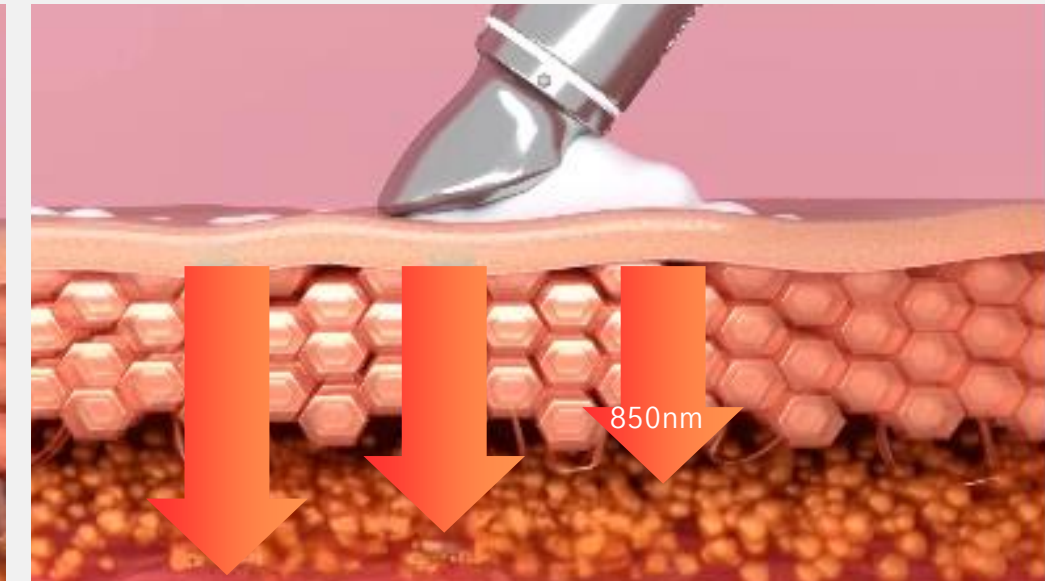


## **NEAR-INFRARED (850nm)**

Warmth-Induced  
Regeneration Enhancement



## **MICRO-VIBRATION**

Lifting Effect with Approximately  
12,000 Micro-Vibrations Per Minute

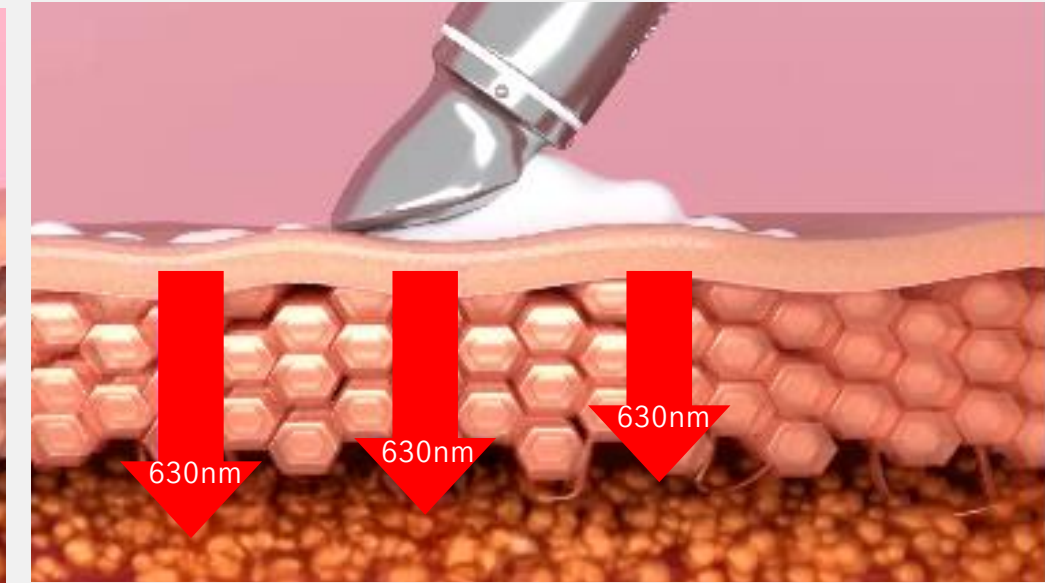

## **RED LED (630~660 nm)**

Light Exposure to Aid  
Collagen Production



# SMART SPECIAL CARE

Galvanic, LED, Micro-Vibration,  
Near-Infrared emitted upon skin contact

Hold the chrome handle  
and touch it to the skin  
**the LED light will turn on and  
start operating immediately.**

If you pump the cream cylinder at the  
bottom of the product one or two times,  
cream will be dispensed from the head.

## Diamond Kneading

Optimally designed  
for skin massage  
with 42 surfaces





# ALL- AROUND SKINCARE

Simple application to any area  
instantly provides comprehensive,  
all-around skincare.




**EYE AREA**



**FOREHEAD**



**NASOLABIAL  
FOLDS**

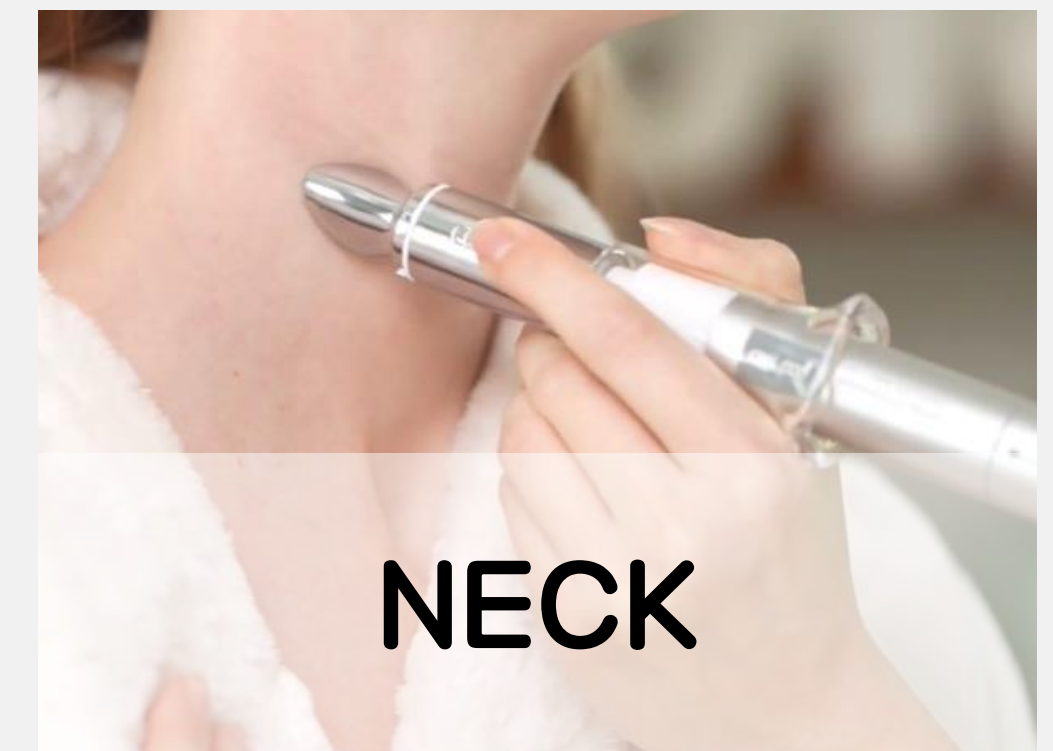

**NECK**

# GALVANIC WRINKLELESS CREAM

Total solution eye cream addressing concerns such as under-eye wrinkles and smile lines

Whitening  Wrinkle Improvement  Moisturizing  Vitality  Soothing

Containing a complex of 10 vitamins, energizes the skin

Niacinamide, panthenol, ascorbic acid peptide, menadione, Biotin, linoleic acid, cyanocobalamin, tocopheryl, Sodium Ascorbate Phosphate, Ascorbic Acid Glucoside

Enriched with Idebenone, it smoothens and revitalizes the skin

Through research from the American Academy of Dermatology and other papers, it is known to have antioxidant effects up to 4 times that of Vitamin C and 10 times that of Coenzyme Q10.

Revitalize damaged and lackluster skin with **glutathione**

Contains **patented ingredients that excel in inhibiting melanin production and aiding in skin damage improvement and prevention**, resulting in resilient and radiant skin

\*Patented Ingredients\*


Korean Angelica Root Extract (Patent No. 10-1578461)

European Strawberry Leaf Extract (EP1761271B1, Synrise Patent)







# EASY TO USE

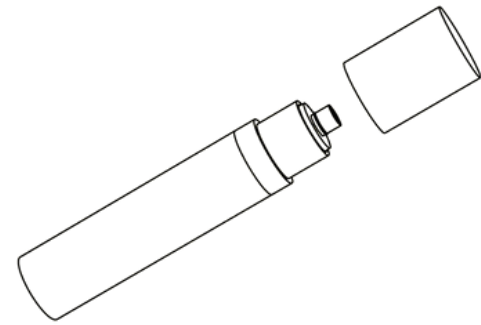

1. Remove the lid of the cream cylinder.

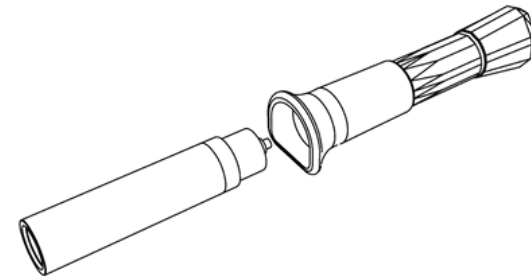

2. Attach the cream cylinder to the main device. Press down until you hear a click

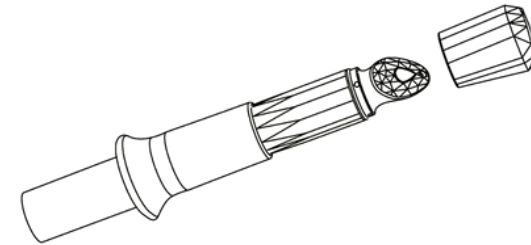

3. Open the lid of the main device head and proceed with use.

# HOW TO REPLACE THE BATTERY

Easy battery replacement with LR44 batteries for excellent portability, no need for recharging.



1. Hold the main device and rotate the cream container (transparent) to detach it from the device.



2. Slowly lower the main body downwards. (Be careful as the battery may fall out)



3. Replace the battery. Check the polarity markings (+) (-) and insert the new battery correctly.




4. After replacing, slowly rotate the cream cylinder connection part (transparent) on the main body to secure it in place.



# CLINICAL TRIALS AND CERTIFICATIONS

## Results of skin absorption test in human subjects

Improvement of skin absorption by 248.584% when using Wrinkle-less Perfector in conjunction with cream compared to cream alone.

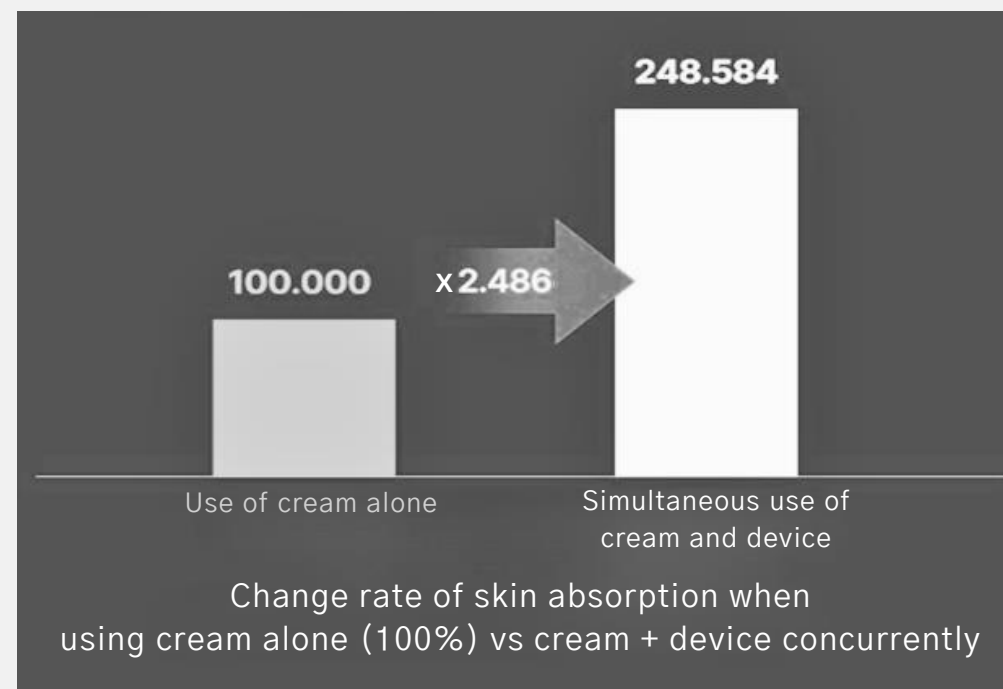

Testing Institution: P&K Skin Clinical Research Center Co., Ltd.

Testing Period: November 01, 2021 ~ November 09, 2021

Number of Participants: 21

Test Area: Entire face

Efficient delivery of active ingredients via microcurrents and gentle vibrations exceeding 5,000 times per minute.

**Abstract:** This study aimed at investigating the effect of electrical current profile upon the iontophoretic transport of (i) ascorbic acid (AA) and (ii) ellagic acid (EA), into porcine skin in vitro, and the impact of the physicochemical properties of both actives on their mechanism of transport when formulated in cosmetic compositions. The experiments were performed using a proprietary iontophoretic device containing a roller to apply the formulation. Three current profiles were tested: (i) galvanic direct current (DC), (ii) square unipolar pulse current (SPC), and (iii) galvanic direct current (DC) + pulse current (PC). The skin samples were collected at different sampling points, extracted and analyzed by HPLC. **Results suggested that the DC + PC mode for only 5 min was able to significantly increase the delivery of AA from o/w cosmetic compositions.** The use of this current profile might improve the skin penetration of AA due to electromigration and passive diffusion, the latter being facilitated by the physical enhancement method. The SPC mode significantly improved the passage of EA in its neutral form from cosmetic o/w formulations by electroosmosis. Tailoring specific electrical current modes considering the ionization state of active ingredients would allow the design of short and personalized cosmetic treatments that significantly improve the penetration efficiency of the active ingredients and possibly reduce the doses applied.

**Keywords:** iontophoresis; topical penetration; ascorbic acid; ellagic acid; pulsed current; galvanic current; square unipolar pulse current

\*Published in the SCI-grade journal *Pharmaceutics*, efficacy data regarding absorption capabilities.

## KC Compliance Assessment Certification

Certification of Compliance Registration for Broadcasting and Communication Equipment

방송통신기자재등의 적합등록 필증 Registration of Broadcasting and Communication Equipments	
상호 또는 성명 Trade Name or Registrant	주식회사 에코텔
기자재명칭(제품명칭) Equipment Name	아이 갈바닉
기기부호/추가 기기부호 Equipment code /Additional Equipment code	UBD
기본모델명 Base Model Number	551-002-1010
파생모델명 Series Model Number	
등록번호 Registration No.	R-R-Ec W-551002-1010
제조사/제조국가 Manufacturer/Country of Origin	주식회사 에코텔 / 한국
등록연월일 Date of Registration	2021-09-15
기타 Others	
<p>위 기자재는 「전파법」 제58조의2 제3항에 따라 등록되었음을 증명합니다. It is verified that foregoing equipment has been registered under the Clause 3, Article 58-2 of Radio Waves Act.</p> <p style="text-align: right;">2021년(Year) 09월(Month) 15일(Day)</p> <p style="text-align: center;">국립전파연구원 국립전파연구원장 Director General of National Radio Research Agency</p> <p style="text-align: center; font-size: small;">※ 적합등록 방송통신기자재는 반드시 "적합성 평가표시" 를 부착하여 유통하여야 합니다. 의반시 과태료 처분 및 등록이 취소될 수 있습니다.</p>	



# WRINKLELESS PERFECTOR & GALVANIC WRINKLELESS CREAM REBORN UNDER YOUR BRAND

**Manufactured as OEM/ODM and distributed  
and sold under your own brand**

Customization options are available for printing and engraving your brand logo, changing product and package colors, and even customizing the cream content within the cream cylinder (also available without cream content) to align with your brand.

For more details, please contact the person in charge at the Ecowell Sales Department via phone, email, or inquire through our website.

EMAIL: [sales@ecowell.co.kr](mailto:sales@ecowell.co.kr)

Homepage: [ecowell.co.kr](http://ecowell.co.kr)



# PRODUCT SPECIFICATIONS

<b>Model Name</b>	5651-002-1010
<b>KC Certification Number</b>	CTK2109-EK-065
<b>Equipment Name</b>	Wrinkless Perfector & Galvanic Wrinkless Cream
<b>Rated Voltage/Power</b>	3.0 V x 2 ea.
<b>Package Size (L x D) / Weight</b>	156 mm x 178 mm x 37 mm / 250 g
<b>Material</b>	ABS, Chrome Plating
<b>Manufacturer / Country</b>	Ecowell Co., Ltd. / South Korea
<b>Cosmetics Responsible</b>	Ecowell Co., Ltd.
<b>Battery</b>	LR44 / 2 EA
<b>Continuous Time (Actual Usage Lifespan)</b>	80 Mins (120 Mins)
<b>Motor</b>	Cylindrical vibration motor
<b>Function</b>	Galvanic therapy, Red LED, Near-infrared, Micro-vibration massage
<b>ON / OFF</b>	Biometric switch





# FAQ

## **Q. Where is the switch located?**

A. This product operates automatically when the main body is held and the galvanic head part is in contact with the skin, without a separate switch.

## **Q. The product keeps running without stopping.**

A. After using the product, please clean any serum or cream residue on the device (especially the head part). You can hold the device with one hand using a tissue or cloth and gently wipe with the other hand. If the same symptom persists even after completely removing the residual cream, please contact the customer service center.



### **[Caution during Use]**

Clean the device body to prevent cream from connecting the galvanic head and the main body.



### **[Close the protective cap]**

Align the narrow part of the protective cap downwards when closing.