SMOOTH SKIN
pure
ice

Powerful · Intelligent · *Skin Friendly*

IPL HAIR REMOVAL DEVICE

USER GUIDE

English..... 4

Français36

Deutsch.....69

Italiano 102

Español 134

Portuguese.....167

Welcome to **SMOOTHSKIN PURE ICE!**

Congratulations! You've opened up a world of unparalleled beauty, smoothness and convenience. Come and join in the conversation at **www.smoothskin.com** where you can find expert advice and support throughout your treatment.

CONTENTS


1. THE SMOOTHSKIN PURE ICE SYSTEM.....	5	10. ICE COOLING INDICATORS	25
2. HOW DOES SMOOTHSKIN PURE ICE WORK?.....	6	11. CHOOSING THE RIGHT INTENSITY LEVEL	26
3. IMPORTANT SAFETY INFORMATION	7	12. CLEANING, MAINTENANCE, STORAGE.....	27
4. POTENTIAL RISKS, SIDE EFFECTS AND SKIN REACTIONS	13	13. TRAVELLING WITH SMOOTHSKIN PURE ICE.....	28
5. WHAT IS THE INTENDED USE OF SMOOTHSKIN PURE ICE	15	14. WARRANTY	28
6. YOUR SUITABILITY	15	15. TROUBLESHOOTING.....	29
7. HOW TO USE SMOOTHSKIN PURE ICE	17	16. LABELS AND SYMBOLS	31
8. THE SKIN TONE SENSOR	23	17. TECHNICAL SPECIFICATIONS.....	32
9. SKIN CONTACT SENSOR	24		

1. THE SMOOTHSKIN PURE ICE SYSTEM

PACKAGE CONTENTS

Inside the packaging of **SmoothSkin Pure Ice** you will find the following:

- Handset and Base Unit
- Power Cable
- User Manual
- Device Storage Bag

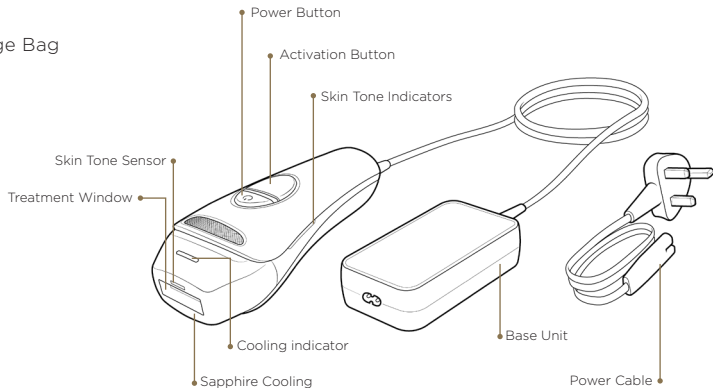2. HOW DOES SMOOTHSKIN PURE ICE WORK?

SmoothSkin Pure Ice is designed to help break the cycle of hair growth. Light energy is transferred through the skin's surface and is absorbed by melanin present in the hair shaft. The absorbed light energy is converted to heat energy (below the surface of the skin), which reduces the hairs growing back.

Treated hairs naturally fall out over the course of a few days to 1-2 weeks.

WHAT TO EXPECT

Immediately post treatment, you should not see any significant side effects from treatment (Refer to section 4 for more information on side effects).

In the first few weeks following the initial treatments, you will still see some hairs growing. These are likely to be hairs that were missed during treatment i.e. they were lying dormant and not in their growth phase (Anagen) during treatment; when IPL is most effective.

Around 6 weeks into the 12 week treatment programme you should see a reduction in hair growth. However, many hairs may still not have been treated in their growth phase. It is important to continue the weekly treatments.

After the 12 week programme, you should see a significant reduction in hairs within the treated area. Any remaining hairs should be finer and lighter in colour. Continued monthly treatments, or treatments as required, should maintain the reduction in unwanted hair.

3. IMPORTANT SAFETY INFORMATION

CONTRAINDICATIONS

DO NOT USE THE DEVICE...

- if you have very dark skin (Darker than skin tone v). The Skin Tone sensors located near the output window of the device will determine the tone of your skin. If the Skin Tone sensors determine your skin is too dark for safe use, operation will be prevented.
- if your skin has been artificially or naturally tanned recently or is burned from sun overexposure. Your skin may be extra sensitive following sun exposure and particularly susceptible to burning, blistering, discoloring, or scarring if you use the device. Wait 1 week until sunburn or tan disappears before using the device. Treating tanned skin with IPL may result in permanent Hyper (darkening) or Hypo (lightening) changes to the treated skin.
- directly on the nipples, genitals or in the anus. These areas may have a darker skin color and/or greater hair density. Using the device in these areas may cause discomfort/pain or injure (burn, discolor or scar) your skin.
- if you have a history of skin cancer or pre-cancerous lesions (e.g. nevi or a large number of moles).
- on dark brown or black spots such as birthmarks, moles or warts in the area you wish to treat. Using the device may injure your skin or make existing conditions worse. You may experience side effects such as burns, blisters and skin color changes or scarring.
- on areas where you have tattoos or permanent make up. The pigment will absorb the light which may cause side effects such as burns, blisters and skin color changes or scarring.

- if you have a known sensitivity to sunlight (photosensitivity) or are taking medication that makes the skin more sensitive and causes photosensitivity, e.g. Retin A, Accutane and/or other topical retinoids. Always check the instruction leaflet that comes with your medicinal product, to see if photosensitivity is listed as a side effect. Use of the Device on photosensitive skin could result in skin damage such as swelling or blisters. Seek medical advice prior to application of the appliance as use may harm the skin.
- if you have any skin condition in the treatment area, including psoriasis, vitiligo, eczema, acne, herpes simplex or active infections or wounds. Using the device may injure your skin or make existing conditions worse.
- on the scalp, lips, nostrils and ears.
- if you are under the age of 18 as this device has not been tested on these individuals.
- if you are pregnant, lactating, as this device has not been tested with these individuals.
- if you have a Chronic Medical Condition.
- if you have varicose veins in the area you wish to treat.
- if you have had a skin peel treatment on the area you wish to treat in the last 30 days as the skin can be more sensitive to the light which may cause pain or burns.

If you are unsure whether this device is safe for you to use, please consult your physician or dermatologist.

If you take any kind of medication regularly and over a long period of time, please consult your physician for any potential impact on skin sensitivity to light.

GENERAL WARNINGS

DO NOT...

- leave any visible hair on the skin in the area to be treated, as this could lead to increased pain, redness of the skin or potential minor burning of the skin.

- expose treated areas to the sun. Wait at least 7 days after treatment before exposing treated skin to direct sunlight. Your skin may be extra sensitive following IPL treatment and particularly susceptible to sunburn. Apply sunscreen (SPF 15 or greater) to the treated area or cover the treated area with suitable clothing.

The Device should not be used over any fillers or botox as the skin can be more sensitive to the light which may cause pain or burns.

This appliance is not intended for the treatment of medical conditions. Medical advice should be sought to address, for example, moles, skin rash, itchy skin, skin fungus or infection, skin bumps, or skin tags. Using the device may injure your skin or make existing conditions worse.

The Device should be kept out of the reach of children under 18 years of age. Children should not use, play with, maintain, or clean the device. There are a number of potential risks – exposure to the light output causing burns or eye damage, electric shock, or strangulation from the cables.

The Device should be kept out of reach of pets and pests which could damage the device. Potential risks include reduced function and electric shock.

The Device is not intended for use by persons with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge. Unless they have been given supervision or instruction concerning use of the equipment to avoid the risk of fire or burns.

The Device Sapphire Block is a critical component that may get warm during use.

The Device is not to be used on skin parts which have an increased likelihood of skin cancer, open wounds or rashes, or swollen, red, irritated, infected, or inflamed areas or skin eruptions.

WARNING: Damage of the skin may occur after prolonged or repeated surface application on one site. Follow the instructions provided and avoid over-usage which included excessive passes, stacking of pulses and/or increased frequency of use.

EYE SAFETY

WARNING: Possible eye injury (Potentially leading to loss of vision) or skin injury if the instructions provided are not followed.

The device emits flashes of Intense Pulsed Light. Direct exposure is potentially harmful to your eyes. Take care to follow the safety precautions below.

The device includes Skin Contact sensors that only allow the device to flash when in good contact with skin.

DO NOT:

- attempt to override the sensors.
- look directly into the treatment window of the handset when the system has electrical power and is switched on
- attempt to activate (flash) the device towards the eyes. It is recommended to look away from the handset when flashing onto skin.
- treat areas around the eye (eyelids, eyebrows, or eyelashes). Misuse can lead to eye damage, protect eye from exposure.

ONLY attempt to activate (flash) the system when the front of the handset is in good contact with the area you wish to treat.

If at any time the device activates when both sensors are **NOT** in contact with the skin, **STOP** use immediately and contact your retailer.

PRECAUTIONS

Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage treated on the face and neck.

The Device is intended for a single user.

DO NOT...

- use the Device on any area where you may want hair to grow back. The treatment is intended for permanent reduction in hair regrowth.
- use the Device on man's face, jaw, or neck. The hair in a man's beard area is too dense and using the Device in these areas may result in skin injury.
- use on areas that have recently received local anaesthesia. The skin may be more susceptible to damage such as redness, swelling or burns.
- use any other components, other than those provided for use with the Device. This may damage the unit and cause it to stop working.
- use if the Device becomes uncomfortably hot to touch. This may indicate that the unit is damaged.
- use the Device if the treatment window is not clean and free from debris. Any debris on or around the treatment window during operation may damage the device, increase the pain of treatment, or result in redness/swelling of the skin. To clean the device, follow the instructions provided.
- subject the equipment to heavy shocks (e.g. drop) to avoid risk of damage to the lamp. If the lamp is damaged, the device will not flash.
- attempt to activate (flash) the Device onto any surface other than skin. This could damage the surface or cause a fire.
- exceed the usage regimen specified. This will not improve the results and could result in increased skin sensitivity or redness or swelling of the skin.

The Device is not recommended for use on red, grey, white, or blonde hair as it is not effective and unsuitable on these hair colors.

ELECTRICAL AND FIRE SAFETY

DO NOT...

- try to open the device, as this may expose you to dangerous electrical components which could lead to electric shock. The device is not intended to be modified or serviced.
- use near bathtubs, showers, basins, or other vessels containing water. Water and electricity are a dangerous combination! To avoid risk of electric shock, **DO NOT** let the device become wet. If the Device becomes wet, it should not be used. To clean the device, follow the instructions provided.
- block the cooling vents. The device may overheat and pose an electric shock risk or fire hazard.
- operate the device if the mains plug is not easily accessible. You should be able to disconnect at any time.
- use if the Sapphire Block is cracked or missing or not in correct place.
- attempt to use this device other than for removing unwanted hair.
- use the device on a surface other than skin. Do not press the treatment button when the treatment window is not placed in full contact with the skin.

Regularly inspect the device housing and Sapphire Block for damage. **DO NOT** use if the Device is damaged; for example, cracked Handset, cable damage, cracked or broken glass on the handset, etc. There could be a risk of electric shock or overexposure to light resulting in skin damage.

WARNING: If the wall socket used to power the equipment has poor connections, the plug of the equipment may become hot. Make sure you plug the equipment into a properly installed wall socket to avoid the risk of fire and burns.

WARNING: Do not leave the equipment unattended when it is switched on to avoid the risk of fire or burns.

4. POTENTIAL RISKS, SIDE EFFECTS AND SKIN REACTIONS

Use of SmoothSkin Pure Ice could result in side effects.

The table below lists the known skin reactions that may occur following treatment with SmoothSkin Pure Ice.

In case of adverse events please contact your national competent authority and the device manufacturer.

POSSIBLE SIDE EFFECT	HOW TO ASSESS & REACT
Mild pain / discomfort in the area being treated.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed, and the pain should diminish with continued use.
Warm feeling or tingling sensations during treatment, which typically disappear after a few seconds to a minute and decrease with continued use.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed.
Itchiness in the treatment area.	This is quite common for IPL treatments, and should subside after a short period. You can keep on using the device as instructed. Do not scratch the area.

POSSIBLE SIDE EFFECT	HOW TO ASSESS & REACT
Skin redness during or after treatment which disappears within minutes.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed once skin redness has disappeared.
Skin redness which does not disappear within minutes and lasts 24-48 hours.	Stop using the device immediately and consult your physician before using it again.
Pain or discomfort that is intense during treatment or persists after a treatment.	Stop using the device immediately and consult your physician before using it again.

IN VERY RARE CASES	HOW TO ASSESS AND REACT
Swelling and redness that does not disappear within two to three days.	Stop using the device immediately and consult your physician before using it again.
Temporary changes in skin colour (lightening or darkening).	If your skin color changes, stop using the device immediately and consult your physician.
Blistering, scarring or burning of the skin.	Stop using the device immediately and consult your physician before using it again. Apply a cold pack to the area. Treat with antiseptic or burn cream.

5. WHAT IS THE INTENDED USE OF SMOOTHSKIN PURE ICE

INTENDED USE/PURPOSE

This device has a non-medical purpose and is indicated for the removal or permanent reduction in unwanted hair.

INDICATIONS FOR USE

The SmoothSkin Pure Ice Hair Removal Device is indicated for the removal of unwanted hair. The SmoothSkin Pure Ice is also indicated for the permanent reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs regrowing when measured at 6, 9 and 12 months after the completion of a treatment regime.

6. YOUR SUITABILITY

BODY AREAS

SmoothSkin Pure Ice is suitable for use on body hair and female facial hair below the cheek line.

SKIN COLOR

SmoothSkin Pure Ice is suitable for use on light, medium and dark skin tones, up to and including skin tone 5.

SmoothSkin Pure Ice incorporates a skin tone sensor that measures the colour of your skin and will prevent use if your skin tone is too dark.

Dark Skin tones are rich in melanin which could absorb too much energy and potentially lead to skin tone damage.

HAIR COLOR

SmoothSkin Pure Ice is suitable for use on naturally black or brown hair. It may not be as effective on white, grey, blonde or red hair.

7. HOW TO USE SMOOTHSKIN PURE ICE

WARNING: Follow each step in the treatment process ensuring you complete all aspects of each step before progressing to the next.

Get the most out of your device with our "how to use" videos - visit www.smoothskin.com/user-guides

STEP 1: PREPARE THE AREA TO BE TREATED

Remove all visible hair by shaving. Make sure there is no hair remaining above the surface of the skin as this may cause injury. This also prevents debris from covering the front of your device.

Cleanse the treatment area and pat dry and thoroughly. **DON'T** use any gels or creams with this device!

STEP 2: POWER THE DEVICE

Plug the power supply into an electrical outlet. Switch on Pure Ice by pressing the power button. You will see a sweep of white lights as the device switches on and the fan will start up.

STEP 3: DO A PATCH TEST

Before your first treatment we RECOMMEND you test your skin for reaction to SmoothSkin Pure Ice.

The patch test area should be approximately 3cm x 2cm in size (equivalent to 2 flashes applied to the skin side by side). Follow Step 5 to treat the area.


Wait 24 hours following the patch test, to ensure your skin is suitable for treatment and there is no adverse reaction to the light energy. If there is no reaction after 24 hours, you may treat the area around the patch test. The 'patch tested area' should **NOT** be re-treated for at least 1 week.

STEP 4: START YOUR TREATMENT WITH SMOOTHSKIN PURE ICE

- A. Turn on the device by pressing the power button. Please refer back to step 1 of this section, for more details.
- B. Press the handset firmly against the area to be treated ensuring the Skin Tone Sensor and the Skin Contact Sensors are in full contact with skin.
If the skin tone sensor does not detect a valid skin tone and / or the skin contact sensors are not in full contact with skin, then the handset cannot be activated.
When the skin tone sensor measures a valid skin tone and the skin contact sensors are in full contact with skin the Skin Tone Indicators will illuminate WHITE to indicate the device is ready to flash.

Refer to section 8 and 9 of this manual for further details.

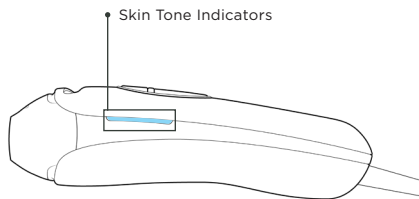CORRECT

Indicator On

Skin Detection

INCORRECT

Indicator Off

No Skin Detection


C. Select the right treatment method:

Glide: recommended for larger areas such as legs. Place the handset on the skin, **press and hold** the activation button continuously and glide the handset along the skin between flashes. Try to move the handset at a steady speed that gives full coverage without overlap or missed areas.

Stamp: for small areas, like underarms and face. Place the handset on the skin, **press and release** the activation button. Remove from the skin before carefully re-applying to the next treatment area

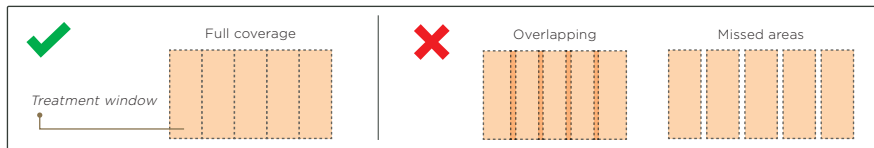HOW TO TREAT THE PUBIC AREA (FEMALES ONLY)

Before treating the pubic area, make sure you have done a skin test. Remove any visible hair, lotions or creams from the area to be treated. You might want to confirm full hair removal using a mirror.



- Stand, sit or crouch down to reach all areas to be treated.
- For pubic area treatment, we recommend using the stamping method.

Select Intensity Level 1 by pressing the Power button until you only see one indicator light:

1. Place the treatment window firmly on your skin, ensuring all skin sensors have skin contact. Stretch your skin with your fingers to improve skin contact.
2. The Skin Contact Indicators show white LEDs.
3. Press and release the treatment button. The device will emit a flash.
4. Lift off the device and place it onto the next adjacent spot to be treated. Repeat steps 2 to 4 and continue spot by spot until you have completed the entire area.

NOTE: Since your skin might have a darker color or pigmentation in the pubic area, the device might not flash even if it has good skin contact.

- If redness or a burning sensation occurs and lasts for several minutes, apply a cool cloth/pad to soothe the area.
- If any pigment changes in the skin, e.g., lightening or darkening, are observed, stop using the device. Seek medical advice to determine the cause which may be due to an underlying medical condition.

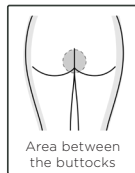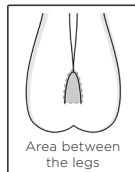WARNING: Make sure not to cover the vents on the device whilst treating, as this can lead to the device overheating.

POWER MODE Timeout

When the unit is left “inactive” i.e. no buttons are pushed for 10 minutes then the unit will revert to STANDBY mode, and the skin sensor lights will switch off. To return to READY mode for treatment, press the power button.

EMERGENCY STOP Function: In case of an emergency/abnormal device operation, remove from skin, unplug the device immediately.

STEP 5: YOUR SMOOTHSKIN PURE ICE TREATMENT SCHEDULE

SmoothSkin Pure Ice should be used once a week until the desired results are achieved.

STEP 6: SWITCHING OFF SMOOTHSKIN PURE ICE


Switch off SmoothSkin Pure Ice by pressing the power button until the indicator lights and fan switch off. Unplug the device from the electrical socket.

Refer to Section 12 for details of how best to look after your SmoothSkin Pure Ice.

8. THE SKIN TONE SENSOR

The skin tone sensor measures the colour of your skin and sets the correct level of output power needed. If your skin is too dark for treatment with SmoothSkin Pure Ice, the skin tone sensor will detect this and prevent the unit from operating. RED lights will be shown in the skin tone indicator if your skin is too dark.

The status of the skin tone sensor can be seen on the skin tone indicator located on both sides of the handset. If a valid skin tone is detected, the skin tone indicators will illuminate WHITE. The number of lights illuminated will be dependent upon your exact skin tone.

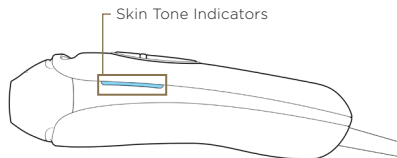- 1 Red light** = Skin is too dark, no treatment possible
- Less lights** = Darker skin and lower energy setting
- More lights** = Lighter skin and higher energy setting


9. SKIN CONTACT SENSOR

SmoothSkin Pure Ice is fitted with skin contact sensors, that check your skin contact before every flash to ensure no misfires.

The sensors must be in full contact with the treatment area in order for the device to flash. If the device is not in full contact with the treatment area, the device will not fire.


When there is not good skin contact no lights will be shown in the Skin Tone Indicator.

Try adjusting the positioning of the device to ensure full contact is achieved, particularly on more awkward areas like arms and shins.

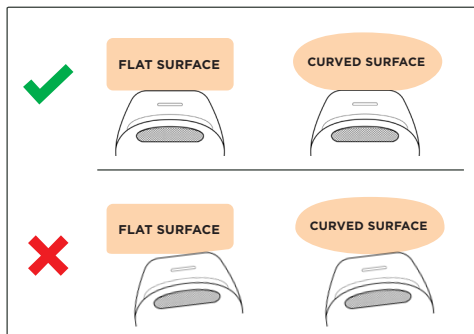10. ICE COOLING INDICATORS


When the device is on, the cooling indicators will be ON indicating the status of the cooling:

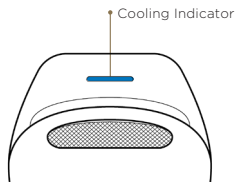Blue Light

INDICATORS ARE BLUE -


The device is ready to use.

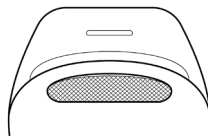White Light

INDICATORS ARE WHITE -

The device is warming up. If it's uncomfortable, stop using the device and wait a few minutes for the cooling indicators to turn blue.

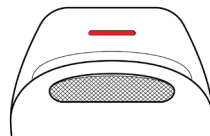Maximum >
Flashing RED

INDICATORS ARE RED -

The device is too warm to use. Keep the device plugged in and wait for a few minutes for it to cool. The cooling indicators will turn blue when it is ready to use.

11. CHOOSING THE RIGHT INTENSITY LEVEL

Push the power button to reduce the light intensity. You can toggle through the 3 following settings:

Intensity level 3 - is the highest intensity setting of the device. It will deliver the best results, as it is using the maximum power output possible for your unique skin tone. This level works particularly well on stubborn areas where the hair tends to be denser.

Intensity level 2 - Default level when the device is switched on, has a lower light intensity than level 3. Recommended for first time IPL users.

Intensity level 1 - is the lowest but fastest intensity level. Recommended for sensitive areas like the pubic area, and for quick maintenance treatments after your initial treatment regime to keep hair regrowth at bay.

When in any of these modes, you can also choose whether to use Glide Method or Stamp Method. Refer to section 7 for more information on the treatment methods.

12. CLEANING, MAINTENANCE, STORAGE

After treatment, always switch off SmoothSkin Pure Ice by unplugging the electrical outlet.

After use, the handset, skin tone sensor and treatment window should be inspected for damage and cleaned with a lint free cloth.

WARNING: The Sapphire window can get hot. Wait 2 minutes after use to allow it to cool down before cleaning.

NEVER use abrasive materials or chemical cleaning fluids such as alcohol or IPA, as these can damage the device and skin sensors leading to a reduction in performance.

WARNING: SmoothSkin Pure Ice is not waterproof. Although a slightly damp cloth can be used to clean the handset, **NEVER** expose the device to excessive amounts of liquid as this could damage it or create a safety hazard.

To ensure maximum performance of your SmoothSkin Pure Ice, it is particularly important to keep the front of the handset and reflectors free from dirt, hair or other debris.

If necessary, the handset can be carefully cleaned by using a slightly damp cloth or cotton bud.

Store SmoothSkin Pure Ice in a cool, dry place. Make sure that the treatment window and skin tones sensor on the handset are protected from damage.

Regularly check the device (including cords) for visible signs of damage. In case of damage or cracks, stop using the device - and visit the website for support details.

To prevent strangulation or tripping hazard, always unplug the device from power socket when not in use.

13. TRAVELLING WITH SMOOTHSKIN PURE ICE

The power supply unit of SmoothSkin Pure Ice is universal and will work with all commonly encountered mains supplies. Providing you have a suitable adaptor for the mains power cable plug, SmoothSkin Pure Ice will operate as normal.

14. WARRANTY

We grant a 2 year warranty on the product starting on the date of purchase. Within the warranty period, we will check faults in materials or workmanship free of charge and replace the complete device, at our discretion.

This warranty extends to every country where this device is supplied by CyDen or its appointed distributor.

This warranty does not cover damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The warranty becomes void if repairs are undertaken by unauthorised persons and if original SmoothSkin Pure Ice parts are not used. To obtain service within the warranty period, please visit our website

www.smoothskin.com

15. TROUBLESHOOTING

PROBLEM	SOLUTION
The device does not start, or the skin tone indicators on the side of the handset do not light up when placed against your skin.	• Make sure the power supply is plugged in. • Make sure the handset is in READY mode (press the Standby/Ready button once). • Ensure the skin sensors are pressed against the skin. • Ensure the skin tone sensor is pressed against the skin.
Device indicates error indication during use	Switching the device off and back on may clear the error. If so, it is safe to continue using the device. If the error persists, then stop using the device and refer to the website.
Whilst pressed against your skin, the first skin tone indicator on the side of the handset is illuminated RED, and SmoothSkin Pure Ice will not work.	The area you are trying to treat is too dark for treatment. SmoothSkin Pure Ice includes an intelligent skin tone sensor that automatically detects your skin tone and selects the optimal treatment setting. Due to the nature of Intense Pulsed Light technology, SmoothSkin Pure Ice is not for everyone and cannot treat very dark skin.
The skin tone indicators on the side of the handset are on or are flashing in an unfamiliar pattern including a red LED indicator light.	SmoothSkin Pure Ice will display any error messages through combinations of lights the on the skin tone indicators on the sides of the handset. This indicates that there is a fault with the device. Stop using the device and refer to the website.

PROBLEM	SOLUTION
SmoothSkin Pure Ice handset, Sapphire block, power supply or cable are damaged, broken, cracked or appear defective.	DO NOT USE. If you are in any doubt about the safety of SmoothSkin Pure Ice or suspect it is damaged in any way, it must not be used, and you should refer to the website.
The Cooling Indicators flash red, and the device will not flash	• The Sapphire block is too hot for use. Keep the device plugged in and wait a few minutes for it to cool. The Cooling indicators will turn blue when it is ready to use. • Ensure the vents on the underside and top of the device are not blocked. • Ensure the room temperature is within the recommended operating range. See section 17 for details.
The Sapphire Block feels hot during use.	There may be something wrong with the device. Do not continue to use if the Sapphire block feels hot.

16. LABELS AND SYMBOLS

	Type BF Applied Part		UK Conformity Marking
	Class II Equipment		Conformity mark of the EU with Notified Body number
	Follow Instructions for use		FCC EMC compliance for North America
	Keep Dry		MET Mark United States & Canada
	WEEE waste disposal mark of the European Union		Australia Regulatory Compliance Mark
	Warning Optical Radiation from device May Cause Eye Injury		China RoHS Environmental Protection Use Period
	Importer		Emergo Europe

FCC Statement:

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

17. TECHNICAL SPECIFICATIONS

SmoothSkin Pure Ice is a filtered broadband Intense Pulsed Light system with the following technical specifications:

- **Repetition Rate:** Manually pulsed every 0.30 to 3.5 seconds (dependent on treatment setting).
- **Pulse:** Partial Free discharge FWHM = 2.5 ± 0.5 ms
- **Max Optical Output:** 2 - 6J/cm²
- **Wavelength:** 510 - 1100nm
- **Power Line Input:** 100-240V-, 50-60Hz, 2A-1A
- **Operating temperature:** Between 5°C and 40°C
- **Operating humidity:** Up to 93%R.H. non condensing
- **Operating Pressures:** 700hPa to 1060hPa
- **Handset Treatment area (spot size):** 3 cm² (30mm x 10mm)
- **The Homogeneity of the treatment area:** Does not vary by more than +/- 20%
- **Pulse Profile:** Decaying pulse waveform with range of pulse lengths (2.5-5.5ms)

The area directly around the treatment window is classified as the Applied Part, and can reach a maximum temperature of 60°C.

Under continued use the temperature of the device may reach a maximum of 60 °C. After treatment, remove the device from the skin.

Ensure SmoothSkin Pure Ice is used as directed in Section 7.

TRANSPORT CONDITIONS:

- **Temperature:** -25°C to +70°C
- **Humidity:** Up to 93%R.H. non condensing
- **Pressures:** 500hPa to 1060hPa

If the device has been stored at temperatures outside of the stated operating range, then please allow at least 30 minutes before use to allow the device to acclimatise.

SmoothSkin Pure Ice has unlimited flashes and is intended for a single user.

Minimum expected service life is 10 years when following the recommended regimen.

This device has no clinical benefit.

CE marking certifies that this appliance conforms to the following EU Directives:

- Low Voltage Directive 2014/35/EU.
- Electromagnetic Compatibility Directive 2014/30/EU.
- Energy related Product Directive 2009 / 125 / EC.

Smoothskin Pure Ice complies with the latest EMC Standards and Amendments applicable to Household and Medical Devices.

These are: EN 55014-1:2021, EN 55014-2:2021, EN 61000-3-2:2019+A1:2021, EN 61000-3-3:2013+AMD1:2019+AMD2:2021, IEC 60601-1-2:2014 (Ed 4.0), FCC CFR 47 Part 15.107 & 15.109, ICES-003 Issue 7.

SmoothSkin Pure Ice complies with the latest Safety Standards and Amendments applicable to Household and Medical Devices.

These are: IEC 60335-1:2010 + Ams, EN 62233, IEC 60335-2-113:2016 IEC 60601-1:2005 (Ed 3.1) +A2:2012, IEC 60601-2-57:2011, IEC 60601-2-83:2019, IEC 60601-1-6:2010 + Amd1:2013, IEC 60601-1-11:2015, IEC 62366-1:2015, IEC 62304:2006 + Ams, CAN/CSA - C22.2 No 60601-1:2008 and 2014, ANSI/AAMI ES 60601-1:2005 + Ams.


SmoothSkin Pure Ice is classified as Risk Group Exempt for all Hazards.

For further details contact CyDen Ltd.

DISPOSAL (END OF LIFE)

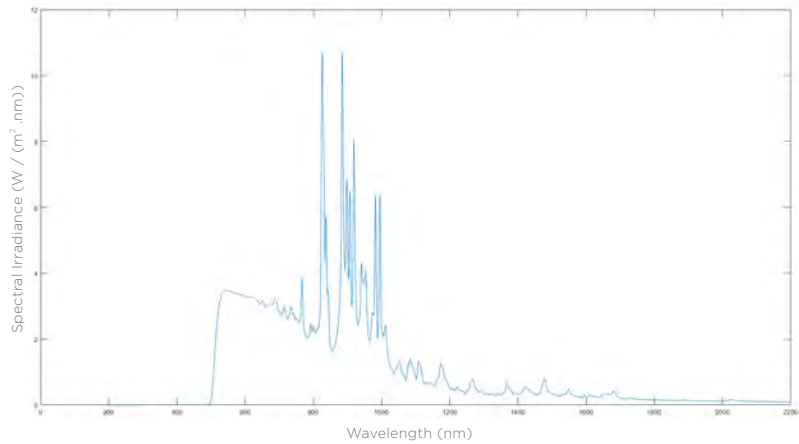In order to minimize hazards to health and the environment and ensure that materials can be recycled, this product should be disposed of at a separate collection facility for waste electrical and electronic equipment. The WEEE symbol marked on the product is there to remind you.

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Para mais informações, contacte a CyDen Ltd.
CyDen Limited, Block A, Bay Studios Business Park, Fabian Way, Swansea SA1 8QB UK
Cyden Ltd. Factory,