

VANITY PLANET

Forever Young


INSTANT HOT & COLD
REJUVENATING WAND

USER GUIDE

Features

- Boosts absorption of serums, creams & moisturizers
- Helps reduce puffiness & redness
- Tightens skin's surface to lessen the appearance of fine lines
- Sonic vibrations enhance both hot & cold mode
- Teardrop-shaped head perfect for under-eye areas

What's the point of facial care if it isn't absorbed by your skin? There isn't one. Which is why Forever Young uses hot and cold settings along with sonic vibrations to make the most of your skincare routine. First, hot mode instantly heats to 107°F (42°C) to help open up your pores and boost the absorption of your favorite serums, creams and moisturizers. Then, change to cold mode for an immediate 43°F (6°C) to shrink pores, lock in serums/moisturizers and lessen the appearance of fine lines. Even better, it helps reduce red spots and puffiness all while soothing your skin. Forever Young features sonic vibrations to enhance either hot or cold mode while rejuvenating your complexion (and giving you a therapeutic massage). So, lock in all that good moisture, tackle those under-eye bags and bring out the most beautiful version of you.

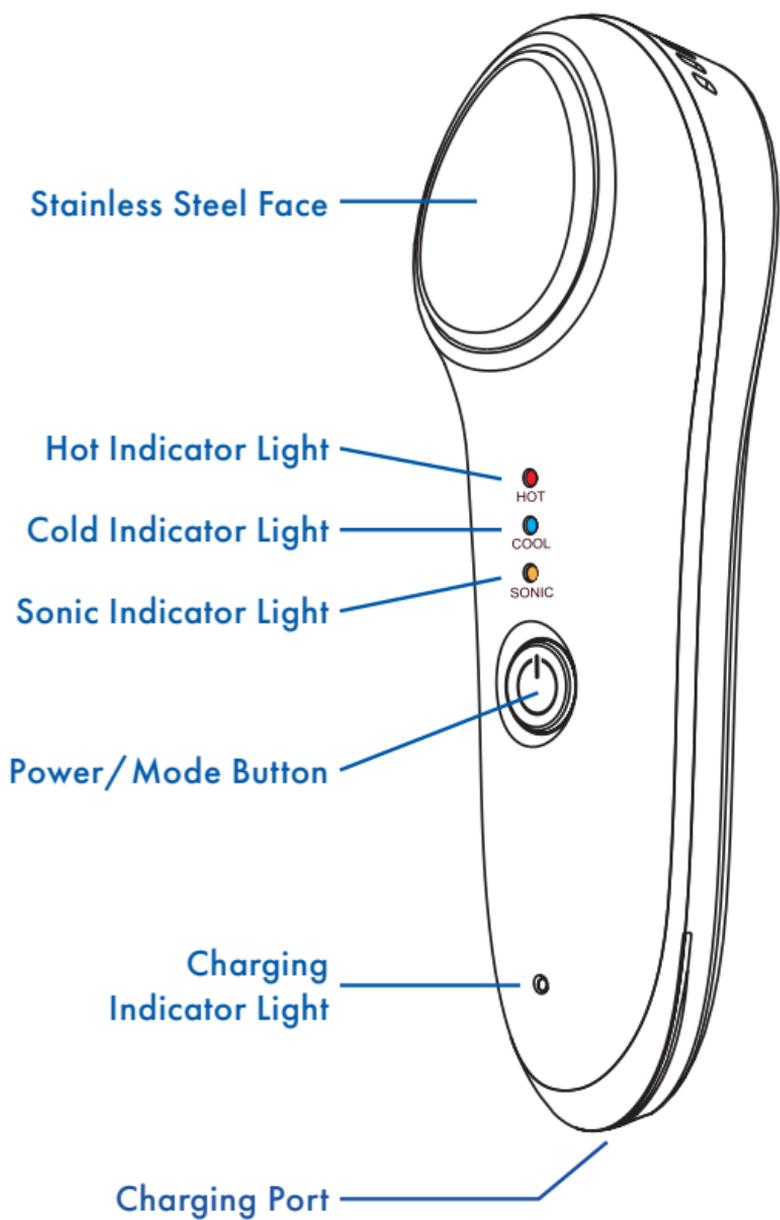

Charging

Your Forever Young takes three hours to fully charge. After it has been fully charged, you can use it for up to one hour before needing to recharge the device. To charge your Forever Young, simply plug the USB cable into the USB port at the bottom of the handle and connect it to a power source. While charging, the indicator light will be red. Once it is fully charged, the light will turn blue.

Using Forever

TURNING YOUR DEVICE ON/OFF & SWITCHING MODES

1. Hold the "On/Off" button for three seconds to activate the hot mode
2. To add sonic vibrations while using the hot mode, simply press the "On/Off" button a second time
3. To shift to the cold mode, press the "On/Off" button a third time
4. Press the "On/Off" button a fourth time to add sonic vibrations while using the cold mode
5. To turn off your Forever Young, hold down the "On/Off" button for three seconds. The device will automatically shut off after two and a half minutes on each mode

Young

INSTRUCTIONS

1. After cleansing your face, apply your favorite serum or moisturizer. Then, use your Forever Young on hot mode (with or without vibrations) to open up your pores and help your skin absorb more of the skincare product.
2. Once your product has been fully absorbed, switch to the cold mode (with or without vibrations) to help shrink your pores and tighten your skin. Now your face is ready and prepped for sleep or for a full face of makeup.
3. If you want to combat puffiness, redness and dark circles, simply use your Forever Young on cold mode (with or without vibrations) to soothe problem areas.
4. Enhance either hot or cold mode with sonic vibrations to make the most of your skincare routine.



say hello
to true
beauty.

Taking Care

To keep your Vanity Planet Forever Young in optimum condition, avoid using abrasive cleaners or chemicals to clean any part of this device. Only clean when the device is turned off and cooled to room temperature.

CLEANING THE HANDLE

To remove any residue or buildup, clean the handle at least once a week by taking a slightly damp cloth and wiping down the body of the device when it is turned off. Do not use soap.

CLEANING THE STAINLESS STEEL FACE

Wipe down the stainless steel face with a slightly damp cloth while the device is turned off and cooled. Make sure to carefully clean around the edges of the face to remove any skincare product residue.

Safety

To ensure you experience all the benefits of your Vanity Planet Forever Young, we recommend the following:

- Do not submerge the handle in water as it is not water resistant
- Do not use with harsh exfoliating cleansers or physical exfoliants such as nuts, shells, seeds or other jagged exfoliating particles
- Make sure that the stainless steel face is properly cleaned between uses to avoid product build up and residue
- Do not use if you are pregnant
- Keep out of reach of children

NOTE: If you are under a doctor's care for any skin condition or if you think you may have a skin condition, please consult your physician before use.

One-Year Warranty

Visit vanityplanet.com/warranty
to register your product and
activate your warranty.

This Vanity Planet product is warranted to be free of manufacturer's defects in materials or workmanship for one year from the date of purchase. Damage or wear resulting from an accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair are not covered under this warranty.

Should this product require warranty service (or replacement at our discretion), please contact client services. Proof of purchase is required. Products returned without prior authorization from Vanity Planet will not be accepted and the sender will not receive a refund or replacement.

vanityplanet.com
info@vanityplanet.com
1.888.870.2132

There are no express warranties except as listed above. This warranty gives you specific legal rights and you may have other rights which vary from state to state.

Please do not return this product to the retailer. Contact us directly for warranty service.

VID