**EFFECTIVE
LIP & EYE**

N Ū F A C E[®]

USER MANUAL

ENSURE YOUR BEST RESULTS!

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NuFACE® EFFECTIVE LIP & EYE TRINITY+ ATTACHMENT

Note: Refer to NuFACE Trinity+ User Manual for more device information.
Remarque: Reportez-vous au manuel d'utilisation de NuFACE Trinity+ pour plus d'informations sur l'appareil.
Nota: consultar con el Manual del Usuario de NuFACE Trinity+ para más información sobre el dispositivo.
Hinweis: Weitere Geräteinformationen finden Sie im NuFACE Trinity+ Benutzerhandbuch.
참고: 자세한 기기 정보는 NuFACE 트리니티+ 사용 설명서를 참조하십시오.

GETTING TO KNOW YOUR ATTACHMENT SE FAMILIARISER AVEC L'APPAREIL CONOCIENDO SU ACCESORIO IHR GERÄT / 기기 알아보기

1) Precise Wands (Electrodes)

Baguettes précises (électrodes)
Barras de precisión (Electrodos)
Behandlungsstäbe (Elektroden)
정확한 완드(전극)

2) Magnetic Effective Lip & Eye Attachment

Accessoire magnétique efficace pour lèvres et yeux
Accesorio magnético efectivo para ojos y labios
Effektiver Lippen- und Augen-Magnetaufsatz
마그네틱 이펙티브 립 앤 아이 부가장치

Magnetic Base (not shown)

Base magnétique (non illustrée)
Base magnética (no se muestra)
Magnetfuß (nicht abgebildet)
마그네틱 베이스(표시되지 않음)

3) Microcurrent Intensity Booster

Amplificateur d'intensité de microcourant
Elevador de la intensidad de la microcorriente
Mikrostrom-Intensitätsbooster
미세전류 강도 부스터

4) Dual Function Power Button ON/OFF + Intensity Control Up/Down

Bouton d'alimentation à double fonction: March/Arrêt + Contrôle d'intensité haut/bas
Botón de encendido doble función: ENCENDIDO/APAGADO + Control de intensidad Arriba/Abajo
Power-Knopf mit doppelter Funktion: Ein-/Ausschalter und Intensitätsregelung auf/ab
전원 온/오프 버튼

INTENDED USE

The NuFACE Trinity + Effective Lip & Eye attachment is intended to be used in conjunction with the NuFACE Trinity+ device for facial stimulation and is indicated for over-the-counter cosmetic use.

IMPORTANT INFORMATION

- Use the Trinity+ Effective Lip & Eye attachment as described within this User Manual.
- The NuFACE Trinity+ device User Manual may also be referenced.
- The risks and dangers of using the attachment in any way that does not follow this User Manual are unknown.
- Your Effective Lip & Eye Attachment has a unique serial number located on the underside of the attachment.
- Register your new device at www.NuFACE.com.

CONTRAINDICATIONS FOR USE

- **DO NOT** use this device if you have a history of epilepsy or seizures.
- **DO NOT** use if you have had recent surgery or other surgical procedures.
- **DO NOT** use if you have any active cancer or over any suspicious or cancerous lesion.
- **DO NOT** use if you are pregnant. The safety of electrical stimulation during pregnancy has not been established.
- **DO NOT** use if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor.
- If you use a body-worn medical device, such as an insulin pump, consult your healthcare provider.
- **DO NOT** use on infants and children under the age of 18.
- Consult your healthcare provider if any of the contraindications above apply to you.

⚠️ WARNINGS

- If you are under the care of a physician, consult your physician before using this device.
- This device is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Apply stimulation only to normal, intact, clean healthy skin.
- **DO NOT** use if areas of your skin are swollen, infected, or inflamed areas or skin eruptions (e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc).
- The product may not be safe for use on the neck in individuals with a known history of carotid artery disease, stroke, or transient ischemic attack (TIA), carotid stenosis, unstable blood pressure, easy fainting, thyroid disease, or prior surgical procedures in the neck area.
- **DO NOT** operate the Effective Lip & Eye Attachment near electronic monitoring equipment (e.g., heart monitors, heart tracing) or in close proximity (e.g., 1 m / 3ft) to shortwave, microwave, or electromagnetic therapy equipment as it may cause the medical equipment to malfunction. Such equipment may produce instability in the device output and cause burns at the stimulation site.
- **DO NOT** disassemble the device as this may cause damage, malfunction, electrical shock, fire, or injury.
- Keep the device and related electrical accessories away from direct contact with water (including baths, showers, and sinks).
- **DO NOT** use across your chest as the introduction of electrical current into the chest may cause rhythmic disturbances to your heart, which could be lethal.

⚠️ WARNINGS (CONTINUED)

- Keep the device out of the reach of children and pets.
- **DO NOT** drop or insert any object(s) into openings within the device.

Consult your healthcare provider before using the device:

- If you have any pre-existing or potential health conditions that may be impacted by use of the device.
- If you have painful areas.
- If you suspect or have a neuromuscular disorder.
- If you have a spinal cord injury.
- If you have any underlying skin conditions.
- If you have undiagnosed or known pain syndromes.
- If you have a tendency to excessively bleed following an injury.
- If you experience constant headaches and/or painful sensations.
- If you have suspected or diagnosed heart disease.

DO NOT apply stimulation/use the device:

- When there is a tendency to excessively bleed following acute injury or breaks in bones.
- When the device is on the charging cradle.
- On the eye area (circular muscle within the orbital rim).
- On the temples, the breast or groin area.
- Over painful areas of your skin or in areas that lack normal sensation.
- While bathing or showering.
- At extreme temperatures.
- If you suspect the device is damaged or not functioning properly.
- While driving or operating machinery.

Discontinue use of the device and contact your healthcare provider if you experience any of the following:

- Persistent pain or discomfort.
- Discomfort or any difficulty breathing or swallowing.
- Negative or harmful reaction(s) or injury.
- Symptoms of hypertension (dizziness or lightheadedness).
- Constant headaches and/or painful sensations.
- Skin irritation or burning.
- Unpleasant sensations on skin during use.
- Consistent flashing lights when without the use of your device.

If you experience any negative reactions while using the device and attachments, stop use immediately and contact your healthcare provider.

PRECAUTIONS

- The Effective Lip & Eye attachment is designed for cosmetic use only.
- The device and attachments should be operated, transported, and stored in the manner specified in this manual. Failure to use and maintain the device and corresponding attachments in accordance with the instructions outlined in this manual may result in product damage or physical injury.

PRECAUTIONS (CONTINUED)

- The long-term effects of transcutaneous electrical stimulators are unknown.
- Apply stimulation only to normal, intact, clean, healthy skin.
- Avoid applying excessive pressure when using the device.
- Skin irritation or hypersensitivity may occur due to the electrical stimulation and / or use of the NuFACE® Microcurrent Activator.
- A microcurrent sensation may occur while using the device. Using a lower intensity and / or additional activator may be needed to reduce the sensation.
- A perception of flashing lights may occur if device is used near or around the eye region.
- **DO NOT** use the device if you see any signs of damage or leaks from the device.
- **DO NOT** remove the battery of the device.
- **DO NOT** store the Trinity+ device in the sun or on a hot surface.
- **DO NOT** place wands from the attachment directly on or into the nose.
- **DO NOT** place wands from the attachment directly on or into the tear duct or eyelid.
- Activator must be reapplied if skin becomes dry during treatment.
- Failure to use and maintain the device and attachments in accordance with the instructions outlined in this manual will void the warranty.
- The use of accessories, transducers, and cables other than those sold by NuFACE may result in increased emissions, decreased life of the device, or injury.
- Children should be supervised to ensure that they **DO NOT** play with the device.

The NuFACE® Trinity+ Effective Lip & Eye attachment is intended for cosmetic use on the face as described. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty wands, expired or non-NuFACE conductive alternative, or any other improper application is neither the responsibility of NuFACE nor its affiliates.

The NuFACE® Trinity+ Effective Lip & Eye attachment has been designed to provide freedom from unacceptable risk when used under normal conditions and reasonably foreseeable hazards. In the event of complete loss of treatment function, the Trinity+ Effective Lip & Eye attachment does not pose a risk to the User.

NUFACE MICROCURRENT ACTIVATOR

The essential first step to your microcurrent treatment to ensure optimal results.

- **Always apply Microcurrent Activator when using the Trinity+ device EFFECTIVE LIP & EYE ATTACHMENT.**
- Refer to packaging for use by date.

TREATMENT AREAS AND TECHNIQUE

The instructions below are designed to provide users with the most effective results when followed correctly.

TREATMENT

The Trinity+ Effective Lip & Eye attachment is designed to stimulate smaller surface areas of the face, such as the lip and eye area.

TECHNIQUES

Using light to medium pressure, glide or hold both wands over each treatment area position.

IMPORTANT: Always ensure that both wands are in contact with the skin during treatment. Always ensure that you are correctly holding the device in either a vertical or horizontal position, as required per the specified treatment area.

- Refer to the Quick Start section in this manual for proper positioning of the wands and techniques for conducting treatment.
- For more advanced treatment options and holds, download the NuFACE App, or visit our website at www.NuFACE.com.

BEFORE USING YOUR NEW ATTACHMENT

Follow these important steps before using your new attachment.

Follow steps 1-5 of the NuFACE Trinity+ Device (sold separately) Manual "**BEFORE USING YOUR NEW DEVICE**" section.

1) **READ** the NuFACE Trinity+ Device User Manual (sold separately) Reference the section in the Trinity+ Device User Manual titled "BEFORE USING YOUR NEW DEVICE" before changing attachments.

2) **CONNECT YOUR TRINITY+ EFFECTIVE LIP & EYE ATTACHMENT**

The NuFACE Trinity+ device is designed with magnetic attachments.

- Remove the Facial Trainer attachment.
- Connect the Effective Lip & Eye attachment onto the NuFACE Trinity+ device.

The Trinity+ device will NOT turn on if the attachment is not securely in place.

3) **TEST YOUR SKIN SENSITIVITY**

IMPORTANT: Test your skin to make sure you are not sensitive to the NuFACE Microcurrent Activator.

- Rub a small amount of Microcurrent Activator on your inner forearm and leave for **15-30 minutes** before washing it off.
- If an allergic reaction (rash, redness, etc.) develops, immediately contact NuFACE for recommendations on alternative NuFACE formulas. **DO NOT continue application** of the Microcurrent Activator with the Trinity+ device until a resolution such as an alternate Microcurrent Activator has been identified.

HOW TO USE YOUR TRINITY+ EFFECTIVE LIP & EYE ATTACHMENT

Get ready to Tone, Lift, and Contour.

The instructions below are designed to provide the most effective results when followed correctly.

1) CLEANSE YOUR FACE

Cleanse your face and neck with an oil-free cleanser. Men should shave before treatment — facial hair can interfere with the conductivity of the wands.

Note: Any cream, lotion, self-tanner, and/or oil-based products on the skin can block the microcurrent.

2) APPLY ACTIVATOR

IMPORTANT: Always apply a NuFACE Microcurrent Activator when using your Trinity+ device.

The instructions below are designed to provide the most effective results when followed correctly.

- Apply a mask-like layer of the NuFACE Microcurrent Activator (sold separately) section-by-section as you treat.
- Reapply activator if skin becomes dry during treatment in order to continue safely using device.

3) PREPARE YOUR DEVICE AND ATTACHMENT

- Remove your device from the charging cradle.
- Ensure your Effective Lip & Eye attachment is securely in place and turn it ON by pressing the ON/OFF button.
- An ascending audible beep is emitted, the energizing hum will begin, and the lights will illuminate indicating that your device is ON.
- Adjust the intensity level to your comfort by pressing the ON/OFF button on the front of your device.
- As you increase or decrease intensity, an audible beep is emitted, and the number of illuminated lights will change.

INTENSITY LEVEL	INDICATOR LIGHT	
Minimum Intensity	First Light Illuminated	●○○
Medium Intensity	First & Second Light Illuminated	●●○
Maximum Intensity	All (three) 3 Lights Illuminated	●●●

DEVICE AUDIBLE SOUND	DEVICE INDICATION
Ascending Beeps & Hums	Powered ON
Descending Beeps	Powered OFF
Single Beep	Intensity Adjusted
Dual-Tone Beep	End of Glide Interval (Tone occurs only during treatment)
Low-Tone Beep	Battery Needs To Be Charged (Tone occurs only during treatment)

HOW TO USE YOUR TRINITY+ EFFECTIVE LIP & EYE ATTACHMENT (CONTINUED)

4) BEGIN TREATMENT

You can also refer to the Quick Start section for basic glides treatment for the lip, eye, and brow area.

- Place the Effective Lip & Eye Attachment in direct contact with the treatment area.
- **Always keep both wands in contact with the skin during treatment to keep microcurrent flowing.**
- Using light to medium pressure, begin treatment.
- The device emits a dual tone beep **every seven (7) seconds** to indicate when to move to the next treatment area.

IMPORTANT: If you feel the electrical current and / or experience a tingling sensation on your skin, do the following:

- Add more NuFACE Microcurrent Activator to the area you are treating.
The device will NOT contract the muscles when used as instructed.
- Lower the intensity on the device.
Lowering the intensity may also reduce or eliminate this sensation.


QUICK START

Following steps 1–4 under 'HOW TO USE YOUR TRINITY+ EFFECTIVE LIP & EYE ATTACHMENT' section:


- Glide or hold both wands over each treatment area position as shown for **one (1) beep seven (7) seconds**.
- Complete each sequence **two (2) times** beginning with the Lip Area on the right side of your face, then **three (3) times** on the Eye and Brow area.
- Then repeat the order by completing the same sequence on the left side of your face.

BEGIN: RIGHT SIDE OF FACE		
ORDER	TREATMENT AREA	SEQUENCE
A	Lip	x2
B	Lip	x2
C	Lip	x2
D	Eye	x3
E	Eye	x3
F	Brow	x3
G	Brow	x3


FINISH: LEFT SIDE OF FACE		
ORDER	TREATMENT AREA	SEQUENCE
A	Lip	x2
B	Lip	x2
C	Lip	x2
D	Eye	x3
E	Eye	x3
F	Brow	x3
G	Brow	x3

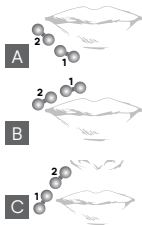Lip Area

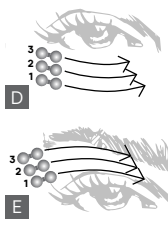Eye Area

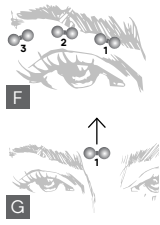Brow Area

HOW TO USE YOUR TRINITY+ EFFECTIVE LIP & EYE ATTACHMENT (CONTINUED)

5) FINISH TREATMENT

- After treatment is complete, massage remaining NuFACE Microcurrent Activator into skin.
- Turn **OFF** the Trinity+ device after treatment.
- A descending audible beep is emitted notifying you that your device is turning OFF.

The Trinity+ device will automatically shut off after 20 minutes of use.

Go to the section "CARING FOR YOUR EFFECTIVE LIP & EYE ATTACHMENT" for instructions on cleaning and maintaining your device.

If you need to clean your device, go to section titled "CLEANING" located in this manual.

TREATMENT FREQUENCY

- Use on targeted areas (3 minutes for the lips, 5 minutes for the eyes).
- For the first 60 days, use 5 days per week.
- After 60 days, use 2-3 times per week for maintenance.
- Allow 24 hours between treatment sessions.

For more advanced treatment options and holds, download the NuFACE App, or visit our website at www.NuFACE.com.

RESULTS

Results may vary based on skin type and skin condition. Best results are achieved with regular use.

NUFACE APP

Download the NuFACE APP for additional how-to instructions, and customized treatment options.

- You can download it from the App store or the Google Play Store.
- To pair your device, in the App, click on 'Add Device' and follow the on-screen instructions to connect. The NuFACE app is available in many countries around the world.

CARING FOR YOUR EFFECTIVE LIP & EYE ATTACHMENT

Cleaning:

- Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device and attachments after every treatment.
- **DO NOT submerge your device or any of the attachments in water at any time.**

Storage:

- The Trinity+ device and accessories are intended to be stored and used at room temperature. Use caution prior to using the device if it has been stored or not used for an extended period of time.

Disposal:

IMPORTANT: If disposing the accessories, follow local regulations for proper disposal.

TROUBLESHOOTING

Device Not Turning On:

1) Check to see if your Trinity+ device is fully charged.

- If the device does not turn ON, the battery may be drained.
- A full charge cycle of **ten (10) hours** may be required for the device.
- Reference the 'NuFACE Trinity+ Device User Manual' for additional troubleshooting.

2) Check to see if the attachment is completely secured to the Trinity+ device.

IMPORTANT: The Trinity+ device will NOT turn on if the attachment is not securely in place.

DO NOT attempt to troubleshoot by disassembling the product as this will void the warranty and may cause injury.

For additional troubleshooting or information specific to the device, reference the NuFACE® Trinity+ Device User Manual.

CUSTOMER SUPPORT

If you are still experiencing issues with your device, contact NuFACE Client Services at info@nuface.com or call **888-360-9171** (US Only) between 8am PST to 5pm PST.

WARRANTY

NuFACE® warrants that all NuFACE® Trinity+ devices and attachments are free from manufacture defects in materials and workmanship for a period of two (2) years from the date of your original purchase. Our obligation under this Warranty is limited to replacement of the product which will not extend the period of your warranty. This Warranty does not cover damage covered by abuse, improper care, application of harmful chemicals, product that may have been altered in any way, loss of the device and product purchased from an unauthorized seller.

NO OTHER WARRANTY EXPRESSED OR IMPLIED, WHETHER OF FITNESS FOR A PARTICULAR PURPOSE OR OF MERCHANTABILITY OR OF ANY OTHER KIND, SHALL EXIST IN RESPECT TO SUCH PRODUCT.

For Devices Purchased in the US: For more information about your warranty or to initiate a warranty claim and receive a return authorization number, contact our Client Services Team via our website at www.NuFACE.com or by calling 888-360-9171. You must show a receipt from an authorized seller as proof of purchase.

For Devices Purchased Outside of the US, International Returns, and International Warranties:

For all product purchases and returns outside of the United States, warranty and returns are per your place of purchase: spa, retailer, or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

Conforms to EN IEC/UL: Device - 60529, 60335-1, 60335-2-23, 60601-1, 60601-1-11, 60601-1-2, 60601-1-6, 60601-2-10, 62133-2, 62366-1

SYMBOLS

One or more of the following markings may appear on the Trinity+ Device, User Manual or packaging.

	European Conformity		Comply with European Restriction of Hazardous Substances Directive 2002/95/EC
	European Representative		Product must be disposed of in accordance with the WEEE directive
	Part Number		FCC-approved equipment authorization mark
	Operator's Manual; Operating Instructions		Class II double insulated electrical device
	Refer to Instruction Manual/Booklet		
	Manufacturer		
	Detachable Power Supply Part		
	Power On/Off		
	Prevention of access to hazardous parts with fingers and protects the equipment inside the enclosure against ingress of solid foreign objects having a diameter of 12.5 mm and greater. Protected against vertically falling water drops when enclosure tilted up to 15°.		
	Type BF applied part, which is equipment that provides a particular degree of protection against electric shock. Particularly regarding allowable leakage current and of the protective earth connection.		

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All Rights Reserved. / Tous droits réservés. / Alle Rechte vorbehalten. /
NuFACE® and Trinity® are registered trademarks of the Carol Cole Company DBA NuFACE.
Medical Device Safety Service GmbH