


Advanced Light Technology
Red, Infrared & Blue Light Therapy

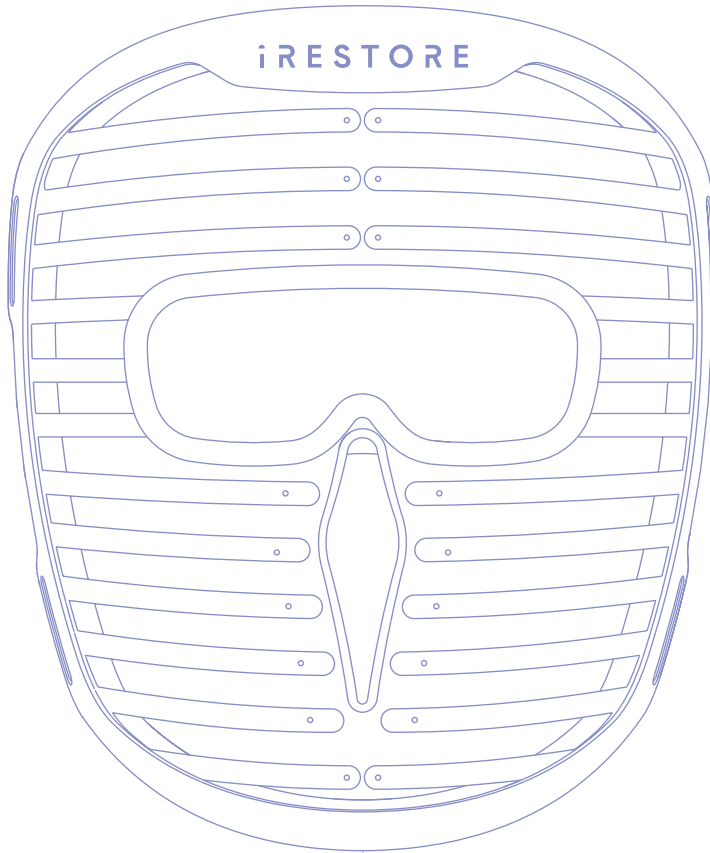

iRESTORE

User Manual

ILLUMINA FACE MASK

iRESTORE

We're on a
different
wavelength.

2

iRESTORE



iRESTORE



Start here!

Scan the QR code for step-by-step videos on setting up and using your Illumina Face Mask. Discover how to seamlessly integrate this into your routine for best results. Let's begin our journey together!



- ◆ Understand the three different modes of your device
- ◆ Best practices for preparing for a device treatment
- ◆ Instructions on how to clean and maintain your device
- ◆ Troubleshooting tips to address any issues that may arise during the use of your device

Table of Contents

6	Product Description
6	Warnings and Precautions
9	Product Diagram
10	How to Use
12	Product Specifications
13	Warranty

PRODUCT DESCRIPTION

The Illumina Face Mask is a low-level light therapy device designed to improve the appearance of skin and help skin look younger and healthier. The Illumina Face Mask utilizes 360 light-emitting diodes (LEDs) with wavelengths of 415nm, 635nm and 830nm.

WARNINGS AND PRECAUTIONS

Do not use the Illumina Face Mask if you have a medical condition that makes your skin sensitive to sun, light or if you are taking photosensitizing medications.



Consult with a doctor before using if you have any skin diseases, cancers or tumors on the face.



Do not use the Illumina Face Mask if you have sensitivity to light or if you are allergic to this specific wavelength of light (415nm, 635nm and 830nm). Skin sensitivity or an allergic reaction can result in an itchy sensation or sunburn-like redness. If these occur, seek medical care from a physician.




The Illumina Face Mask should not be used by or on children under 18 years of age.

Discontinue the use of the Illumina Face Mask if the itching or tingling occurs and lasts more than one (1) hour after completion of a session. If the problem continues for more than one (1) hour after use, discontinue usage and consult with a physician.



Do not pull on the cord while handling, untangling or straightening it. The pressure may break the cord, which could result in a malfunction or electrical hazard.

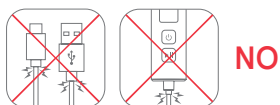


Do not get the Illumina Face Mask or charging cord wet. This may result in an electric shock.



Do not use the mask while driving or operating heavy machinery.

Do not insert any objects into any opening of the mask.


Do not use the device if it is overheating or you suspect that it is malfunctioning.

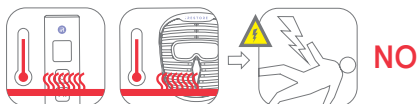


Discontinue use if this device or charging cord is not working properly or appears damaged in any way.

WARNINGS AND PRECAUTIONS

Before using the Illumina Face Mask, ensure your face is towel-dried, and do not use it when your face is wet from showering, swimming or bathing. Only use the RENEW Calming Facial Spray prior to mask sessions.

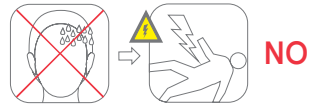


Use the USB-C Charging Cable provided with the Illumina Face Mask or other USB-C Charging Cable that complies with relevant safety regulations.

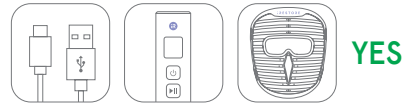

Keep the device out of reach of children and pets.

Always unplug the Illumina Face Mask from the power source prior to cleaning the device. Failure to do so may result in an electrical shock that may cause harm or injury.

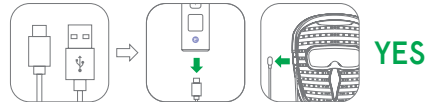

It is important to only use the Illumina Face Mask according to the instructions provided. Using it for any other purpose is not recommended and will void the warranty.

Never attempt to repair the Illumina Face Mask on your own. Please contact your distributor or Freedom Laser Therapy to request a repair or replacement.

Avoid dropping the Illumina Face Mask, as this may cause the device to malfunction and become damaged. If you drop the device or suspect device malfunction (e.g., device exterior cracks, LEDs do not light up or the device does not turn on), contact Freedom Laser Therapy to request a repair or a replacement.



This Illumina Face Mask is equipped with a lithium-ion battery and it is recommended to always carry it in your carry-on luggage while traveling. When it is time to dispose of your device, please take it to specific recycling or household hazardous waste collection points due to the nature of the battery.

DISCLAIMER

This device is not intended to diagnose, cure or prevent any disease. Freedom Laser Therapy makes no claims, representations or warranties regarding the ability of this product to cure any physical, skin or mental condition. A qualified health professional should be consulted with regard to any condition requiring medical attention.


PACKAGE INCLUDES (ACCESSORIES LIST):

Accessories	Quantity
The Illumina Face Mask	1
Controller with Battery	1
USB-C Charging Cable	1
Adjustable Straps	2
Storage Pouch	1
Microfiber Cloth	1


The Illumina Face Mask components are shown below:



FRONT

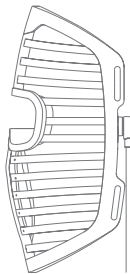

SIDE

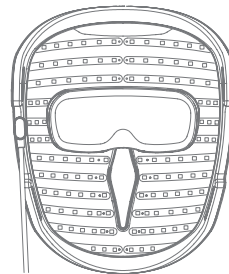


INSIDE

OPERATION MODES

Mode	Face Mask LED Light	Indicator Light
P1	Red (635nm) and infrared (830nm) turn on	Red light
P2	Red (635nm), infrared (830nm) and blue (415nm) turn on	Purple light
P3	Blue (415nm) turns on	Blue light

PRODUCT DIAGRAM

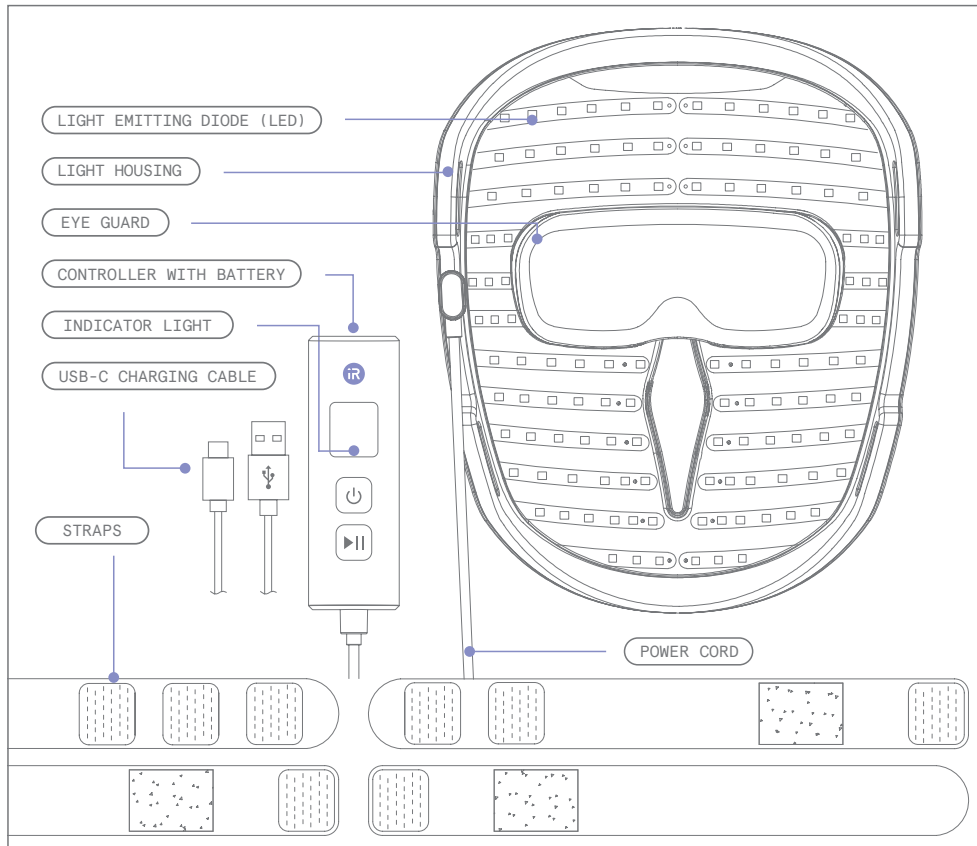


	Description	Function
1	Light Housing	Configured with 360 LEDs to provide a low-level light source to the face
2	Eye Guard	Padding around the light housing to prevent light from shining into eyes and provide better comfort
3	Light Emitting Diode (LED)	360 LEDs (415nm, 635nm and 830nm) are the light source
4	Controller with Battery	External system control to operate On/Off, indicator light and Start/Pause
5	Power Cord	Connect mask and controller
6	Straps	Adjustable elastic bands to customize your fit*
7	USB-C Charging Cable	Use a USB power adapter (adapter not included) to charge the mask

* See labels on straps for installation: upper left (UL), upper right (UR), lower left (LL) and lower right (LR).

HOW TO USE THE ILLUMINA FACE MASK

For best results, please use the Illumina Face Mask for 10 minutes, 3-5 times in a week.

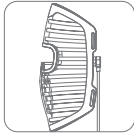


Set Up

Connect the power cord to the mask. Attach the two parts of the upper strap to the mask (labeled UL for Upper Left and UR for Upper Right). If desired, attach the lower straps as well (labeled LL for Lower Left and LR for Lower Right). The lower straps are optional and are designed to allow the user to pull the mask closer to the face.



Preparation

Cleanse your face with RENEW Hydrating Foam Cleanser to remove all makeup, skincare, dirt, and oils. Follow with RENEW Calming Facial Spray. Place the mask comfortably on your face and lightly tighten the straps. Ensure the mask is snug but not too tight. Adjust the upper straps upward to relieve pressure.

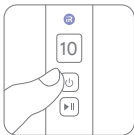

Step 1

Press and hold ⏻ (Power) button to turn on the mask. The mask enters P1 mode automatically. Press the ⏻ (Power) button again to enter P2 mode. Press the ⏻ (Power) button a third time to enter P3 mode.




Step 2

Once a mode is selected, press ▶|| (Start/Pause) button to begin a 10-minute session, and the indicator light will turn on. Once the session is complete, the indicator light will turn off, and the mask will turn off automatically in 30 seconds. Note: If the ▶|| (Start/Pause) button is not pressed within 30 seconds once a mode is selected, the mask will turn off automatically.



Step 3

If you need to pause the session at any time, press the ▶|| (Start/Pause) button, and the indicator light will flash. Press the ▶|| (Start/Pause) button to resume the session, and the indicator light will turn on. If the ▶|| (Start/Pause) button is not pressed within 30 seconds during a paused session, the mask will turn off automatically.

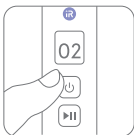

Step 4

Hold the ⏻ (Power) button for several seconds if you need to turn off the mask. The indicator light will turn off immediately; then, the mask will turn off in 30 seconds.

IMPORTANT: If you want to change mode in a session, press the ▶|| (Start/Pause) button to pause a session, then follow steps 1-4, listed above.

HOW TO RECHARGE THE ILLUMINA FACE MASK

- Plug the charging cord into the controller and connect its USB to a USB power adapter (adapter not included).
- Plug your adapter into an AC outlet. The display will turn on, showing power status indicating battery percentage (25, 50 and 75).
- The display will turn off when the device is fully charged.
- Unplug the USB power adapter from the AC outlet and remove the charging cord.

IMPORTANT:

- If the indicator light starts to flash during a session and the battery power is low, please recharge the device before next use.
- The buttons will not work when charging to prevent electrical hazards.
- If the red indicator light is quickly flashing while charging and the built-in battery's temperature is abnormal, the device will stop charging automatically. When the battery temperature returns to normal, the device will resume charging.

CLEANING INSTRUCTIONS

- The Illumina Face Mask should be unplugged prior to any cleaning.
- Use included microfiber cleaning cloth to clean the mask.
- For stains that cannot be removed with a dry cloth, you may use a lightly dampened cloth with water to clean the mask.
- Do not use any harsh solvents to clean the mask; please only use the iRESTORE Device Cleaning Solution (not included).
- Use iRESTORE's alcohol-based Device Cleaning Solution to help kill germs and bacteria that could potentially be harmful if left on the device during use, ensuring a hygienic experience every time.

STORAGE INSTRUCTIONS

- After each use, store the mask in a dry location at room temperature, away from sunlight, windows, sinks, bathtubs, kitchens and bathrooms.


TROUBLESHOOTING

- My Illumina Face Mask is not powering on.
The battery power is too low for the mask to work; please charge it before use.
- My Illumina Face Mask is not lighting up, but my controller is on.
Contact Freedom Laser Therapy to discuss the possibility of receiving a replacement mask.

PRODUCT SPECIFICATIONS

Item Number.....	ILLUMINA FACE MASK ID-610
Color.....	WHITE
Wavelength of LED.....	415NM±10NM; 635NM±10NM; 830NM±10NM
Net Weight.....	0.40 LBS
Size.....	8.54 X 7.17 X 3.87 IN.
Number of LED	360 (3 PER LED BULB) 120 LED BULBS
Maximum Operation Time.....	10 MINUTES
Battery.....	18650/3350MAH
Charging Power.....	USB DC 5V, 2A
Charging Current.....	1.7A
Charging Time.....	2 HOURS
Operating Environment Temperature.....	32°F-104°F (0°C-40°C)
Operating Environment Humidity.....	15% TO 90% RH
Operating Environment Pressure.....	700 HPA TO 1060 HPA
Storage & Transportation Temperature.....	-4°F TO 113°F (-20°C TO 45°C)
Storage & Transportation Humidity.....	10% TO 90% RH

Labels on device

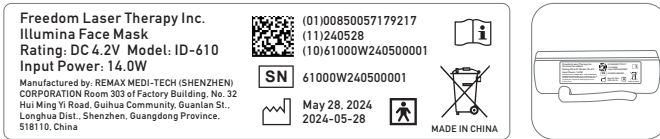

MANUFACTURED BY:
 REMAX MEDI-TECH (SHENZHEN) CORPORATION
 ROOM 303 OF FACTORY BUILDING, NO. 32,
 HUI MING YI ROAD, GUIHUA COMMUNITY,
 GUANLAN STREET, LONGHUA DISTRICT, SHENZHEN,
 GUANGDONG PROVINCE, 518110, CHINA
 +86-755-28012955

WARRANTY

Illumina Face Mask has a 12-month warranty against defects in material or workmanship from Freedom Laser Therapy. Items damaged or abused due to user error or failure to follow instructions provided are not covered under the warranty and will be charged at standard repair rates.

Within 30 days of purchase, validate your warranty by visiting www.irestorelaser.com/warranty

ADDITIONAL TERMS – LIMITED PRODUCT WARRANTY

Please keep this warranty in a safe place.

During the limited product warranty period, Freedom Laser Therapy will repair or replace, at its sole discretion, any part or parts determined to be defective. In repairing or replacing the equipment, Freedom Laser Therapy may replace it with a substitute device or repair the equipment with substitute parts of like kind, quality and functionality. Substitute equipment or parts may have been remanufactured or may contain non-original manufactured parts.


This limited warranty does not provide coverage for equipment which has been subjected to abuse, misuse, exposure to moisture or water, unauthorized modification or repair, damage due to external causes or an Act of God, normal wear and tear which does not affect the original manufacturer's intended use.

This limited warranty is non-transferable and does not apply to equipment on which the serial number has been altered, defaced or is missing.

If a covered product failure occurs, please call Freedom Laser Therapy toll-free at 1-888-373-7425.

THIS WARRANTY AND THE REMEDIES SET FORTH ABOVE ARE EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES, REMEDIES AND CONDITIONS, WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED. ANY IMPLIED WARRANTIES, INCLUDING WITHOUT LIMITATION, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE OR USE, SHALL BE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY; OTHERWISE, THE REPLACEMENT OR REPAIR AS PROVIDED UNDER THIS LIMITED WARRANTY IS THE EXCLUSIVE REMEDY. TO THE FULL EXTENT PERMITTED BY LAW, IN NO EVENT SHALL FREEDOM LASER THERAPY BE LIABLE, WHETHER IN CONTRACT OR TORT, INCLUDING NEGLIGENCE, FOR ANY INDIRECT, INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES OF ANY KIND OR LOSS OF REVENUE OR PROFITS, LOSS OF BUSINESS, LOSS OF INFORMATION OR DATA, SOFTWARE OR APPLICATIONS OR OTHER FINANCIAL LOSS ARISING OUT OF OR IN CONNECTION WITH THE ABILITY OR INABILITY TO USE THE EQUIPMENT OR REPAIR THE EQUIPMENT.



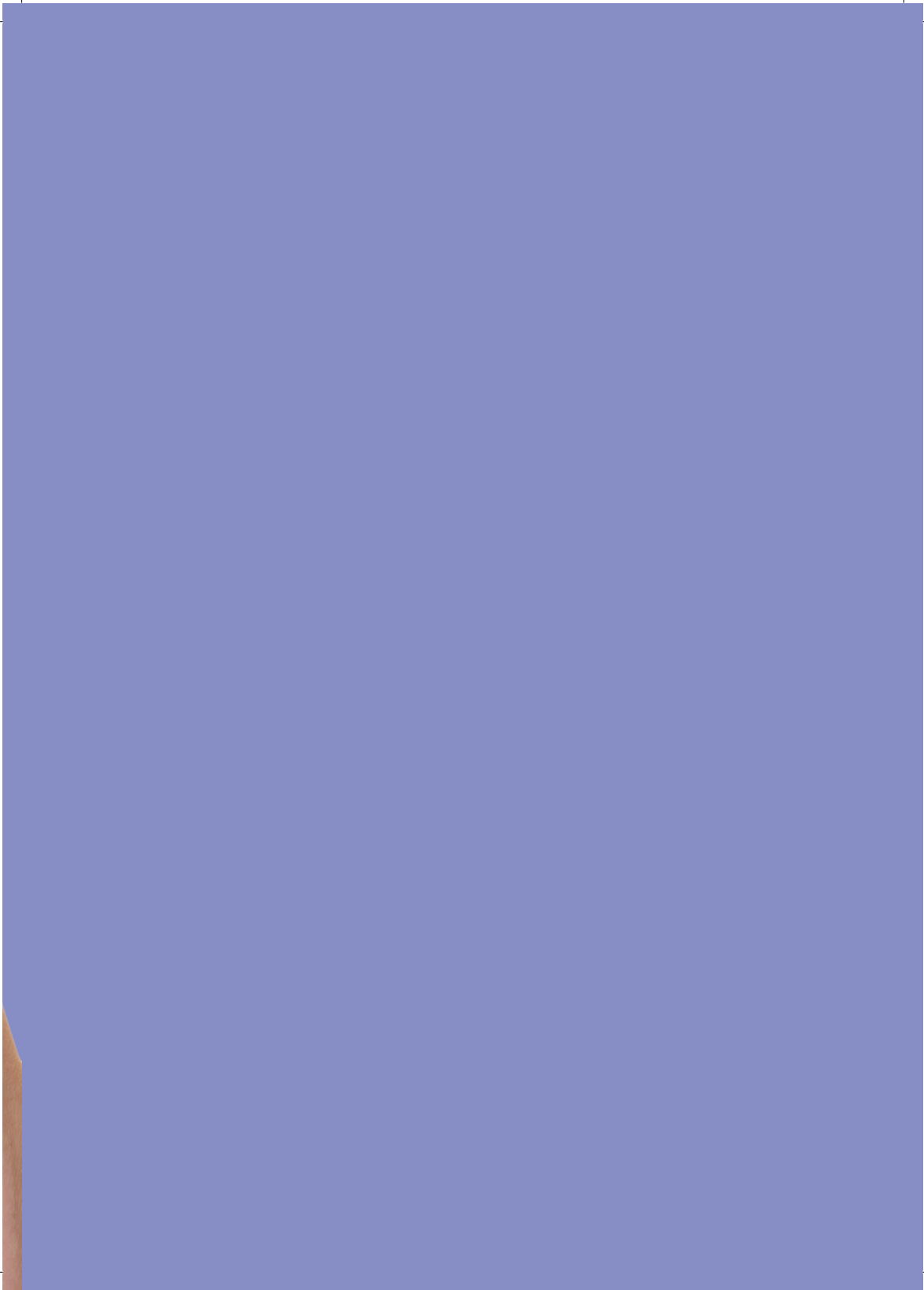

i R E S T O R E

FREEDOM LASER THERAPY, INC.
16782 VON KARMAN AVE, UNIT 15
IRVINE, CA 92606 USA
TOLL FREE 1-888-373-7425
WWW.IRESTORELASER.COM

©2024 FREEDOM LASER THERAPY, INC.
PATENT PENDING
VERSION: A/1