



CMYK 75 0 5 0



CMYK 24 44 0 0



CMYK 46 0 90 0



CMYK 23 16 13 46



CMYK 0 0 0 0



CMYK 5 0 19 0



CMYK 14 1 24 0



CMYK 21 2 15 0



CMYK 10 0 1 0



### User Program

Take a peek in Little Martin's Drawer™ and you will find exclusive products that will both surprise you and enrich your family's life together!

Participate in Little Martin's User Program to become a member of our product testing team and receive samples of pre-release products.

### Customer Care

Because our goal is to continuously create innovative and high quality baby products, we welcome you to visit our website or email us and let us know what you think.

### Our Mission

Little Martin's Drawer™ provides a trendy lineup of innovative baby products that makes parenting tasks easier and babies smile.

### Our Values

Love of patience.  
Love of protection.  
Love of companionship.

Little Martin's Drawer™ is a registered trademark of Ebzns LLC. Copyright© 2017 Ebzns Holdings LLC. All rights reserved.

[support@littlemartinsdrawer.com](mailto:support@littlemartinsdrawer.com)  
[www.littlemartinsdrawer.com](http://www.littlemartinsdrawer.com)



In USA imported for: Little Martin's Drawer,  
235 Montgomery Street, Suite 350,  
San Francisco, CA 94104

Made in China



keep  
baby  
healthy

## Little Martin's Electric Toothbrush

the reliable way to  
clean baby's teeth



LED flash light



16,000 strokes  
per min



PET  
brush head



smart timer


Thank you for choosing Little Martin's Drawer. We advise that you carefully read this instruction manual to ensure safe and correct use.

- Please keep this manual in a safe place.
- If giving as a gift, please be sure to include this manual.


### Precautions

- Adult supervision required at all times during use.
- Keep device and small parts out of reach of children.
- Assembly, cleaning and maintenance of the device is required by an adult.
- Do not open, repair, or alter the product without authorization.
- Do not soak product in water for long period of time.
- Do not place product in high temperatures or direct sunlight.
- Do not use if product is damaged.
- Do not use product for other uses other than intended.

### Charging

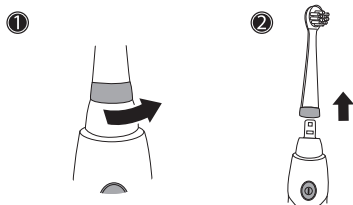


To detach: Turn brush head anti-clockwise and pull upwards.

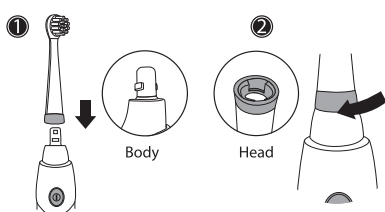

To attach: Align brush head and gently push downwards, then turn clockwise.

Note: Notch on neck should align with inner notch of brush head to lock in place.


### Battery Replacement

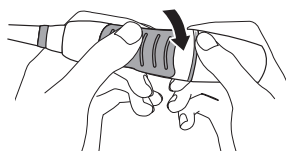

Turn battery cap anti-clockwise to remove. Insert AAA battery according to polarity signs. Align line marker on battery cap with line marker on brush handle, then turn clockwise to replace battery cap.

### Operating Instructions

1. Wet bristles and apply appropriate amount of toothpaste.
2. Divide oral area into up, down, left, and right sections.
3. Gently place toothbrush bristles against teeth at a slight angle towards the gum line.
4. Press the power button to turn on, then move brush head across teeth and gums with small upward and downward movements.
5. Move to the next section when the brush pulses (at 30 seconds) until all four sections of the mouth have been cleaned.
6. Brush will automatically power off after 2 minutes, or press the power button to turn off.

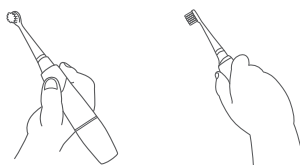

[When operating toothbrush, grip as if holding a pencil or a rock]

### Specifications

Material:	ABS, TPE, nylon
Speed:	16,000 strokes per minute
Battery Type:	AAA
Working Time:	120 days (twice per day, 2 mins per use)
Product Size:	164 x 25 x 25 mm
Weight:	31 grams



CMYK 0 56 18 0



CMYK 24 44 0 0



CMYK 46 0 90 0



CMYK 23 16 13 46



CMYK 0 0 0 0



CMYK 5 0 19 0



CMYK 14 1 24 0



CMYK 21 2 15 0



CMYK 0 8 1 0



### User Program

Take a peek in Little Martin's Drawer™ and you will find exclusive products that will both surprise you and enrich your family's life together!

Participate in Little Martin's User Program to become a member of our product testing team and receive samples of pre-release products.

### Customer Care

Because our goal is to continuously create innovative and high quality baby products, we welcome you to visit our website or email us and let us know what you think.

### Our Mission

Little Martin's Drawer™ provides a trendy lineup of innovative baby products that makes parenting tasks easier and babies smile.

### Our Values

Love of patience.  
Love of protection.  
Love of companionship.

Little Martin's Drawer™ is a registered trademark of Ebzns LLC. Copyright© 2017 Ebzns Holdings LLC. All rights reserved.

support@littlemartinsdrawer.com  
www.littlemartinsdrawer.com



In USA imported for: Little Martin's Drawer,  
235 Montgomery Street, Suite 350,  
San Francisco, CA 94104

Made in China



8 54878 00741 4

keep  
baby  
healthy

## Little Martin's Electric Toothbrush

the reliable way to  
clean baby's teeth



LED flash light



16,000 strokes  
per min



PET  
brush head



smart timer


Thank you for choosing Little Martin's Drawer. We advise that you carefully read this instruction manual to ensure safe and correct use.

- Please keep this manual in a safe place.
- If giving as a gift, please be sure to include this manual.


### Precautions

- Adult supervision required at all times during use.
- Keep device and small parts out of reach of children.
- Assembly, cleaning and maintenance of the device is required by an adult.
- Do not open, repair, or alter the product without authorization.
- Do not soak product in water for long period of time.
- Do not place product in high temperatures or direct sunlight.
- Do not use if product is damaged.
- Do not use product for other uses other than intended.

### Charging

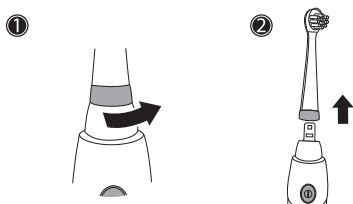


To detach: Turn brush head anti-clockwise and pull upwards.

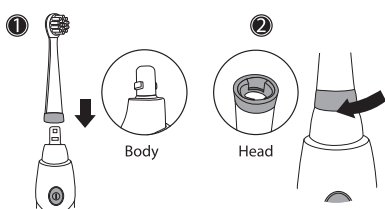

To attach: Align brush head and gently push downwards, then turn clockwise.

Note: Notch on neck should align with inner notch of brush head to lock in place.


### Battery Replacement

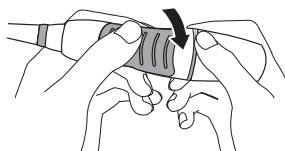

Turn battery cap anti-clockwise to remove. Insert AAA battery according to polarity signs. Align line marker on battery cap with line marker on brush handle, then turn clockwise to replace battery cap.

### Operating Instructions

1. Wet bristles and apply appropriate amount of toothpaste.
2. Divide oral area into up, down, left, and right sections.
3. Gently place toothbrush bristles against teeth at a slight angle towards the gum line.
4. Press the power button to turn on, then move brush head across teeth and gums with small upward and downward movements.
5. Move to the next section when the brush pulses (at 30 seconds) until all four sections of the mouth have been cleaned.
6. Brush will automatically power off after 2 minutes, or press the power button to turn off.

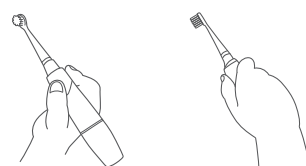

[When operating toothbrush, grip as if holding a pencil or a rock]

### Specifications

Material:	ABS, TPE, nylon
Speed:	16,000 strokes per minute
Battery Type:	AAA
Working Time:	120 days (twice per day, 2 mins per use)
Product Size:	164 x 25 x 25 mm
Weight:	31 grams