



User Manual

Mini Stand

 JOOVV.COM

 [JOOVVSOCIAL](https://www.instagram.com/JOOVVSOCIAL)

 [JOOVVSOCIAL](https://www.facebook.com/JOOVVSOCIAL)



Corporate Office

Joovv, Inc.
555 N. El Camino Real, Suite 342
San Clemente, CA 92672

Customer Service

Email: help@joovv.com
Web: www.joovv.com
Setup: www.joovv.com/setup

Important Safety Instructions

CAUTION: Do not use this product for any purpose not explicitly specified by Joovv.

Improper installation may cause property damage or personal injury. If you do not understand these directions, or have doubts about the safety of the installation, contact Customer Service at help@joovv.com. Joovv is not responsible for damage or injury caused by incorrect installation or use. This product is not designed to hold devices other than the Joovv Mini.

This product contains small items that could be choking hazards if swallowed. Keep these items away from young children.

Parts List

1 x Mini Stand	1 x M4 Allen Wrench
4 x M4x6mm Allen Screws	1 x Instruction Manual



If you'd prefer instructional setup videos, please visit: www.joovv.com/setup

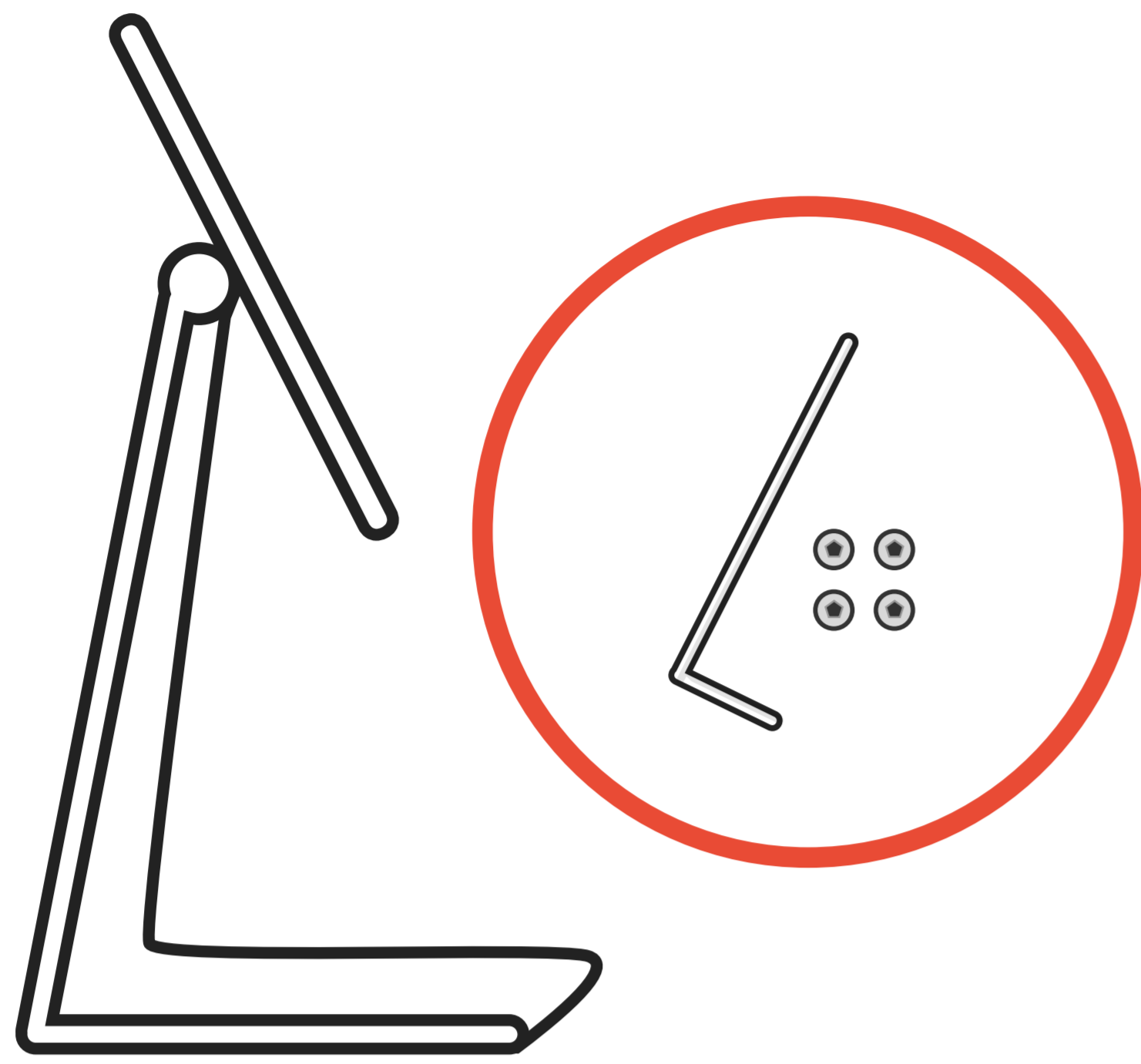


FIG. A

STEP 1

Lay out all hardware included on a padded surface, such as carpet or a blanket, as shown in **Figure A**.

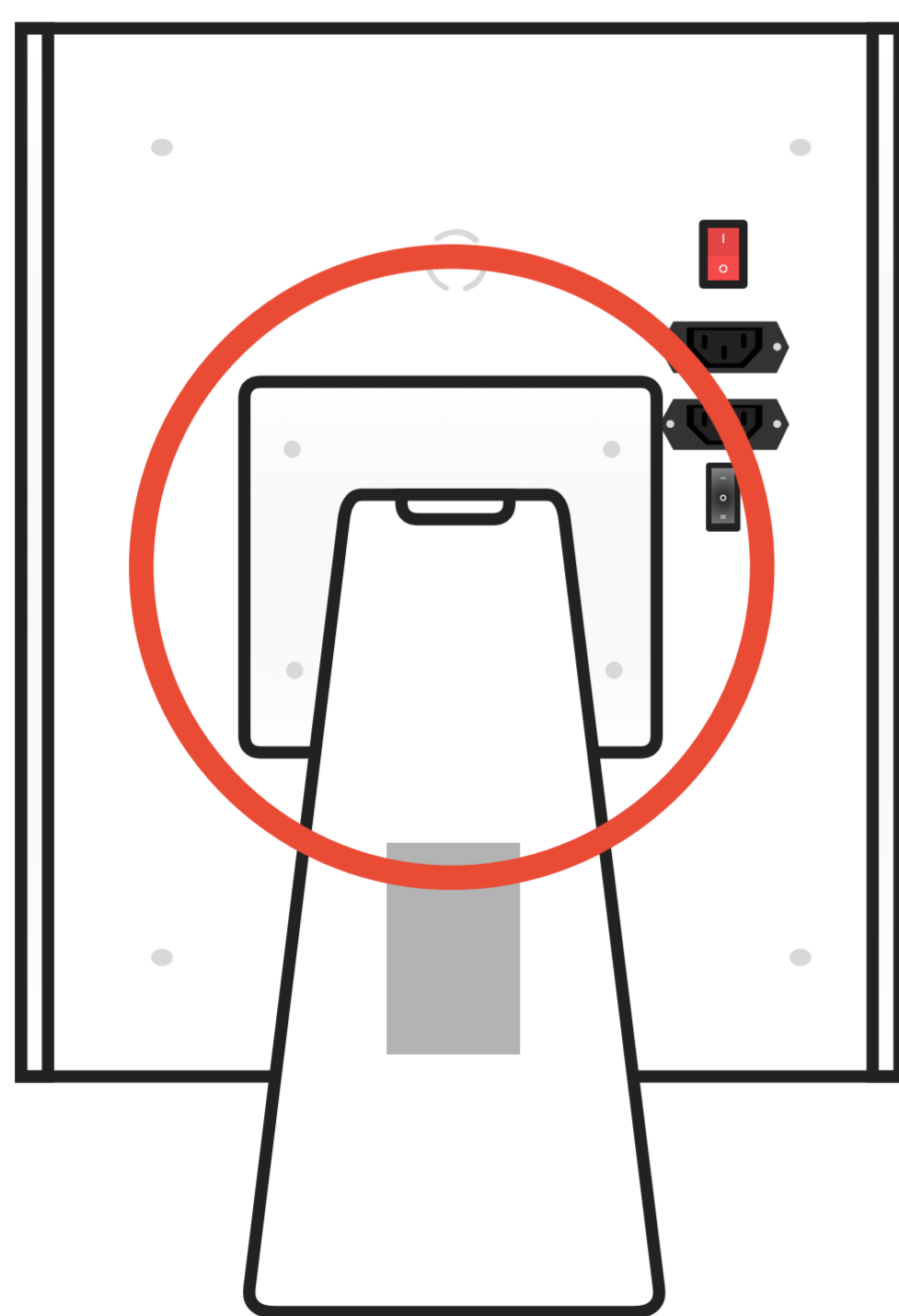


FIG. B

STEP 2

Carefully lay the Joovv Mini face down on a padded surface. Align the faceplate on the Mini Stand with the designated screw holes on the back of the Joovv Mini, as shown in **Figure B**.

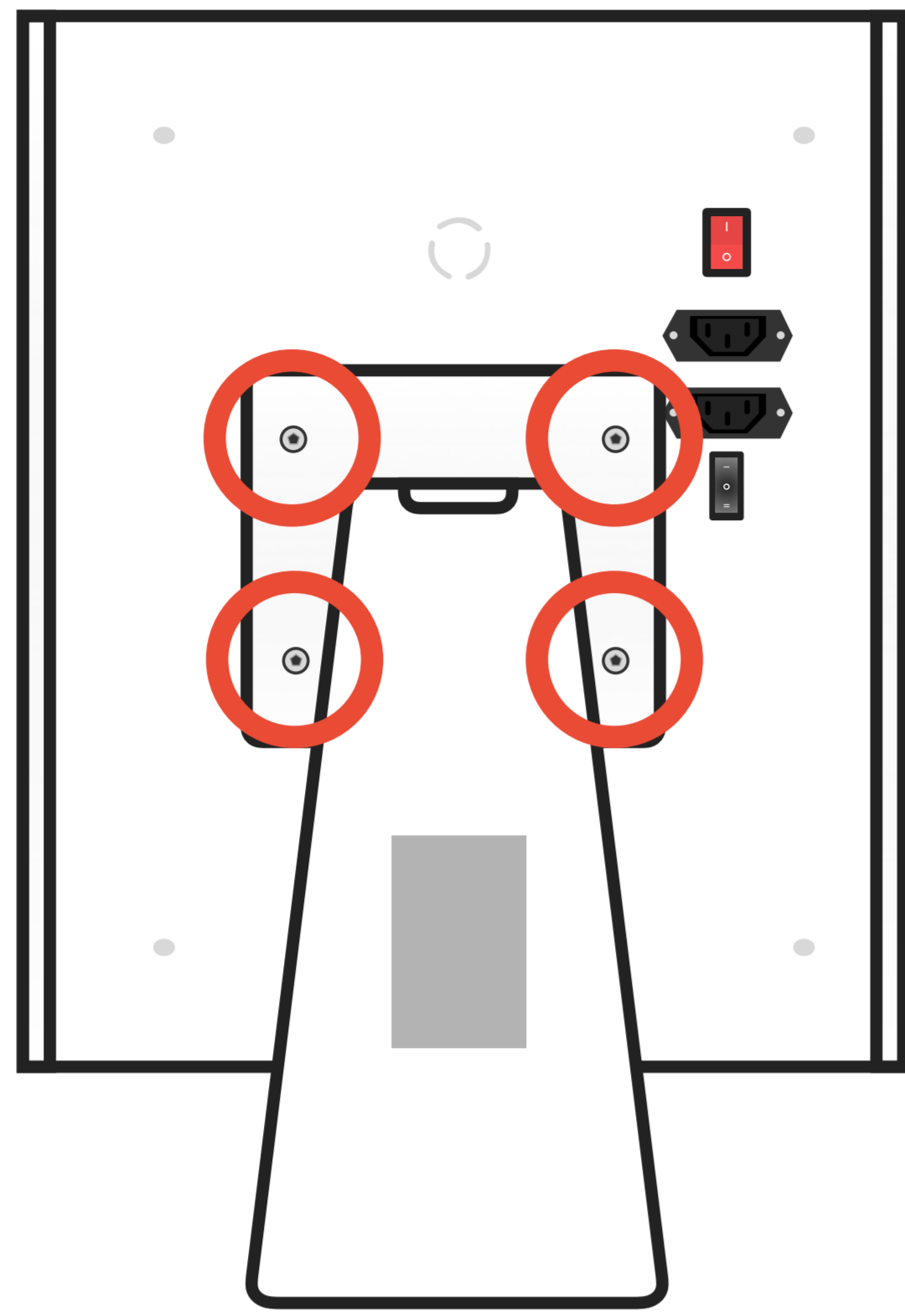


FIG. C

STEP 3

Using the included Allen wrench, screw in the four (4) provided screws in order to attach the Mini Stand faceplate to the Joovv Mini, as shown in **Figure C**.



FIG. D

STEP 4

Set the Mini Stand on any stable, flat surface and adjust the angle of the Joovv Mini to target the desired area, as shown in **Figure D**. Then, plug in your Joovv Mini and you're ready to start Joovvin'!*



*Please refer to your Joovv Mini manual for full operating instructions.

Warranty

Joovv, Inc. warrants that its devices will function and perform within the specifications for two years from the date of delivery. The warranty covers all parts and labor during the warranty period with some exclusions. If service is required, the buyer is responsible to return the device to Joovv within the 2-year warranty period. The warranty shall not apply to equipment that has been abused, modified, damaged, or fitted with substitute parts. Joovv shall not be liable for any consequential loss or inconvenience caused by any failure of its products. If the device is repaired or replaced, Joovv reserves the right to provide parts in different colors from the original device.

To obtain warranty service, contact Joovv Customer Service at help@joovv.com or visit www.joovv.com/warranty to request an authorization number.



Corporate Office

Joovv, Inc.
555 N. El Camino Real, Suite 342
San Clemente, CA 92672

Customer Service

Email: help@joovv.com
Web: www.joovv.com
Setup: www.joovv.com/setup