


Lifesmart

Massage Chair

LC3200-S
Manual

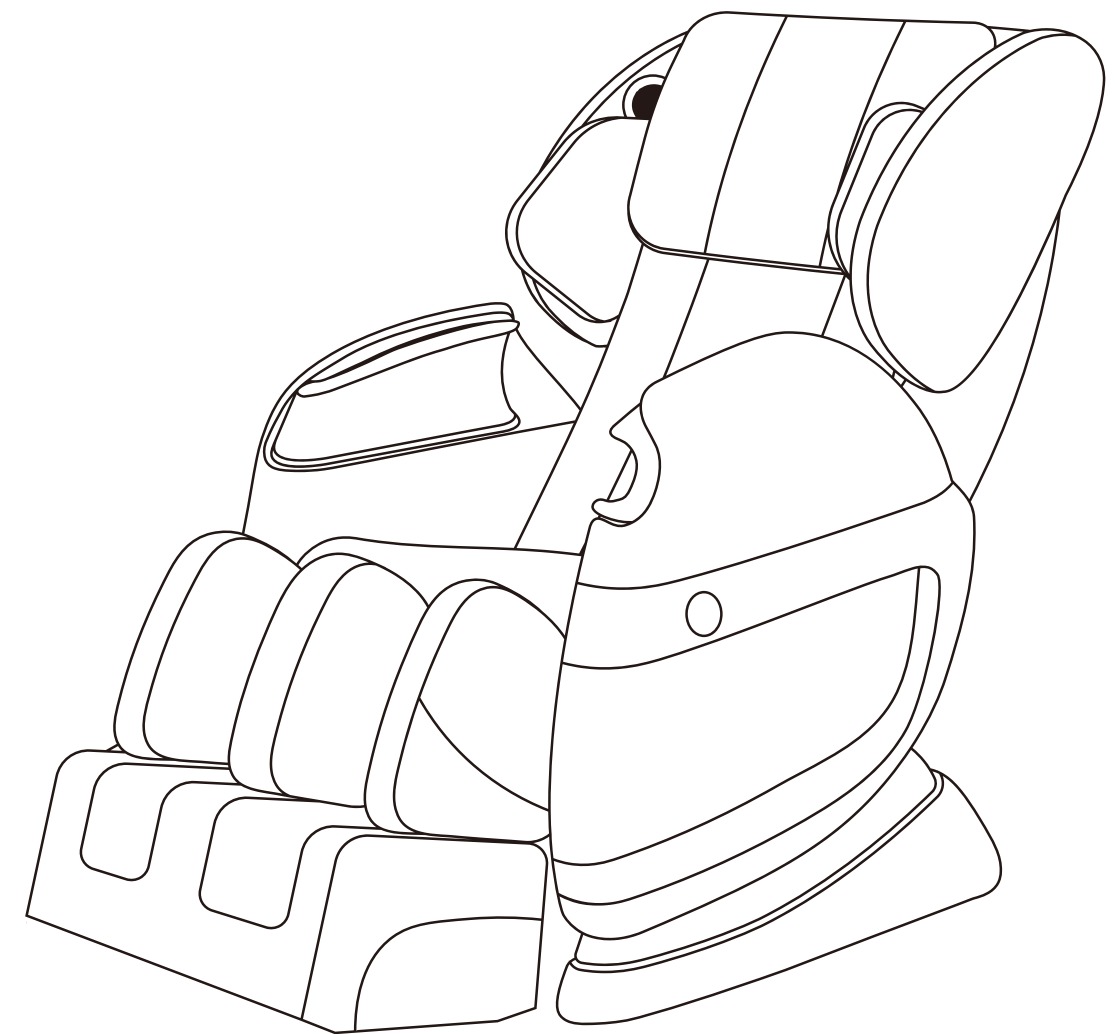

Questions, problems, missing parts? Contact technical support department at **657-341-0362, 9:00 a.m. - 5 p.m., PST, Monday - Friday** or visit us online at www.esmartpro.com
Please have your model and serial number ready



M_LC3200-S_2019
Printed in China

Contents

Thank you for purchasing eSmart massage chair.
Please read the owner’s manual (operating manual) before use.

Safety Precautions.....1-4

Grounding or double-insulation5

Components & Descriptions.....6-7

Before Use.....8

Bluetooth Speakers.....10

Adjusting Massage Chair.....13

Start Massage.....14

After Use.....14

Troubleshooting.....15

Safety Check.....15

Storage and Care.....16

Specification.....17

- This product is intended for residential use only.
- Use this product only as described in this manual. Any other use is not recommended by the manufacurer and may cause fire, electric shock or injury to persons.
- Remember to keep the manual in an easily accessible place.

Safety Precautions

IMPORTANT SAFETY INSTRUCTIONS

This massage chair is intended for household only .When using an electrical appliance, basic precautions should always be followed, including the following:
Read all instructions before using this massage chair.

DANGER-To reduce the risk of electric shock:

- 1)Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- 2)Never use pins or other metallic fasteners with this appliance.
- 3)Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
- 4)Keep Dry - Do not operate in a wet or moist condition.

WARNING-To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet. when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons. Keep children away from extended foot support.
- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.
- 10)Do not use outdoors.
- 11)Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 12)To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13)Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
- 14)Do Not Crush - Avoid sharp folds.
- 15)Do not use massager in close proximity to loose clothing or jewelry.
- 16)Keep long hair away from massager while in use.
- 17)Connect this appliance to a properly grounded outlet only. See Grounding Instructions.
- 18)Remove control box when not in use, and store out of reach of children.
- 19)Burns will result from improper use, read instructions before using, do not wet-do not use pins, never remove cover.

IMPORTANT SAFETY INSTRUCTIONS

● Safety Precautions

- To ensure safe and correct use of the massage chair, do not operate it without reading this operation manual carefully.
- Do not use the massage chair when the upholstery or seat is removed. This may cause injury or malfunction.
- Do not use the massage chair in combination with other therapeutic equipment or electric blanket, etc., as this may result in ineffectiveness or injury.
- Do not use the massage chair an hour before or after eating as this may have adverse effects on the user.
- Do not allow children or pets to play around the massage chair, e.g., behind the backrest, under the seat or leg rest of the massage chair. This may result in injury.
- Do not rest or place heavy objects on the armrests, leg rest or backrest as this may result in malfunction or injury.
- Do not use the massage chair when your body is wet.
- Do not operate the massage chair with wet hands.
- It is not recommended to use the massage for more than 15 minutes.
- Do not perform continuous massage on the same spot of your body for more than 5 minutes at a time as this may result in excessive stimulation and may have adverse effects.
- Do not insert your hand or foot along the paths of the massage rollers during use as this may cause injury.
- Stop operation immediately if you experience any discomfort.
- This product is not intended for self-treatment of conditions that should be managed by a qualified health care provider.
- People on medication or with medical conditions, please consult your doctor before use. Please do not use the product if you are not feeling well.
- Remove all scarves, neckties, necklaces and jewelry before using the massager.

People with Medical Conditions

Consult your doctor before use if you:

- are under medical rest as ordered by a doctor.

have spinal disorders, an abnormal spinal condition or have suffered a spinal injury.

have back problems.

have diabetes, osteoporosis or sensory impairment.

have joint dysfunction.

have a pacemaker or other electronic medical devices.

may be pregnant.

have phlebitis or thrombosis.

have an increased risk for blood clots.

had recent surgery.

have surgical pins, screws, or anything mechanical implanted in your body.

SAVE THESE INSTRUCTIONS

Safety Instruction

● WARNING

The power plug must be fully inserted into and appropriate outlet. Improper connection of the equipment-grounding conductor could result in an electric shock.

Check with a qualified electrician and servicemen if you are in doubt as to whether the product is properly grounded. Do not modify the power plug provided by the product. If the outlet is not suitable. please ensure that a proper outlet is installed by qualified Electrician.

When using the massage chair, basic precaution should always be followed to reduce the risk of electric shock, an injury or burns to the person.

Read the following details before usage:

1. Do not use the massage chair outdoors and avoid places with direct sunlight.
2. Do not use the massage chair in places where the humidity is high .
3. Do not use the massage chair near any heating device or fire .
4. Do not use the massage chair when smoking.
5. Recommended to unplug the power supply after use.
6. Please ensure that a proper outlet is installed by qualified electrician, Improper connection of the equipment-grounding conductor could result in an electric shock.

⚠ Caution

- The massage chair is to be used indoors only.
- If you use the massage chair in a cold room, do not increase the room temperature abruptly. It is recommended that you increase the temperature gradually to a normal level.
- If the massage chair is stored in a cold place and is being brought into a warm environment, it is advisable to wait one hour before use. This is because its performance may not be at optimum level, as there may be water droplets on the mechanical parts due to condensation. Using the massage chair under such conditions may result in malfunction.

User-Maintenance Instruction

User-maintenance instructions shall include:

- a) Instructions for cleaning and user maintenance operations are commended by the manufacturer, such as lubrication or nonlubrication , and a statement to the user that any other servicing should be performed by an authorized service representative or that the appliance has no user serviceable parts.
- b) Instructions for an appliance employing an automatically reset thermal limiter that shuts off the entire appliance shall inform the user what to expect if the thermal limiter operates
- c) Specific instructions for the proper method of storage of the cord, the total appliance, and the like, when the appliance is not in use, and for care of the cord while in use (for example, for a hand-supported product, untwisting)
- d) In the case of an appliance intended to be used with water additives, conditioners, or other solutions with or without water, specific instructions regarding the proper liquid or additive to use and the amount to be used in conjunction with the appliance

Safety Precautions

⚠ Caution

- The massage chair is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the massage chair by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the massage chair.
- Avoid using the massage chair under the influence of alcohol.
- Do not apply massage rollers directly on the head, elbow or knee joints, torso or abdomen.
- If your legs or feet slip off the leg rest during massage, do not force them back in as this may result in injury.
- Remove the power plug from the electrical outlet immediately should liquid be accidentally spilled onto the massage chair. Do not use wet hands to connect or disconnect the power plug from the electrical outlet.
- Do not tamper with the power cord or use the massage chair with a damaged cord as this may result in malfunction and electric shock.

● Instances When the massage chair Should Not Be Used

In the following events, please switch off the main power switch immediately and disconnect the power cord from the electrical outlet:

- If water is accidentally spilled onto the massage chair. This may result in electric shock or malfunction.
- If the fabric of the backrest is torn and the internal components are exposed.
- If pain or discomfort is felt during massage, stop operation immediately and consult your doctor.
- If you detect any malfunction or any other abnormal condition during operation.
- If there is a power failure. Injury may occur when power is restored unexpectedly.
- If there is lightning.

● Assembly and Repair of the massage chair

- Do not disassemble the backrest cover of the massage chair. Touching the internal components may result in malfunction or electric shock.
- The massage chair must only be professionally serviced or repaired. Do not attempt to disassemble or repair the massage chair yourself.

● Things to Note About the Power Plug and Cord

- Check that the Alternating Current (AC) voltage corresponds with the specifications indicated for the massage chair.
 - Do not connect or disconnect the power plug from the electrical outlet with wet hands. This may result in malfunction or electric shock.
 - When disconnecting the power plug, pull it by holding the plug, not the cord.
 - Although operation will be stopped automatically by the Auto-Timer function, always remember to switch off the main power switch after use.
 - Do not place the power cord under the massage chair or any other heavy objects.
 - Do not wind the power cord around the massage chair as this may damage the cord and result in a fire or electric shock.
 - Do not operate the massage chair with a damaged cord or extension cord.
 - Do not use the massage chair if the electrical outlet is loose.
- If the power cord is damaged, it should be replaced by a qualified technician.

Grounding or double-insulation.

This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

DANGER Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product-if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

b) For a grounded, cord-connected product rated less than 15 A and intended for use on a nominal 120-v supply circuit, the instructions in either (1) or (2):


1) This product is for use on a nominal 120-v circuit, and has a grounding plug that looks like the plug illustrated in sketch A in Figure 86.1. A temporary adapter that looks like the adapter illustrated in sketches b and c may be used to connect this plug to a 2-pole receptacle as shown in sketch B if a properly grounded outlet is not available. The temporary adapter should be used only until a properly grounded connected outlet (sketch A) can be installed by a qualified electrician. The green colored rigid ear, lug, or the like extending from the adapter must be held in place by a metal screw.

2) This product is for use on a nominal 120-v circuit and has a grounding plug that looks like the plug illustrated in sketch a in Figure 86.1. Make sure that the product is connected to an outlet having the same configuration as the plug. No adapter should be used with this product.

C) For all other grounded, cord-connected products:


This product is for use on a circuit having a nominal rating more than 120 V (or This product is rated more than 15 a and is for use on a circuit having a nominal rating of 120v") and is factory-equipped with a specific electric cord and plug to permit connection to a proper electric circuit. Make sure that the product is connected to an outlet having the same configuration as the plug. No adapter should be used with this product. If the product must be reconnected for use on a different type of electric circuit the reconnection should be made by qualified service personnel.

d) For a permanently connected product

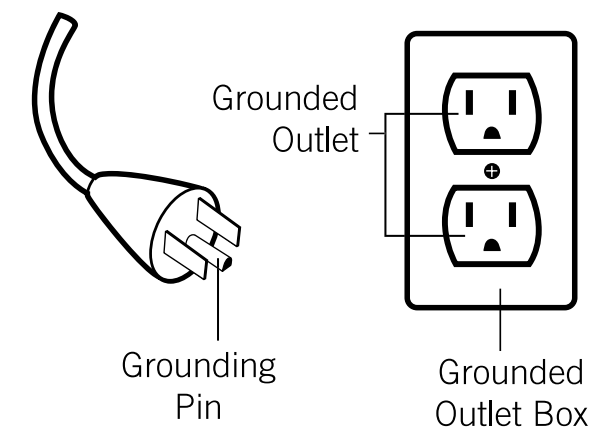

Components & Descriptions

● Chair

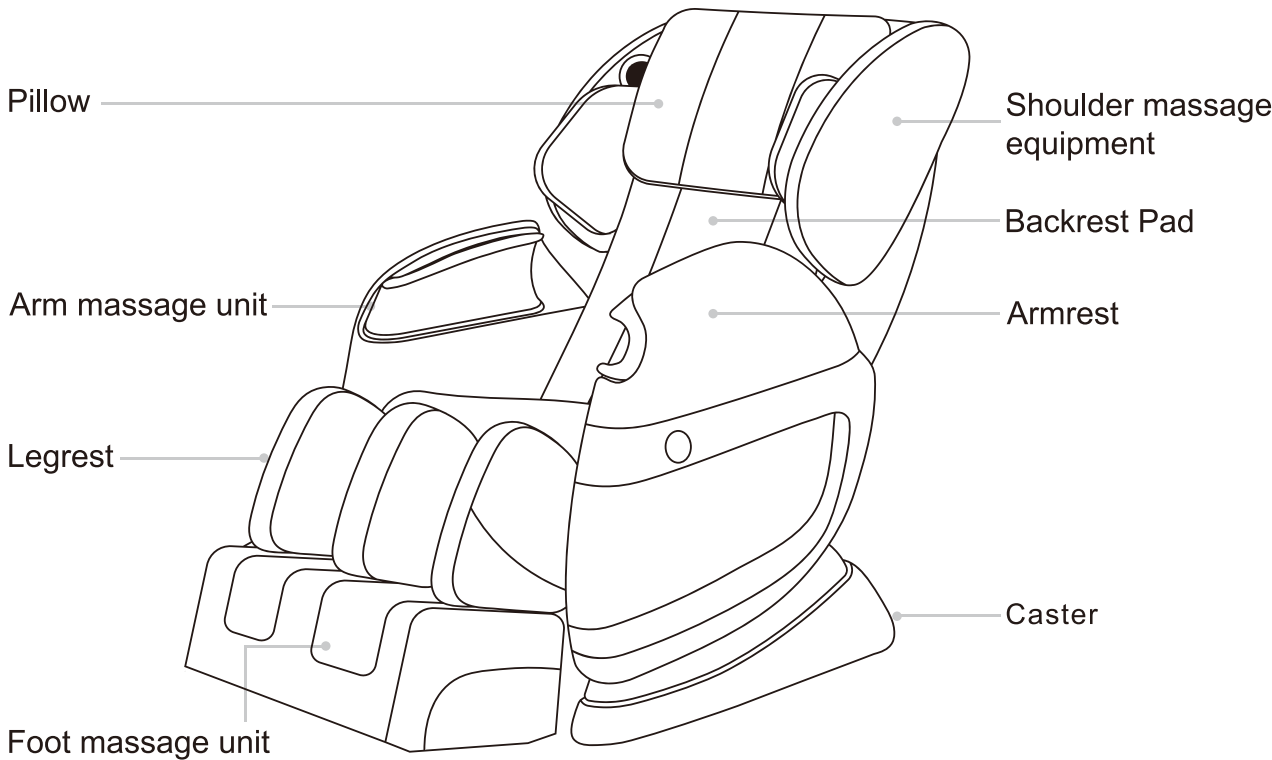


● Power cord and power socket



Attention1:Before inserting power cord,please ensure main power switch is off.
Attention2:Do not disassemble or replace the power cord when it is damaged.A retailer or a qualified professional shall be needed.



● Main power switch

Turn on the main power switch, massage chair into standby mode

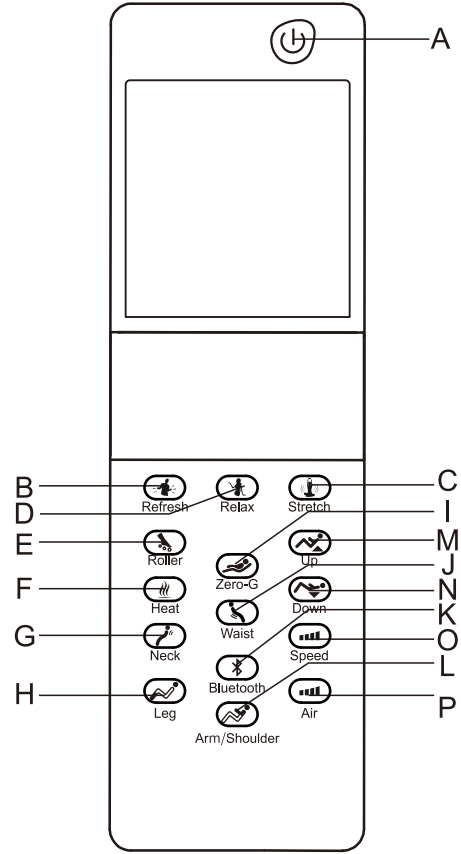● Controller

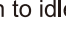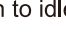

A	Power	Press the button to light on the LCD screen. Press any massage function button to start massage. Press the button again to stop all massage function, the massage rollers will return to its original position
B	Refresh function	Press the button to begin Refresh function
C	Stretch	Press this button to begin Stretch massage
D	Relax	Press this button to begin Relax massage
E	Roller	Press the button to activate the rolling function at the bottom of the feet.Press again to stop the function.
F	Heat	Press the button to activate the heating on the back,Press again to stop the heating function.
G	Neck massage	Press the button to start massage on the neck. Press the button again to stop the neck massage.
H	Leg air massage	Press the button to start rhythmic air pressure massage on legs. Press the button to stop the function.
I	Zero Gravity	Press this button to activate the Zero-Gravity function. The icon () will be lighted. Press the button to return to idle position. The icon () will not be lighted.
J	Waist massage	Press this button to start roller massage on your waist. Press the button again to stop the function.
K	Bluetooth	Press the button to power Bluetooth ON/OFF.
L	Arm/Shoulder	Press the button to start Arm and Shoulder massage with rhythmic air pressure.Press again to stop the function.
M	Up	Angle Adjustment control button for massage chair. Press and hold“up”button until the chair is at the desired angle.
N	Down	Press and hold“down”button until the chair is at the desired angle.
O	Speed	The speed Adjustment button for Kneading speed. Each time the button is pressed, the intensity will increase by 1 level. The maximum level is 3. The speed can be adjusted in the cycle.
P	Air	The Air Intensity Adjustment button for Air-pressure Each time the button is pressed, the intensity will increase by 1 level. The maximum level is 3.


 WARNING	● Please do not spray water or other liquid on the chair. ● Do not operate the chair with wet hand.
 CAUTION	● Do not sit on the controller, do not drop or step on the controller. Do not pull over the wire. ● Do not press the button with sharp object or it will result in damage.

Before use

- When the temperature is too low to run the machine, please increase the ambient temperature gradually, do not increase the temperature in sudden. Please move the chair to the warm place and keep on waiting for 1 hour before starting.
→ Or it will result in moisture of parts and abnormal operation.


How to assemble the shoulder

- Unzipping the shoulder device, put the iron part into the slots on the end of the shaft, let the hole "A" match with the slots "B"
- Connect the air hose "E" to the air hose "F" exposed on the backrest.
- Use the screw (M6 X 16) to fix on the shaft of the backrest.
- Zip the cover
- Repeat above steps to assembly the other shoulder.

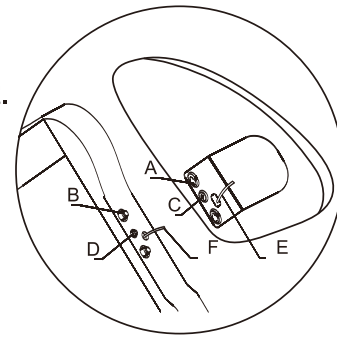

How to assemble backrest pad and seat cushion

- Connect the wire of the heater plug from the backrest to underneath of the seat (Fig 1)
Connect the cable of remote controller with the cable under seat.(Fig 1)
- Connect the backrest and the seat by zipper. (Fig 2)
- Attached backrest pad with zipper. (Fig 3)
- Attached the pillow pad using the velcro fasten (fig 4)

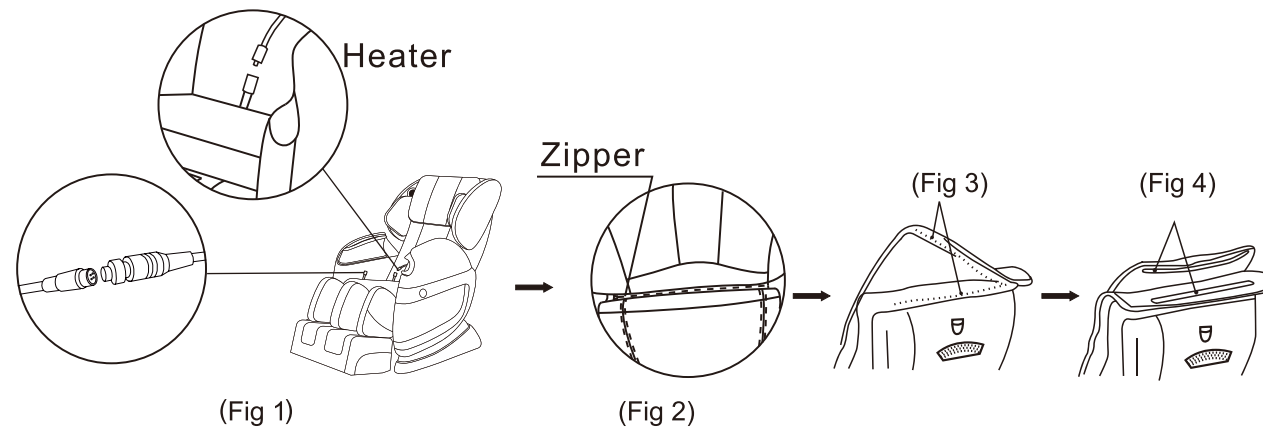


Where to place the chair

Please make sure that the massage chair has about 19.685in in front of it, leaving about 19.685in behind the chair to allow free reclining of the backrest.




- Please do not use in a humid place such as sauna or bathroom.
- Place on a levels surface before use.
- Please do not use in the following situations: in a utility room or dusty place or in high temperatures (over 40°C).
- Do not place it in the sun or near heating appliances.
- Please keep cigarettes, ashtrays or flammable goods away from the massage chair.

- To prevent damage to the floor and reduce noise, please lay carpet under the massage chair (do not use an electric blanket).
- Please be careful placing this machines, as it is heavy.
- To avoid excessive noise, make sure the massage chair is a minimum of 1 meter away from the TV, radio, and other AV machine.

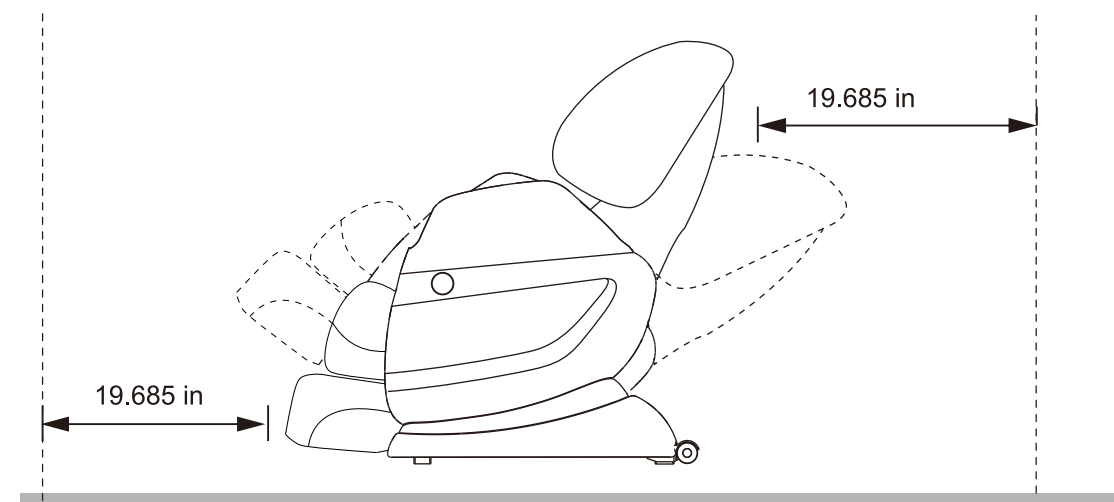

- Please ensure there is a chair and accessories inside the packing box

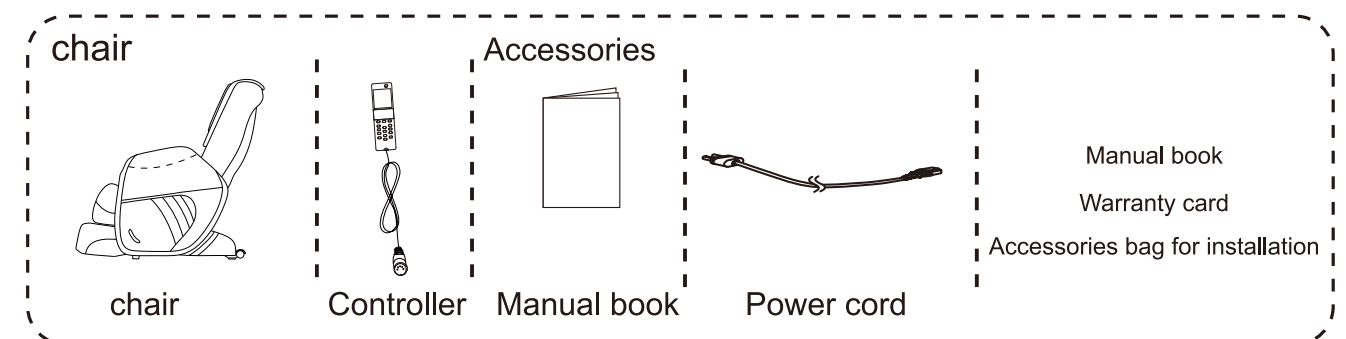


How to move

When moving, please turn off main power switch, then unplug the power plug.

Using caster to move

Please put the back of the chair to the bottom and then hit the rear wheel, then move

Lift up the calvesrest then recline the chair to about 45 degrees to easily move the chair with a slight pull or push

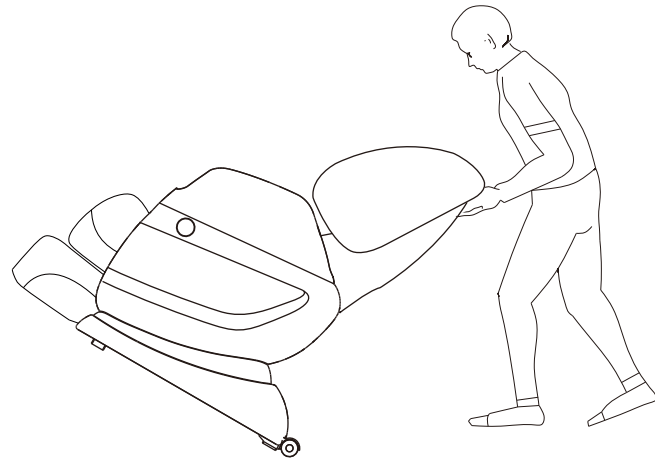


- The chair cannot be moved if the passage is less than 29.528in in width.
- Contact the dealer in case you cannot move the chair normally.
- Due to the heavy weight of the chair, make sure that it will not damage your waist or other parts of your body during moving.
- Controller is connected to the chair, make sure it does not fall during moving.
- Do not move the chair with person loaded.
- Please move the chair after unplugging.
---> Otherwise it may cause damages to the power plug, causing fire disaster or electric.

Bluetooth

This massage chair has two Bluetooth speakers built in to play your music via a compatible Bluetooth enabled music device.


1. Press Bluetooth button on the remote control.
2. Using your Bluetooth music player, locate the massage chair name "LC3200-S"
3. Connect and enjoy your music.

How to open the backrest

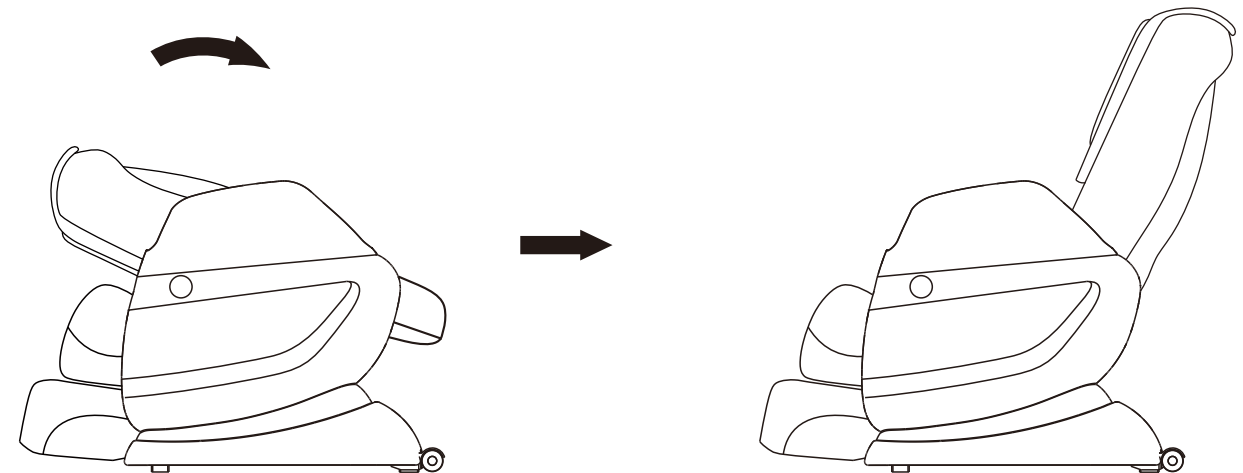

Turn the backrest to back until a "click" sound is heard.


Fold backrest

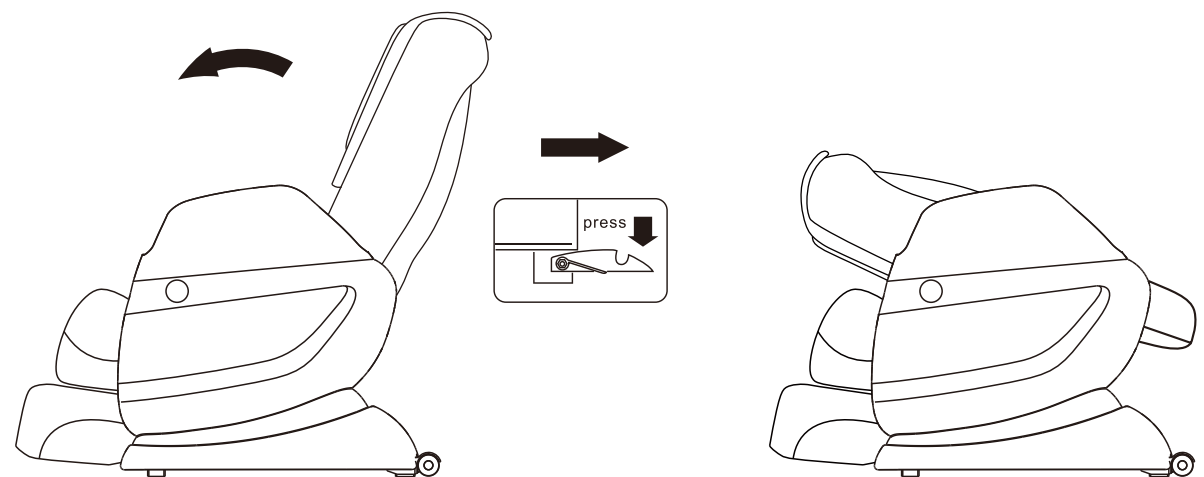


Release lock lever as figured to lay to backrest onto the seat cushion.

How to switch on the power

① Insert power cord into outlet

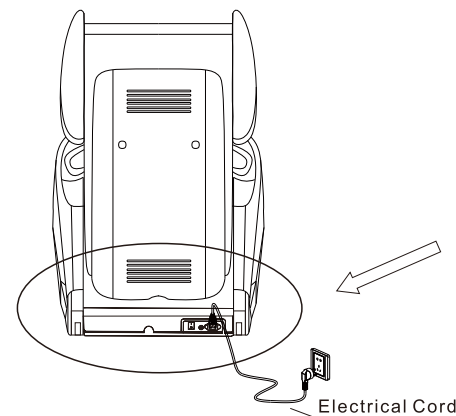

② Turn on main power switch on the back of the chair

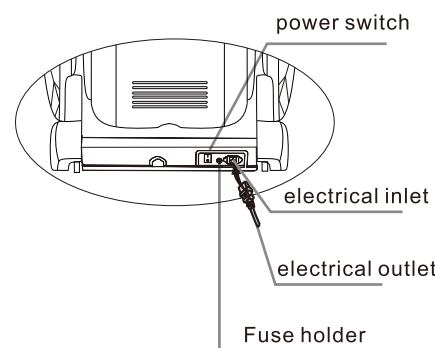


※ If the rollers get stuck, turn off the switch and press the ON/OFF button, the roller will be back to its initial position.

! WARNING	<ul style="list-style-type: none"> Do not connect power cord to an external transformer
! CAUTION	<ul style="list-style-type: none"> Ensure the main power switch is off before plugging the power socket or controller. Please insert power plug and controller plug fully Read this manual book carefully before using, set the chair according to the requires, ensure normal operation

Connecting the power


Check the power cord, power plug, socket on the chair, main power switch, controller are unbroken for the first time use or use after long time storage. Then plug power cord.

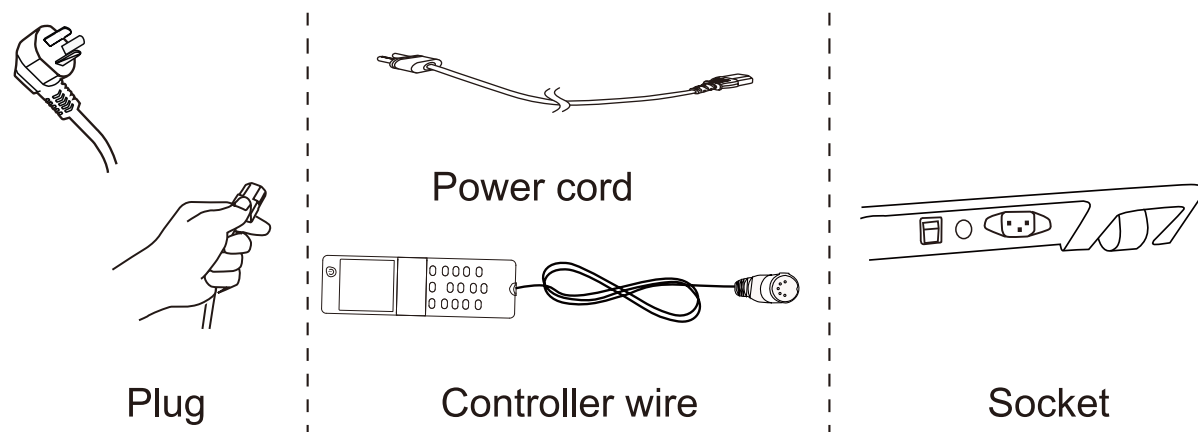


Confirm the fabric

Make sure if the fabric is damaged according to the following procedures before use.

① Open zipper and back cushion from one side

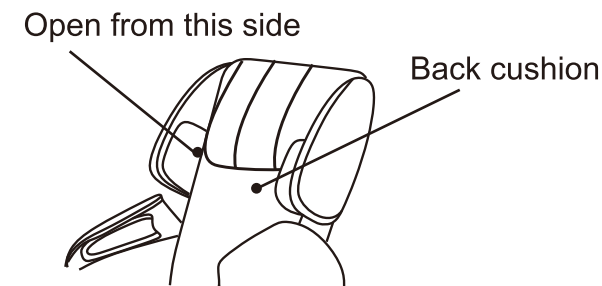

② Put your hand between cushion and backrest to see if the fabric is damaged.

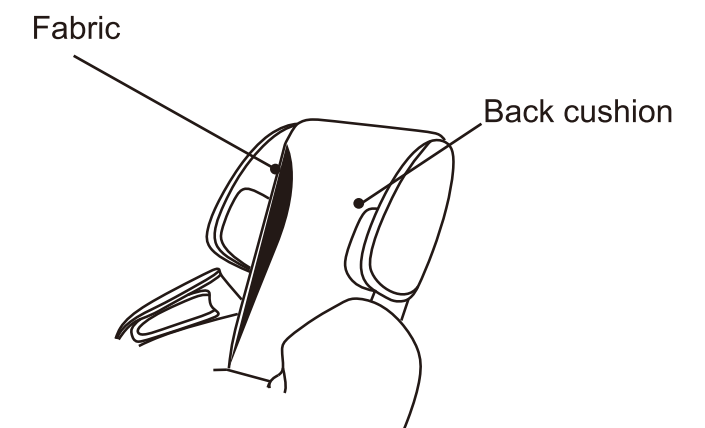

! CAUTION	<p>Before using the chair, confirm that whether there is damage in the material of the armchair or of other parts. In case of damage, even the smallest tearing, in the chair, please stop using it immediately, turn off the power switch, unplug the chair, and contact the customer service center. Never use the armchair with dismantled cushion or damaged material, otherwise there will be serious safety hazard.</p>
-----------	---

Adjusting massage chair


Sitting in the chair

Adjust the backrest to your desired position by pressing and holding on the "UP" "DOWN" buttons.

! CAUTION	<ul style="list-style-type: none"> Do not sit in the chair when legrest is rising, or after legrest rises, it may result in accident or reclining chair, causing personal injury. Raise the chair before using. Enjoy the massage from light massage course. Adjust the angle of your reclining angle according to your body figure. Do not stand on the legrest, it will break the legrest or even the whole chair.
! WARNING	<p>When adjusting angle, please make sure there is no person (Especially children) or pet around. Besides, there must be at least 11.811in space at the front or back of the chair, and there is no foreign matter in the way when rotating.</p>

Start massage

1. Connecting the power,turn on main power switch

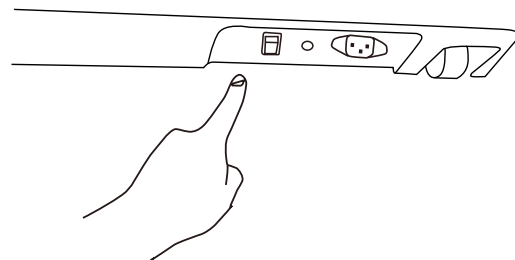

● In the Auto Mode

You may change from one massage program to another at any time without having to wait for it to complete. The twinkled button is the current massage mode

● In the Manual Mode

You may change from one massage program to another at any time without having to wait for it to complete.

After use

Rise the chair up and lower the footrest completely.

Press ON/OFF key to turn off the chair,wait until chair gets back to its initial status.

After airbags deflate completely,rollers get back to its initial position,sole massage stops, heating stops,then step out of the legrest.Do not stand in the legrest.Turn off main power switch,unplug power plug.

! CAUTION

After you use the chair, turn off the main power switch, and unplug the power plug from the socket.

→Otherwise, accident or injuries may be caused by children or pets.

After you use the chair, be sure to raise back completely and lower the footrest completely.

After you use the chair, please clean it and place it in dry places.

When unplugging the power cord, please do not pull it directly; instead, you should unplug it pinching the power plug.

Please be sure to put the power cord away in prevention of tripping over.

Troubleshooting

If there is still a problem after trying the solutions given in this troubleshooting guide, stop using the chair immediately to prevent accident or injury, and unplug the power plug from the wall outlet. Contact the vendor or your nearest service center for assistance.

Problem	Causes and troubleshooting
The massage chair can not operate even though it has been connected to the power supply.	<ul style="list-style-type: none"> ● Check whether the plug is well connected and working well ● Check whether the main switch on the massage chair is turned to position “I”
During operation, the controller is working, but the chair can not work	<ul style="list-style-type: none"> ● The massage chair could be overheated after prolonged use. It will automatically cool down to resume its normal condition. ● Please turn off the power and allow the chair to cool down for 30 or 50 minutes
During operation, both the controller and the massage can not work	<ul style="list-style-type: none"> ● If the set-time in the Timer is complete, please press the timer function again to resume operation ● Switch OFF the power and switch ON again to restart the chair

Safety check

Regular Maintenance Will Preserve the Life of Your Chair

Have you experienced any of these problems?

- The chair does not start with the Main POWER switch set to ON.
- The power cord has been cracked or scratched.
- The power cord or plug becomes abnormally hot.
- The chair produces a burnt smell, abnormal sound or vibration.
- The auto timer function does not work in some cases.
- Other failures.

If you find anything listed on the left, stop using and unplug the power cord from the outlet to prevent failure or an accident. Be sure to contact the vendor or our service center for maintenance.

To use safely

The chair fabric (roller fabric) is a consumable part, and its service life is 3 years (conditions of use: two 15-minute sessions a day, every day). For safe use, periodic inspection by the vendor is recommended before the end of its service life.

*Durable years vary depending on conditions of use.

* After about five years of use, a professional inspection is recommended even if there are no symptoms or problems. For inspection fees, contact the vendor.

Storage and care

Storage

- After cleaning the dirt and dust, place the chair in a dry place
- If the chair is left unused for a while, please cover the chair to prevent from dust.
- Pets sometimes chew on power cord, so try to keep the pets away from the chair even when it is stored.

CAUTION Do not install the chair where it can come into contact with direct heat, such as direct sunlight or a radiator.

Care and Maintenance

CAUTION ● When cleaning the chair, be sure to unplug the power plug from the wall outlet first. Do not plug or unplug the power plug with wet hands.

[Caring for the Back Pad and Covers]

Soak a piece of soft cloth in lukewarm water in which a small amount of mild detergent has been added. Then thoroughly wring it out before cleaning off any spots or grime.

CAUTION ● The fabric color can change if the chair is in contact with vinyl covers or wall coverings for a prolonged period of time.

[Caring for Other Plastic Parts]

Soak a piece of soft cloth in lukewarm water in which a small amount of mild detergent has been added, thoroughly wring it out, and clean off any spots or grime. Then wipe with a dry cloth so that no detergent residue remains.

CAUTION ● Do not use alcohol-based or other household cleaning products or abrasives because they can scratch or crack the finish or cause color to fade.

[Caring for the Remote Control]

Wipe dirt off the remote control with a dry cloth.

CAUTION ● Never use a moistened or wet cloth. May result in a mechanical failure.

Product specification

Type	LC3200-S
Power	AC120V~ 60Hz
Rated power	120W
Rated time	15 minutes
Dimension	Approx 29.33in*54.17in(72.28in)*44.49 (30.51) in (W*L*H) The length and height figures in () represent the dimensions when fully reclined.
Weight	About 149.91lb
Outer material	PU leather
Accessories	Main power cord, Controller, manual book, warranty card, accessories bag for installation

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this device.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

-Reorient or relocate the receiving antenna.

-Increase the separation between the equipment and receiver.

-Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

-Consult the dealer or an experienced radio/TV technician for help.

To satisfy FCC / IC RF exposure requirements, a separation distance of 20 cm or more should be maintained between the antenna of this device and persons during device operation.

To ensure compliance, operations at closer than this distance is not recommended.