
MAGNITONE
LONDON

XO LIGHTS
OUT

Tri Colour LED
Micro-Sonic Silicone
Cleansing Brush

INSTRUCTION
MANUAL

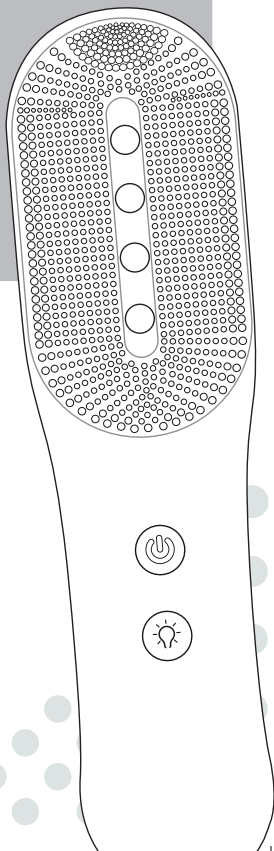Professional product.
Read the instructions of use carefully.
Keep the box and user manual for future reference.

CONTENTS:

1 x Micro-Sonic Cleansing Brush
1 x Additional Charging Socket Cover
1 x USB Charging Cable

TECHNICAL SPECIFICATION

Product: MAGNITONE XO LightsOut Tri-Colour LED
Micro-Sonic Silicone Cleansing Brush

Battery: 3.7V
Output Voltage: 5V
Model Number: MLEDX02G

© Magnitone 2022
Registered trademark, owned by:
HSNF Ltd., London, N10 2JE, U.K.
HSNF IE Ltd., Dublin 7, D07 P4AX, Ireland

Made in China

CUSTOMER SERVICE
help@magnitone.co.uk
www.magnitone.co.uk

 [magnitonelondon](https://www.instagram.com/magnitonelondon)

 [MagnitoneSkincare](https://www.facebook.com/MagnitoneSkincare)

 [magnitoneofficial](https://www.tiktok.com/magnitoneofficial)

WARNINGS

Trust us, you want to read these

DO NOT USE IF YOU:

- You are under the age of 14.
- Have scarred, broken, inflamed or infected skin.
- Have eczema, dermatitis or another skin condition.
- Have extreme sunburn.
- Are taking medication that can cause skin sensitivity.
- You have a pacemaker.
- Close supervision is necessary when this product is used by, on or near children or physically or mentally challenged persons.
- Do not use on or around the delicate eye area.
- Strictly for external use only.

TECHNICAL - DO NOTS:

- Always disconnect the power from the USB source before use.
- Keep away from water (including baths, showers, and sinks).
- Wipe the device with a dry cloth ONLY, do not wash it with water directly.
- Do not use the device while bathing or showering.
- Do not charge near flowing water.
- Keep out of the reach of children and animals.
- Discontinue using the device immediately if you suspect the device to be malfunctioning, or not operating properly.
- Do not reach for a charging cable that has fallen into water.
- Turn off power at the mains before retrieving the unit or it may result in electric shock and / or death.
- Do not attempt to disassemble or repair any part of the product as this may cause damage, harm, malfunction, electrical shock, fire or injury.
- The battery of the device is not removable by the user. If the device is to be disposed of, do not put it into a fire. Dispose of the device as per local regulations.
- Do not use the device if you see any signs of internal battery leakage.
- Do not misuse the product - used only for its intended purpose.
- Never operate the brush if it has been damaged.
- Do not use near to shortwave or microwave equipment. Other electrical equipment in the vicinity may be affected by the usage of this product.
- LightsOut is intended to be stored and used at room temperature and normal humidity.
- Do not store LightsOut in the sun or on a hot surface or in a wet room such as a bathroom. High temperatures may be damaging.
- The device should be operated, transported and stored in the manner specified in this manual.

EN

GETTING STARTED

YOUR FIRST TIME:

'LightsOut' is a non-invasive LED photon Cleansing Brush, intended to cleanse skin with target specific skin concerns, using red, blue and yellow light therapy.

Please note, this is not a medical device and should not be used to treat any medical conditions. LightsOut is designed to be used only by individuals in good health. If this is not the case for you, please consult your doctor before use.

If you experience any abnormal skin irritation, discontinue use immediately.

LightsOut is intended for cosmetic use on the face as described in the Instructions of Use. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty contacting elements, or any other improper application is neither the responsibility of HSNF Ltd nor its affiliates.

DISPOSAL:

This product is classed as Electrical or Electronic equipment. Please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.



HOW OFTEN:

Suitable for daily use.

EN


PARTS AND ACCESSORIES

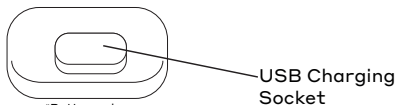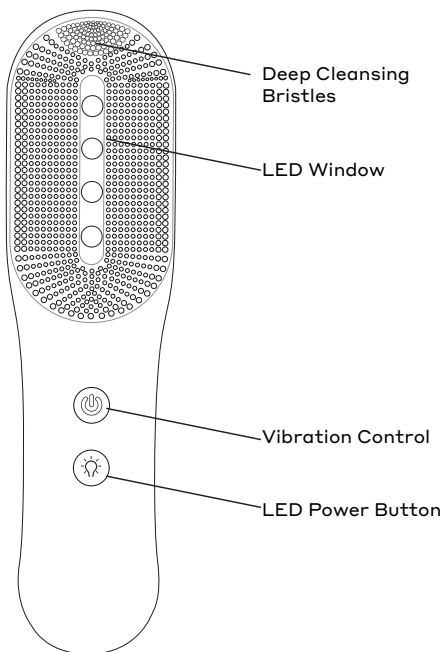*Bottom edge

Always make sure to insert the tab after charging to ensure waterproof seal.

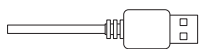1 x USB Charger

LIGHT THERAPIES:

What's the best light treatment for you? The colours of light penetrate the skin at different levels, offering a variety of benefits and treatments:

Blue Light (476-495nm)


Blue light therapy helps to control sebum production. It also removes *Propionibacterium acnes* and improves the environment in the hair follicle. Regular use can treat dark spots and reduce the appearance of scars.

Yellow Light (590-620nm)

Yellow light therapy can be used to fade hyperpigmentation signs. Improves cell oxygenation. Regular use leaves the skin smooth and glowing.


Red Light (620-750nm)

Red light therapy promotes skin natural renewal. Improves the appearance of fine lines and wrinkles. Enhances skin radiance and glow.

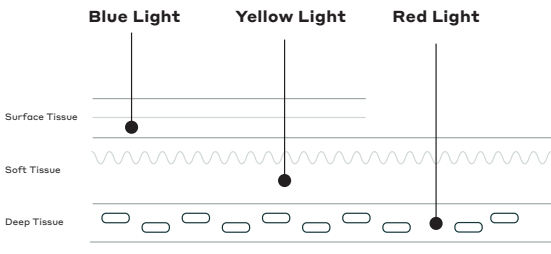HOW TO USE


GETTING STARTED:

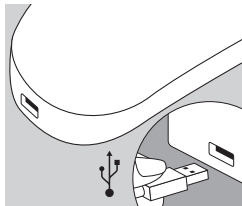
1) Make sure your device is fully charged! To charge your brush, pull out the small rubber cover on the bottom of the device, and connect the USB charger. ALWAYS make sure the cover is secured back onto the device once charging is complete. A full charge will take 3 hours.

2) Before you begin, clean your face to remove makeup, surface dirt and pollution. Have you heard of our amazing WipeOut range?! Both our cleansing SWIPES and Cloths remove make-up with just water! Check them out at magnitone.co.uk.

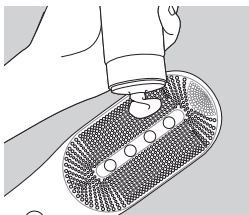
3) Once your skin is clear of makeup and surface dirt, apply your cleanser onto the brush head or directly onto your skin. XO LightsOut is compatible with any cleanser.

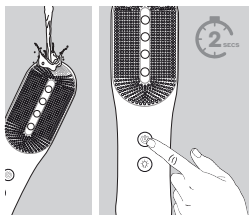4) Wet the head of the brush slightly, and press the top button down for 2 seconds to activate the micro-sonic vibrations. Continue to press the button, to select your preferred speed, the faster the speed the deeper the vibrations!

HOW TO USE

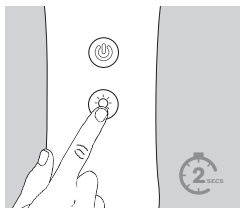
5) Press the lower button down for 2 seconds to activate the LED mode, continue to press the button to select your desired LED treatment. Both modes can operate at the same time.

6) Don't scrub! Let the micro-sonic vibrations do all the work! Glide the brush gently over your face and neck (avoiding the eye-area) in small slow circular motions. **The device will automatically switch-off after 2 minutes.** To switch off **before 2 minutes** press and hold the operating buttons.

*Some devices do not have 2 minute auto-shut off timer feature.

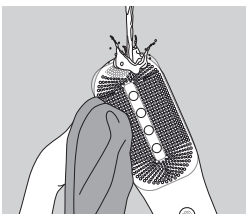7) Now rinse the brush and pat dry with a dry, clean cloth.

8) Complete the routine and apply your usual creams and potions.

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