# MAGNITONE

LONDON

Since 2013, it has been our mission to energise and condition skin naturally; no knives, needles or baffling science! We create quality devices at affordable prices, making **#NoFilterSkincare** accessible to all.

MAGNITONE  
LONDON

## BAREFACED3

### Vibra-Sonic™ Facial Cleanse + Massage Brush

Read the instructions of use carefully. Keep the box and user manual for future reference.

#### THIS MANUAL INCLUDES

- TECHNICAL SPECIFICATION
- WHAT'S IN THE BOX
- REGISTER YOUR WARRANTY
- WARNINGS
- RECYCLING + DISPOSAL
- CLEANING + STORAGE
- CHARGING YOUR DEVICE
- PARTS + ACCESSORIES
- HOW TO GO BAREFACED!
- YOUR BRUSH HEADS

#### TECHNICAL SPECIFICATION

**Product:** MAGNITONE BareFaced3 Vibra-Sonic Cleanse + Massage Brush

**Product Code:** ML05G

**Input:** 5V  200mA 

**Battery:** Ni-MH 2.4V 600mAh

**Waterproof Rating:** IP7

© MAGNITONE 2023

Registered trademark, owned by:

HSNF Ltd., London, N10 2JE, UK. HSNF IE Ltd., Dublin 7, D07 P4AX, Ireland.

Made / Fabriqué / Prodotto / Hecho / Hergestellt / Wyprodukowano /  
Vervaardigd in China.

## WHAT'S IN THE BOX

- 1 x BareFaced3 Vibra-Sonic Cleanse + Massage Brush
- 1 x SiliSoft Silicone Cleansing Head
- 1 x ToneltUp Massage Head
- 1 x Bathroom Caddy
- 1 x USB Charging Cable

Scan the  
QR code

## REGISTER YOUR WARRANTY

Your device is covered by a 12 month warranty from your date of purchase. You must register your device to validate your warranty.

Scan the QR code to register your warranty and unlock an extra 12 months cover. That's 2 years in total!

This warranty does not cover fault or defect caused by:

- Faults arising from failure to observe operating and cleaning instructions
- Careless handling or misuse
- Transit damage
- Alterations, repairs or adjustments undertaken by anyone other than our qualified product engineers.

## CUSTOMER SERVICE

help@magnitone.co.uk  
www.magnitone.co.uk

## WARNINGS

BareFaced3 is a deep cleansing and massaging facial brush. Powered by Vibra-Sonic™ technology; BareFaced3 helps to deep cleanse by removing dirt and impurities, and stimulating blood flow. Some may experience breakouts or other adverse skin changes after the first few days using your brush but do not panic! This is completely normal and is your skin's way of adjusting to a new and deeper method of cleansing. Skin should re-balance itself after 2 weeks.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

## CONTRAINDICATIONS

### DO NOT USE IF:

- You have broken skin, irritated skin, varicose veins, rashes, spots and moles.
- You have eczema, dermatitis or other skin conditions.
- You have extreme sunburn.
- You are taking medication that can cause skin sensitivity.

## SAFETY CAUTIONS

- DO **NOT** use on or near the delicate eye area.
- DO **NOT** drop or misuse the product.
- Do **NOT** leave the charger plugged in once the brush is fully charged as this may damage the appliance.
- DO **NOT** disassemble or repair any part of the appliance.
- Children shall not play with the appliance.

## RECYCLING + DISPOSAL

Your MAGNITONE BareFaced3 must not be disposed of with other household or commercial waste at the end of its working life. When you have no further use for your product, please dispose of it according to your local authority's recycling procedure.

## CLEANING + STORAGE

Make sure your device is disconnected from the USB charging cable before cleaning. Rinse the brush under running water after every use, and let the brush air dry or wipe down with a dry cloth.



## CHARGING YOUR DEVICE

Your BareFaced3 may have a little bit of charge when you take it out of the box but, to maximise battery life, you should give your brush a full charge before you use it.


To charge, simply attach the magnetic end of the USB charger to the bottom of your device (A) and plug in the USB to a computer or USB charging port/plug.

The lights will flash during charging and will turn solid once fully charged. (B) If the lights turn red when the brush is in use, this means the battery is low and you need to charge it.

A full charge should last up to 1 month, depending how often you use it.
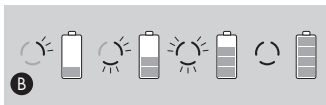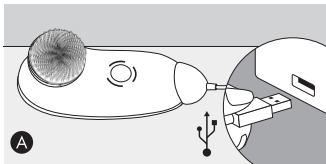
## PARTS + ACCESSORIES
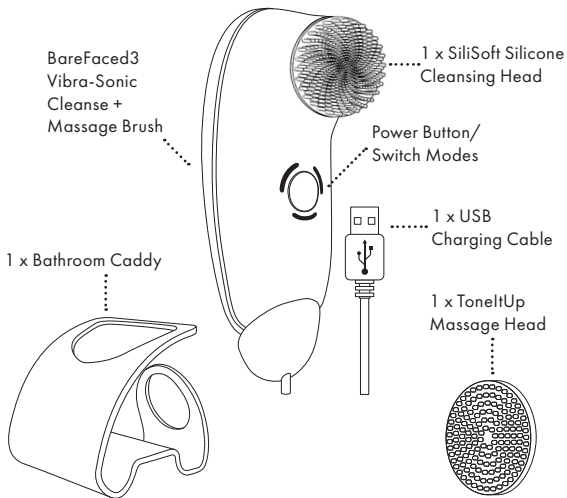
## LIGHT INDICATOR

Gentle Cleanse

Deep Cleanse




PulseLift Massage

## HOW TO GO BAREFACED!



### TO CLEANSE:

1. First cleanse first - remove all makeup. Try WipeOut Microfibre Cleansing Cloths to remove makeup with just warm water. (C)
2. Apply your preferred cleanser either to your face or directly onto the SiliSoft brush head. (D)
3. Choose your mode. Press once for a Gentle Cleanse or twice for a Deep Cleanse.
4. Hold the BareFaced3 to your skin and use small circular motions to help 'wobble' out any dirt.
5. BareFaced3 will run for one minute, with a beep every 20 seconds to remind you to move on to the next area. Start with the cheeks and chin, then move to the T-zone and finish with the neck. (E)
6. When you're done, rinse your BareFaced3 to remove any cleanser (F) and place it in the bathroom caddy to drip dry. (I) Rinse your face with warm water to remove any cleanser (or use your WipeOut cloth) and carry on with the rest of your skincare routine.
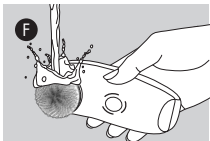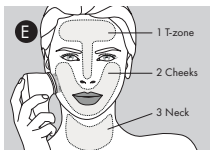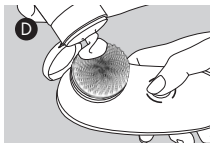
## TO MASSAGE:


1. When you reach the step in your skincare routine to apply your serums or moisturiser, give them an extra boost with your BareFaced3 in PulseLift massage mode. Attach the ToneltUp Massage Brush Head to your BareFaced3. (G)
2. Apply your preferred serum or moisturiser either to your skin or directly on the brush head. (D)
3. Press the button 3 times to enter PulseLift massage mode.
4. BareFaced3 will run for one minute, with a beep every 20 seconds to remind you to move on to the next area. Start with the cheeks and chin, then move on to the T-zone and finish with the neck. (H)
5. When you're finished, rinse the ToneltUp Head to remove any serum/moisturiser (F) and place your BareFaced3 in the bathroom caddy to drip dry. (I)
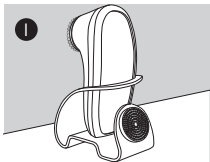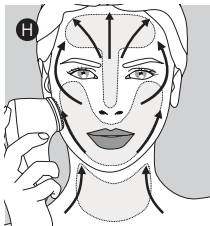
## YOUR BRUSH HEADS


### NEW SILISOFT™ SILICONE CLEANSE HEAD

Suitable for all skin types, the super soft silicone bristles and unique sweep-swirl design give a deep yet gentle cleanse. Unlike traditional plastic bristles, the high-grade silicone means SiliSoft will keep its shape and won't need replacing. Better for your face, for your wallet and for the planet!

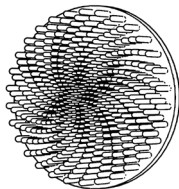
### TONEITUP MASSAGE HEAD

When you reach the step in your skincare routine to apply your serums or moisturiser, give them some extra oomph with your BareFaced3 in PulseLift massage mode. The short, firm nodules of the ToneItUp head will help to push the products into the skin, while boosting micro-circulation.
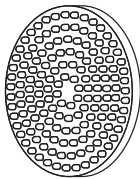
## HOW TO CHANGE THE BRUSH HEADS

1. To switch between brush heads, hold the BareFaced3 firmly in one hand.
2. Then, pull the brush head gently at an angle (not straight) and it'll pop off. (J)
3. Fix the desired head on by pushing until it clicks into place.
4. Store the 'other' head in your Bathroom Caddy until you need it next. (K)
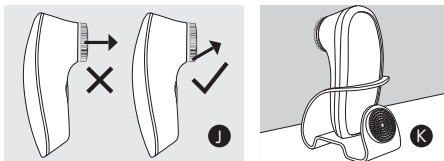
## HOW TO USE THE POWER TRAVEL LOCK

Avoid the awkward moment when your bag starts vibrating as you go through airport security! To activate the Power Travel Lock, press and hold the power button for 5 seconds to turn it on. The red light will flash 3 times when it is active. Repeat to unlock.

## SHARING IS CARING!

We love seeing your pictures and videos. Take a photo, make a video, tag us on social and we'll give it a share!

magnitoneofficial

magnitonelondon

magnitoneskincare

youtube.com/c/MAGNITONELondon

Make sure to tag us in and use **#NoFilterSkincare**.

## OUR PROMISE

We believe that all skin is beautiful.

Lines, blemishes, freckles, pimples - they're what make us unique and make us human. We want to represent real people and real skin, so we promise that all of our imagery is 100% un-retouched.

Our devices and skincare are simply tools to help you feel confident in the gorgeous skin you're in.

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