

# MAGNITONE

LONDON

Since 2013, it has been our mission to energise and condition skin naturally; no knives, needles or baffling science! We create quality devices at affordable prices, making **#NoFilterSkincare** accessible to all.

**MAGNITONE**  
LONDON

## **FACEROCKET®**

### 5-in-1 Facial Firming + Toning Device

Read the instructions of use carefully. Keep the box and user manual for future reference.

#### **THIS MANUAL INCLUDES**

- TECHNICAL SPECIFICATION
- WHAT'S IN THE BOX
- REGISTER YOUR WARRANTY
- WARNINGS
- CLEANING + STORAGE
- RECYCLING + DISPOSAL
- CHARGING YOUR DEVICE
- PARTS + ACCESSORIES
- TECHNOLOGIES
- HOW TO USE FACEROCKET®
- LIFT OFF SUPERCONDUCTIVE GEL

#### **TECHNICAL SPECIFICATION**

**Product:** MAGNITONE FaceRocket®  
5-in-1 Facial Firming + Toning Device

**Product Code:** MFR01G

**Input:** 5.0V  1000mA 

**Battery:** Li-Ion 3.7V 2000mAh



#### **© MAGNITONE 2023**

Registered trademark, owned by:  
HSNF Ltd., London, N10 2JE, UK. HSNF IE Ltd., Dublin 7, D07 P4AX, Ireland.

Made / Fabriqué / Prodotto / Hecho / Hergestellt / Wyprodukowano /  
Vervaardigd in China.

#### **WHAT'S IN THE BOX**


- 1 x FaceRocket® 5-in-1 Facial Firming + Toning Device
- 1 x LiftOff Superconductive Prepare + Glide Gel 60ml
- 1 x USB Charging Cable

#### **REGISTER YOUR WARRANTY**

Your device is covered by a 12 month warranty from your date of purchase.

You must register your device to validate your warranty.

Scan the QR code to register your warranty and unlock an extra 12 months cover. That's 2 years in total!

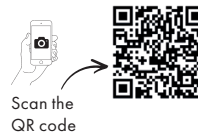

This warranty does not cover fault or defect caused by:

- Faults arising from failure to observe operating and cleaning instructions.
- Careless handling or misuse.
- Transit damage.
- Alterations, repairs or adjustments undertaken by anyone other than our qualified product engineers.

#### **CUSTOMER SERVICE**

help@magnitone.co.uk  
www.magnitone.co.uk

### WARNINGS

Get your skin looking and feeling out of this world with FaceRocket® - the ultimate facial firming and toning device. Clinically proven to improve skin firmness by 29%, dermal density by 35% and reduce wrinkles by 10%, FaceRocket® delivers real results in just 4 weeks.\*

FaceRocket® is intended for facial stimulation and is for cosmetic use only. This is not a medical device and should not be used to treat any medical conditions. Any consequence resulting from misuse, application to other parts of the body, connection to improper voltage sources, dirty contacting elements, or any other improper application or use is neither the responsibility of HSNF Ltd nor its affiliates.

Some users may find the sensation quite strong, especially in comparison to other microcurrent products they have used. A tingling sensation (like a small electric shock) and some involuntary muscle twitching may occur while using MAGNITONE FaceRocket®. Lowering the mode setting may reduce or eliminate this sensation. If this is uncomfortable, we suggest using a lower strength setting. Normally this sensation subsides after the first few uses. If you experience any pain or discomfort, stop using the device immediately.

*\*Independent study of 32 individuals aged 33-50. Average result based on regular use over 4 weeks.*



- For use on the face only.
- Adverse effects may occur if the instructions are not followed.
- Prolonged or repeated application on one site may potentially damage the skin.
- **DO NOT** use over the eyes (the circular muscle within the orbital rim), on or around the thyroid glands (midline of the neck, in-front of and either side of the wind-pipe) on the temples (above the ear), breasts or groin area.
- **DO NOT** apply stimulation across your chest as the introduction of an electrical current into the chest may cause rhythmic disturbances to your heart, which can be lethal.
- **DO NOT** use the appliance in conjunction with flammable based cosmetic products, unless they are recommended by the manufacturer.
- When using on the jaw, care should be taken to minimise treatment directly over the midline of the neck, specifically avoiding the thyroid.
- There may be a perception of flashing lights during the use of your FaceRocket®, resulting from the stimulation of the eye area. If you notice this condition constantly without the use of your device, discontinue use and seek medical advice.
- Some users may experience a metallic taste in their mouth during the use of the device, resulting from stimulation and varying sensitivity around the mouth. If you notice this condition constantly without the use of your device, discontinue use and seek medical advice.
- Some users may experience a metallic taste in their mouth during the use of the device, resulting from stimulation and varying sensitivity around the mouth. If you notice this condition constantly without the use of your device, discontinue use and seek medical advice.

**CONSULT YOUR DOCTOR BEFORE USE IF YOU HAVE ANY UNDERLYING SKIN CONDITIONS.**

## CONTRAINDICATIONS

### DO NOT USE IF:

- You have areas where the integrity of the skin is damaged (such as burns, blisters, lesions, incision, wounds, active skin disease or infection); are suffering from skin disease or inflammation of the skin, have sensitive skin or an allergic constitution.
- You have a history of epilepsy or seizures.
- You have had recent surgery or other surgical procedures.
- You have any electronic implanted device such as cardiac pacemaker.
- You use a body-worn medical device, such as an insulin pump.
- You have any active cancer or have recently recovered from or had treatment for any suspicious or cancerous lesions.
- You are pregnant.
- You are under the care of a doctor for any ongoing medical condition, consult them before using this device.
- You have a known history of carotid artery disease, stroke or transient ischemic attack (TIA), carotid stenosis, unstable blood pressure, easy fainting, thyroid disease, or prior surgical procedures in the neck area.
- You have a tendency for bleeding following acute injury or breaks in bones.
- You are wearing implants not compatible with its use.
- This appliance shall not be used by or on persons up to and including 18 years of age.

## SAFETY CAUTIONS

- **DO NOT** immerse the device in water.
- **DO NOT** drop or misuse the product.
- **DO NOT** handle the device or charger with wet hands.
- **DO NOT** operate the device if it has been damaged.
- **DO NOT** disassemble or attempt to repair the product.
- **DO NOT** place the product in direct sunlight or expose the battery to high temperatures, or heat-generating equipment. Other electrical equipment in the vicinity may be affected by the usage of this product.
- **DO NOT** bypass or remove the inherent safety device of the appliance.
- Make sure your device is disconnected from the USB charging cable before use.
- The appliance **MUST** only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Strictly for external use.

## CLEANING + STORAGE

Cleaning and user maintenance shall not be made by children without supervision. Make sure your device is disconnected from the USB charging cable before cleaning. Use a slightly damp towel to clean the device and dry with a soft cloth.

## RECYCLING + DISPOSAL

Your MAGNITONE FaceRocket® must not be disposed of with other household or commercial waste at the end of its working life. When you have no further use for your product, please dispose of it according to your local authority's recycling procedure.




## CHARGING YOUR DEVICE

Before your first use, charge your FaceRocket®. Plug into any USB Socket on a charge plug or computer to charge your device. When the device is fully charged, the lights will go green.

How long the FaceRocket® battery charge lasts completely depends on how long you use it for, on what mode and on what level (for example, radiofrequency uses a lot more power than cryotherapy). If you are using your device daily, you will likely need to charge it at least once a week.

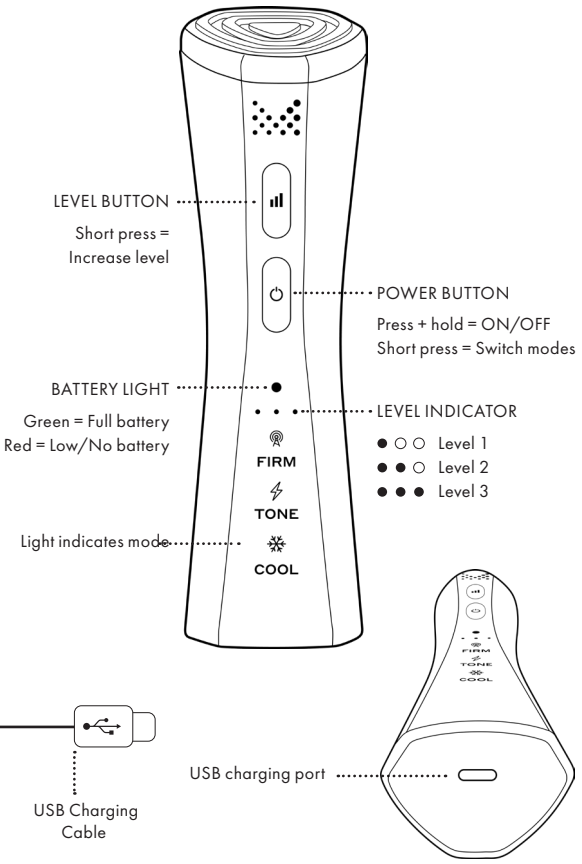


PARTS + ACCESSORIES



TECHNOLOGIES

FaceRocket® uses a combination of technologies, all designed for a firmer, brighter and tighter look + feel.

<p><b>RADIOFREQUENCY</b></p>	<p>Radio waves heat below the skin's surface to stimulate collagen and elastin production, resulting in firmer skin.</p> <p>42°C</p>
<p><b>MICROCURRENT</b></p>	<p>Small pulses of electricity stimulate the facial muscles to contract, strengthening and toning them, helping to lift the face.</p> <p>300 - 500 mA</p>
<p><b>RED + GREEN LED LIGHT THERAPY</b></p>	<p>Red light waves penetrate below the skin's surface to stimulate collagen + elastin production.</p> <p>Green light waves penetrate below the skin's surface encourage blood flow, helping skin to appear more radiant.</p> <p>Red: 630 - 680 nm / Green: 475 - 545 nm</p>
<p><b>CRYOTHERAPY</b></p>	<p>Cool and soothe the skin, reducing puffiness and encouraging lymphatic drainage.</p> <p>4°C</p>
<p><b>SONIC MASSAGE</b></p>	<p>Gentle vibrations stimulate below the skin's surface, encouraging blood flow and lymphatic drainage.</p>

MODE	TECHNOLOGIES	LEVELS
 <b>FIRM</b>	<ul style="list-style-type: none"> <li>• Radiofrequency</li> <li>• Red LED</li> </ul>	n/a
 <b>TONE</b>	<ul style="list-style-type: none"> <li>• Microcurrent</li> <li>• Red LED</li> </ul>	<ol style="list-style-type: none"> <li>1. ~ 300 mA</li> <li>2. ~ 400 mA</li> <li>3. ~ 500 mA</li> </ol>
 <b>COOL</b>	<ul style="list-style-type: none"> <li>• Cryotherapy</li> <li>• Green LED</li> <li>• Sonic massage</li> </ul>	<ol style="list-style-type: none"> <li>1. No massage</li> <li>2. Massage 1</li> <li>3. Massage 2</li> </ol>


### HOW TO USE FACEROCKET®

#### BEFORE YOU BEGIN

Make sure you have read and understood all warnings, cautions and instructions before you begin.

Take your 'Before' photo! Remember to take a before photo of your cleansed face, without makeup, in a well-lit room. Take photos every week in the same place and position to track your progress. You may not notice changes daily, but regular selfies will help you to see the difference!

Share your journey with us on social - use the hashtag **#FaceRocket**

-  [magnitoneofficial](#)
-  [magnitonelondon](#)
-  [magnitoneskincare](#)
-  [youtube.com/c/MAGNITONELondon](#)


















#### HOW OFTEN

For best results, follow our 28 Day Face Programme! This daily calendar will guide you through your first 28 days of treatments. Use the calendar on the next page or scan the QR code to download and print your own.



















When your 28 days have finished, maintain the results by using FIRM Mode at least once a week and TONE mode 1-3 times a week. COOL mode can be used daily, as required.

# FACEROCKET®

<b>0</b>  <b>BEFORE YOU BEGIN</b>	<b>1</b>  ● ○ ○	<b>2</b>  ● ○ ○	<b>3</b>  ● ○ ○
<b>KEY</b>  TAKE A SELFIE  <b>FIRM mode</b> 5-10 minutes  <b>COOL mode</b> 3-5 minutes  <b>COOL mode</b> 1-3 minutes ● ○ ○ Level 1 ● ● ○ Level 2 ● ● ● Level 3	<b>8</b>  ● ● ○	<b>9</b>  ● ● ○	<b>10</b>  ● ● ●
	<b>15</b>  ● ● ●	<b>16</b>  ● ● ●	<b>17</b>  ● ● ●
	<b>22</b>  ● ● ●	<b>23</b>  ● ● ●	<b>24</b>  ● ● ●


# 28 DAY FACE PROGRAMME

<b>4</b>  ● ○ ○	<b>5</b>  ● ○ ○	<b>6</b>  ● ○ ○	<b>7</b>  ● ○ ○
<b>11</b>  ● ● ○	<b>12</b>  ● ● ○	<b>13</b>  ● ● ●	<b>14</b>  ● ● ○
<b>18</b>  ● ● ●	<b>19</b>  ● ● ●	<b>20</b>  ● ● ●	<b>21</b>  ● ● ●
<b>25</b>  ● ● ●	<b>26</b>  ● ● ●	<b>27</b>  ● ● ●	<b>28</b>  ● ● ●

## FIRM


### RADIOFREQUENCY + RED LED LIGHT THERAPY

Firm + tighten the skin with Radiofrequency + Red LED, working together below the surface to stimulate collagen + elastin production. Do not use FIRM Mode more than once every 3 days.

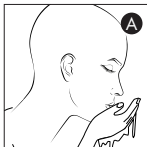

1. Before you begin, make sure your skin is fully cleansed - free of makeup, surface oils and daily grime (A). (Men should shave before treatment as facial hair can interfere with conductivity).

*IMPORTANT: Any cream, lotion, self-tan, or other oil-based products on the skin can block the radiofrequency.*

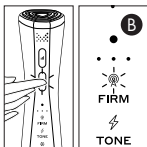


*We recommend a full double cleanse; try our WipeOut Amazing Microfibre Cleansing Cloths to get rid of makeup and surface dirt, then supercharge your second cleanse with one of our award-winning cleansing brushes!*




2. Press and hold the power button to turn on FaceRocket®. It will automatically start in FIRM Mode (B).

3. While the device warms up, apply a generous layer of Superconductive Gel to your skin (C). You can apply to the whole face or section by section, and re-apply throughout the treatment as necessary.


*IMPORTANT: You must use a Superconductive Gel to help the radiofrequency reach below the skin's surface. Without it, the treatment will not work.*

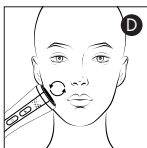

4. Starting at the chin, place the head of FaceRocket® flush to your skin and gently move in small circular motions (D). After 5 seconds you will hear a beep, this means you need to move on to the next section of your face. Always keep FaceRocket® moving in FIRM Mode. Do NOT hold in one place.

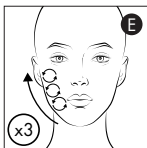


5. Work along the jawline up to the ear on one side of your face (E), then repeat this section of your face 2 more times. So you're doing 3 x 'repetitions' in total.

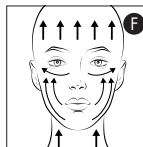

6. Using this same technique move on to the next section of your face and continue to work around your face (F):

- Corner of mouth to ear.
- Under eye/cheekbone to ear (Do NOT use on the eye lid / delicate eye area).
- Forehead.
- Under chin.


Apply more gel as necessary.

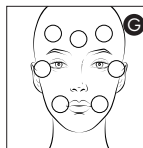

7. Once you have completed your whole face, you can go back and work on any specific problem areas such as smile line, frown lines or crows feet (G). To do this, gently move FaceRocket® in small circular motions for 15 seconds (3 x 5 second beeps).



8. When you're done, remove any excess gel with a clean, damp cloth (H) and complete your skincare routine as normal.

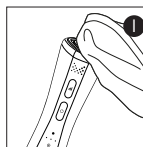

9. Wipe clean the head of your FaceRocket® using a damp cloth (I). Do NOT rinse with water.






## TONE

### MICROCURRENT + RED LED LIGHT THERAPY

Lift, tone and contour with microcurrent. Small pulses of electricity stimulate facial muscles, strengthening + toning them.

In TONE Mode, DO NOT use FaceRocket® over the eyes (the circular muscle within the orbital rim), on or around the thyroid glands (midline of the neck, in-front of and either side of the wind-pipe) on the temples (above the ear), breasts or groin area. Make sure you have read and understood all WARNINGS before you begin.

TONE Mode can be used daily. For best results, use 2-3 times a week.

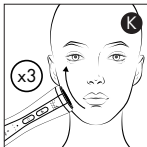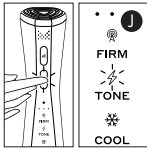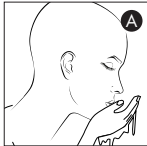

1. Before you begin, make sure your skin is fully cleansed - free of makeup, surface oils and daily grime (A). (Men should shave before treatment as facial hair can interfere with conductivity).




*IMPORTANT: Any cream, lotion, self-tan, or other oil-based products on the skin can block the radiofrequency.*

*We recommend a full double cleanse; try our WipeOut Amazing Microfibre Cleansing Cloths to get rid of makeup and surface dirt, then supercharge your second cleanse with one of our award-winning cleansing brushes!*



2. Press and hold the power button to turn on FaceRocket®. Press the power button again to select TONE Mode (J).  
FaceRocket® in TONE mode will start at Level 1 Microcurrent intensity. Always start at Level 1 and work your way up after a few treatments. Press the Level button to increase the intensity.
3. Apply a generous layer of Superconductive Gel to your skin (C). You can apply to the whole face or section by section, and re-apply throughout the treatment as necessary.

*IMPORTANT: You must use a Superconductive Gel to help the microcurrent reach below the skin's surface. Without it, the treatment will not work.*

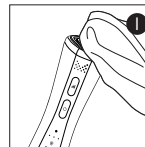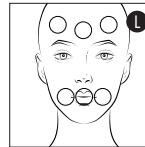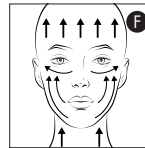


4. Starting at the chin, place the head of FaceRocket® flush to your skin and gently glide along the jaw, up to your ear (K), until you hear a beep. Repeat this section of your face 2 more times. So you're doing 3 x 'repetitions' in total.
5. Using this same technique move on to the next section of your face and continue to work around your face (F):
  - Corner of mouth to ear.
  - Under eye/cheekbone to ear (Do NOT use on the eye lid / delicate eye area).
  - Forehead.
  - Under chin.

Apply more gel as necessary.
6. Once you have completed your whole face, you can go back and work on any specific problem areas such as lips, smile lines or frown lines (L). To do this, hold FaceRocket® in one place for up to 15 seconds (3 x 5 second beeps).
7. When you're done, remove any excess gel with a clean, damp cloth (H) and complete your skincare routine as normal.
8. Wipe clean the head of your FaceRocket® using a damp cloth (I). Do NOT rinse with water.

COOL CRYOTHERAPY + GREEN LED LIGHT THERAPY + SONIC MASSAGE  


Cooling Cryotherapy soothes and reduces puffiness while Green LED and Sonic Massage stimulate blood flow, help lymphatic drainage and improve radiance.

COOL Mode can be used daily, as often as you like or need. It can be used in conjunction with the other 2 modes. For best results, use following the FIRM or TONE treatments.

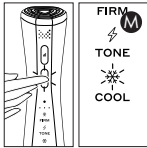


1. Make sure your skin is clean and makeup free.
2. Press and hold the power button to turn on FaceRocket®. Press the power button again twice to select COOL Mode (M). Press the Level button to choose your massage setting:

- a. Cryotherapy + Green LED Light Therapy only.
- b. Cryotherapy + Green LED Light Therapy + Sonic Massage 1.
- c. Cryotherapy + Green LED Light Therapy + Sonic Massage 2.

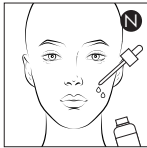


3. Cooling mode does not require pre-application of any product however, for best results, apply a generous layer of a nourishing serum (N) such as Hyaluronic Acid or Vitamin C (not included) to help FaceRocket® to glide and to hydrate your skin. This is optional and Cooling Mode can be used without any product at all.

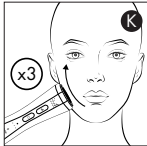

4. Starting at the chin, place the head of FaceRocket® flush to your skin and gently glide along the jaw, up to your ear (K), until you hear a beep.

*Optional: Repeat this section of your face 2 more times. So you're doing 3 x 'repetitions' in total.*

5. If you have any areas you want to focus on, you can hold the device in place for up to 15 seconds - still or in small circular motions.
6. When you're done, any serum should be absorbed. Complete your skincare routine with SPF and moisturiser.
7. Wipe clean the head of your FaceRocket® using a damp cloth (L). Do NOT rinse with water.

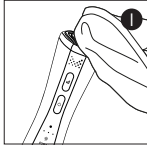MAGNITONE  
LONDONLIFT OFF SUPERCONDUCTIVE  
PREPARE + GLIDE GEL

Effortlessly glide and ensure the proper delivery of microcurrent pulses deeper to the skin with this superconductive, water-based gel.

## HOW TO USE

Test your skin for sensitivity to LiftOff Prepare + Glide Gel. Rub a pea-sized amount of gel on the underside of your forearm and leave for 15-20 minutes before washing it off. If a reaction (rash or redness) develops, immediately contact MAGNITONE for recommendations on alternative products. You should not begin treatment with MAGNITONE FaceRocket® until a suitable gel is found.

Ensure your skin is clean and dry. Apply a generous layer of the gel to all your treatment areas, immediately before using your FaceRocket® device. Wipe or rinse off after use.

## INGREDIENTS

Aqua (Water), Carbomer, Glycerin, Methyl Paraben, Propyl Paraben, Propylene Glycol, Sodium Hydroxide

## CAUTIONS

For external use only. Not to be used on sore or damaged skin. In the unlikely event of rash or irritation, discontinue use. Avoid contact with eyes.





## STORAGE

Store in a cool dry place away from sunlight, with the lid firmly closed.



### SHARING IS CARING!

We love seeing your pictures and videos. Take a photo, make a video, tag us on social and we'll give it a share!

-  [magnitoneofficial](#)
-  [magnitonelondon](#)
-  [magnitoneskincare](#)
-  [youtube.com/c/MAGNITONELondon](#)

Make sure to tag us in and use **#NoFilterSkincare**.

### OUR PROMISE

We believe that all skin is beautiful.

Lines, blemishes, freckles, pimples - they're what make us unique and make us human. We want to represent real people and real skin, so we promise that all of our imagery is 100% un-retouched.

Our devices and skincare are simply tools to help you feel confident in the gorgeous skin you're in.

**xoxo**

Look out for the stamp

