MAGNITONE

LONDON

Since 2013, it has been our mission to energise and condition skin naturally; no knives, needles or baffling science! We create quality devices at affordable prices, making **#NoFilterSkincare** accessible to all.

MAGNITONE
LONDON

GETLIT

Tri Colour LED Face Mask

Read the instructions of use carefully. Keep the box and user manual for future reference.

THIS MANUAL INCLUDES

- TECHNICAL SPECIFICATION
- WHAT'S IN THE BOX
- REGISTER YOUR WARRANTY
- WARNINGS
- CLEANING + STORAGE
- RECYCLING + DISPOSAL
- GETTING STARTED
- PARTS + ACCESSORIES
- HOW TO GETLIT!

TECHNICAL SPECIFICATION

Product: MAGNITONE GETLIT Tri Colour LED Face Mask

Product Code: MLEDGL01W

Input: 5V  350mA 

Battery: Li-Ion 3.7V 155mAh

© **MAGNITONE 2021**

Registered trademark, owned by:

HSNF Ltd., London, N10 2JE, UK. HSNF IE Ltd., Dublin 7, D07 P4AX, Ireland.

Made / Fabriqué / Prodotto / Hecho / Hergestellt / Wyprodukowano /
Vervaardigd in China.

WHAT'S IN THE BOX

1 x GetLit Tri Colour LED Face Mask

1 x USB Charging Cable

REGISTER YOUR WARRANTY

Your device is covered by a 12 month warranty from your date of purchase. You must register your device to validate your warranty.

Scan the QR code to register your warranty and unlock an extra 12 months cover. That's 2 years in total!

This warranty does not cover fault or defect caused by:

- Faults arising from failure to observe operating and cleaning instructions
- Careless handling or misuse
- Transit damage
- Alterations, repairs or adjustments undertaken by anyone other than our qualified product engineers.

Scan the
QR code

CUSTOMER SERVICE

help@magnitone.co.uk

www.magnitone.co.uk

WARNINGS

'GetLit' is a non-invasive LED photon Face Mask, intended to target specific skin concerns, using red, blue and yellow light therapy.

Please note, this is not a medical device and should not be used to treat any medical conditions. GetLit is designed to be used only by individuals in good health. If this is not the case, please consult your doctor before use.

GetLit is intended for cosmetic use on the face as described. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, or any other improper application is neither the responsibility of HSNF Ltd nor its affiliates.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. The appliance **MUST** only be supplied at safety extra low voltage corresponding to the marking on the appliance.

CONTRAINDICATIONS

DO NOT USE IF:

- You have scarred, broken, inflamed or infected skin.
- You have eczema, dermatitis or any other skin condition.
- You have sunburn.
- You are taking medication that can enhance skin sensitivity (e.g. Certain tablets or creams for acne).

SAFETY CAUTIONS

- Do **NOT** use the mask for any other purpose other than cosmetic facial treatment.
- Do **NOT** use this product in the bathroom, near basins or in close proximity with any other vessels containing water.
- Do **NOT** submerge the device in water.
- Do **NOT** use the device while bathing or showering.
- Do **NOT** place any heavy objects on top of the mask.
- Do **NOT** press down, or fold the mask.
- Do **NOT** drop or insert any non-specified object on the device.
- Do **NOT** disassemble the device as this may cause damage, harm, malfunction, electrical shock, fire or injury. There are no user-serviceable parts inside. **Note:** Disassembly voids your warranty.
- Do **NOT** use the device if you see any signs of internal battery leakage.
- The device is intended to be stored and used at room temperature and normal humidity. Do not store GetLit over the sun or on a hot surface. High temperatures may be damaging.
- The device should be operated, transported and stored in the manner specified in this user manual.
- Before use – make sure that the USB wire is disconnected from the power source.
- If the USB cable becomes damaged, it should be replaced.
- If the device begins to overheat, turn the device **OFF** immediately.
- Discontinue use immediately if you suspect the device to be malfunctioning, or not operating properly.

CLEANING + STORAGE

Cleaning and user maintenance shall not be made by children without supervision. Always make sure your USB cable is unplugged before cleaning/use. Use a soft, dry cloth to clean the mask and safely store away, ideally in the original box.

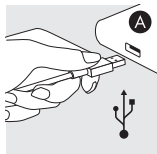
RECYCLING + DISPOSAL

Your MAGNITONE GetLit must not be disposed of with other household or commercial waste at the end of its working life. When you have no further use for your product, please dispose of it according to your local authority's recycling procedure.

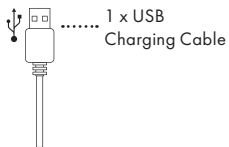
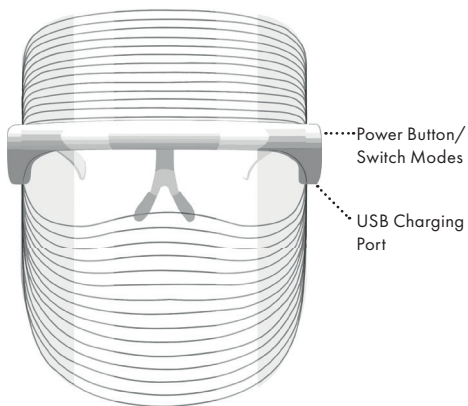
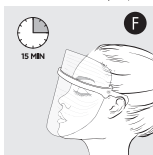
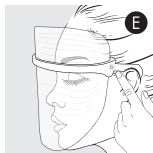
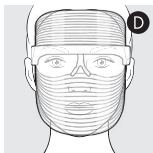
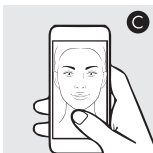
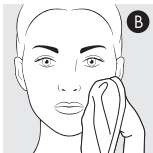
GETTING STARTED

Before you begin, charge your GetLit LED Mask for 2 hours. Plug into any USB socket, on a charge plug or computer, to charge your device (A). Whilst on charge, the internal port will glow red, once fully charged the light will turn blue.

PARTS + ACCESSORIES

HOW TO GETLIT!

1. Make sure your face is clean and free of makeup (B).
2. Take a selfie! Regular selfies are the best way to keep track of any change to your skin (C).
3. Place the mask over your face, making sure that the frames are securely positioned behind your ears (D).
4. Press the Power button down for 3 seconds to turn the device on. The red LED light will automatically activate. Press the Power button to switch between the red, yellow and blue LED light settings (E).
5. Relax! Just chill and let the LED work its magic. The mask will shut off automatically after 15 minutes. If you wish to shut off the mask earlier, simply press and hold the Power button to turn off (F).
6. When you're done, complete your skincare routine with your usual serums, creams or moisturiser.

WHAT'S THE RIGHT LIGHT FOR YOU?

 RED LIGHT 630nm~680nm	An easy rhyme to remember it - 'red for wrinkles'! Red LED light has been shown to improve skin firmness and reduce the appearance of fine lines and wrinkles by stimulating collagen + elastin production.
 YELLOW LIGHT 600nm~605nm	'Yellow for glow'. Yellow LED light has been shown to stimulate blood flow to the skin. More blood flow = increased delivery of oxygen + nutrients to the skin, leaving it feeling more radiant and glowing.
 BLUE LIGHT 420nm~460nm	And finally, 'blue for blemishes'! Blue LED light is widely used to target spots and blemishes. The blue LED light has been shown to have antibacterial properties and make oil producing follicles less active.

**Rather watch a video than read?
Us too!**

Scan the QR
to watch

SHARING IS CARING!

We love seeing your pictures and videos. Take a photo, make a video, tag us on social and we'll give it a share!

magnitoneofficial

magnitonelondon

magnitoneskincare

youtube.com/c/MAGNITONELondon

Make sure to tag us in and use **#NoFilterSkincare**.

OUR PROMISE

We believe that all skin is beautiful.

Lines, blemishes, freckles, pimples - they're what make us unique and make us human. We want to represent real people and real skin, so we promise that all of our imagery is 100% un-retouched.

Our devices and skincare are simply tools to help you feel confident in the gorgeous skin you're in.

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Batterie: Li-Ion 3.7V 155mAh

© MAGNITONE 2021

Registered trademark, owned by:
HSNF Ltd., London, N10 2JE, UK.
HSNF IE Ltd., Dublin 7, D07 P4AX, Ireland.

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