**MAGNITONE**  
LONDON

## LIFTOFF

### Microcurrent Facial Toning + Lifting Device

Read the instructions of use carefully. Keep the box and user manual for future reference.

#### THIS MANUAL INCLUDES

- TECHNICAL SPECIFICATION
- WHAT'S IN THE BOX
- REGISTER YOUR WARRANTY
- WARNINGS
- CLEANING + STORAGE
- RECYCLING + DISPOSAL
- PARTS + ACCESSORIES
- HOW TO LIFTOFF!
- LIFT OFF SUPERCONDUCTIVE PREPARE + GLIDE GEL

#### TECHNICAL SPECIFICATION

**Product:** MAGNITONE LIFT OFF Microcurrent Facial Toning + Lifting Device

**Product Code:** MLF01P

**Input:** 5.0V  250mA 

**Battery:** Li-Ion 3.7V 250mAh

UK  
CA CE

© **MAGNITONE 2021**

Registered trademark, owned by:

HSNF Ltd., London, N10 2JE, UK. HSNF IE Ltd., Dublin 7, D07 P4AX, Ireland.

Made / Fabriqué / Prodotto / Hecho / Hergestellt / Wyprodukowano /  
Vervaardigd in China.

## WHAT'S IN THE BOX

- 1 x LiftOff Microcurrent Facial Toning Device
- 1 x Superconductive Prepare + Glide Gel 60ml €
- 1 x USB Charging Cable

## REGISTER YOUR WARRANTY

Your device is covered by a 12 month warranty from your date of purchase. You must register your device to validate your warranty.

Scan the QR code to register your warranty and unlock an extra 12 months cover. That's 2 years in total!

This warranty does not cover fault or defect caused by:

- Faults arising from failure to observe operating and cleaning instructions
- Careless handling or misuse
- Transit damage
- Alterations, repairs or adjustments undertaken by anyone other than our qualified product engineers.

## CUSTOMER SERVICE

help@magnitone.co.uk  
www.magnitone.co.uk

## WARNINGS (1/2)

Are you ready for LiftOff? Lift, tone + contour your face in just 5 minutes with this handheld microcurrent facial toning device. Essentially a workout for your face, similar to facial massage, small pulses of electricity penetrate deep into the skin, keeping it looking plump, lifted and youthful. Used in expensive salon treatments for years, this handy device will get you professional results at home.

LiftOff is intended for facial stimulation and is for cosmetic use only. This is not a medical device and should not be used to treat any medical conditions. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty contacting elements, or any other improper application is neither the responsibility of HSNF Ltd nor its affiliates.

Some users may find the sensation quite strong, especially in comparison to other microcurrent products they have used. A tingling sensation (like a small electric shock) and some involuntary muscle twitching may occur while using MAGNITONE LiftOff. Lowering the mode setting may reduce or eliminate this sensation. If this is uncomfortable, we suggest using a lower strength setting. Normally this sensation subsides after the first few uses. If you experience any pain or discomfort, stop using the device immediately.

- For use on the face only.
- Adverse effects may occur if the instructions are not followed;
- Prolonged or repeated application on one site may potentially damage the skin;
- **DO NOT** use over the eyes (the circular muscle within the orbital rim), on or around the thyroid glands (midline of the neck, in-front of and either side of the wind-pipe) on the temples (above the ear), breasts or groin area.

**WARNINGS (2/2)**

- **DO NOT** apply stimulation across your chest as the introduction of an electrical current into the chest may cause rhythmic disturbances to your heart, which can be lethal.
- **DO NOT** use the appliance in conjunction with flammable based cosmetic products, unless they are recommended by the manufacturer;
- When using on the jawl, care should be taken to minimise treatment directly over the midline of the neck, specifically avoiding the thyroid.
- There may be a perception of flashing lights during the use of your LiftOff, resulting from the stimulation of the eye area. If you notice this condition constantly without the use of your device, discontinue use and seek medical advice.
- Some users may experience a metallic taste in their mouth during the use of the device, resulting from stimulation and varying sensitivity around the mouth. If you notice this condition constantly without the use of your device, discontinue use and seek medical advice.

**CONSULT YOUR DOCTOR BEFORE USE IF YOU HAVE ANY UNDERLYING SKIN CONDITIONS.**

## CONTRAINDICATIONS

### DO NOT USE IF:

- You have areas where the integrity of the skin is damaged (such as burns, blisters, lesions, incision, wounds, active skin disease or infection); are suffering from skin disease or inflammation of the skin, have sensitive skin or an allergic constitution.
- You have a history of epilepsy or seizures.
- You have had recent surgery or other surgical procedures.
- You have any electronic implanted device such as cardiac pacemaker.
- You use a body-worn medical device, such as an insulin pump.
- You have any active cancer or have recently recovered from or had treatment for any suspicious or cancerous lesions.
- You are pregnant.
- You are under the care of a doctor for any ongoing medical condition, consult them before using this device.
- You have a known history of carotid artery disease, stroke or transient ischemic attack (TIA), carotid stenosis, unstable blood pressure, easy fainting, thyroid disease, or prior surgical procedures in the neck area.
- You have a tendency for bleeding following acute injury or breaks in bones.
- You are wearing implants not compatible with its use.
- This appliance shall not be used by or on persons up to and including 14 years of age;

## SAFETY CAUTIONS

- **DO NOT** immerse in water.
- **DO NOT** drop or misuse the product.
- **DO NOT** handle the device or charger with wet hands.
- **DO NOT** operate the device if it has been damaged.
- **DO NOT** disassemble or attempt to repair the product.
- **DO NOT** place the product in direct sunlight or expose the battery to high temperatures, or heat-generating equipment. Other electrical equipment in the vicinity may be affected by the usage of this product.
- **DO NOT** bypass or remove the inherent safety device of the appliance.
- Make sure your device is disconnected from the USB charging cable before use.
- The appliance **MUST** only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Strictly for external use.

## CLEANING + STORAGE


Cleaning and user maintenance shall not be made by children without supervision. Make sure your device is disconnected from the USB charging cable before cleaning. Use a slightly damp towel to clean the device and dry with a soft cloth.

## RECYCLING + DISPOSAL

Your MAGNITONE LiftOff must not be disposed of with other household or commercial waste at the end of its working life. When you have no further use for your product, please dispose of it according to your local authority's recycling procedure.

## PARTS + ACCESSORIES
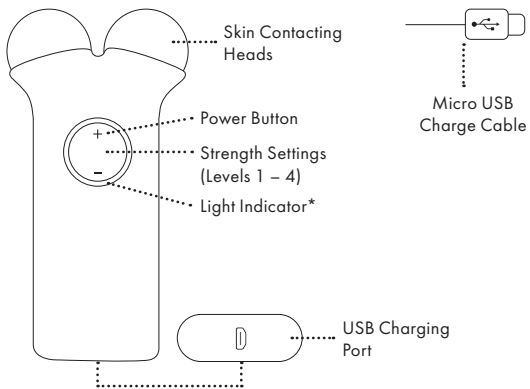
### \*LIGHT INDICATOR

Low Battery  
(Strength setting light  
flashes during usage)

Charging

Battery Full

## AREAS OF USE


LiftOff is designed to stimulate large surface areas of the face, such as the jawline, jowl, cheekbones and forehead. **DO NOT** use over the eyes (the circular muscle within the orbital rim), on or around the thyroid glands (midline of the neck, in-front of and either side of the wind-pipe) on the temples (above the ear), breasts or groin area. Make sure you have read and understood all **WARNINGS** before you begin. Follow the instructions below for best results.

## HOW OFTEN

- Use LiftOff for 3-5 minutes per treatment.
- Always allow at least 24 hours between treatments.
- For best results, use LiftOff a minimum of 4 times per week for the first 2 months.
- To maintain the results, use LiftOff a minimum of once per week ongoing.



**REMEMBER** - LiftOff is like going to the gym for your face! You can't go once and expect a six-pack straight away, you've got to train and maintain.
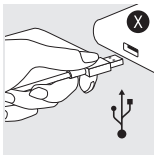
## BEFORE YOUR FIRST TIME

1. Charge LiftOff for 2 hours. Plug into any USB Socket on a charge plug or computer to charge your device (X). When the device is fully charged, the lights will go solid.
2. Test your skin for sensitivity to LiftOff Prepare + Glide Gel (Y). Rub a pea-sized amount of gel on the underside of your forearm and leave for 15-20 minutes before washing it off. If an allergic reaction (rash or redness) develops, immediately contact MAGNITONE for recommendations on alternative formulas. You should not begin treatment with MAGNITONE LiftOff until your reaction to the gel has been resolved. Make sure you have read and understood all of the **WARNINGS**.
3. Take your 'Before' photo! Remember to take a before photo of your cleansed face, without makeup, in a well-lit room (Z). Take photos every few weeks in the same place and position to track your progress. You may not notice changes daily, but regular selfies will help you to see the difference!
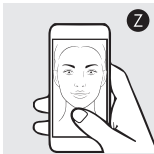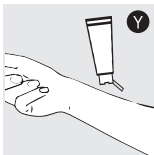
## HOW TO LIFT OFF!

The instructions below are designed to provide the most effective results. Best results are achieved and maintained with regular use. All skin is different, and sensitivity varies. During use, you should feel a tingling pulsation and it may cause your facial muscles to twitch. This is normal! Some people feel no sensation at all.

What you should feel is a pulsed sensation - this is often barely noticeable but it can feel like small electric shocks that can cause the facial muscles to twitch whilst the product is in contact with the skin. However, it shouldn't cause pain or prolonged discomfort or restricted movement after use. If the sensation is painful on the lowest setting, then you may wish to discontinue using it.





Alternatively, try switching from our Super-conductive Gel to a normal moisturiser or serum, or to no gel at all. The oil base in normal moisturisers can help to dampen and block the current, meaning it penetrates less well into the skin and may be helpful for you to build up use from.

1. Before you begin, make sure your skin is fully cleansed - free of makeup, surface oils and daily grime (A). (Men should shave before treatment as facial hair can interfere with conductivity)

**IMPORTANT:** Any cream, lotion, self-tan, or other oil-based products on the skin can block the microcurrent.

\*We recommend a full double cleanse; try our WipeOut Amazing Microfibre Cleansing Cloths to get rid of makeup and surface dirt, then supercharge your second cleanse with one of our award-winning cleansing brushes! Check out the website for more skinspiration.

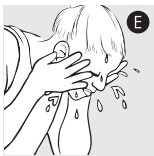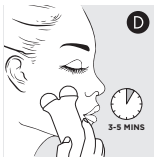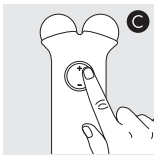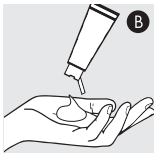- Apply a generous layer of LiftOff Prepare + Glide Superconductive Gel (B). It is best to do this section by section. You can reapply during treatment if your skin becomes dry. If you need to buy more LiftOff Gel, please visit [www.magnitone.co.uk](http://www.magnitone.co.uk)
- Press the '+' button to switch the device into LEVEL 1. Press '+' or '-' to select your desired strength setting, from 1 [Low], to 4 [High]. When you begin using LiftOff, always start at LEVEL 1 and work your way up (C).
- Follow the LiftOff Glide Techniques below (a beep sounds every 5 seconds to help you know how long to hold each glide). You should use the device for around 3-5 minutes in total (D).
- To switch the device off, either press + on LEVEL 4 or - on LEVEL 1. (**NOTE:** The device will automatically switch off after 10 mins).
- After use, rinse off the Glide Gel with water or a WipeOut Microfibre Cleansing Cloth (E) and complete your skincare routine. (**NOTE:** The use of products containing Retinol within 24 hours before/after treatment may cause skin irritation.)
- Use a slightly damp, lint-free cloth to gently wipe clean your device after every use.  
**Never wash LiftOff under a tap!**

## LIFTING TECHNIQUES

Always keep both contacting points touching the skin throughout each lift. For each lift, press down gently at the bottom and very slowly glide the device upwards until you hear the beep. Then move on. For a full facial, complete each glide listed below 3 times.


### JAW / JOWL

Place contacting points slightly off centre, taking care to avoid the midline of the neck. Glide toward the back of the neck until you hear the beep (F).

### CHEEKS

Place contacting points below the naso-labial fold ("Smile lines") or by the side of the nose and glide upwards toward the ear, over the cheekbone avoiding the eye and temple area, until you hear the beep (G).
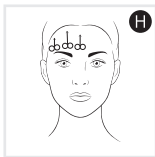
### FOREHEAD

Place contacting points slightly above the eyebrows and glide towards the hairline until you hear the beep (H).

Watch the demo video on YouTube.

**MAGNITONE**  
LONDON

## LIFT OFF SUPERCONDUCTIVE PREPARE + GLIDE GEL

### PRODUCT INFORMATION

Effortlessly glide and ensure the proper delivery of microcurrent pulses deeper to the skin with this superconductive, water-based gel. For use with MAGNITONE LiftOff Microcurrent Facial Toning Device only.

### HOW TO USE

Ensure your skin is clean and dry. Apply a generous layer of the gel to all your treatment areas, immediately before using your LiftOff device. Wipe or rinse off after use.

### CAUTIONS

For external use only. Not to be used on sore or damaged skin. In the unlikely event of rash or irritation, discontinue use. Avoid contact with eyes.

### STORAGE

Store in a cool dry place away from sunlight, with the lid firmly closed.

### INGREDIENTS

Aqua (Water), Carbomer, Glycerin, Methyl Paraben, Propyl Paraben, Propylene Glycol, Sodium Hydroxide.

crueltyfree  
andvegan

## SHARING IS CARING!

We love seeing your pictures and videos. Take a photo, make a video, tag us on social and we'll give it a share!

magnitoneofficial

magnitonelondon

magnitoneskincare

youtube.com/c/MAGNITONELondon

Make sure to tag us in and use **#NoFilterSkincare**.

## OUR PROMISE

We believe that all skin is beautiful.

Lines, blemishes, freckles, pimples - they're what make us unique and make us human. We want to represent real people and real skin, so we promise that all of our imagery is 100% un-retouched.

Our devices and skincare are simply tools to help you feel confident in the gorgeous skin you're in.

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