# MAGNITONE

LONDON

Since 2013, it has been our mission to energise and condition skin naturally; no knives, needles or baffling science! We create quality devices at affordable prices, making **#NoFilterSkincare** accessible to all.

**MAGNITONE**  
LONDON

## SCRUBUP

3-in-1 Ultrasonic + Microcurrent  
Scrubber

Read the instructions of use carefully. Keep the box and user manual for future reference.

### THIS MANUAL INCLUDES

- TECHNICAL SPECIFICATION
- WHAT'S IN THE BOX
- REGISTER YOUR WARRANTY
- WARNINGS
- CLEANING + STORAGE
- RECYCLING + DISPOSAL
- CHARGING YOUR DEVICE
- PARTS + ACCESSORIES
- THE SCIENCE-Y BIT
- HOW TO SCRUB UP!

### TECHNICAL SPECIFICATION

**Product:** MAGNITONE SCRUB UP 3 in 1 Ultrasonic + Microcurrent Scrubber

**Product Code:** MSU01W

**Input:** 5V  1000mA

**Battery:** Li-Ion 3.7V 500mAh

© MAGNITONE 2021

Registered trademark, owned by:  
HSNF Ltd., London, N10 2JE, UK.  
HSNF IE Ltd., Dublin 7, D07 P4AX, Ireland.

Made / Fabriqué / Prodotto / Hecho / Hergestellt /  
Wyprodukowano / Vervaardigd in China.

### WHAT'S IN THE BOX

- 1 x ScrubUp 3-in-1 Ultrasonic + Microcurrent Scrubber
- 1 x USB Charging Cable
- 1 x Instruction Manual

### REGISTER YOUR WARRANTY

Your device is covered by a 12 month warranty from your date of purchase. You must register your device to validate your warranty.

Scan the QR code to register your warranty and unlock an extra 12 months cover. That's 2 years in total!

This warranty does not cover fault or defect caused by:

- Faults arising from failure to observe operating and cleaning instructions
- Careless handling or misuse
- Transit damage
- Alterations, repairs or adjustments undertaken by anyone other than our qualified product engineers.

### CUSTOMER SERVICE

help@magnitone.co.uk  
www.magnitone.co.uk

## WARNINGS

### AREAS OF USE:

ScrubUp is intended for cosmetic use on the face as described. Avoid applying the device to the delicate eye area, the lips and the ears. This is not a medical device and should not be used to treat any medical conditions. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty contacting elements, alternative preparation gel, or any other improper application is neither the responsibility of HSNF Ltd nor its affiliates.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

## CONTRAINDICATIONS

### DO NOT USE IF

- You have eczema, dermatitis or other any skin disorders.
- You are taking medication that can cause skin sensitivity
- You have areas where the integrity of the skin is damaged (such as sunburns, burns, blisters, lesions, incision, wounds, active skin disease or infection); are suffering from skin disease or inflammation of the skin, have sensitive skin or an allergic constitution.
- You have a history of epilepsy or seizures.
- You have had recent surgery or other surgical procedures.
- You have any electronic implanted device such as cardiac pacemaker.
- You use a body-worn medical device, such as an insulin pump.
- You have any active cancer or have recently recovered from or had treatment for any suspicious or cancerous lesions.
- You are pregnant.
- You are under the care of a doctor for any ongoing medical condition, consult them before using this device.
- You have a known history of carotid artery disease, stroke or transient ischemic attack (TIA), carotid stenosis, unstable blood pressure, easy fainting, thyroid disease, or prior surgical procedures in the neck area.
- You have a tendency for bleeding following acute injury or breaks in bones.
- You are wearing implants not compatible with its use.

## SAFETY CAUTIONS

- **DO NOT** use over the eyes (the circular muscle within the orbital rim), on or around the thyroid glands (midline of the neck, in-front of and either side of the wind-pipe) on the temples (above the ear), breasts or groin area.
- **DO NOT** use the appliance in conjunction with flammable based cosmetic products, unless they are recommended by the manufacturer.
- **DO NOT** immerse the device in water.
- **DO NOT** drop or misuse the product.
- **DO NOT** handle the device or charger with wet hands.
- **DO NOT** operate the device if it has been damaged.
- **DO NOT** disassemble or attempt to repair the product.
- **DO NOT** place the product in direct sunlight or expose the battery to high temperatures, or heat-generating equipment. Other electrical equipment in the vicinity may be affected by the usage of this product.
- **DO NOT** bypass or remove the inherent safety device of the appliance.
- Make sure your device is disconnected from the USB charging cable before use.
- The appliance **MUST** only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Adverse effects may occur if the instructions are not followed.
- Prolonged or repeated application on one site may potentially damage the skin.
- If you experience any skin irritation, stop using the device immediately.
- When using on the jawl, care should be taken to minimise treatment directly over the midline of the neck, specifically avoiding the thyroid.

- Some users may experience a metallic taste in their mouth during the use of the device, resulting from stimulation and varying sensitivity around the mouth. If you notice this condition constantly without the use of your device, discontinue use and seek medical advice.
- Close supervision is necessary when this product is used by, on, or near children or physically or mentally challenged persons.
- Children shall not play with the appliance.
- Strictly for external use only.

## CONSULT YOUR DOCTOR BEFORE USE IF YOU HAVE ANY UNDERLYING SKIN CONDITIONS.

## CLEANING + STORAGE

Cleaning and user maintenance shall not be made by children without supervision.

Make sure your device is disconnected from the USB charging cable before cleaning. Use a slightly damp towel to clean the device and dry with a soft cloth.

## RECYCLING + DISPOSAL


Your MAGNITONE ScrubUp must not be disposed of with other household or commercial waste at the end of its working life. When you have no further use for your product, please dispose of it according to your local authority's recycling procedure.


## CHARGING YOUR DEVICE

Charge ScrubUp for 12 hours. Plug into any USB Socket on a charge plug or computer to charge your device. When the device is fully charged, the lights will go solid.

## PARTS + ACCESSORIES

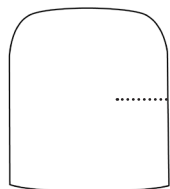..... Dust cover

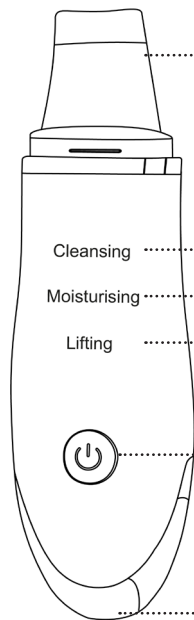..... Stainless steel probe

Cleansing ..... Cleaning mode

Moisturising ..... Moisturising mode

Lifting ..... Lifting mode

..... Power button  
(Function Switch Key)

..... USB Charging port

## THE SCIENCE-Y BIT

ScrubUp uses a combination of technologies for each mode - Cleansing, Moisturising + Lifting.

ULTRASONIC

High frequency sound waves create pressure to push out grime + excess oil from pores.

POSITIVE IONS

Positive ions attract the 'bad' stuff, working to help to cleanse + extract pores.

NEGATIVE IONS

Negative ions aid the 'good' stuff by helping the skin to absorb moisture + nutrients.

MICROCURRENT

Small electrical pulses stimulate the muscles, encouraging blood flow + collagen production.

MODES	CLEANSING	MOISTURISING	LIFTING
TECHNOLOGIES			
USE WITH	WATER	SERUM	SUPER CONDUCTIVE GEL
SPATULA POSITION	EDGE	FLAT	FLAT
HOW OFTEN	1-2 PER WEEK	DAILY	1-4 PER WEEK

## HOW TO SCRUB UP!

### BEFORE YOU BEGIN

Make sure you have read and understood all of the **WARNINGS**.

The instructions below are designed to provide the most effective results. Best results are achieved and maintained with regular use. All skin is different, and sensitivity varies. During use, you may feel a tingling pulsation. This is normal! Some people feel NO sensation at all. What you should feel is a pulsed sensation - this is often barely noticeable but it can feel like small electric shocks that can cause the facial muscles to twitch whilst the product is in contact with the skin. However, it shouldn't cause pain or prolonged discomfort or restricted movement after use. If the sensation is painful, then you may wish to discontinue using it.

### TOP TIP!


Remember to take a before photo of your cleansed face, without makeup, in a well-lit room. Take photos every few weeks in the same place and position to track your progress. You may not notice changes daily, but regular selfies will help you to see the difference, especially in Lifting Mode!

Rather watch than read?  
Watch the How To video  
on YouTube:

## CLEANSING MODE

The ultrasonic waves + positive ion technology create pressure that dislodges the dirt from your pores.  
Recommended use - no more than 1-2 times per week.

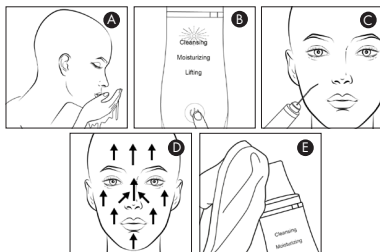
1. Make sure your face is clean and makeup free.
2. Splash water onto your face, **your skin must be wet** while using ScrubUp in Cleansing Mode (A). You may need to continue to wet your face throughout the treatment.
3. Switch on by holding the power button for 3 seconds to activate Cleansing Mode (B).
4. Place the end of the metal spatula to your skin with a light pressure (C) and glide in upwards strokes across your skin, avoiding the delicate eye area (D). **DO NOT** press hard or repeat more than twice over the same area of skin.
5. You may need to wipe the spatula clean throughout the treatment to remove the extracted dirt, oil + grime. Use a damp WipeOut cloth.
6. When you've finished your treatment, rinse your face with warm water and pat dry.
7. Press and hold to switch off your device, wipe clean the spatula using soap and water (E). **DO NOT** rinse or submerge in water.
8. Finish your treatment by applying your favourite serum or moisturiser. Help it to infuse using Moisturising Mode.

## MOISTURISING MODE

The ultrasonic waves + negative ion technology help the skin to absorb and infuse the product. Can be used daily.

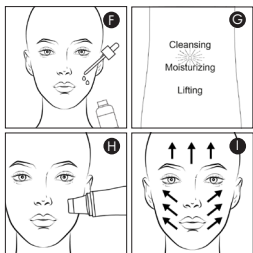
1. After cleansing, generously apply your favourite nourishing or hydrating skin serum or moisturiser (\*not included) (F). **DO NOT** use exfoliating serums, acids or activated products as this may cause irritation.
2. Switch the device on by holding the power button for 3 seconds and then press to activate Moisturising Mode (G).
3. Using the flat side of the metal spatula (H), glide upwards across your skin using a gentle pressure (I) and watch your serums sink in!
4. Once you've covered your whole face, press and hold to switch off your device, wipe clean the spatula using a damp cloth with soap and water (E). **DO NOT** rinse or submerge in water.

## LIFTING MODE

Small pulses of microcurrent stimulate the muscles, encouraging blood flow + collagen production. Recommended use - 3-4 times per week for 4 weeks (to see results), then 1-2 times per week ongoing (to maintain).

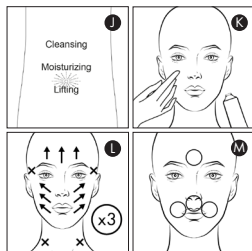1. Before you begin, make sure your skin is fully cleansed - free of makeup, surface oils and daily grime. (Men should shave before treatment as facial hair can interfere with conductivity) **IMPORTANT:** Any cream, lotion, self-tan, or other oil-based products on the skin can block the microcurrent.
2. Switch the device on by holding the power button for 3 seconds and then press twice to activate Lifting Mode (J).
3. Apply a generous layer of LiftOff Superconductive Gel (\*not included) (K). It is best to do this section by section. You can reapply during treatment if your skin becomes dry.  
\*To buy LiftOff Superconductive Gel, please visit [www.magnitone.co.uk](http://www.magnitone.co.uk)
4. Starting at the chin, place the flat side of the spatula flush to your skin (H) and gently glide along the jawline, up to the ear, on one side of your face. Repeat 2 more times. 3 x 'reps' in total.
5. Using this same motion, move on to the next section of your face and continue to work clockwise/anti-clockwise around your face (L). **DO NOT** use

on - the delicate eye area, temples or sides of throat (glands). Apply more gel as necessary.

- Once you have completed your whole face, you can go back and work on any specific 'problem' areas such as lips, smile lines or frown lines (M). To do this, hold ScrubUp in one place for 15 seconds.
- Press and hold to switch off your device, wipe clean the spatula using a damp cloth with soap and water. **DO NOT** rinse or submerge in water.
- When you're done, remove any excess gel from your face with a clean, damp WipeOut cloth (E) and complete your skincare routine. Finish your treatment by applying your favourite serum or moisturiser. Help it to infuse using Moisturising Mode!

## SHARING IS CARING!

We love seeing your pictures and videos. Take a photo, make a video, tag us on social and we'll give it a share!

 [magnitoneofficial](#)

 [magnitonelondon](#)

 [magnitoneskincare](#)

 [youtube.com/c/MAGNITONELondon](https://youtube.com/c/MAGNITONELondon)

Make sure to tag us in and use **#NoFilterSkincare**.

## OUR PROMISE

We believe that all skin is beautiful.

Lines, blemishes, freckles, pimples - they're what make us unique and make us human. We want to represent real people and real skin, so we promise that all of our imagery is 100% un-retouched.

Our devices and skincare are simply tools to help you feel confident in the gorgeous skin you're in.

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