

ToothShower®

Sonic Toothbrush



User Manual



Congratulations on purchasing your ToothShower Sonic Toothbrush! The ToothShower Sonic Toothbrush has a patented, high-efficiency sonic motor. The bristles are designed to flare and massage at high speeds. The vibration frequencies offer up to 48,000 strokes per minute. The five modes and unique vibrations enable you to easily reach your gum line for improved plaque removal and cleaner teeth.

Toothbrushing alone is not enough to clean all the surfaces of your teeth. The American Dental Association (ADA) recommends brushing your teeth twice a day and cleaning between the teeth (interdental) once a day. ToothShower Sonic Toothbrush and ToothShower Water Flosser Suite work in combination to fulfill the ADA daily recommendations for brushing and interdental cleaning.

Food and bacteria get trapped under this loose gum tissue (sulcus).

There is a space under the gum tissue.



If the food and bacteria are not removed from within this space, the gum tissue will become inflamed and bleed.

This tissue (free gingiva) surrounding the teeth is not attached to the bone underneath.

Table of Contents



Warnings and Precautions	3
.....	
Appliance Description	4
.....	
Operating Instructions	5
.....	
Brushing Modes	8
.....	
Toothbrushing Instructions	9
.....	
Care and Cleaning	12
.....	
Limited One Year Warranty	13

This product is intended to be used as directed. It uses electricity, and as such, could cause electrocution if not used as directed. Failure to use this product as intended and failure to follow the instructions, warnings, and precautions stated below could lead to additional physical injuries including, but not limited to, injuries to mouth and gums, and injuries to other body parts.

Warnings and Precautions



- Store the charger at a safe distance from water, water source and heated surfaces. Do not use the charger outdoors.
- Unplug immediately if the product has fallen into water and prior to reaching for it.
- Never use a charger with a damaged electrical cord or plug, and never force it into an outlet. Electrical cord cannot be replaced, discard if damaged in any way.
- Do not clean any parts of this product in the dishwasher.
- Do not use any attachments other than ToothShower Sonic Toothbrush recommended parts.
- Do not plug into an outlet with a voltage other than 100 to 240 volts. Voltage converters do not warrant voltage compatibility.
- Closely instruct and supervise children and individuals with special needs in the proper use of this appliance. Do not attempt to use the appliance as a toy.
- Discontinue use if any parts of the product appear damaged in any way.
- This product is designed for cleaning teeth, gums and tongue. Do not use it for any other purpose.
- If you have had oral or gum surgery in the previous 3 months, consult your dental professional prior to using this product.
- Initially your gums may bleed when using this product, consult a dental professional if there is excessive bleeding or if bleeding continues after two weeks of use.
- This product contains no user-serviceable parts. Refer to warranty and support if the product no longer works properly or needs repair.
- Stop using the appliance and contact a dental professional if you experience any discomfort or pain.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.



Contents

- (1) ToothShower[®] Sonic Toothbrush Waterproof Handle
- (2) Dupont Bristled Toothbrush Heads
- (1) Charging Base

Soft Dupont Bristle Toothbrush Head

Toothbrush Waterproof Handle

On/Off and Brushing Mode Button

Working Mode Indicators:

Clean
Whiten


Polish

Massage

Sensitive

Power and Charger Indicator

Charging Base

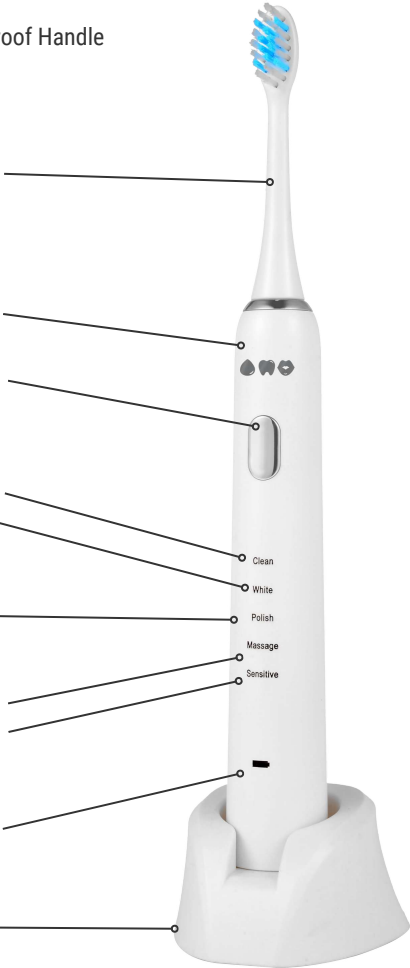



1 Charge toothbrush waterproof handle

- Place the charging base onto a dry and safe location. Plug the charger into an electrical outlet. Place the toothbrush handle onto the ToothShower charging base. The charging base indicator should be lit.
- A full charge will take at least 24 hours and enables approximately 6 weeks of brushing time, depending upon how often the toothbrush is used. The battery indicator light will blink when the charge is getting low.



Charger for model RLT 222

- The constant blue LED light indicates charging is in progress. Once fully charged, the light turns off automatically.
- Unplug the charger from the electrical outlet to save energy or leave it plugged in with the toothbrush handle on the charging base. There is no risk of overcharging the battery.



Charger for model RLT 202

Important: Do not clean the charging base by immersing it under running or still water. See cleaning instructions on page 12.



2

Attach the brush head onto the handle


- Firmly press the included brush head down onto the metal shaft of the handle until it stops. The brush head bristles should be facing the front of the handle.



3

Wet the toothbrush head bristles

- Wetting toothbrush bristles prior to applying the tooth paste softens the toothbrush bristles and rinses off debris. Only use ToothShower Sonic Toothbrush heads.

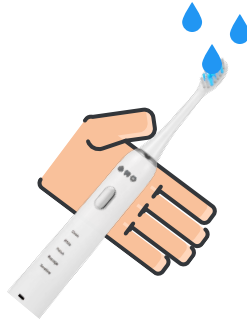



4 Apply a small amount of toothpaste

- You can use a toothpaste of your choice. A fluoridated toothpaste is recommended to help fight decay.



5 Place brush into your mouth before turning on

- To avoid splashing, guide the tooth brush head to your teeth before switching on the appliance. Press the power On/Off button to turn on the ToothShower Sonic Toothbrush.



Massaging Angled Bristles: Your ToothShower Sonic Toothbrush is equipped with the ability to have the perimeter bristles of the toothbrush head flare in unique angles to clean and massage the line where the tooth and gum meet. People often miss the gum line, where the majority of plaque accumulates.



Switching between brushing modes

- Press the On/Off and brushing mode button once, this will automatically start the Clean mode. By pressing the button again you can switch to the next mode. Each mode will run for 2 minutes before stopping.

Five brushing mode choices for varying oral care needs:

Clean - up to 31,000 strokes per minute, along with unique bristles at the perimeter of the brush head which will flare slightly to **clean** the teeth and gum line.

White - up to 41,000 strokes per minute, along with unique bristles at the perimeter of the brush head which will flare wider to **whiten** more surface area.

Polish - up to 48,000 strokes per minute, along with unique bristles at the perimeter of the brush head which will flare more vigorously for a **polish** effect.

Massage - up to 31,000 strokes per minute, along with unique bristles at the perimeter of the brush head which will flare wider and add a **massage** motion for a gentle stimulation to the gum line.

Sensitive - up to 31,000 strokes per minute, along with bristles at the perimeter of the brush head which will flare in a gentler motion for **sensitive** areas.

Brush for the full 2 minutes, switching quadrants of the mouth at the indicated 30 second delay intervals. The device will turn off after the Sensitive mode or by holding the power On/Off and brushing mode button down for 5 seconds.



Toothbrushing Instructions



Start with the outside of the teeth, then the insides, ending with the chewing surfaces. Do not scrub or press hard, simply let the brush do all of the work.


1 Brush at a 45 degree angle

Place the brush at a 45-degree angle to the front teeth surfaces. Bristles must contact both the tooth surface and gum line. This ensures you are brushing into the gum line to remove food and bacteria in the sulcus (space under the loose gum tissue).



2 Brush in a circular motion

Move the brush in a small, circular motion. The circular motion ensures you are reaching into the gum line and down along the tooth surface.

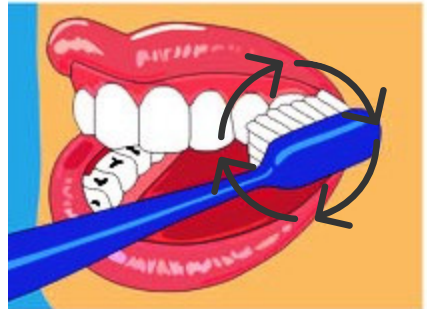

Toothbrushing Instructions



Guide the toothbrush head slowly from tooth to tooth. Hold the brush head in place for a few seconds before moving on to the next tooth. Brush along the gum line as well as the teeth.

3 Brush the tongue side of teeth

Clean the inside surfaces of the back teeth by moving the brush in a small back and forth motion. Make sure to brush all the way back behind the last molar. Do not let your tongue keep you from reaching these areas.



4 Brush behind the upper front teeth

Clean the inside surfaces of the front upper teeth by tilting the brush vertically using small up and down strokes. Be sure to have the bristles reach up into the gum line behind these teeth.



Toothbrushing Instructions



You can divide your mouth into four quadrants and brush all the teeth surfaces of each quadrant for 30 seconds. Each brushing mode offers a vibration frequency delay at 30 second intervals to know when to switch quadrants.

5 Brush behind the lower front teeth

Clean the inside surfaces of the front lower teeth by tilting the brush vertically using small up and down strokes. This is a difficult area to keep clean. Spend more time here and make sure to reach into the gum line with the bristles.



6 Brush the chewing surfaces of the teeth

Move the toothbrush in a back and forth motion to clean the chewing surfaces of the teeth.





1. After brushing, rinse your brush head

- Thoroughly rinse any remaining toothpaste off of the toothbrush head bristles. Once a week, remove the brush head and rinse the entire brush head with warm water. It is recommended to replace the toothbrush head every 3 months for effective brushing.

2. Wipe the handle dry

- Wipe the handle dry prior to placing it onto the charging base.

3. Clean charging base

- Unplug the charger before you clean it. Wipe the surface of the charger with a damp cloth.

4. Storage of toothbrush

- If you are not going to use your ToothShower Sonic Toothbrush for an extended period of time, unplug the charger from the electrical outlet. Clean it and store it in a cool, dry place away from direct sunlight.

5. Disposal of toothbrush

- This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in your household waste at the end of its useful life. Disposal can take place at an appropriate lithium ion battery collection point provided by your county or municipality.

Limited One Year Warranty



ToothShower, LLC warrants that this product is free from defects in materials and workmanship for 1 year from the date of purchase to the original purchaser only. Purchaser will be required to submit the original purchase receipt as proof of purchase date and the entire product to support a warranty claim. The product cannot have been abused, misused, altered or damaged after purchase and must have been used only according to ToothShower, LLC instructions with accessories or consumable parts approved by ToothShower, LLC. ToothShower, LLC will replace at its discretion any part of the product, which in its opinion is defective. Installation is the responsibility of the consumer and is not covered by the warranty. Consumable parts such as tips and brush heads are not covered by this limited warranty. This warranty gives you specific legal rights and you may have other rights that vary, depending on the law in your place of residence.

BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter. Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim. Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Pennsylvania and any such arbitration shall be subject to the Federal Arbitration Act and Pennsylvania law.

WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY.

Further, by purchasing, installing, or using this product you agree that arbitration must be on an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of ToothShower products or the general public. Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

ToothShower®

105 Airport Road, Pottstown, PA 19464

Contact us at: info@toothshower.com
www.toothshower.com

© 2021 ToothShower, LLC
All rights reserved
Color and images may not reflect actual product

ToothShower® and ®

Registered trademarks of ToothShower, LLC