




Black Electric Mobile Overbed Table

Instruction Manual

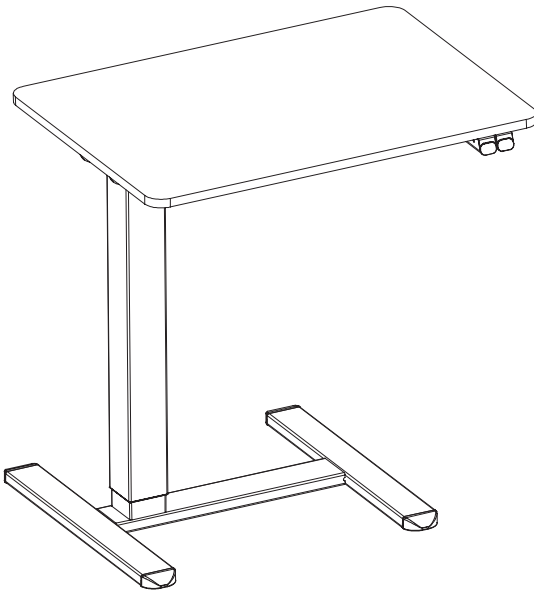

SKU: CART-E1TB



Scan the QR code with your mobile device or follow the link for helpful videos and specifications related to this product.

<https://vivo-us.com/products/cart-e1tb>

GET IN TOUCH | Monday-Friday from 7:00am-7:00pm CST



help@vivo-us.com



www.vivo-us.com
Chat live with an agent!



309-278-5303



WARNING!

If you do not understand these directions, or if you have any doubts about the safety of the installation, please call a qualified technician. Check carefully to make sure there are no missing or defective parts. Improper installation may cause damage or serious injury. Do not use this product for any purpose that is not explicitly specified in this manual. Do not exceed weight capacity. We cannot be liable for damage or injury caused by improper mounting, incorrect assembly or inappropriate use.



WARNING: CHOKING HAZARD

SMALL PARTS - NOT FOR CHILDREN UNDER 3 YEARS. ADULT SUPERVISION IS REQUIRED.




TIPOVER WARNING

SERIOUS OR FATAL CRUSHING INJURIES CAN OCCUR FROM TIPOVER. TO HELP PREVENT TIPOVER:


- NEVER ALLOW CHILDREN TO CLIMB, STAND, HANG, OR PLAY ON ANY PART OF CART.
- USE TIPOVER RESTRAINT OR ANCHOR CART TO WALL

USE OF TIPOVER RESTRAINTS MAY ONLY REDUCE, BUT NOT ELIMINATE RISK OF TIPOVER.


PACKAGE CONTENTS

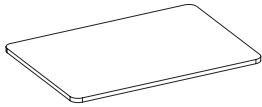


A (x1)
Desktop

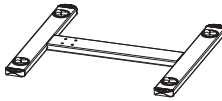


B (x1)
Base

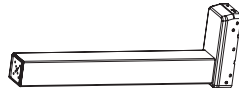


C (x1)
Adjustable
Leg

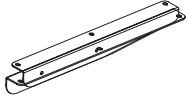

D (x1)
Side Bracket
Set

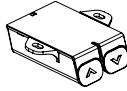

E (x1)
Control
Panel

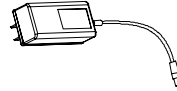

F (x1)
AC Power
Adapter



G (x4)
Cable Clips



S-A (x5)
M6x35mm
Screw



S-B (x5)
M6x11mm
Screw



S-C (x9)
M5x12mm
Screw



T-A (x1)
5mm Allen
Wrench

NOTE: SOME HARDWARE MAY NOT BE USED DURING ASSEMBLY.

DO NOT EXCEED WEIGHT CAPACITY.

Failure to do so may result in serious injury.



TOOLS NEEDED




Phillips
Screwdriver

ASSEMBLY STEPS


STEP 1

Install Base (B) to the underside of Adjustable Leg (C) using M6x35mm Screws (S-A) and 5mm Allen Wrench (T-A).

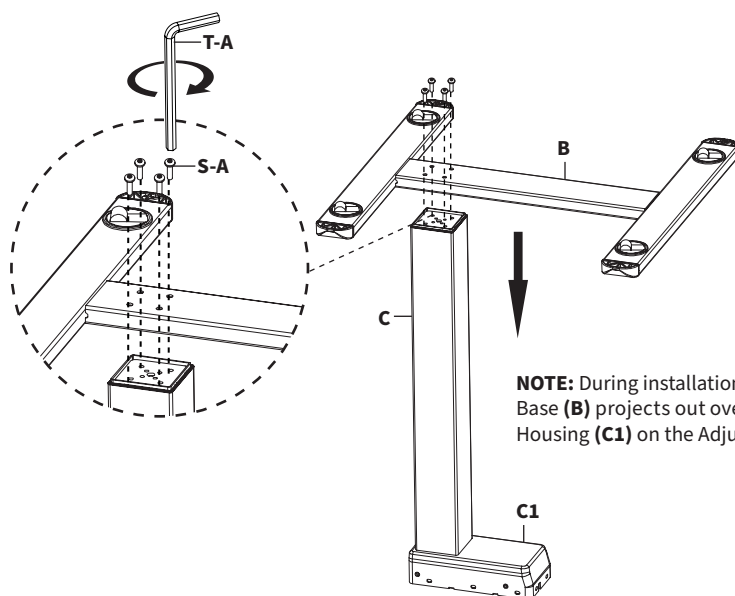


STEP 2

Install Side Bracket Set (D) to Adjustable Leg (C) using M6x11mm Screws (S-B) and 5mm Allen Wrench (T-A).

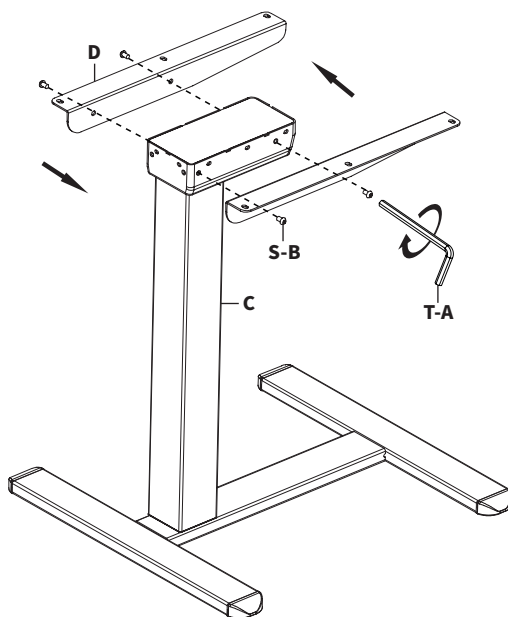

STEP 3


Install Desktop (A) to Side Bracket Set (D) using M5x12mm Screws (S-C) and a Phillips screwdriver.

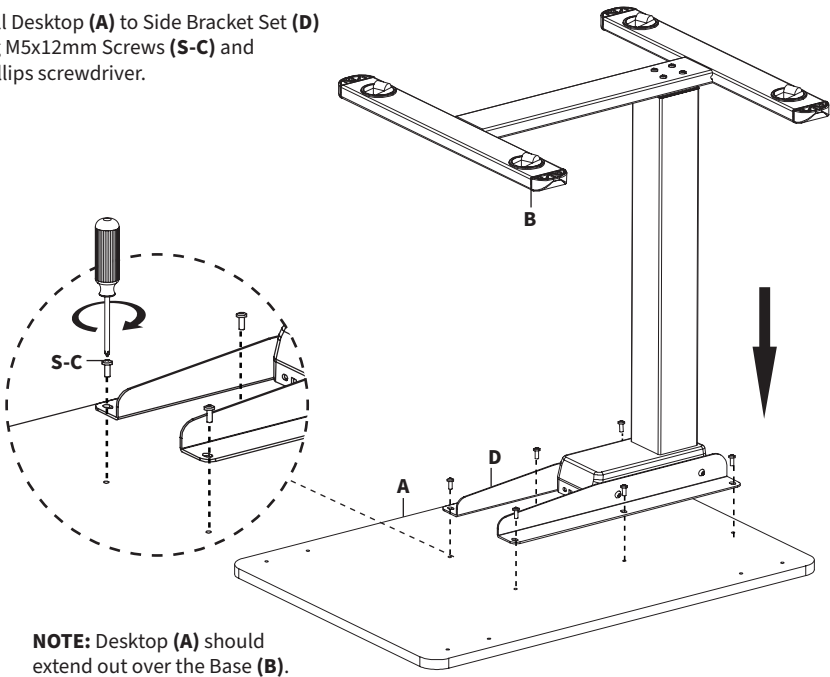

NOTE: Desktop (A) should extend out over the Base (B).


STEP 4

The Control Panel (E) can be installed on any of the four corners of the Desktop (A). Select the mounting location best for your use case and install Control Panel (E) to Desktop (A) using M5x12mm Screws (S-C) and a Phillips screwdriver.

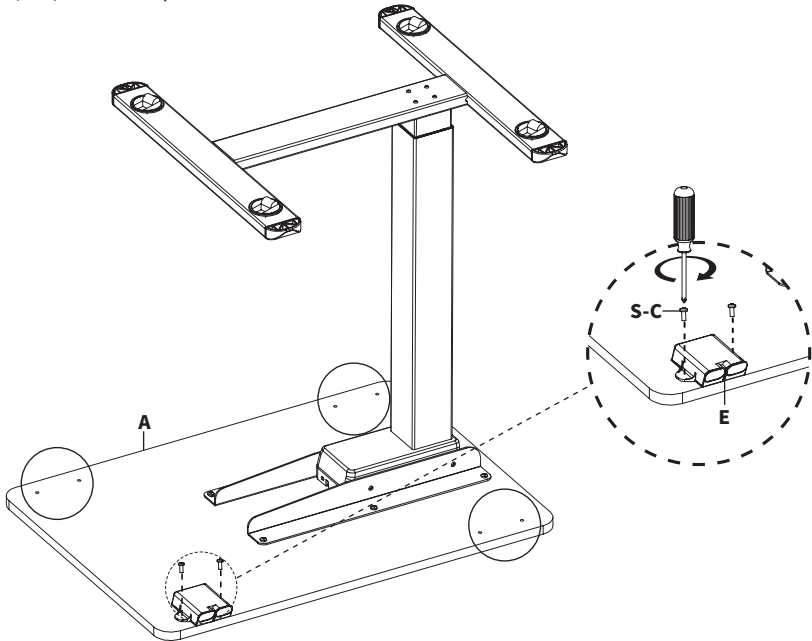


STEP 5

Plug the Motor Cable (**C1**) from Adjustable Leg (**C**) and AC Power Adapter (**F**) into Control Panel (**E**). Remove protective film from Cable Clips (**G**) and apply them to the underside of Desktop (**A**) to aid in cable management.

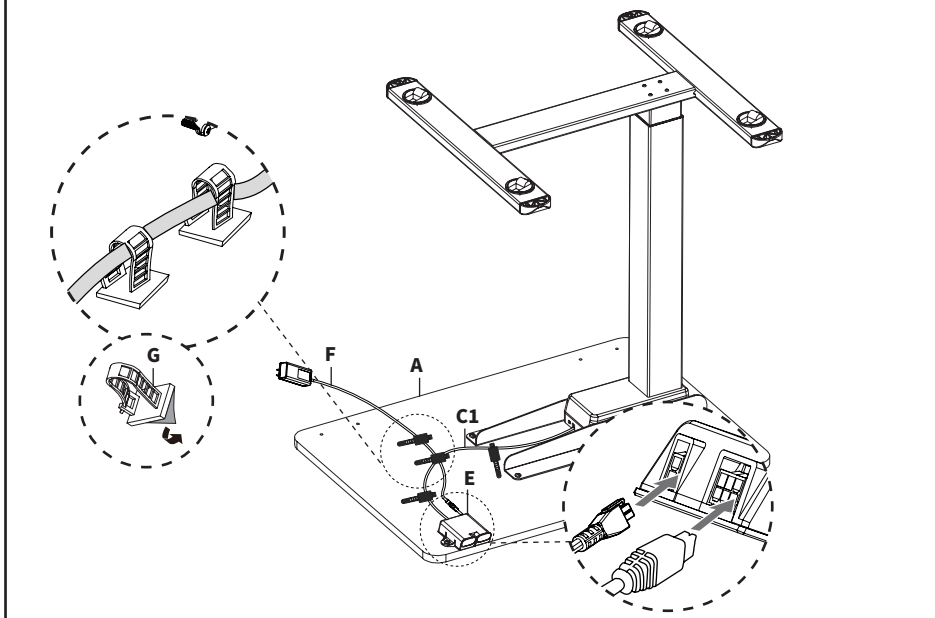


STEP 6

Plug AC Power Adapter (**F**) into a compatible outlet and follow the operating instructions on next page before use.

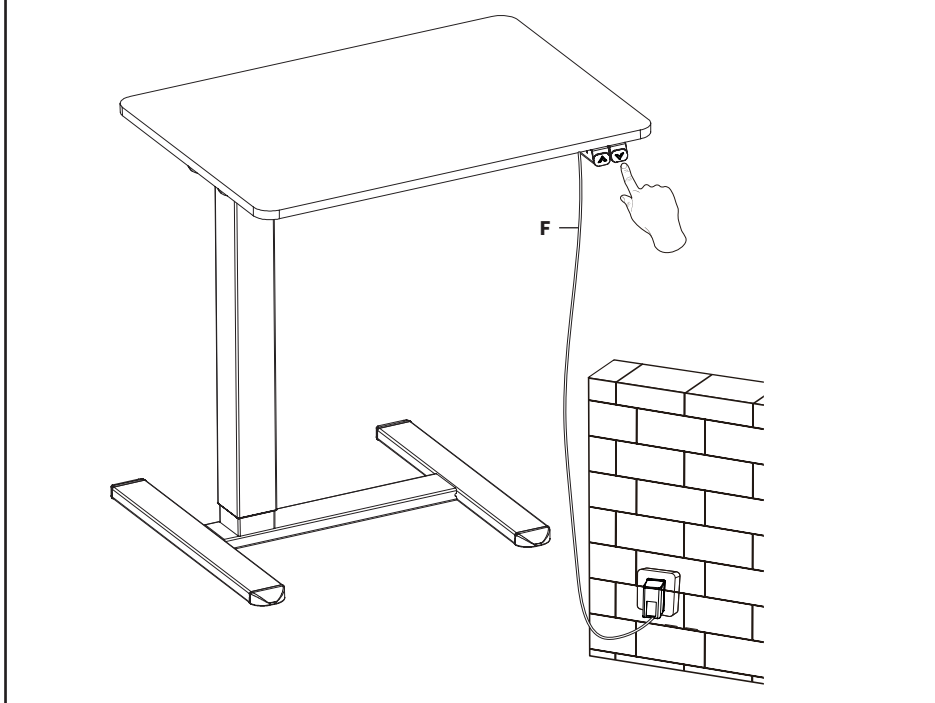

CONTROL PANEL

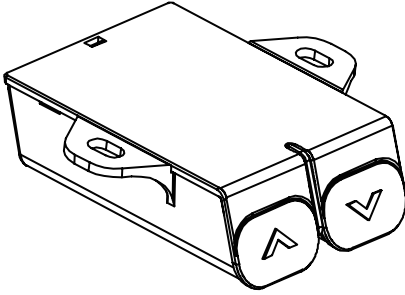

**Up Arrow:**


Press and hold to lift the desktop; release button to stop.

**Down Arrow:**

Press and hold to lower the desktop; release button to stop.

RESET PROCEDURE

Press and hold the  until the desk lowers completely; continue hold for 5 seconds. Release  when desk rebounds 1cm. Reset is now complete.

If the reset needs to be halted for any reason, simply release the . Perform reset again.



Open **Monday - Friday 7:00am - 7:00pm** CST,
our dedicated support team can offer immediate assistance with rapid response times. If any parts are received damaged or defective, please contact us. We are happy to replace parts to ensure you have a fully functioning product.



309-278-5303

AVG. RESOLUTION TIME (*within office hrs*): **5M 4S**



www.vivo-us.com
Chat live with an agent!

AVG. RESOLUTION TIME (*within office hrs*): **< 15 M**



help@vivo-us.com

AVG. RESPONSE TIME (*within office hrs*): **1HR 8M**

- 23% within < 15m
- 38% within < 30m
- 61% within < 1hr
- 83% within < 2hr
- 92% within < 3hr