
verilux®

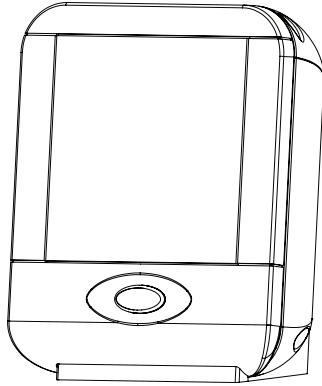To reduce the risk of damage during shipping, the bulb for your new HappyLight Liberty is packaged in protective bubble wrap. **You must remove the top cover, install the bulb and lens, and reinstall the top cover prior to using your lamp. See the Bulb/Lens Installation/Replacement section.**

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READ ALL INSTRUCTIONS BEFORE USING

Important Safeguards

⚠ DANGER:

- To avoid electrocution, do not operate this device near water.

⚠ WARNING:

- Do not use with power supply voltage other than 120 VAC.
- Do not defeat the safety purpose of the polarized plug. A polarized plug has two blades with one wider than the other. The wider blade is provided for your safety. If the provided plug does not fit in your outlet, consult an electrician for replacement of the obsolete outlet.
- Close supervision is recommended when these lamps are used by or near children, or others who may not understand the need for these precautions.
- To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the bulb to cool.
- Do not block vents. Vents are necessary for proper function of the lamp. Blocking them may cause the lamp to overheat.
- Do not shorten or cut the power cord.
- Always unplug lamp before performing bulb replacement or installation and allow sufficient time for the bulb to cool.
- Light bulbs get hot quickly! Do not touch the bulb while the lamp is on.
- This bulb contains a trace amount of mercury. For safe handling, cleanup, disposal and recycling in accordance with applicable laws in US see www.epa.gov/cfl and in Canada, please visit www.ec.gc.ca/mercure-mercury/ and www.lightrecycle.ca
- This product contains mercury, a chemical known to the State of California to cause birth defects or other reproductive harm.
- This product contains magnets. If you have a pacemaker, implantable defibrillator, insulin pump, any other type of electro-mechanical device, or brain shunts that contain magnets, keep sufficient distance (about six inches) between the device and the lamp.
- Some people may experience light sensitivity to fluorescent lighting. Discontinue use if you experience any discomfort. Follow directions for pretest on page four prior to use.
- Certain medications (e.g., anti-depressants) may induce light sensitivity which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.
- If you have a pre-existing eye condition, consult with your health care provider before using this product.

Important Safeguards — continued

▲ CAUTION:

- This lamp is for indoor use only.
- Avoid placing this lamp in areas that are exposed to direct sunlight or close to heat-radiating products such as heaters. Do not place the lamp on the top of appliances that radiate heat.
- Protect the power cord from being walked on or pinched, particularly at plugs, convenience receptacles and the point where the power cord exits the lamp.
- Do not operate this lamp if it has been damaged in any way. For example, the power supply cord or plug has been damaged, liquid has been spilled on or objects have fallen onto the lamp, the lamp has been exposed to rain or moisture, the lamp does not operate normally or has been dropped.
- This lamp may cause interference with radios, cordless telephones or devices that use a wireless remote control, such as televisions. If interference occurs, move the lamp away from the device, plug the lamp or the device into a different outlet or move the lamp out of the line of sight of the remote control receiver.
- Top of lamp is secured with strong magnets. Avoid lifting or carrying lamp by the top.
- Avoid placing the lamp in areas that are dusty, humid/moist, lack ventilation or are subject to constant vibration.
- Do not dismantle other than to replace the bulb. There are no user serviceable parts in these lamps.
- This lamp is not for use with light dimmers, timers, motion detectors, voltage transformers or extension cords.
- Do not leave the lamp unattended during use.
- Do NOT use lubricating agent in assembly. Please see troubleshooting section if you are having difficulty assembling your lamp.
- Use only the type of bulb and maximum wattage indicated by Verilux.
- Never cover the lamp or place anything on top of it when it is in operation.
- Do not use solvents or cleaners containing abrasives, or ammonia-based cleaners on this product.
- After cleaning the lamp, all moisture should be properly wiped away and dried before restoring power.

SAVE THESE INSTRUCTIONS

† This device complies with Part 18 of the FCC Rules.

This Class B digital apparatus complies with Canadian ICES-005.

3

Edit#: 002	Title: VT10 Bilingual Manual
Date: 01/11/19	Version: Rev22

Features

Your new Verilux HappyLight Energy Lamp is a remarkable product with many features and benefits. It:

- delivers Natural Spectrum light to your home, office, dorm room or apartment.
- cues your body's own energy enhancers to help you regain focus and concentration.
- has up to 10,000 LUX* of bright Natural Spectrum light to lift spirits, sharpen concentration and improve energy levels.
- is designed with a compact shape and small footprint to make it easy to place in a convenient location.
- includes an energy-efficient, 26-watt bulb that lasts up to 10,000 hours and provides flicker-free illumination.

* 10,000 lux at 10cm/4" from lamp. 5,000 lux at 15cm/6" to 20cm/8" from lamp.

Device Pretest

To help determine if you have sensitivity to light treatment, we recommend that you perform this pretest upon first use of this device:

Step 1 - Turn on the lamp and position it about 45cm/18 inches from your face. (Do not look directly into the lamp.)

Step 2 - Operate the lamp for approximately 10 minutes and turn it off. If the light has caused any discomfort, stop using it and consult a doctor. If you have no discomfort, you may consider gradually decreasing the distance from the lamp and increasing the use time, which results in more lamp intensity and greater effect. For your initial sessions, you should monitor the degree of any discomfort and accordingly adjust the time used and lamp location.

▲WARNING: Some people may experience light sensitivity to fluorescent lighting. Discontinue use if you experience any discomfort. Prior to using, follow the "Device Pretest" directions, shown above.

▲WARNING: Certain medications (e.g., anti-depressants) may induce light sensitivity which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.

▲WARNING: If you have a pre-existing eye condition, contact your health care provider prior to use.

Minimum Recommended Program

30 minutes
to
1 hour*/day

Researchers at the National Institute of Health demonstrated that white light, in a range of 2,500 to 10,000 LUX (a measurement of light intensity), helped to reverse the symptoms of “Winter Blues.”

Exposure strategies affect people differently. We all have different reactions to our environment. It is important to establish a distance and schedule that is comfortable for you. We recommend you use the light regularly and at about the same time each day. It is not necessary to sit in front of the lamp continuously. You can divide the time into several mini-sessions. Continue to use the HappyLight daily to feel an improvement in your mood and activity level.

When you are feeling sluggish and lethargic, spend more time with the light, sit closer to it or increase the light intensity. If you feel a sort of edginess, as though you have consumed too much coffee, it may be time to turn off the unit or reduce the intensity of the light.

* Recommended exposure time is a function of output and distance. Individual results may vary. Also, refer to the warnings at the bottom of the Device Pretest section of this manual.

▲WARNING: Certain medications (e.g., anti-depressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity, or hypomania. Consult with your healthcare provider before using this product.

▲WARNING: If you have a pre-existing eye condition, or are at risk for developing Age Related Macular Degeneration or other eye conditions consult with your healthcare provider before using this product.

Edit#: 002	Title: VT10 Bilingual Manual
Date: 01/11/19	Version: Rev22

Operation

Use consists of sitting close to the HappyLight with your eyes open. While staying oriented toward the light, feel free to engage in other activities such as reading, writing, hobbies and crafts, sewing, needlepoint, applying make-up or working at your computer or desk.

The HappyLight should be positioned as close as 4 inches, but within 24 inches of your face and should be slightly off-center. (See diagram.) *Note: Do not look directly into the lamp.*

The light from the HappyLight must be directed at your eyes; therefore, your eyes must be open to achieve the full benefit. *Note: Tinted glasses reduce the amount of light reaching your eyes.*

When you are feeling sluggish and lethargic, spend more time with the light, sit closer to it or adjust the settings to increase the intensity. If you experience edginess or are over-stimulated, move the light farther away or reduce session times.

▲WARNING: Do not block vents. Vents are necessary for proper function of the lamp. Blocking them may cause the lamp to overheat.

▲CAUTION: Avoid placing this lamp in areas that are exposed to direct sunlight or close to heat-radiating products such as heaters. Do not place the lamp on top of appliances that radiate heat.

Bulb/Lens Installation/Replacement

Your HappyLight uses a bulb manufactured only by Verilux. To order a replacement, please call Verilux toll-free, at 1-800-454-4408.


Initial Bulb Installation	Bulb Replacement	Lens Replacement	Instruction Steps
1.	1.	1.	Be certain the lamp is unplugged prior to beginning.
2.	2.	2.	Remove the top cover (A) by lifting firmly. The top is secured by two (strong) magnets (B) that will release when lifted properly.
X	3.	3.	Slide the lens (C) out of its channel by pulling it gently upward.
3.	X	X	Carefully remove the bulb from the bubble wrap.
X	4.	X	Remove the old bulb by twisting it counterclockwise until the bulb pins align with the openings in the bulb mount and lifting upward. (See illustration 2.)
4.	5.	X	Insert the new bulb into the bulb mount and rotate it clockwise until it clicks in place. (Either pin in either opening is fine. See illustration 3.)
5.	X	X	Remove the lens from the plastic bag that also held this instruction manual.
6.	6.	4.	(Re)install the lens and top cover, reversing the procedure shown in "Lens Replacement" steps two and three above.
7.	7.	5.	Dispose of any old parts and packaging according to local regulations.

▲WARNING: Always unplug lamp before performing bulb replacement or installation and allow sufficient time for the bulb to cool.


▲WARNING: This bulb contains a trace amount of mercury. For safe handling, cleanup, disposal and recycling in accordance with applicable laws in US see www.epa.gov/cfl and in Canada, please visit www.ec.gc.ca/mercure-mercury/ and www.lightrecycle.ca

Edit#: 002	Title: VT10 Bilingual Manual
Date: 01/11/19	Version: Rev22


Bulb/Lens Installation/Replacement — continued

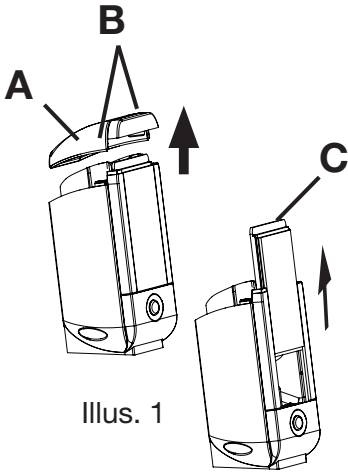Illus. 1

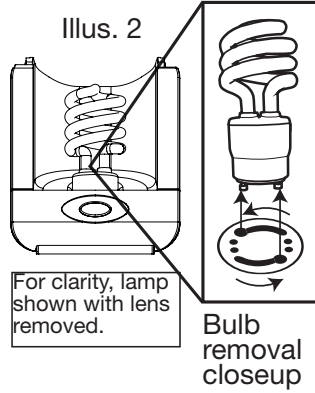Bulb removal closeup

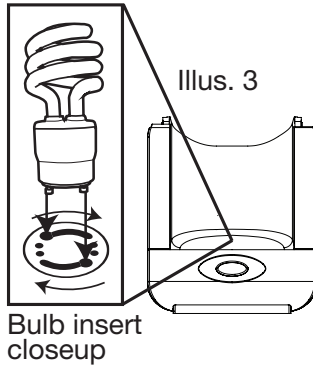Illus. 3

Care and Cleaning

Your lamp is made from high-quality materials that will last for many years with minimum care. You may want to periodically clean the lamp using a mild non-abrasive cleaner and soft cloth. When cleaning, make sure you have turned off and unplugged the unit and allowed sufficient time for the unit to cool.

Troubleshooting

Before requesting service on your Verilux lamp, please:

- Make sure all power connectors, cords and plugs are inserted fully and securely.
- Be sure the bulb is securely seated in the socket. (See Bulb/Lens Installation/Replacement section.)
- Make sure there is power to the wall outlet, or try another outlet – one that is not connected to a dimmer.

PROBLEM	CHECK	SOLUTION
Light Won't Come On	Has the plug come loose?	Insert plug securely into wall outlet.
	Is the bulb old?	Check bulb for blackening/darkening at bottom of spiral where it joins the base (an indicator of bulb life) and replace bulb.
	Is the bulb connected securely into socket?	Be sure the bulb is securely seated in the socket. Repeat the bulb installation process. [WARNING] Always unplug lamp before performing bulb replacement or installation.
Blackening/Darkening at Ends of Bulb	Is the bulb new?	Sometimes when a new bulb is turned on for the first time the tubing will darken, but this will fade after a few minutes.
	Is the bulb old?	With older bulbs, darkening of the tubing indicates end of bulb life. (This doesn't affect operation except to reduce the amount of light available, since the phosphor in that area has become inactive.) Bulb replacement is recommended.
Light Flickers	Has the plug come loose?	Insert plug securely into wall outlet.
	Is the bulb connected securely to the socket?	Be sure the bulb is securely seated into the socket. Repeat the bulb installation process. [WARNING] Always unplug lamp before performing bulb replacement or installation.
	Is the lamp connected to a dimmer switch?	Connect lamp plug only to wall outlet without dimmer switch.

Technical Specifications

HappyLight Compact Energy Lamp

Environment: Operating Temperature: 10 to 30°C (50 to 86°F)

Storage Temperature: -20 to 50°C (-4 to 122°F)

Relative Humidity: 20 to 90% non-condensing

Bulb: Replacement bulb: 26-watt 10,000-hour bulb. Replace with Verilux brand only, Model number CFS26GU24VLX by calling 1-800-454-4408 or visiting verilux.com.

Power: Input Power: 120 VAC, 60Hz, 0.50 A

ETL Listed in Canada and the USA

Edit#: 002	Title: VT10 Bilingual Manual
Date: 01/11/19	Version: Rev22

One-Year Limited Warranty

ATTENTION! ONCE OPENED, PLEASE DO NOT RETURN THIS PRODUCT TO THE STORE WHERE IT WAS PURCHASED FOR REPAIR OR REPLACEMENT!

Many questions can be answered by visiting www.verilux.com, or you may call our Customer Service Department at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Verilux, Inc., 340 Mad River Park, Waitsfield, VT 05673

Verilux warrants this product to be free from defects in material and workmanship for a period of one year from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux, Inc., will, at its option, repair or replace defective parts of this product at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance or delivery fees. This warranty does not cover damage, defect or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

Return Authorization is required for all returns. To obtain a Return Authorization, please contact the Verilux Customer Service Department at 1-800-786-6850.

If, during the first three years of ownership, this product fails to operate properly, it should be returned as specified at www.verilux.com/warrantyreplacement or as instructed by a Verilux customer service representative at 1-800-786-6850.

Note: *Verilux recommends using a quality surge suppressor on all electronic equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges. Suppressors may be purchased at electronics stores.*

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

Please visit our website at: www.verilux.com or call our Customer Service Department at 1-800-786-6850

Representatives are available Monday — Friday, 9am — 5pm EST

verilux[®]


340 Mad River Park, Waitsfield, VT 05673

Register warranty online at:
verilux.com/warranty

Printed in China for: Verilux, Inc.

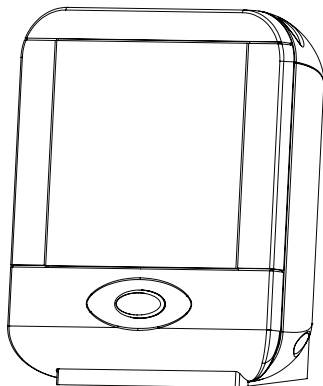© Copyright 2019 Verilux, Inc. All rights reserved.

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© Copyright 2019 Verilux, Inc. Tous droits réservés.
La garantie limitée est assurée par Verilux, Inc., 340 Mad River Park, Waitsfield, VT 05673
Imprimé en Chine pour : Verilux, Inc.