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
**verilux®**

340 Mad River Park, Waitsfield, VT 05673 USA

**verilux**<sup>®</sup>  
The Healthy Lighting Company<sup>®</sup> Est. 1956
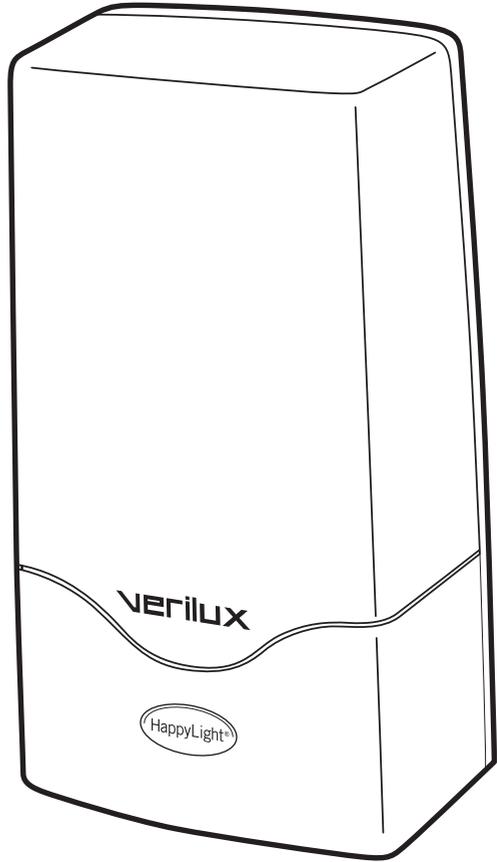
# HappyLight<sup>®</sup>

Compact Energy Lamp

**Model: VT01C**

**SAVE THESE INSTRUCTIONS**

Dear Customer,

Thank you for purchasing the Verilux® HappyLight® Energy Lamp. You now own an innovative and quality manufactured product, backed by a one-year limited warranty.

This product brings Natural Spectrum® light into your home or office. Many other healthy lighting products are available through our-catalog and online. Visit us on the web at [www.verilux.com](http://www.verilux.com) to request a free catalog and learn more about all our quality Verilux products, or call our toll-free number at 1-800-454-4408.

As a Verilux customer, your satisfaction means everything to us. We look forward to serving you now and in the future.

Have a Bright Day!

Nicholas Harmon  
President, Verilux, Inc.

## Contents

Important Safeguards . . . . .	4
Features. . . . .	5
Device Pretest . . . . .	6
Operation. . . . .	6
Selecting a Location . . . . .	6
Instructions for Use . . . . .	6
Recommended Program . . . . .	7
Bulb Removal and Replacement. . . . .	8
Re-installing Bulb and Plastic Diffuser . . . . .	9
Troubleshooting. . . . .	10
Care and Cleaning. . . . .	10
Technical Specifications . . . . .	11
Warranty . . . . .	12

## READ ALL INSTRUCTIONS BEFORE USING

### Important Safeguards

---

#### **▲ DANGER:**

- To avoid electrocution, do not operate this device near water.

---

#### **▲ WARNING:**

- Do not use with power supply voltage other than 120 VAC.
- Do not defeat the safety purpose of the polarized plug. A polarized plug has two blades with one wider than the other. The wider blade is provided for your safety. If the provided plug does not fit in your outlet, consult an electrician for replacement of the obsolete outlet.
- Close supervision is recommended when these lamps are used by or near children, or others who may not understand the need for these precautions.
- Do not operate lamp in close proximity to flammable or combustible vapors, such as aerosol spray products, or where oxygen is being administered.
- To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp, and allowed sufficient time for the bulb to cool.
- If this lamp has been damaged in any way, for example, the power supply cord or plug has been damaged, liquid has been spilled or objects have fallen onto the lamp, the lamp has been exposed to rain or moisture, the lamp does not operate normally or has been dropped, it should not be operated. There are no user serviceable parts in this lamp.
- Do not block vents. Vents are necessary for proper function of the electronic ballast; blocking them may cause overheating of the electronics and increase risk of fire.
- This bulb contains a trace amount of mercury. For more information on clean-up and disposal, see [www.epa.gov/cfl](http://www.epa.gov/cfl).
- This product contains a chemical known to the State of California to cause birth defects or other reproductive harm.
- Some people may experience light sensitivity to fluorescent lighting. Discontinue use if you experience any discomfort. Follow directions for pretest on page 6 prior to use.
- Certain medications may induce light sensitivity (e.g. anti-depressants) which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.
- If you have a pre-existing eye condition, consult with your health care provider before using this product.

---

#### **▲ CAUTION:**

- These lamps are for indoor use only.
- Avoid placing this lamp in areas that are exposed to direct sunlight or close to heat-radiating products such as heaters. Do not place the lamp on the top of appliances that radiate heat.
- Protect the power cord from being walked on or pinched, particularly at plugs, convenience receptacles and the point where the power cord exits the lamp.

---

**▲ CAUTION: — continued**

---

- This lamp may cause interference with radios, cordless telephones or devices that use a wireless remote control, such as televisions. If interference occurs, move the lamp away from the device, plug the lamp or the device into a different outlet or move the lamp out of the line of sight of the remote control receiver.
- Unplug lamp during lightning storms or when unused for long periods.
- Avoid placing the lamp in areas that are dusty, humid/moist, lack ventilation or are subject to constant vibration.
- This lamp is not for use with light dimmers, timers, motion detectors, voltage transformers or extension cords.
- Do not leave the lamp unattended during use.
- Do NOT use lubricating agent in assembly. See troubleshooting section if you are having difficulty assembling your lamp.
- Light bulbs get hot quickly! Do not touch the bulb while the lamp is on.
- Use only the type of bulb and maximum wattage indicated by Verilux.
- Never cover the lamp or place anything on top of it when it is in operation.
- Always unplug lamp before performing bulb replacement.
- Do not use solvents or cleaners containing abrasives, or ammonia-based cleaners on this product.
- After cleaning the lamp, wipe away all moisture and properly dry the lamp before restoring power.

## SAVE THESE INSTRUCTIONS

### Features

Your new Verilux HappyLight Energy Lamp is a remarkable product with many features and benefits:

- Delivers Natural Spectrum light to your home, office, dorm room or apartment with no dangerous ultraviolet (UV) wavelengths.
- Has up to 5,000/2,500 LUX\* of bright Natural Spectrum light per unit to lift spirits, sharpen concentration and improve energy levels.
- Cues your body's own energy enhancers to help you regain focus and concentration.
- Is designed with a compact shape and small footprint to make it easy to place in a convenient location.
- Includes an energy efficient, cool burning 26-watt, 10,000-hour bulb.
- Equipped with a practical switch on the cord for more comfort.
- Includes best-in-class power system with electronic ballast and Flicker Elimination Technology® to assure no-buzz operation and instant-on, flicker-free illumination.

\*5,000 Lux (scotopic measurement) / 2,500 Lux (photopic measurement) at 15cm to 20cm/6" to 8"

## Device Pretest

To determine if you are sensitive to light treatment, perform this pretest:

Step 1 - Place the lamp about 45 cm/18" from your face and turn it on. **(Do not look directly into the lamp.)**

Step 2 - Operate lamp for approximately 10 minutes and turn it off. If you feel any discomfort, discontinue use and consult a doctor. If you feel no discomfort, for subsequent use, you may want to gradually decrease the distance between the lamp and your face and/or lengthen the session time, which will result in increased lamp intensity and greater effect. If you develop any discomfort, increase the distance from your face and reduce the session time until you have found the time/distance combination that works for you.

**⚠WARNING:** Discontinue use if you feel any sensitivity to fluorescent lighting.

**⚠WARNING:** Certain medications may induce light sensitivity (e.g. anti-depressants) which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.

**⚠WARNING:** If you have a pre-existing eye condition, contact your health care provider prior to use.

## Operation

### Selecting a Location

- Place your Verilux HappyLight Energy Lamp on a table or other flat surface.
- For best performance, place the HappyLight as close as 15cm/6", but within 60cm/24" from your face.

**⚠WARNING:** Do not block vents. Vents are necessary for proper function of the electronic ballast. Blocking them may cause overheating of the electronics and increase risk of fire.


**⚠CAUTION :** Avoid placing this lamp in areas that are exposed to direct sunlight or close to heat-radiating products such as heaters. Do not place the lamp on the top of appliances that radiate heat.

### Instructions for Use

Use consists of sitting close to the HappyLight Energy Lamp with your eyes open. While staying oriented toward the light, feel free to engage in other activities such as reading, writing, hobbies and crafts, sewing, needlepoint, applying make-up or working at your computer or desk.

## Operation — continued


1. Turn on the lamp and allow a few minutes for the bulb to achieve full brightness.
2. The HappyLight Energy Lamp should be positioned as close as 15cm/6", but within 60cm/24" from your face, and should be slightly offset from center (See diagram.) *Note: Do not look directly into the lamp.*

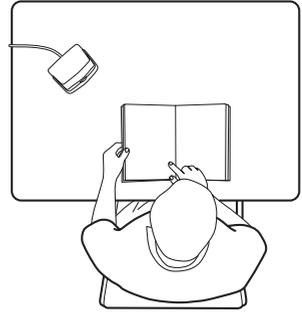The light from the HappyLight Energy Lamp must be directed at your eyes; therefore, your eyes must be open to achieve the full benefit. *Note: Tinted glasses reduce the amount of light reaching your eyes.*

When you are feeling sluggish and lethargic, spend more time with the light, sit closer to it or add more lighting devices to increase the intensity. If you experience edginess or are over stimulated, move the light farther away or reduce session time.

## Minimum Recommended Program

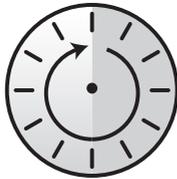30 minutes  
to  
1 hour\*

Researchers at the National Institute of Health demonstrated that white light, in a range of 2,500 to 10,000 LUX (a measurement of light intensity), helped alleviate the “Winter Blues.”

Different exposure strategies work differently for different people. We each have distinct reactions to our environment. For light therapy, it is important to establish a distance and schedule that is comfortable for you. We recommend you use the light regularly and at about the same time each day. It is not necessary to sit in front of the lamp continuously. You can divide the time into several mini-sessions. Continue to use the HappyLight daily to feel an improvement in your mood and activity level.

When you are feeling sluggish and lethargic, spend more time with the light, sit closer to it or add more lighting devices to increase the intensity. If you experience edginess or are over stimulated, move the light farther away or reduce use time.

\* Recommended exposure time is a function of output and distance. Individual results may vary.

## Operation — continued

**▲WARNING:** Discontinue use if you feel any sensitivity to fluorescent lighting. (See page 6.)

**▲WARNING:** Certain medications may induce light sensitivity (e.g. anti-depressants) which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.


**▲WARNING:** If you have a pre-existing eye condition, contact your health care provider prior to use.

## Bulb Removal and Replacement

Your HappyLight uses a bulb manufactured only by Verilux. To order a replacement bulb, please call Verilux, toll-free, at 1-800-454-4408.

### To Replace the Bulb:


1. Unplug the lamp from power outlet.
2. To access the bulb, remove the plastic diffuser by pressing the tabs (D) and lifting out the diffuser.
3. Firmly hold the bulb base (B) (do not hold the glass portion of the bulb) and turn bulb in a slight counter-clockwise direction until it releases from the socket.
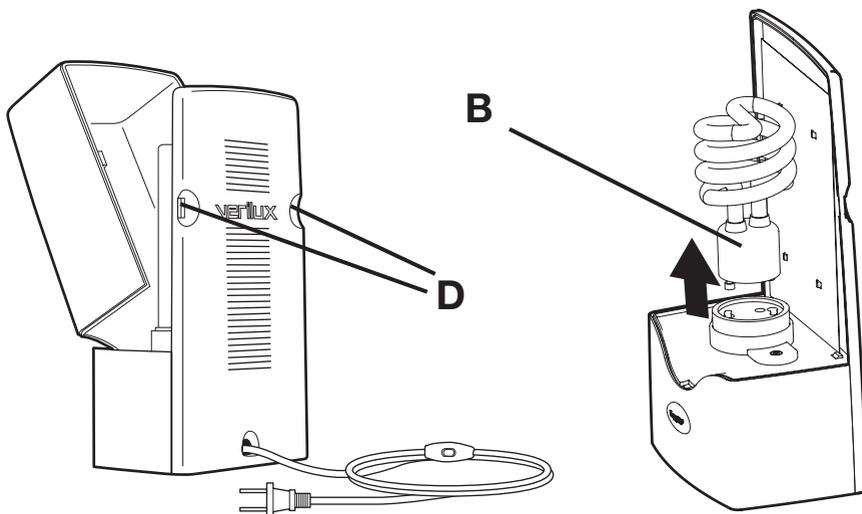
## Bulb Removal and Replacement — continued

### Re-installing Bulb and Plastic Diffuser

1. To install the bulb, insert bulb into lamp by firmly holding bulb base and pushing down and rotating slightly in a clockwise direction to seat properly.
2. Attach the diffuser to the lamp body by first inserting the tab located at the bottom of the diffuser (A) into the opening of the lamp base (B).
3. Insert one diffuser corner at a time into the side track and gently down into place into the bottom track.
4. Snap side tab into place on side slot. Entire side of diffuser should be seated into the track. Repeat for other side of diffuser. **Note:** Alignment of diffuser corner into side track.
5. Firmly push diffuser down until **Tab A** snaps into place. Align top of diffuser in top track as necessary.
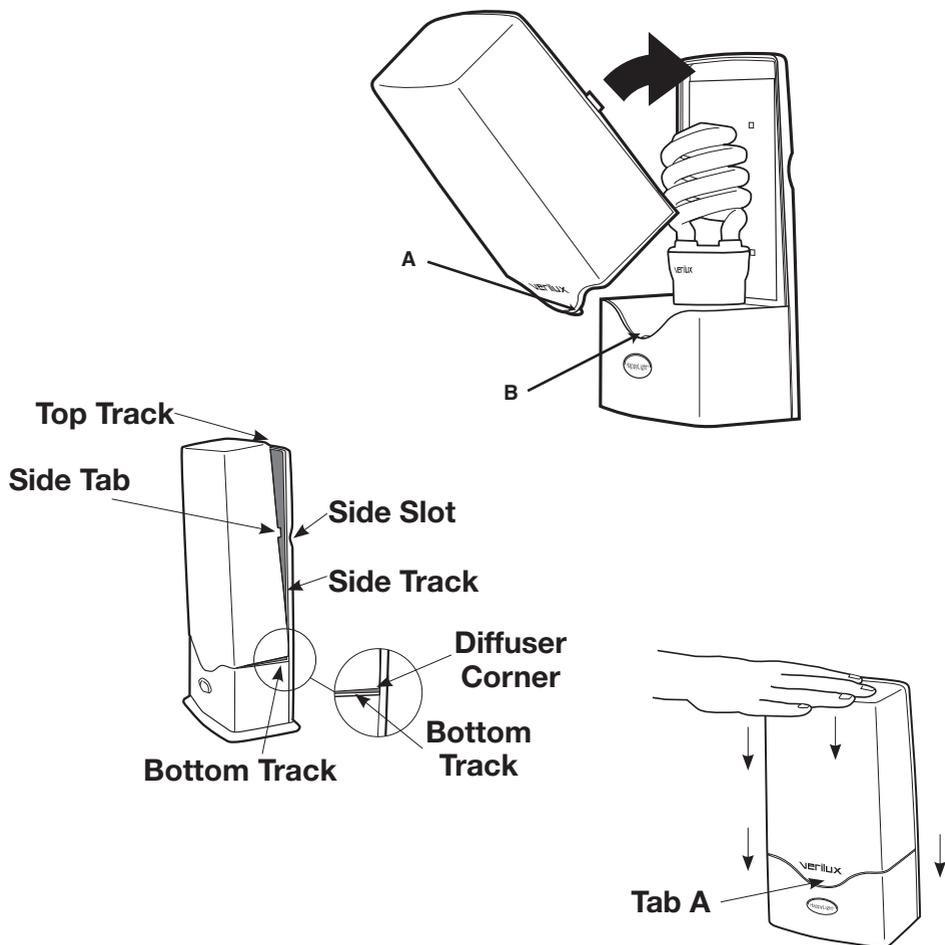
## Troubleshooting

Before requesting service on your Verilux lamp, please read the following:

- Make sure all power connectors, cords and plugs are inserted fully and securely.
- Be sure the bulb is securely seated in the socket. (See pages 8 and 9.)
- Make sure there is power to the wall outlet or try another outlet – one that is not connected to a dimmer.

PROBLEM	CHECK	SOLUTION
Light Won't Come On	Has the plug come loose?	Insert plug securely into wall outlet.
	Is the bulb old?	Check bulb for blackening/darkening at bottom of spiral where it joins the base (an indicator of bulb life) and replace bulb.
	Is the bulb connected securely into socket?	Be sure the bulb is securely seated in the socket. Repeat the bulb installation process. <b>⚠ CAUTION:</b> Always unplug lamp before performing bulb replacement or installation.
Blackening/Darkening at Ends of Bulb	Is the bulb new?	Sometimes when a new bulb is turned on for the first time the tubing will darken, but this will fade after a few minutes.
	Is the bulb old?	With older bulbs, darkening of the tubing indicates end of bulb life. (This doesn't affect operation except to reduce the amount of light available, since the phosphor in that area has become inactive.) Bulb replacement is recommended.
Light Flickers	Has the plug come loose?	Insert plug securely into wall outlet.
	Is the bulb connected securely to the socket?	Be sure the bulb is securely seated into the socket. Repeat the bulb installation process. <b>⚠ CAUTION:</b> Always unplug lamp before performing bulb replacement or installation.
	Is the lamp connected to a dimmer switch?	Connect lamp plug only to wall outlet without dimmer switch.

## Care and Cleaning

Your Verilux lamp is made from high-quality materials that will last for many years with a minimum of care.

**⚠ WARNING:** To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp, and allow sufficient time for the bulb to cool.

**⚠ CAUTION:** Clean only with a soft dry cloth – do not use solvents or cleaners containing abrasives or ammonia.

**⚠ CAUTION:** After cleaning the lamp, wipe away all moisture and properly dry the unit before restoring power.

## **Technical Specifications**

### **HappyLight Energy Lamp**

#### **Environment**

Operating Temperature: 15°C to 35°C; 56°F to 95°F

Storage Temperature: 0°C to 40°C; 32°F to 104°F

Relative Humidity: 5% to 95% non-condensing

#### **Bulb**

Replacement bulb: 26-watt 10,000-hour bulb (replace with Verilux brand only, Model number CFS26GU24VLX)

#### **Power**

Input Power: AC 120V~60Hz

**ETL Listed in Canada and the USA**

## One-Year Limited Warranty

**ATTENTION!** ONCE OPENED, PLEASE **DO NOT** RETURN THIS PRODUCT TO THE STORE WHERE IT WAS PURCHASED FOR REPAIR OR REPLACEMENT!

Many questions can be answered by visiting [www.verilux.com](http://www.verilux.com), or you may call our Customer Service Department at 800-786-6850 during normal business hours.

This limited warranty is provided by Verilux, Inc., 340 Mad River Park, Waitsfield, VT 05673 USA

Verilux warrants this product to be free from defects in material and workmanship for a period of one year from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux, Inc., will, at its option, repair or replace defective parts of this product, at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance or delivery fees. This warranty does not cover damage, defect or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

**Return Authorization is required for all returns.** To obtain a Return Authorization, please contact the Verilux Customer Service Department at 800-786-6850.

If, during the first year of ownership, this product fails to operate properly, it should be returned as specified at [www.verilux.com/warrantyreplacement](http://www.verilux.com/warrantyreplacement) or as instructed by a Verilux customer service representative at 800-786-6850.

**Note:** *Verilux recommends using a quality surge suppressor on all electronics equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges and may be purchased at electronics stores.*

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

Please visit our website at [www.verilux.com](http://www.verilux.com) or call our Customer Service Department at


**1-800-786-6850**

Representatives are available Monday — Friday, 9:00a.m. to 5:00p.m. EST

**verilux**<sup>®</sup>

340 Mad River Park, Waitsfield, VT 05673 USA

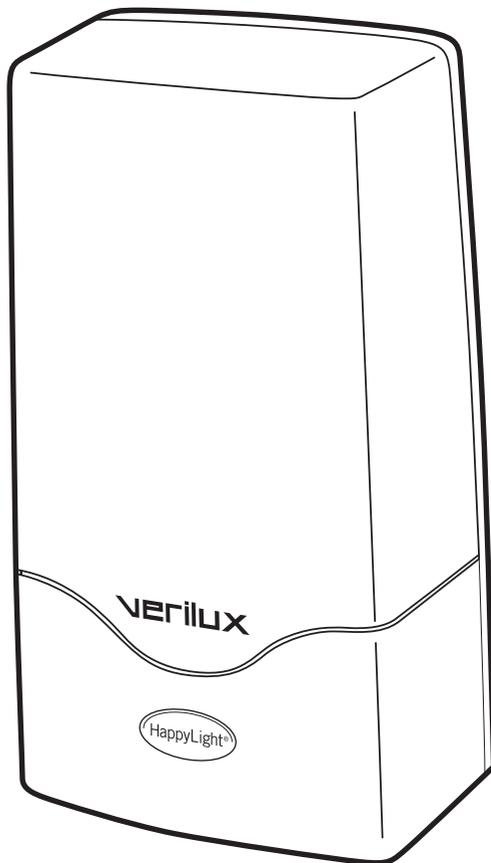**verilux**<sup>®</sup>  
The Healthy Lighting Company<sup>®</sup> Est. 1956

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