




VITALITY CHAIR

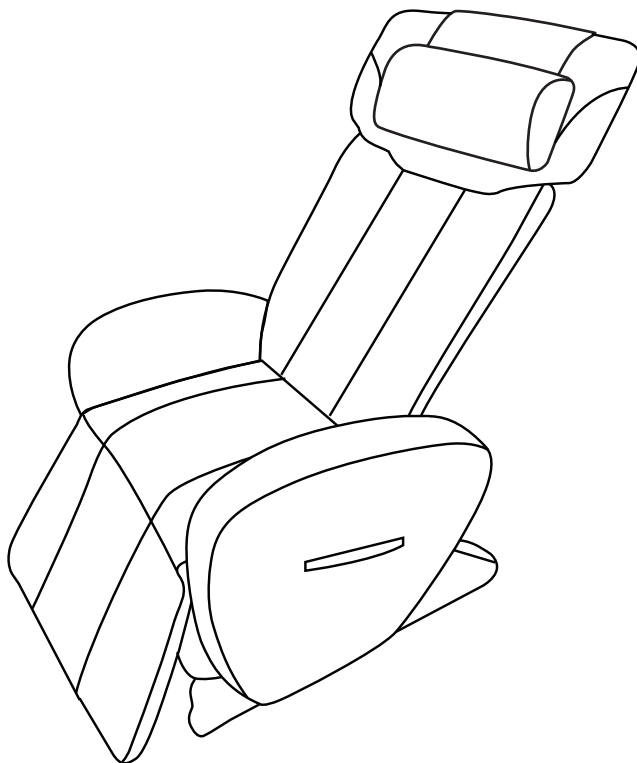

MASSAGE LOUNGER

Before operating this unit, please read these instructions completely.

* Due to continuous quality improvement, the specifications and part names in this instruction manual may be different than the actual product

Table of Contents

Thank you for purchasing the Vitality Chair.
Before using this chair for the first time, please read through all sections of this User Manual carefully, paying special attention to the Safety Instructions.

PLEASE RETAIN THIS MANUAL FOR FUTURE REFERENCE.

IMPORTANT SAFETY INSTRUCTIONS	1
PARTS AND ACCESSORIES	3
SETUP	4
Headrest	6
Backrest	7
Pillow	9
Power Switch	10
Transport	11
Grounding	12
REMOTE CONTROL	13
OPERATION	14
Auto Programs	14
Bluetooth Connection	15
Manual Options	16
FEATURES	18
PRODUCT STORAGE AND MAINTENANCE	19
LIMITED WARRANTY	20

Important Safety Precautions

- ▶ Basic precautions should always be followed when using an electrical unit, including the following.
- ▶ Please ensure that you read all instructions before using the VITALITY CHAIR.

⚠ DANGER – To reduce the risk of electric shock:

Always unplug this unit from the electrical outlet immediately after using and before cleaning.

⚠ WARNING – To reduce the risk of burns, fire, electric shock, or injury:


1. Always insert the power plug all the way into the outlet to reduce the risk of short-circuit and fire.
2. The unit should never be left unattended when plugged in. Always unplug the unit when not in use.
3. The unit is not to be used by people who cannot move or communicate on their own.
4. Only use the unit for the purpose described in these instructions.
5. Do not stand on the unit.
6. Always return the seat to the upright position and make sure the legrest is completely retracted.
7. Make sure that nothing is obstructing the legrest while retracting. Keep children away from the legrest.
8. This unit should not be used by children. Do not allow children to play on or around the unit.
9. Never operate the unit if it has a damaged power cord or plug. Return the unit to the nearest authorized service center if it is not working properly, if it has been dropped or damaged, or if it has been immersed in water.
10. Keep power cords away from heated surfaces. Do not carry the unit by the power cord.
11. Do not use the unit on top of heating appliances. Do not operate the unit under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury.
12. Do not use the unit in bathrooms or other damp or humid places as this may cause electric shock or cause the unit to malfunction. Do not spill water or any liquids onto the controller.
13. Do not use the unit outdoors.
14. Do not use or operate where aerosol spray products are being used or where oxygen is being administered.

Important Safety Precautions


15. When unplugging the unit, turn off all controls before removing the plug from power outlet. Do not use the unit with a transformer as this may lead to malfunction or electric shock.
16. Connect this unit to a properly grounded outlet only. See Grounding Instructions.
17. Persons with any of the conditions listed below or who are undergoing medical treatment or suffer from medical conditions should consult with their physicians before using the unit:
 - Pregnant women, illness, and/or in poor physical shape.
 - Pacemaker, defibrillator, or other personal medical devices prone to electronic interference.
 - Are prohibited from receiving massage by a physician due to thrombosis or aneurysm, acute varicose veins, or other circulatory disorders.
 - Heart disease.
 - Back, neck, shoulder, or hip pain.
 - Irregular curvatures of the spine.
 - Loss of heat sensation.
18. Do not use the unit against bare skin. Exposing the skin directly to the massager may irritate the skin.
19. Do not use the unit while wearing any hair accessories. Do not use the unit for massaging the head, abdomen, neck, elbows or knees.
20. Do not use the unit after drinking alcohol.
21. When the backrest is in a reclined position, do not sit on it or sit with your legs resting on the headrest. Do not sit on the legrest as the unit may fall over.
22. Do not allow children or pets to play on or around the unit, especially during operation.
23. Do not unplug the unit during operation. Unplug the unit immediately if there is a power outage.
24. Keep hands and feet clear of mechanism. Only the occupant should operate it.
25. Empty pockets of loose change and small items.

Parts and Accessories


Front View

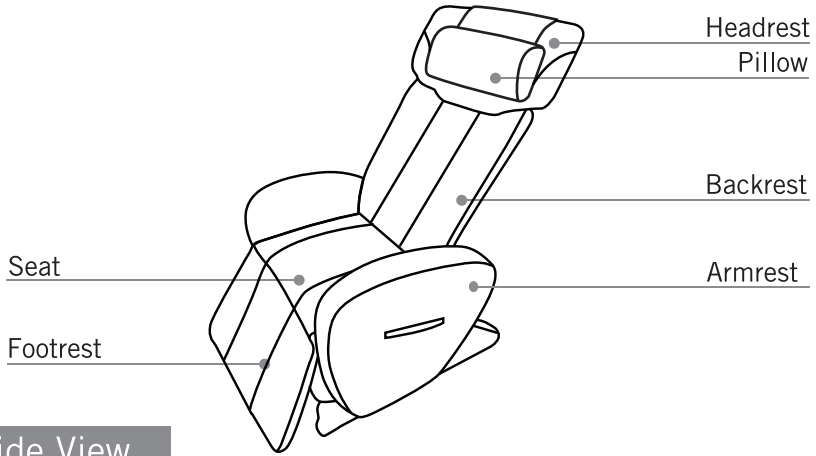


Side View

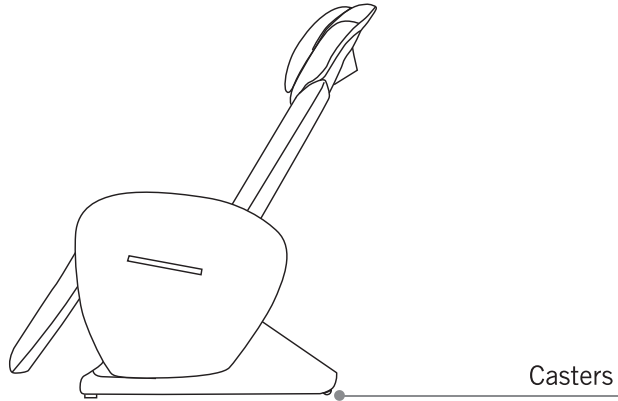


Accessories

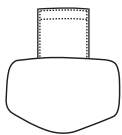

Pillow

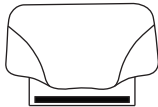

Headrest

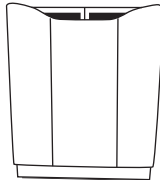

Backrest



Power Cord




Instruction Manual

STEP 1:

1. Remove all of the smaller miscellaneous parts from the package carton.
2. Removing the body of the chair from the carton requires two or more people, as shown here:



GRIP POINT
(armrests on both sides)

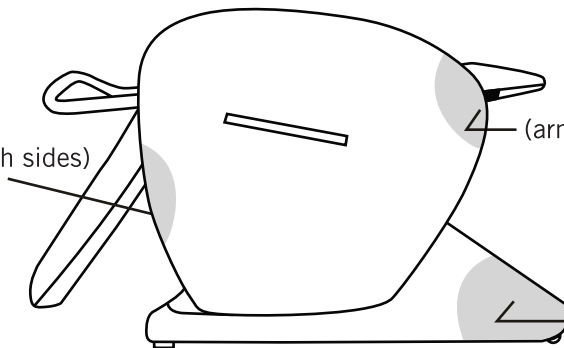


GRIP POINT
(armrests on both sides)

GRIP POINT
(casters)


STEP 2:

1. Set the chair down onto a flat, level surface at your desired location.
2. Raise the backrest to an approx. 115-120° angle, until you hear a “click” sound. Then lock the backrest into position, as shown below:

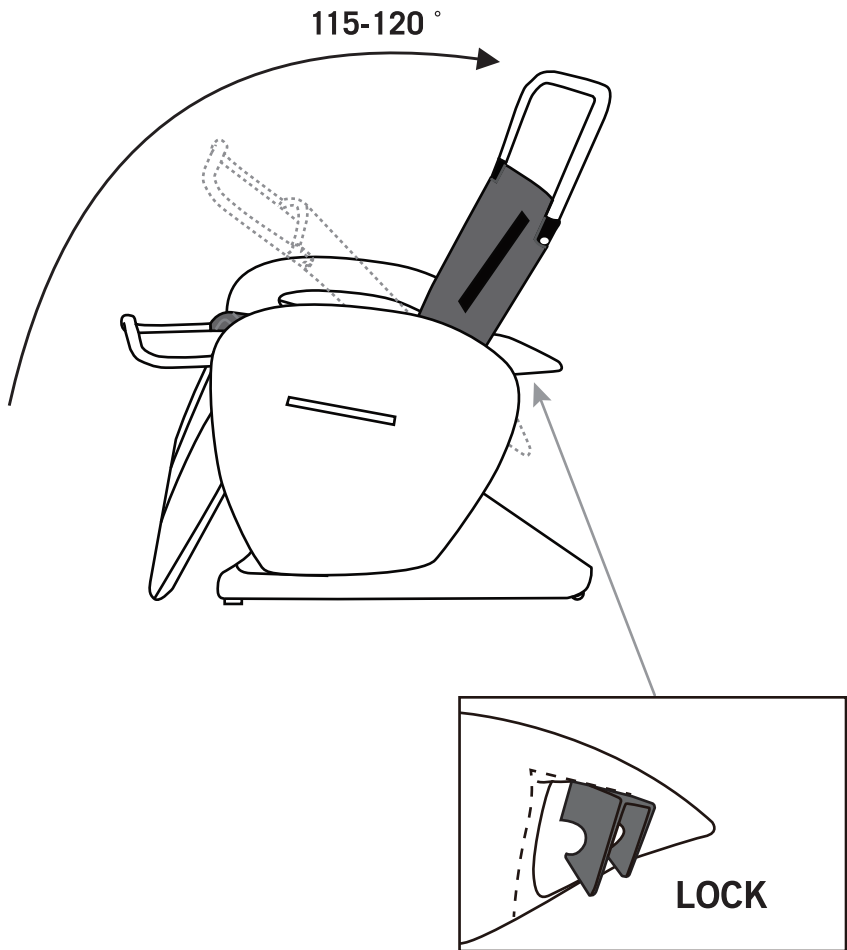


STEP 3:

1. Place the speaker cable in front of the metal frame.

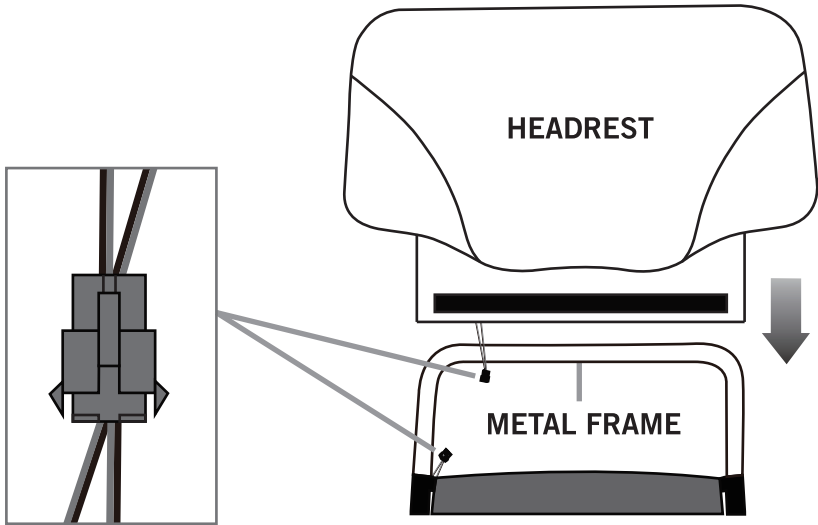


2. Plug both ends of the speaker cable into the connector.

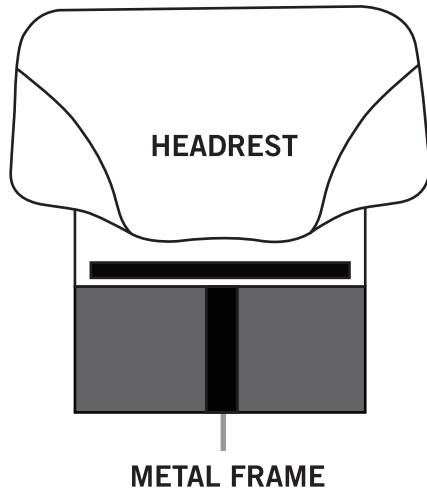


2. Install the headrest through the metal frame.


STEP 4:

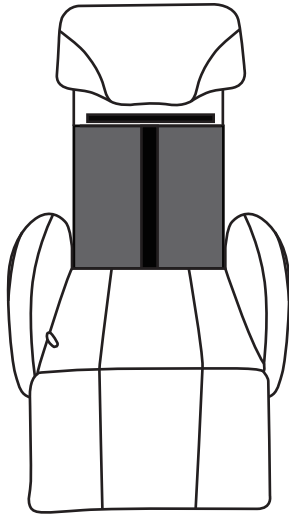

1. Unzip the backrest and place it over the metal frame.

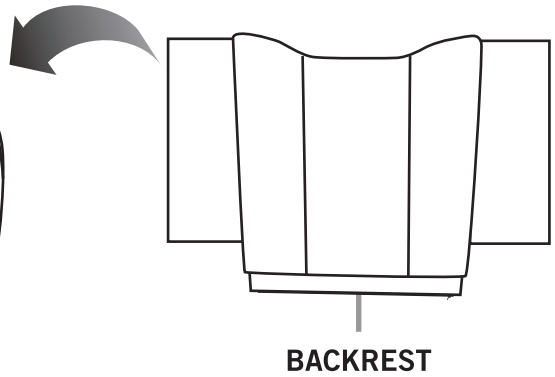


2. Zip up the backrest to the seat.
(See picture below)

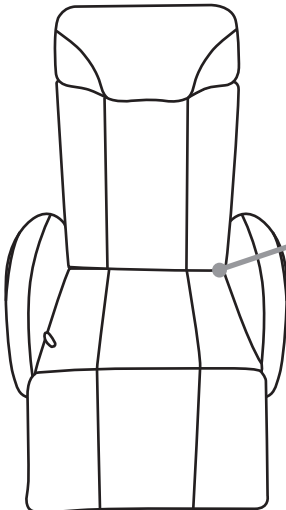


STEP 5:

1. Zip the backrest



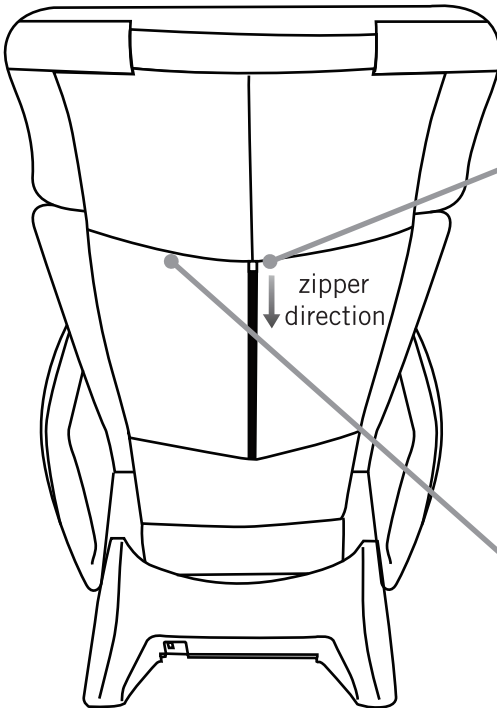


1. Adjust the backrest to the back of the headrest.

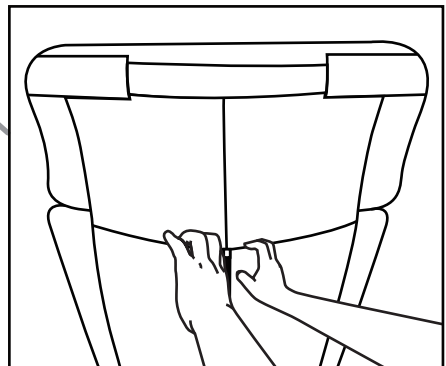


STEP 6:

1. Insert the pillow strap through the gap of the chair's headrest.

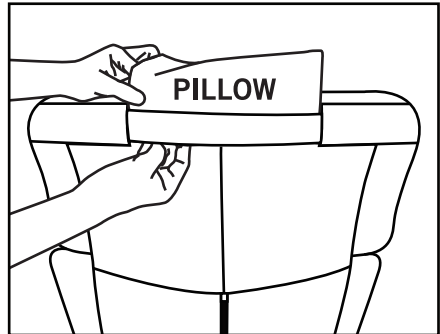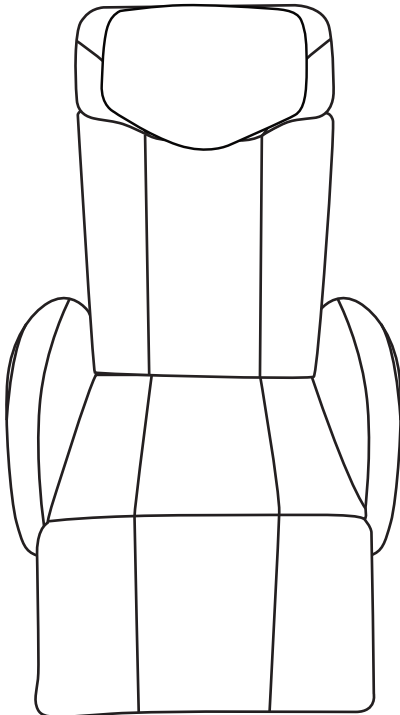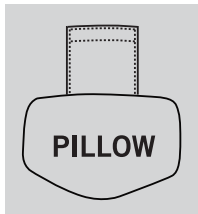

2. Pull the strap to adjust the pillow to suit the height of your head.

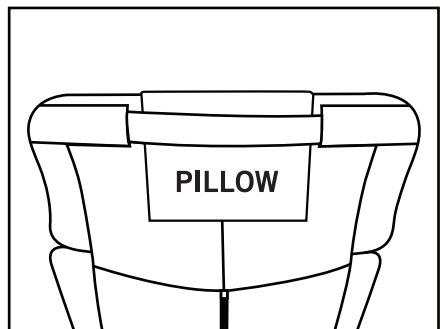


Setup - Power Switch

STEP 7:

1. Inspect the power cord to make sure that no part of it is broken or damaged.


STEP 8:

1. Connect one end of the power cord to the connector at the lower left corner of the chair's back end. Connect the other end of the cord to a standard electrical outlet.
2. Press the power switch to turn on power.

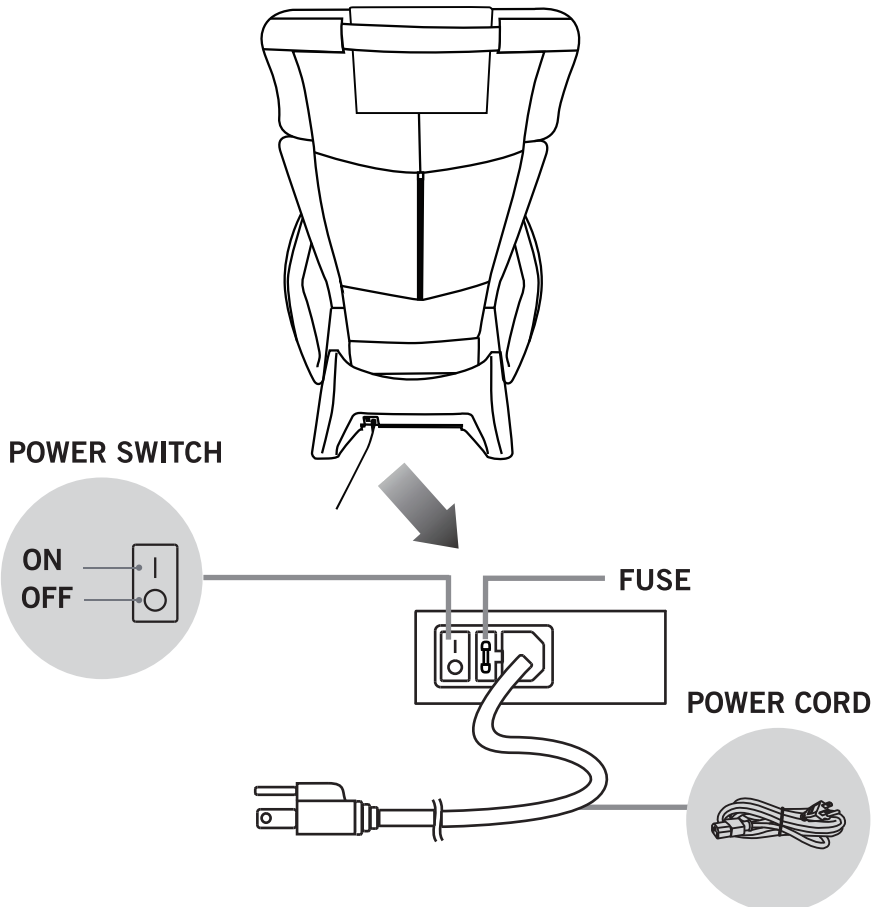

MOVE METHOD:

Transporting the chair requires two people working together. Set the chair in its upright position, with one person on each side grasping an armrest. Carefully move the chair to your desired location.

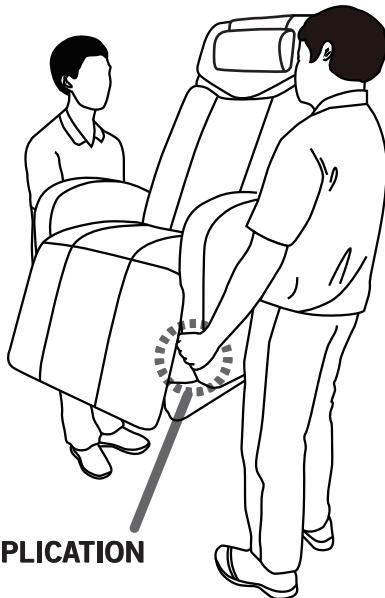


WARNING

Always use the method described above to move the chair. Pushing or dragging the chair will damage the bottom of the unit, and may cause malfunction.

GROUNDING INSTRUCTIONS

This product must be grounded in the event of an electrical malfunction. Grounding provides a low resistance path for the electric current to reduce the risk of electric shock. The cord is equipped with a built-in equipment-grounding conductor and a 3-prong grounding plug. This plug must be used with an appropriate outlet installed and grounded in accordance with all local codes and ordinances.

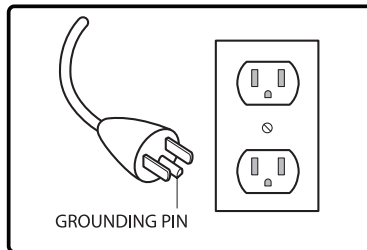

DANGER

- ▶ Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product – if it will not fit the outlet, have a proper outlet installed by a qualified electrician.
- ▶ This product is for use on a nominal 120-volt circuit and has a grounding plug that looks like the plug illustrated. Make sure that the product is connected to an outlet having the same configuration as the plug. No adapter should be used with this product.

GROUNDING

This massage chair belongs to Class I of the IEC Protection Class. Be sure to use a three-core power supply socket. To avoid electric shock, leakage, or other electric risk, ensure the ground wire connected to the power supply has obtained a stable ground connection.

Remote Control

1 POWER ON / OFF

2 RELAX PROGRAM

3 VITA RELAX

4 REFRESH PROGRAM

5 VIGOR PROGRAM

6 POWER NAP PROGRAM

7 MUSIC SYNC. PROGRAM

8 BLUETOOTH

9 VOLUME UP

10 VOLUME DOWN

11 RECLINE UP

12 RECLINE DOWN

13 AUTO PROGRAM 2 ~ 7

14 SHAKE UP

15 SHAKE DOWN

SPECIAL PROGRAM


16 BACK RELAX

17 HIP RELAX

18 LEG RELAX

19 ZERO GRAVITY

20 PAUSE

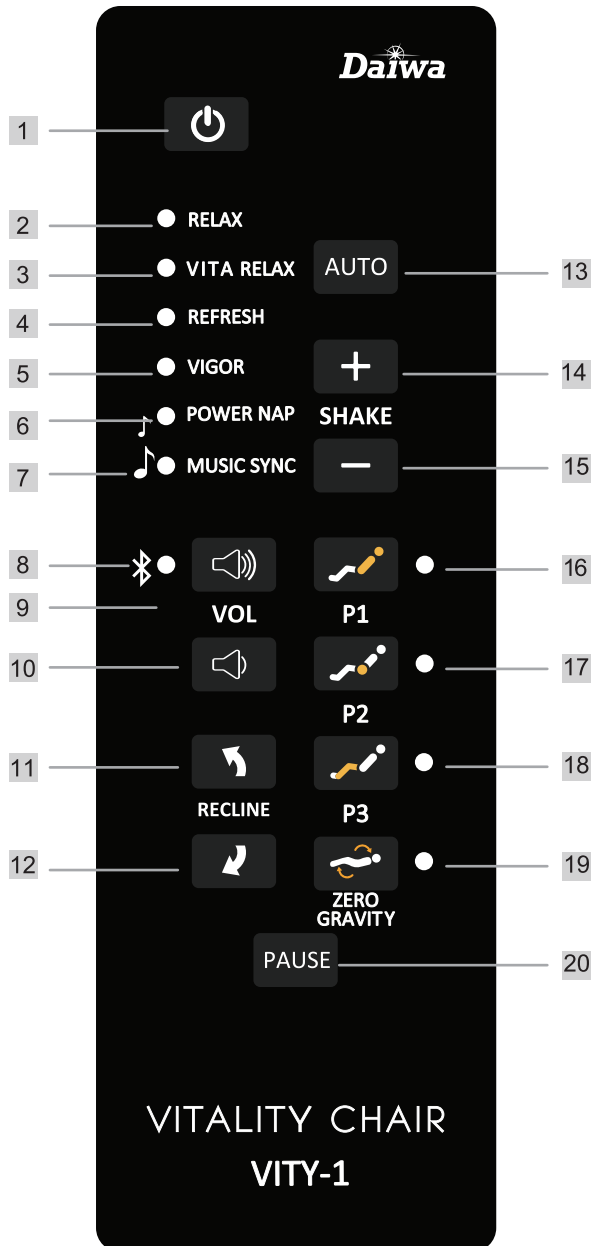

POWER BUTTON

1. Press once to turn on power.



- The chair will enter its Standby mode and await further instruction. LED lights will flash in sequence.
- You may then select your desired massage options by pressing the corresponding buttons.

2. Press again to shut off power.

- The back of the chair will automatically return to its upright position.

AUTO BUTTON

1. Press once to enter the Automatic Programs mode. All programs are 15 minutes.



- LED lights for each option will illuminate in sequence. Press AUTO button until you arrive at your desired program option.

RELAX

Induces a feeling of calm, useful for unwinding after a stressful day. Gentle shaking intensity, alternates between two comfortable recline positions.

VITA-RELAX

Relaxing and energizing. Light shaking intensity. Reclining positions and shake speed are automatically varied throughout.

REFRESH

Stimulating and reenergizing. Strong shaking intensity that leaves you with a slight tingling feeling. Reclining positions and shake speed are automatically varied throughout.

VIGOR

Massaging and stimulates muscle contractions, useful after a strenuous workout. Powerful shaking intensity that leaves you with a tingling feeling. Reclining positions and shake speed are automatically varied throughout.




POWER NAP

Calming, encourages napping. Very gentle shaking intensity creates a rhythmic rocking motion.

MUSIC SYNC

Motions based on the rhythm of the music connected through Bluetooth. Range of gentle to strong shaking intensity in sync with your chosen songs.

PAIRING YOUR BLUETOOTH DEVICE:

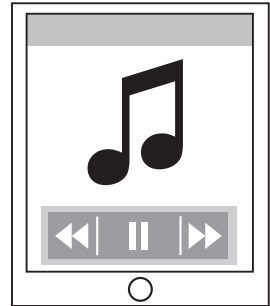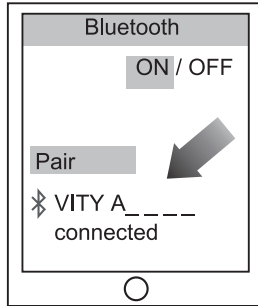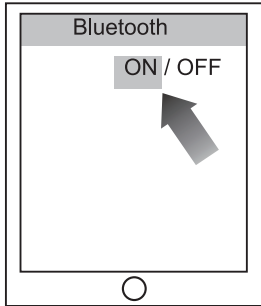

1. Turn on Bluetooth on your smart device.
2. Your Vitality Chair's unique bluetooth device ID will be displayed as VITY A_ _ _ _ , followed by the first 4 digits of the chair's serial number.
3. For the best listening experience, adjust the volume to the highest setting on your smart device.
4. Press [VOL + or -] on the remote control to adjust the volume that comes out the speakers.
5. While in a manual program, press [SHAKE + or -] on the remote control to adjust your preferred shake intensity.



NOTES:

1. If any outside electromagnetic interference is present, it may require 2-3 attempts before pairing.
2. This chair can be paired with only one Bluetooth device at a time.
3. If your smart device has an existing pairing before connecting, delete this pairing before reconnecting with the Vitality Chair.
4. This chair also provides a standard headphone jack to connect with a smartphone, MP3, or other devices for playing music.

VOLUME ADJUSTMENT

10 Volume levels are available while you enjoy your Bluetooth device.



VOL



1. Press the [VOL +] button to increase volume.
 - You will hear two beeps when the highest volume is reached.
2. Press the [VOL -] button to decrease volume.
 - You will hear two beeps when the lowest volume is reached.

MANUAL OPTIONS

This chair features powerful shaking massage action. Each session lasts 15 minutes, during which 14 levels of intensities can be adjusted.




1. Press [SHAKE +] to increase shaking intensity during operation.
2. Two beeps will sound when the highest intensity is reached.
3. Press [SHAKE -] to decrease shaking intensity during operation.
4. Two beeps will sound when the lowest intensity is reached.

NOTE: The load detecting function will significantly decrease its shaking intensity if you leave the chair during operation. The chair will automatically stop all vibration after 5 seconds. If another program is not selected within 3 minutes, the chair will turn itself off.


ADJUSTING BACKREST POSITION

Raising the backrest



1. Press and hold  to raise the backrest. Release this button when the chair reaches your desired position.
2. As the backrest rises, the footrest simultaneously lowers.
3. Two beeps will sound when the chair reaches its farthest original upright position.
4. To quickly set the chair to its full upright position, press this button twice in a row.

Reclining the backrest

1. Press and hold  to lower the backrest. Release this button when the chair reaches your desired position.
2. As the backrest lowers, the footrest simultaneously rises.
3. Two beeps will sound when the chair reaches its farthest reclined position.
4. To quickly set the chair to its furthest reclined position, press this button twice in a row.

MANUAL PROGRAMS



In order to help evenly distribute the users weight across the chair, the Vitality Chair features a zero-gravity reclined position. While in this position, you may select any shake speed. This position may be used with or without shaking massage.



P1

Targets your upper and lower back. This program runs for 10 minutes. Reclining positions and shake speed are automatically varied throughout.



P2

Targets your lower back, hips, and other muscles below the waist. This program runs for 10 minutes. Reclining positions and shake speed are automatically varied throughout.



P3

Targets your buttocks and thighs. This program runs for 10 minutes. Reclining positions and shake speed are automatically varied throughout.

PAUSING OPERATION

PAUSE

1. Press once to temporarily stop all functions. The five program indicator lights on the remote control will flash until further instructions.
2. Press this button again to resume your previous operation.

NOTES: MAXIMUM USER WEIGHT ON CHAIR IS 265 LBS.

1. All operation will stop automatically if excessive weight is loaded onto the chair for more than 15 seconds. A buzzer will sound every other second until the load is removed.
2. If any malfunction occurs in the motor, operation will cease after 3 seconds, and a beeping sound will be heard to alert the user. At this point, please turn off power and contact an authorized dealer or service agent for assistance.

SPECIAL FEATURES:

- ▶ Zero gravity
- ▶ Surround sound speakers
- ▶ Bluetooth
- ▶ Headphone jack
- ▶ Remote control pocket holder
- ▶ Extra padded pillow

AUTOMATIC PROGRAMS:

- ▶ Power Nap
- ▶ Relax
- ▶ Refresh
- ▶ Vita Relax
- ▶ Vigor

BENEFITS OF VITALITY CHAIR:

The Vitality Chair has a variety of settings that can relax you with gentle rocking motions or it can revitalize you with a vibrating massage. It is a form of Whole Body Vibration (WBV, also known as Biomechanical Stimulation), which has been used for over a decade by fitness trainers as a form of exercise and massage. Explore some of the benefits below.

Relaxing and Stimulating

At its most gentle, the combined vibration and shaking effect feels like a soothing rocking motion that is immediately relaxing. Whole-body vibration can increase mood-boosting brain chemicals and decrease stress hormones.

Relief from Joint Pain/Pressure

The Vitality Chair fully reclines into a zero-gravity position, which relieves pressure on your back and joints. The vibration massage effect can help further relieve painful conditions by creating numbness in afflicted areas. Regularly reclining in this position also helps promote circulation to the legs and feet, which can help to relieve or prevent varicose veins.

Fat Burning and Cellulite Reduction

Whole Body Vibration stimulates rapid micro-contractions of the muscles, which burns calories and acts as a kind of passive exercise. Whole body vibration can cause energy expenditure, increase the uptake of oxygen in the blood, and increase flexibility of the muscles. Vibration can also temporarily reduce the appearance of cellulite, by flattening out fat cells, giving skin a smoother appearance.

Reduced Loss of Muscle and Bone Density

Whole-body Vibration can reduce muscle loss and bone density loss in people who cannot move freely. This is also very beneficial for people who do not get enough physical exercise.

Massage

The combined effects of shaking and vibration throughout the body massages muscles, reducing soreness, boosting circulation, and increasing flexibility. If you are out of shape, the Vitality Chair can loosen stiff or tight muscles, making you less prone to injuries.

CLEANING

- Before cleaning, switch off power and remove the power plug from the wall outlet.
- Dip a soft, cotton cloth in neutral detergent and wring dry. Use it to wipe the pillow cushion, back cushion, armrests, legrests, and any plastic parts of the chair.
- Use a piece of wet cloth (without detergent) to re-wipe the pillow cushion and back cushion. Place cushions in a well-ventilated location for drying.
- Use a soft, dry cloth to re-wipe the armrests, legrests, and any plastic parts of the chair.



CAUTION

- Before any maintenance, always disconnect the power plug from the wall outlet.
- To avoid electric shock, never touch the power plug with wet hands.
- Do not use abrasive chemicals, such as benzene or alcohol, to clean this chair. Damage to the chair may result.
- DO NOT apply an iron or any other heated device to the chair.

STORAGE

Before putting into storage, wipe off all accumulated dirt and dust from the surface of the chair. Cover the entire surface of the chair with a dust cover (not included).



CAUTION

Do not store the product under direct sunlight, nor in a location with high temperatures. Doing so may cause the upholstery to fade.

DANGER

Immediately turn off power, disconnect the power plug from the wall outlet, and contact the nearest authorized service center if:

- A burning smell is detected.
- The power cord feels hot.
- Any other abnormal symptoms.
- The cord is plugged in and the power switch flipped on, but no functions operate.

VITALITY CHAIR VITY-1

Consumer Limited Warranty: ONE-YEAR WARRANTY

This warranty guarantees that the product purchased is to be free from defects in material and workmanship under normal usage and local voltage for a period of one year. This warranty expires in **TWO (2) years on parts, and ONE (1) year on labor** after the date of purchase.

U.S. Jaclean, Inc. (referred to as “the Warrantor”) will repair this product with new or rebuilt parts in case of defects in material or workmanship, free of charge, in the USA provided that:

- 1) The product is operated in accordance with the printed instructions supplied.
- 2) The product is used under normal conditions for which designed.
- 3) The product is not subjected to misuse, negligence, or accident.
- 4) The product receives proper care, protection, and maintenance.

This warranty is extended to the original end-user when purchased from an authorized U.S. Jaclean, Inc. distributor. A purchase receipt or other proof of date of the original purchase is required for service and parts replacement under this warranty. This warranty does not cover shipping charges for defective items if a claim is made more than 30 days after the date of purchase.

This warranty only covers failures due to defects in materials and workmanship, and does not cover normal wear or cosmetic damages. This warranty does not cover damages which occur in shipment, or failures which result from: accident, misuse (including, without limitation, failure to comply with instructions and use limitations contained in the Operating Instructions), abuse, neglect, mishandling, faulty installation, set-up adjustments, improper operation or maintenance, alteration, modification, power line surge, lightening, improper line voltage, commercial use such as hotel, office, restaurant, or other business or rental use of product, service by anyone other than from U.S. Jaclean, Inc., or damage that is attributable to “acts of God.”

Continued use of this product after discovering a defect voids this warranty. This warranty shall be governed by the state of California, United States of America.

IMPORTANT: DO NOT return defective product to the store.

Contact U.S. Jaclean, Inc. for repair services.

Please retain this Warranty Certificate, together with the original purchase receipt, and present it to our technician for inspection during service.

For assistance and more information, please contact U.S. Jaclean, Inc.

Online warranty registration at: daiwamassage.com



Distributed by U.S. Jaclean, Inc.
1816 West 135th Street
Gardena, CA 90249

1-310-525-2700
9AM - 5PM PST (Monday - Friday)

daiwamassage.com