



# User Manual

## SEDONA

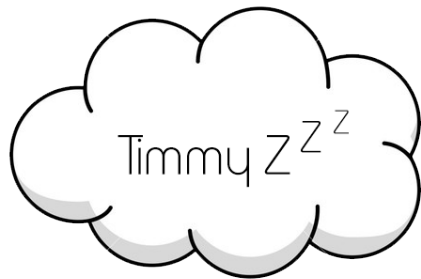

## The SEDONA Remote

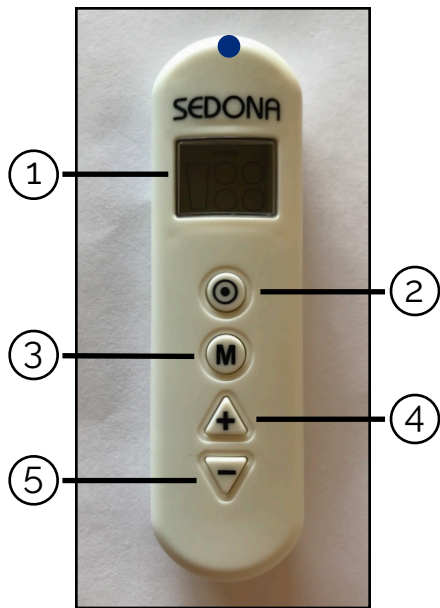

1. Display

2. Power ON and OFF

3. Mode (select program)

4. + Key (increases intensity)

5. - Key (reduces intensity)

# Power Supply

Open the battery cover; place the two AAA batteries into the compartment; close the battery cover. After press the ON/OFF button (if the device doesn't turn on then batteries are empty, please swap the batteries for two new AAA batteries.)



1. Open the battery cover

2. Place the two AAA batteries into the compartment

Another power source is to plug the control remote into an outlet



Plug the cable into the top of the remote

# How to Start Device


Press the ON/OFF button to boot up the unit. The display panel back light will turn on. Then press the M (mode) button to select the desired program. Press "+" or "-" to set the intensity. (Make sure the unit is connected to the pillow otherwise the unit will automatically shut down)

Make sure inductor is connected




Press the ON/OFF button

Press the M button until desired mode

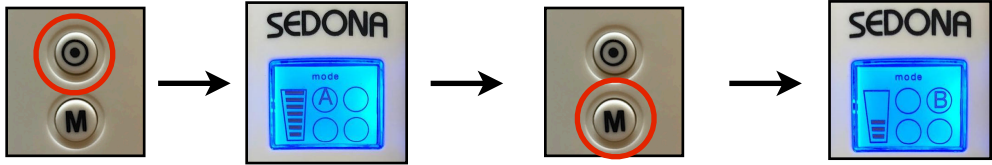

The selected program will start automatically a few seconds after you have reached your desired intensity and run for 30 minutes

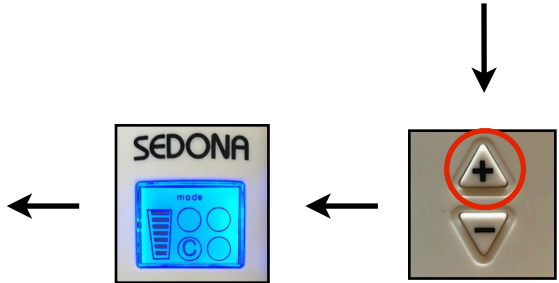


Press the "+" button till desired intensity


To have the program run indefinitely (without stopping)

Press the ON/OFF button

Press the M button and "+" button at the same time

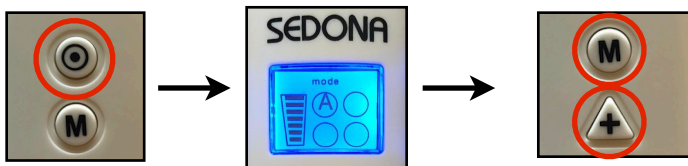

You will notice the lock button in the upper left corner

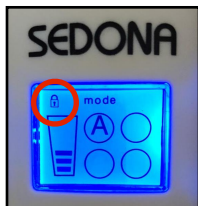

Now select desired program and intensity

Program		Frequency	Intensity	Use
A	Very Low	3Hz	0.6 gauss	Deep Sleep
B	Low	7.8Hz	0.6 gauss	Restorative Sleep/ Schumann frequency
C	Medium	25Hz	0.6 gauss	Awakening/Active Energy/ Concentration
D	High	1,200Hz	0.6 gauss	Pain Relief



SEDONA WELLNESS Disclaimer: The information on this manual is for reference purposes only and is not intended to recommend our pulsed electromagnetic field therapy device system as a drug or as a diagnosis for any illness or disease condition; nor as a product to eliminate disease or other medical condition. The Information has not been evaluated by U.S. Food and Drug Administration. Worldwide, there are no governmental health agencies that recognize a need to supplement natural magnetic fields using pulsed electromagnetic fields. Also, the website, company, employees, practitioners or its distributors make no warranty of any kind, expressed or implied with regard to the Information or how you use it. Sedona Wellness Trading Co. LLC makes no medical claims, real or implied, as to benefit of our device and methods. Our product is not intended to be used to diagnose, treat, cure or prevent any disease. Readers should consult appropriate health professionals relating to their health and well-being. Readers accept all responsibility for self-experimentation.

SEDONA'S "TimmyZzz" IS STRICTLY A SLEEP ENHANCEMENT TOOL

Copyright © 2024 TimmyZzz

