

Please Note:

- The total time required for treatment with the HairMax LaserComb is 8 minutes for Professional 12 and 11 minutes for Lux 9 models.
- It is important to move the LaserComb very slowly. Remember you are bathing the scalp with laser light. If you move it too fast, the treatment will not achieve the optimal results to promote hair growth.
- While there are no negative effects from overuse, it is recommended that you use the LaserComb according to directions.
- The LaserComb has advanced electronics to make sure that the device is working properly. If for some reason you are experiencing difficulties with your device, before sending it back to us for repair, please refer to the Trouble Shooting section of your User Manual.

For further assistance please call +1-561-417-0200.



REGISTER

To validate your warranty
please register online at
www.hairmax.com/warranty
or call Customer Service
toll-free at 1-866-527-3726

CONTACT

Lexington International, LLC
777 Yamato Road, Suite 105
Boca Raton, FL 33431
+1-561-417-0200
HAIRMAX.com

HAIRMAX®
The science of hair growth.

QUICK START GUIDE

Lux 9, Professional 12 LaserComb






The HairMax LaserComb® is a very effective device when used correctly. It is manufactured using precision components and should last 10-15 years under normal usage. Please take time to read the entire User Manual before operating the LaserComb.

Lux 9, Professional 12 LaserComb Models




Use the LaserComb Lux 9 for 11 minutes, and the LaserComb Professional 12 for 8 minutes, 3 times a week. The treatments should be spaced every other day, for example Monday, Wednesday, Friday. Move the LaserComb slowly 1/2 inch every 4 seconds.

It is recommended, but not essential, that the LaserComb be used on a clean scalp after washing. Your hair can be wet, towel dried or dry. The LaserComb should be used prior to applying any topical treatments, gels, sprays or concealing products.

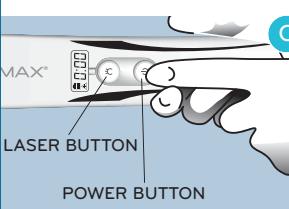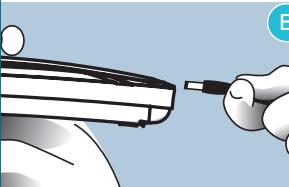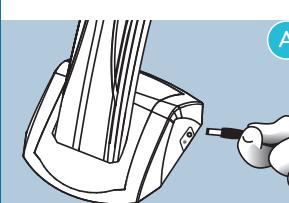

Charge your LaserComb using either the charging cradle (Figure A) or by placing the power cord directly into the handle end of the LaserComb (Figure B). For initial charging, allow 2-3 hours to fully charge your LaserComb. A green light will illuminate on the side of the charging stand indicating that your power cord is operating properly. During charging, the display on the LaserComb will illuminate for 20 seconds and the battery indicator bars will light up indicating the battery level and confirm the charging process. After 20 seconds the display will turn off but the battery will continue to charge.

To turn on the LaserComb press the power button and hold until power button light turns green. Next, press the laser button which will light up red. That means the LaserComb is ready for use (Figure C).

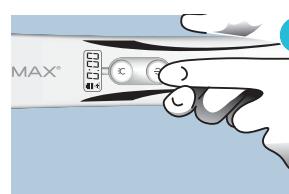

Before beginning treatment, be sure to test the LaserComb for sensitivity (see page 5 of your User Manual).

To start your treatment place the LaserComb flat on your scalp so that both rows of teeth are touching the scalp (Figure D).

Start at your hairline. Move the LaserComb from front to back of scalp in a slow combing motion, stopping every ½ inch (1cm) and pausing for 4 seconds (Figure E). Do not lift the LaserComb off your scalp while combing from front to back.

Next, lift the LaserComb off your scalp and start at a new position at the hair line. Continue in this manner from front to back and then side to side (Figure F) until you cover your entire scalp with multiple passes. When the LaserComb is on it will beep and vibrate every four (4) seconds. This means that it is time to move the LaserComb to the next area of the scalp. Repeat these passes for the total treatment of 8 min. for Professional 12 and 11 min. for Lux 9 models.

If you have long hair, use your fingers to lift your hair out of the way while using the LaserComb to be sure that the laser energy reaches your scalp (Figure G).

After treatment, turn off the LaserComb by pressing and holding the power button (lower button) until green indicator light turns off (Figure H). Store the device in the charging cradle for future use. You may leave the power cord plugged into the charging cradle at all times.