

## Please Note:

- The time required for treatment with the HairMax® LaserBand is as little as 3 minutes.
- Your HairMax turns off automatically after 12 minutes.
- Your HairMax has advanced electronics to ensure the device is working properly. If for some reason you are experiencing difficulties with your device, before sending it back to us for repair, refer to the troubleshooting section of your User Manual on page 14.
- For further assistance, contact HairMax Customer Service at the numbers provided below.

### Corporate Office

Lexington Intl., LLC  
777 Yamato Road, Suite 105  
Boca Raton, Florida 33431 USA

### Customer Service – Warranty Repairs

Lexington Intl., LLC  
1040 Holland Drive, Suite 1  
Boca Raton, Florida 33487 USA

**Direct line** +1 561-417-0200

**Toll Free** 1-866-527-3726

Monday - Friday 9 AM - 5 PM EST

**United States** 1-866-527-3726

**United Kingdom** 0808-134-9868

**Australia** 1-800-937-358

### Register

To validate your warranty, register online at [www.hairmax.com/warranty](http://www.hairmax.com/warranty)


Patent [www.HairMax.com/patents](http://www.HairMax.com/patents)

Made in China

© Copyright 2016 Lexington Intl., LLC.

All rights reserved.

HairMax laser devices are indicated to treat Androgenetic Alopecia, and promote hair growth in males who have Norwood Hamilton Classifications of IIa to V and in females who have Ludwig (Savin) I-4, II-1, II-2, or frontal patterns of hair loss and who both have Fitzpatrick Skin Types I to IV.

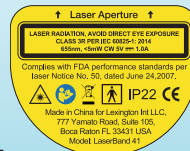

**HAIRMAX**  
The science of hair growth.



HairMax®  
LaserBand 41

QUICK START GUIDE

The HairMax® LaserBand is a highly effective device when used correctly. It is manufactured using precision components and should last 10-15 years with normal use. Please take time to read the User Manual before operating your HairMax.

## LaserBand 41 Model



Use your HairMax for as little as 3 minutes, 3 times a week, ideally every other day (for example, Monday, Wednesday, Friday).

### For maximum benefit:

- Scalp should be freshly washed (hair can be wet, dry or damp)
- Use your HairMax prior to applying styling products (these can build up and impede the laser light)
- Regularly clean the laser area with a damp cloth. Do not use harsh detergents or cleaning products on any part of your HairMax

### Charging Options

#### Charge your HairMax for 3-4 hours before initial use.


Prior to charging your HairMax, make sure your HairMax is turned off.

**OPTION 1** Plug the power cord into the charging stand and the other end into an electrical wall outlet (Figure A). Place your HairMax on the charging stand and wait until the indicator light starts flashing. If the indicator light does not start flashing within 5 seconds, reposition your device on the charging stand.


**OPTION 2** Place the power cord directly into the charging socket on the bottom left side of your HairMax and the other end into an electrical wall outlet (Figure B).

The indicator light will flash during the charging process. A steady light indicates your HairMax is fully charged (Figure C).

**Do not use your HairMax while it is plugged into the power cord.**

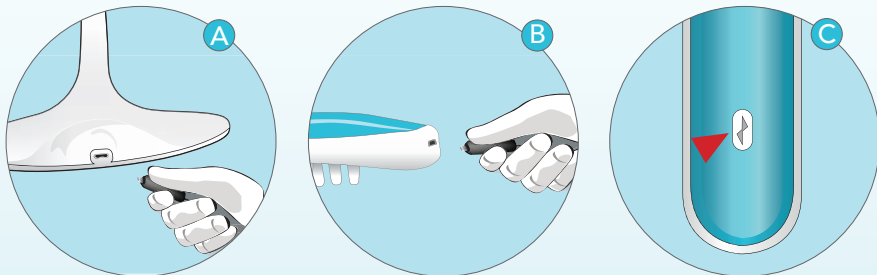

**Prior to treatment**, test your HairMax for sensitivity (Figure D). For more information, see page 5 of the User Manual.

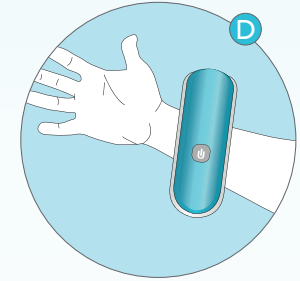


### Treatment Method

Your HairMax will be used in 6 positions on your scalp. Leave it in place for a minimum of 30 seconds for each position. It will vibrate every 30 seconds to let you know when to move it to the next position.

To turn on your HairMax, press and hold the power button until the red lasers turn on (Figure E). Next, place it on your frontal hair line. This is position 1. When you feel the vibration, slide it backwards about 2 inches (5 cm) to position 2. Leave it in position 2 until you feel the vibration again.

Repeat this process for all six positions on your scalp. The entire treatment should be a minimum of 3 minutes.

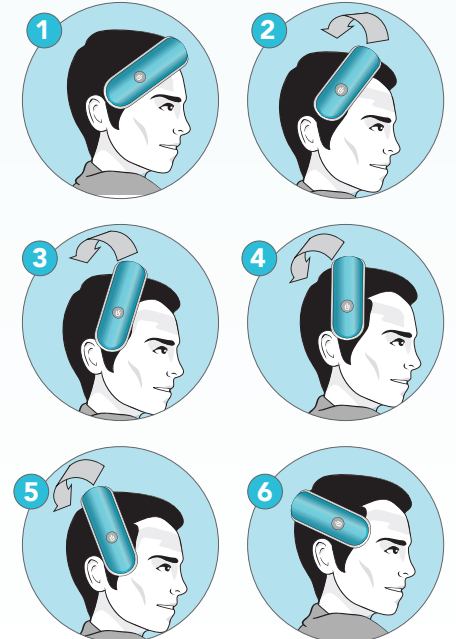

**After treatment**, turn off your HairMax by pressing and holding the power button until the lasers turn off (Figure E). Remove your HairMax from your head and place it on the charging stand for future use.

