



58.6 kg

Bluetooth Smart Body Scale

USER'S MANUAL

glucology™

CONTENTS

1. Getting Started	3
2. Using the Scale	4
3. Features	6
4. Connecting scale to app	8
5. How the app works	9
6. 12 Parameters being measured	12
7. Technical Details	13
8. Support and Safety	14



GETTING STARTED

The Smart Body Scale

Thank you for purchasing the Glucology™ Bluetooth Smart Body Scale. This product provides body measurements that are based on 12 parameters (refer to page 12), using Bioelectrical Impedance Resistance Technology, so that you can remain better informed about your health.

This manual will provide guidance with how to best maintain and use the smart scale alongside the accompanying FitGo app.

The information contained in the app is for general information purposes only and is not meant to substitute professional advice or health-related treatment. You should not use any information gained from the app to make health-related decisions. The provided body data should be used at your discretion and is not intended to diagnose, treat, prevent or cure any disease.

USING THE SCALE

For First Use

Ensure you install 2 X 1.5V AAA batteries by opening the battery door at the bottom of the product.

Place the scale on a hard and flat surface to ensure accurate measurement.

The scale will either auto turn on or press ON to get started.

Connect the "FitGo" app to your new scale (refer to 'Connecting Scale to App' instructions).

Clean the scale with a damp cloth and try to avoid chemical abrasive cleaners. Do not put in dishwasher. This scale is not waterproof.

Store scale in a cool and dry environment.

USING THE SCALE

Reminder

Do not use the scale on people who have body implants such as pacemaker, artificial limbs, contraceptive devices, metal plates or screws. It may cause the device to malfunction or produce an inaccurate result. If in doubt, please consult your physician.

Make sure to have good contact between your bare feet and the metal plates on the scale, otherwise only data on your weight will be displayed.

The appliance is intended for domestic use within the home.

Children should be supervised to ensure that they do not play with the appliance.

FEATURES



LED Invisible Window

Displays weight

High Precision Sensor

Metal electrode plates that send a weak electric current flow to calculate resistance of the body for accurate body data



Smooth tempered glass panel

The thick glass provides durability and strong load bearing capacity



FEATURES

Turning scale on/off

The scale will automatically turn on when the connected app is in use and a person is standing on the scale.

The scale will automatically turn off when it is not in use/left inactive.

Weight Capacity

Note that the scale is designed for maximum capacity of 180kg.

If this limit is exceeded, the scale comes with an overload indicator to protect the machine and help it work longer.

CONNECTING SCALE TO APP

DOWNLOAD APP

Scan the appropriate QR code to download the 'FitGo' app on your device.



iOS App




Android App

HOW IT CONNECTS

Once you open the app, there is no connecting settings required. Ensure your device has its bluetooth connection on.

You can also set your target weight to customise your own health plan.


HOW THE APP WORKS

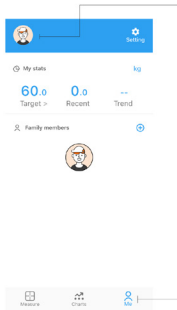

Step 1

Edit your target weight by clicking 'Me' and going to your profile.


You can also add up to 8 users in one account in the 'Family Members' section.

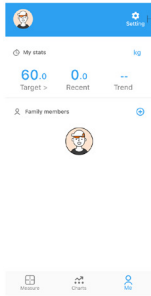

Step 2

On the same page, click 'Setting' then set unit/health index/language as you want.

HOW THE APP WORKS

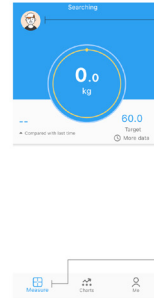

Step 3

Click 'Measure', then choose the user by clicking on the profile picture before proceeding to stand on the scale.



Step 4

Once you're on the scale, 12 kinds of measuring data will be shown on the app.

HOW THE APP WORKS




Step 5

Click 'Charts' to show the history measuring data.

Your previous measurements will be stored automatically.

12 PARAMETERS BEING MEASURED

The following body data will be measured:

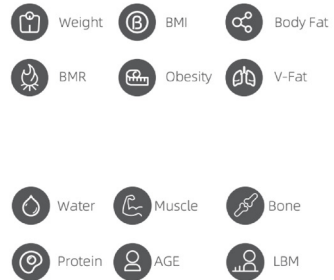

TECHNICAL DETAILS

CAPACITY:	180 kg/397 lb
MINIMAL WEIGHT:	5 kg
BATTERY:	2 x 1.5v AAA
DIMENSION:	280 x 280 x 22mm
DIVISION:	0.1 kg/0.2lb
LED DISPLAY:	75 x 31 mm
APP ACCOUNTS:	Up to 8 users
WARRANTY:	12 months

SUPPORT AND SAFETY

IMPORTANT SAFETY INSTRUCTIONS

Make sure you read these instructions, keep these instructions and follow all instructions.

Do not attempt to repair or service Glucology Smart Kitchen Scale yourself. Refer all servicing to qualified service personnel.

CONTACT

For more information, visit our website:

Australia Customers: www.ibdmedical.com.au
International Customers: www.glucology.store

OTHER CONTACT INFORMATION:

US & Canada Number: 1 800 494 2301
Worldwide Number: +61 2 8398 7520
Customer Support email: cs@glucology.store

Headquarter Address:
31/330 Wattle Street, Ultimo NSW 2007

