KiCA Pro

Massage Device Instructions

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

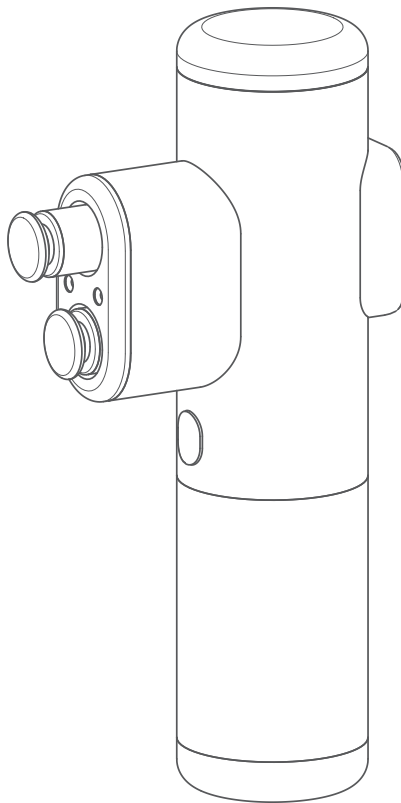
EN V1.1

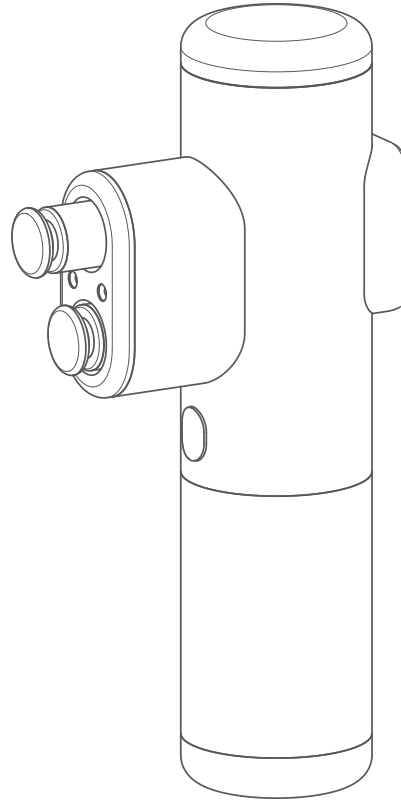Overview

- [1] Flat head (Non-removable)
- [2] Trigger button
- [3] Touch screen
- [4] Power button
- [5] Charging port
- [6] Battery indicator

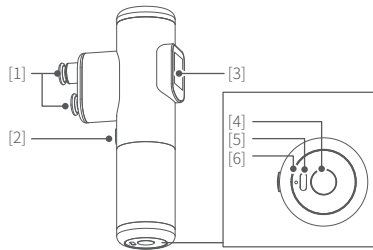Flat Head:

Massage flat major muscle groups on the back.

Gentle Head:

Promote the absorption of essential oil.

Spherical Head:

Relax soft tissues and sensitive muscles.


U-shaped Head:

Massage neck (not cervical vertebrae) and tendon.

Start to Use

1. Massage Head Combination

The device can be used directly with the built-in massage head (two flat head), on it, or combine with another one depending on your needs. See complete combination as below.

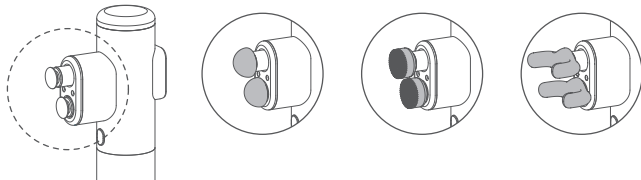
2. Power ON / Start running / Power Off

Power ON: Long press the power button until hear the tone.


Start running:

Method one: Single tap trigger button to start running, single tap again to stop.

Method two: Operate on the touch screen to select the modes and gears.

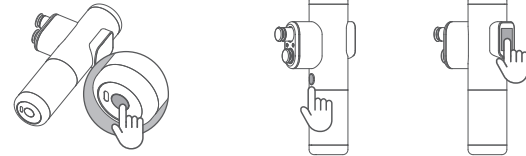Do not touch the middle and the back of the two massage heads while the product is running to avoid injuries.

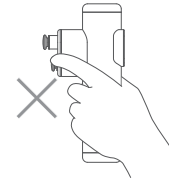Power Off:

Long press the power button until hear the tone to power off.

If there is no operation within 10 minutes, it will power off automatically.

3. Touch screen

Operate on the touch screen to select the modes and gears.

- [1] Battery level
- [2] Current rotation speed
- [3] Massage time for this mode
- [4] Total massage time for current mode
- [5] Current gear (1-10)
- [6] Current mode
- [7] Current pressure


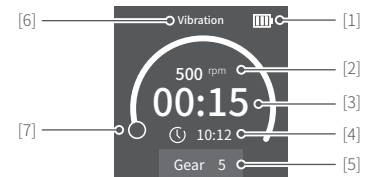
Diagram:



Basic massage modes

Setting

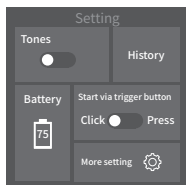Package choosing

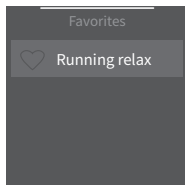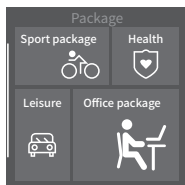Favorites

4. Basic massage modes

Mode	Rated rotate speed (rpm)	Gears
Vibration mode	500-2600	1-10
Kneading mode	100-900	1-10
Percussion mode	50-450	1-10

5. Locked Motor Protection

⚠ Do not press the massage head firmly to block the motor rotation while using, otherwise it will reduce the service life of the product!

This product has locked motor protection function. If the rotation is blocked for 3 seconds, it will stop running automatically after a beep sound. Single tap the trigger button to start again.


6. Intelligent Protection

This product has intelligent protection function. After working for 10 minutes continuously, a beep sound will be issued, and the product will stop running automatically. Single tap the trigger button to start again.


7. Charging

⚠ Do not use it while charging.


The battery indicator is red while charging, and it will turn green after charging completed. Screen shows the charging icon. When the product is on, it will return to the interface before charging, after unplug the power supply while charging.

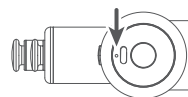Battery indicator

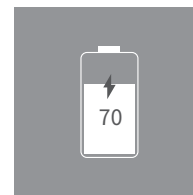Charging

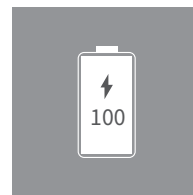Charging completed

8. Specifications

Product Name: KiCA Pro Massage Device

Product Model: A3D

Operating Voltage Range: 10.8V~12.6V

Operating Temperature Range: 0°C ~+40°C

Charging Temperature Range: 0°C ~+40°C

Battery Capacity: 2200mAh

Charging Time: About 4.5 hours

Weight: About 630g

Size: 97.2×46×188.3mm

Theoretical Battery Life: About 20 hours (Test in lab environment under the 1st gear with flat head)

List


- | | |
|------------------|----|
| ① KiCA Pro | ×1 |
| ② Storage Bag | ×1 |
| ③ Spherical Head | ×2 |
| ④ Gentle Head | ×2 |
| ⑤ U-shaped Head | ×2 |
| ⑥ USB-C Cable | ×1 |
| ⑦ Manual | ×1 |

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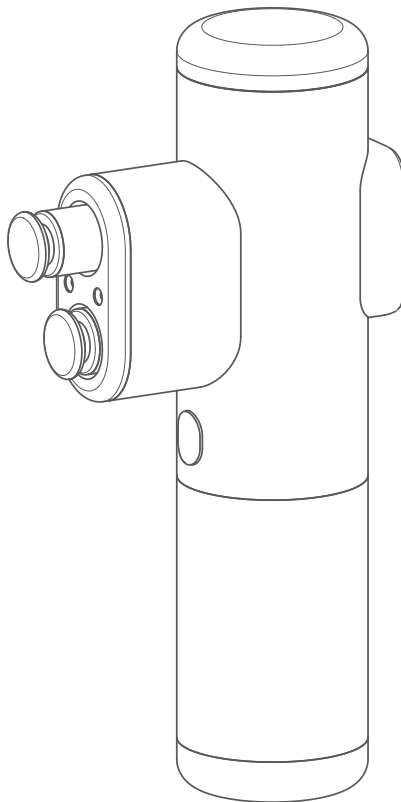Contraindications

1. Cancers
2. Metal bone anastomosis (do not use the device directly in the fracture area)
3. Osteoporosis
4. Vein thrombosis and pulsation
5. Hypertension
6. Pregnant woman
7. Varicose veins
8. Acute inflammation
9. Do not massage the bony area
10. Do not massage surgical site
11. Do not massage the cervical vertebrae and spine
12. Kidney stones (do not use the device near the kidneys)
13. Do not massage one episode for more than 5 minutes

Safety Instruction

1. It is recommend to massage your body after wearing dry clothes or placing a clean and dry towel on massage part, then gently press and move the device.
2. High-speed and high-pressure massage can only be used to soft tissues without feeling pain or maladjustments. Do not massage head and any skeletal sites with high speed and high pressure.
3. Massage the same part of the body with high speed and high pressure for a long time will.
4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the massage device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise,

please operate carefully.

5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

Storage and Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please storage the product in a dry environment.
4. DO NOT overcharge or overuse the battery, otherwise, it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.

Precautions when discarding the product

The lithium ion battery used by the product is a recyclable and renewable valuable resource.

When discarding it, please have it handled by a professional agency that can recycle chargeable batteries. Before discarding the product, power it off and take the lithium ion battery out from it following the steps below. The battery is to be disposed of safely.

1. Disassemble the top housing with help of a professional tool.
2. Pull out the batter terminals from the PCBA.
3. Use a professional cutter to cut off the bottom housing and then take out the battery.

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FCC regulatory conformance:

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

NOTE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE:

The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user's authority to operate the equipment.

Warranty Card / 售后服务卡

Product Model

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Serial Number

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Purchase Date

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Customer Name

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Customer Tel

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Customer Email

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Warranty:

1. Within one year from the date of selling, the product is malfunctioning under normal due to non-artificial reasons.
2. The malfunction of the product is not caused by artificial reasons such as unauthorized disassemble or conversion or addition.
3. The buyer is able to provide the certificate of maintenance service: the warranty card, the legitimate receipts, invoices, or screenshot of purchasing.

The followings cases are not included under the warranty:

1. Unable to provide legitimate receipt and warranty card with buyer's info.
2. The damage caused by human or irresistible factor.

More details about after-sales policy, please refer to after-sales page on website:

<https://www.kica-tech.com/en/support>

* Our company reserves the right of final interpretation on the above-mentioned after sales terms and limitations.

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Facebook

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.....

This document is subject to change without notice.

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The latest user
manual

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KiCA

Manufactured by: Guilin Feiyu Technology Incorporated Company

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Fabricant: Guilin Feiyu Technology Incorporated Company

Fabricado por: Guilin Feiyu Technology Incorporated Company

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Produttore: Guilin Feiyu Technology Co., Ltd