
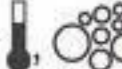
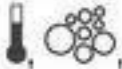
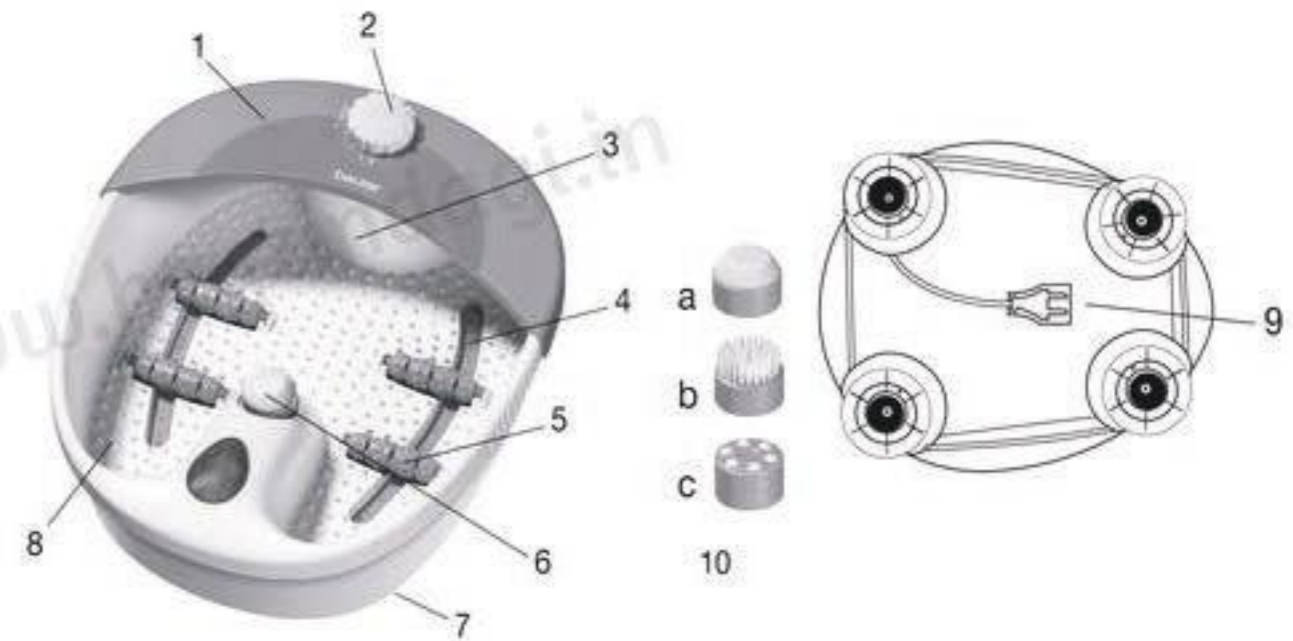


How to Use?

- Fill the foot bath with warm water no higher than the level indicator (approx. 4-5 cm)
- Ensure that the rotary switch is set to ① before the unit is connected to the mains
- Make sure the cable is routed safely so that there is no risk of anyone tripping over it
- Be seated comfortably in a chair
- Place your feet on the footrest and apply slight pressure
- Set the desired function on the rotary switch:
 1. ① Off
 2.  Vibration Massage + Infrared Light,
 3.  Heat + Bubble Massage,
 4.  Heat + Bubble and Vibration Massage + Infrared Light
- Pour the water out over the narrow side of the footbath after use

Massage roller

- Install the rollers, making sure that the roller fixing points downwards, as shown in the illustration given in the instruction manual
- Move your feet back and forth on the four massage rollers integrated into the bottom of the footbath for special massage effect
- Optionally, the massage rollers can be removed by pulling them firmly upwards



1. Splash protection
2. Rotary switch
3. Water level mark
4. Bubble strip
5. Massage attachment
6. Removable pedicure attachments a. Corn and callous remover b. Brush c. Massage attachment
7. Rubber feet (bottom side)
8. Foot rest with massage naps
9. Cord winder (bottom of device)
10. Massage rollers

REMEMBER

- Never connect the foot bath to the power supply when your feet are already in the water. If the device is damaged, this can lead to an electric shock
- Do not use the footbath if you have any changes caused by disease or injury in your foot (e.g. open wounds, warts, fungal infection)
- Do not use the foot bath for longer than 20 minutes to avoid risk of overheating. Allow it to cool for at least 15 minutes before using it again
- This device must not be used by children or persons (including children) with restricted physical or sensory capabilities without supervision

- Do not place the foot bath in the shower, bathtub or wash basin for filling
- Please do not use it for animals or uses beyond what it's been made for
- Switch the device off and disconnect the power supply after every use and before every cleaning
- Never stand up with your feet in the foot bath, as you might slip
- If the device is not working properly check the connectivity thoroughly. Send it to the Beurer Service Centre in India for repairs

P.S.: Consult your doctor before ordering this product especially for those suffering from diabetes, thrombosis, varicose veins, vein inflammations or recovering from foot surgery