



# QUICK START GUIDE

- As per our assessment, the foot warmer can accommodate feet sized up to UK 11 comfortably. Sizes 12 will fit, but struggle for comfort a bit. Kindly order accordingly.
- Please unpack the Foot Warmer and Cord
- You will need a Type G plug point or adaptor for the same. Identify a point close by that does not stretch the cord too much. Will also help in avoiding falls due to tripping
- Plug in the cord, after checking the settings on the controller are in '0' position
- Settle into a comfortable seat and make sure the floor where you are keeping the warmer is dry and cleaned well to avoiding dirtying the fabric and any electrical mishaps. You can comfortably place it on a rug too
- Please take off your shoes and other footwear. While you can use socks, we recommend not using any unless necessary for best comfort. Do not use nylon or synthetic fabrics as they can lead to burns - only cotton
- Switch the Settings as you wish for Heat and Massage.
- If you wish to use it over hours, get initial relief by setting level 2 for both and then bringing them to 1 or zero to maintain the comfort
- Once done, change the settings to 0 for both and allow the warmer to cool down before storing it.

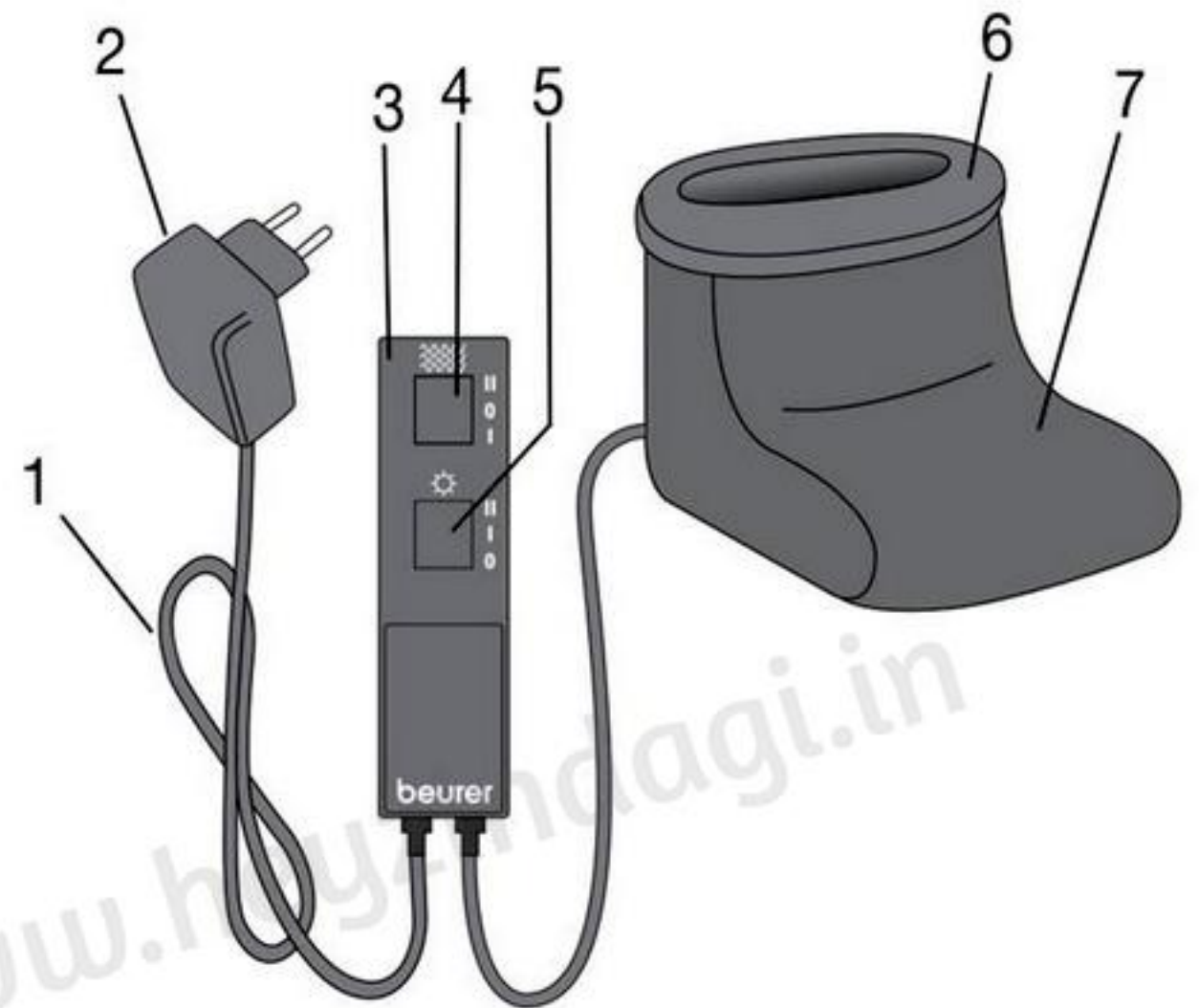

1. Cable
2. Mains Part
3. Controller
4. Sliding switch for massage settings (0-1-2)
5. Sliding switch for heat settings (0-1-2)
6. Fleece lining
7. Foot warmer

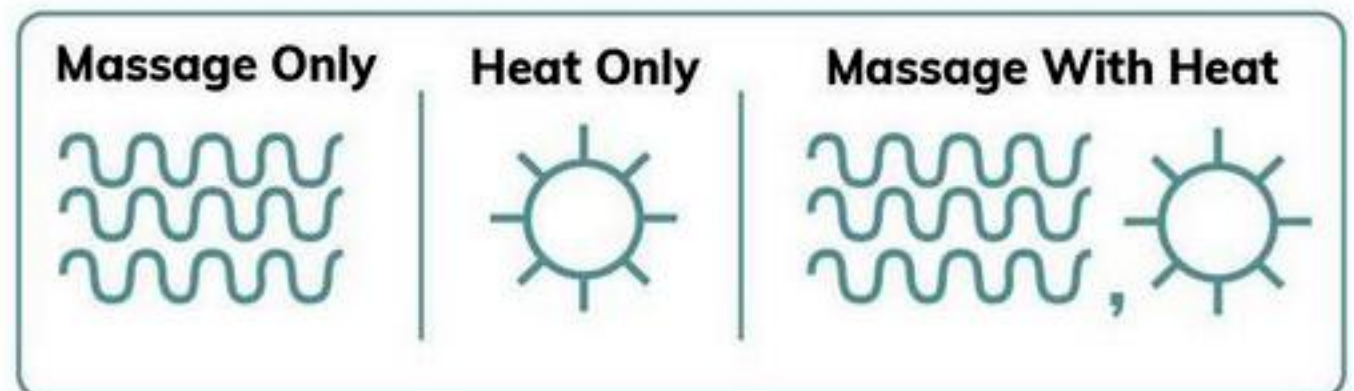

## REMEMBER

- This device must not be used by children or maintained by them without supervision
- It is not designed for use in hospital and institutional environments
- Use and store in dry areas only

- Please do not use it for animals or uses beyond what it's been made for
- If the unit does not switch on, check the connectivity thoroughly. Send it to the Beurer Service Centre ([hyperlink](#)) in India for repairs

**P.S.:** Kindly consult your Physician or Medical Specialist before using this product in case of skin insensitivity or heat intolerance.

To be used under supervision by children and those with limited response to emergencies