







KNOW THE MULTI FUNCTION THERMOMETER

1. Forehead/cover cap
2. Forehead thermometer mode button 
3. Ear thermometer mode/on button  Ear 
4. LED (red) (high temperature alarm)
5. LED (green)
6. LCD
7. Sensor tip with lens (cover cap removed)
8. Memory button 
9. Set button 
10. Battery compartment lid
11. Date
12. Time
13. Measurement mode ("Ear", "Forehead", "Object")
14. Battery warning symbol
15. Measurement display
16. "Memory" symbol and storage space number
17. "Measurement in progress" symbol
18. Unit of measurement ($^{\circ}\text{C}$ or $^{\circ}\text{F}$)

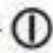

How to Use?

- Press the  button for 1 second to switch the thermometer on
- Make sure the device emits two short beeps signalling a successful self test

Measuring Body Temperature in the ear


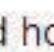




- Remove the cap by gently pushing it up (1) and then pulling it off forwards (2)
- Pull the ear of the patient slightly up and backwards before inserting the sensor tip
- Insert the sensor tip carefully and press the button for 1 second
- Release the button
- Note the measured value of the temperature on the display following a short beep

Measuring Body Temperature in the forehead

- Place the measuring head with the forehead/cover cap fitted on the temples, hold down the Forehead button and move the thermometer smoothly over the forehead to the other temple

- Release the Forehead button
- Note the measured value of the temperature on the display following a short beep

Measuring surface temperature

- Press the  and  buttons simultaneously for 3 seconds to switch to the Object temperature mode. Ensure the  symbol appears in the display
- Press and hold the  button or  Forehead and point the sensor tip towards the object or liquid to be measured maintaining a distance of 3 cm (NOTE: never immerse in liquids)
- Release the  button
- Note the measured value of the temperature on the display following a short beep

For details on setting the device and using other features, please go through the user manual for Multi-Function Thermometer by Beurer Germany, available for download below

REMEMBER

- Do not use the thermometer in taking measurement in the ear in cases of certain acute infectious diseases, inflammatory diseases (e.g. discharging pus or secretion), ear injuries (e.g. eardrum damage) or during post-surgical recovery in ear. Do not use the ear thermometer for infants below 6 months

- Clean the ear properly before measuring as ear wax can affect measurements. Always measure a person's temperature in the same ear as some people produce different readings in their left and right ear
- The thermometer needs to have been in the room in which the measurement is taken for at least 30 minutes before use