Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd.

Data is subject to change without notification. Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

Home Skinovations Ltd. Tavor Building, Shaar Yokneam POB 533, Yokneam 20692, ISRAEL [www.silkn.com](http://www.silkn.com), [contact@silkn.com](mailto:contact@silkn.com)

## Table of Contents

1. What is Flash&Go™	2
1.1. Flash&Go™ device description	2
1.2. Package contents	2
2. Silk'n Flash&Go™ Intended Use	3
3. Safety with Flash&Go™	3
4. Contraindications	3
5. When to avoid using Flash&Go™?	6
6. Precautions – How to use Flash&Go™ Safely	7
7. Reducing the risk of injury	8
8. Risks of using Flash&Go™	9
9. Long Term Hair Removal the Flash&Go™ Way	11
10. How To Use Flash&Go™	12
10.1. Flash&Go™ Setup	12
10.2. Treating with Flash&Go™ for the first time	12
10.3. What to expect when treating with Flash&Go™?	14
10.4. After treating with Flash&Go™	15
11. Flash&Go™ Maintenance	15
11.1. Cleaning Flash&Go™	15
12. Troubleshooting	16
12.1. "My Flash&Go™ does not start."	16
12.2. "A light pulse is not emitted when I press the Pulse trigger."	16
13. Customer Service	17
14. Specifications	17
15. Labeling	54
16. IEC 60601-1-2 Declaration	55

Before using Flash&Go™ for the first time, please read this User Manual in its entirety. Pay particular attention to procedures for using and operating the device and after-use procedures.

We recommend you re-familiarize yourself with this User Manual before each use of Flash&Go™.

Flash&Go™ is a powerful electrical device. As such, it should be used with special attention to safety. Please read all warnings and safety precautions before use, and strictly follow them when using Flash&Go™


## 1. What is Flash&Go™?

**Flash&Go™ is a light-based device for long-term hair removal, designed for home-use.**

### 1.1 Silk'n Flash&Go™ device description

The Silk'n Flash&Go™ HAIR REMOVAL DEVICE has a TREATMENT SURFACE, a CONTROL PAD and a PULSE TRIGGER.

- |                     |   |
|---------------------|---|
| ① Control pad       | ③ Central button                          |
| ② Pulse Trigger     | ④ Energy level indicator lights           |
| ⑤ Skin Color Sensor | ⑦ Ready/Skin tone warning indicator light |
| ⑥ Treatment Surface |   |

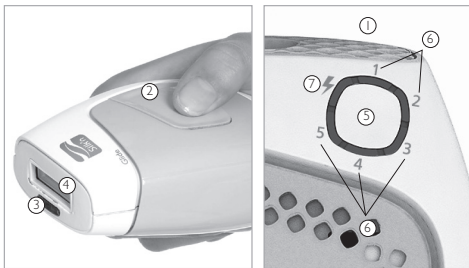On the Hair Removal device, you can find the Flash&Go™ Control Pad including a central Button which is also the Energy Level Setting Button, and POWER, Energy Level, Ready, and System Warning Indicator Lights.

### 1.2 Package contents

Upon opening the Flash&Go™ package, you will find the following parts:

- Flash&Go™ Hair Removal device with TREATMENT SURFACE
- AC cord Adaptor
- This User Manual
- A beauty Case

## 2. Silk'n Flash&Go™ Intended Use

Silk'n Flash&Go™ is intended for removal of unwanted hair. Silk'n Flash&Go™ is also intended to effect long term, or permanent hair reduction. The Flash&Go™ device is also intended for permanent reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs regrowing when measured at 6, 9, and 12 months after the completion of a treatment regime.

## 3. Safety with Flash&Go™

- **HPL™ technology in Flash&Go™ - Superior safety with lower energy level**  
Home Pulsed Light™ technology is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy used in Flash&Go™ reduces its potential to cause harm or complications, and contributes to your overall safety.
- **Flash&Go™ protects your eyes.** The Flash&Go™ Hair Removal Device is intended for use below the cheekbones but not around the eyes. HPL light can be harmful to the eyes even when not pointed directly at them. To protect your eyes from such harm the Flash&Go™ Hair Removal Device has a built-in safety feature for eye protection. It has been designed so that a light pulse cannot be emitted when the treatment surface is facing open air. The Pulse Trigger is enabled only when the treatment surface is in full contact with the skin.
- **Flash&Go™ protects your skin.** Flash&Go™ comes with a Skin Color Sensor enabling use only on lighter skin complexions to ensure skin safety. Furthermore, the 0.9X3cm² Optical Lens through which pulses of light are delivered is recessed inside the device's treatment surface. This enables Flash&Go™ to protect your skin by avoiding direct contact between the optical lens and the skin.

## 4. Contraindications

### Important Safety Information – Read Before Use!

Flash&Go™ is not designed for everyone. Please read and consider the information in the following section before use. For further information and personalized advice, you may also visit [www.silkn.com](http://www.silkn.com) or your local Flash&Go™ domain. **If you are unsure if Flash&Go™ is safe for you to use, please ask your doctor or dermatologist!**

**DO NOT USE Flash&Go™ on a naturally dark skin complexion!**

Flash&Go™ removes unwanted hair by selectively addressing hair pigment. Varied quantities of pigment also exist in the surrounding tissue of skin. The quantity of pigment in a particular person's skin, which is manifested by their skin complexion, determines the degree of risk they are exposed to using Flash&Go™. Treating dark skin with Flash&Go™ can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation)

**DO NOT USE Flash&Go™ on naturally dark skin!**

Hair	Skin						
	Lightest	Light	Medium-Light	Medium	Dark	Darkest	Black
Single hair	energy 4-5	energy 4-5	energy 3-4-5	energy 3-4	×	×	×
Multiple hairs	energy 4-5	energy 4-5	energy 3-4-5	energy 3-4	energy 1-2-3	energy 1-2	×

A unique Skin Color Sensor is embedded in Flash&Go™ to measure the treated skin complexion at the beginning of each session and occasionally during the session. This Skin Color Sensor ensures that a pulse will be emitted only on suitable skin tones.

**DO NOT USE Flash&Go™ on tanned skin or after recent sun exposure!**

Tanned skin, particularly following sun exposure, contains large quantities of the pigment Melanin. This applies to all skin types and complexions, including those which don't seem to tan quickly. The presence of large quantities of Melanin exposes the skin to higher risk when using Flash&Go™. Using Flash&Go™ on skin that has been recently exposed to sunlight can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation).

**DO NOT USE Flash&Go™ on tanned skin or after recent sun exposure! Such use can cause burns or skin injury. Avoid exposure to the sun for 4 weeks before your Flash&Go™ treatment!**

**Never use Silk'n Flash&Go™ around or near the eyes – Flash&Go™ can be used only from the cheekbones down.**

**Do not use the Flash&Go™ Device anywhere on the ears, nipples, or genitals, or around the anus. These areas may have more sensitive skin, a darker skin tone, and/or a greater hair density, and using the Flash&Go™ device there may cause serious skin injury.**

**It is not recommended to use the Flash&Go™ device on man's face, jaw, or neck. The hair in a man's beard area is too dense for this purpose.**

The Flash&Go™ device may not be as effective on naturally white, grey, blond and red body hair. Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage treated on the face and neck.

**Although, Silk'n Flash&Go™ can be used for the treatment of facial hair, extra caution should be used on the face to avoid the eyes.**

## 5. When to avoid using Flash&Go™?

Certain conditions may prevent the use of Flash&Go™ temporarily. **DO NOT USE** Flash&Go™ if any of the following currently apply to you:

- If you are pregnant or nursing (lactating).
- If you were exposed to strong sunlight or an artificial tanning machine during the past 4 weeks.
- If you have a tattoo or permanent makeup on the area to be treated.
- Avoid treating areas with dark brown or black spots, such as birth marks, moles or warts.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such as cold sores in the area to be treated. Wait for the affected area to heal before using Flash&Go™.
- If you have a history of keloidal scar formation, a known sensitivity to light (photosensitivity) or are taking medication that makes the skin more sensitive to light, including non-steroidal anti-inflammatory agents, (e.g., aspirins, ibuprofens, acetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonylraes, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A®, Accutane® and/or topical retinoids.
- If you have abnormal skin conditions caused by diabetes, for example, or other systemic or metabolic diseases.
- If you are currently or have recently been treated with Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A®, topical retinoids or azelaic acid.
- If you have been treated with Accutane® (isotretinoin) within the past 6 months.
- If you have been on a steroid regimen within the past 3 months.
- If you have a history of herpes outbreaks in the area of treatment, unless you have consulted your physician and received preventative treatment before using Flash&Go™.
- If you suffer from epilepsy.
- If you have an active implant, such as a pacemaker, incontinence device, insulin pump, etc.
- If you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- If you have a history of skin cancer or areas of potential skin malignancies.
- If you have received radiation therapy or chemotherapy treatments within the past 3 months.
- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated.

**If you are unsure if Flash&Go™ is safe for you to use, please consult with your doctor or dermatologist.**

## 6. Precautions – How to use Flash&Go™ Safely

### ■ Choose your energy levels CAREFULLY!

Energy level refers to the intensity of the light pulse that is projected on your skin during use, from the lowest level to the highest level. Indicator lights on the control pad illustrate the energy level at which the machine is set. As the energy level increases, so do the results of Flash&Go™ as well as the risk of side effects (see section 8 - "Risks of Using Flash&Go™").

### ■ Always begin your first use of Flash&Go™ at the lowest energy setting!

Only if you experience little or no discomfort during and after use of Flash&Go™ at the lowest energy setting, raise the energy level by one indicator light the next time you use Flash&Go™, and so on for each subsequent hair removal session.

For detailed instructions on energy level setting, see "Energy Level" box in section 10.2 - "Treating with Flash&Go™ for the first time".

### ■ Avoid adverse effects!

Do not treat the same area of skin more than once per hair removal session!

Try to avoid overlapping pulses!

If your skin blisters or burns, STOP USE IMMEDIATELY!

### ■ Avoid complications after use of Flash&Go™!

Do not expose treated areas of skin to the sun!

Sun exposure includes constant unprotected exposure to direct sunlight for over 15 minutes, or constant unprotected exposure to diffuse sunlight for over 1 hour.

To protect recently treated skin when exposed to sunlight, be sure to thoroughly apply sunscreen SPF 30 or higher, for 2 weeks after each hair removal session.

### ■ Always shave the area to be treated and make sure that the skin is clean and dry before using Flash&Go™.

- Cover birthmarks and tattoos before Flash&Go™ application.
- Never look directly at the light coming from the Flash&Go™ treatment surface.
- Do not use Flash&Go™ on nipples and genitals (male or female).
- Do not use Flash&Go™ on any body site where you might later want hair.
- Do not use Flash&Go™ for any purpose other than hair removal.
- Never point the Flash&Go™ treatment surface in an attempt to emit a light pulse into open space. Always make sure that the surface is pointed at, and in full contact with the skin during application.
- Remove the Flash&Go™ treatment surface from the skin if either the skin or the surface is too hot.

- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone to clean the skin before using Flash&Go™.
- Use of Flash&Go™ may cause temporary pigmentation changes (See section 8 - "Risks of Using Flash&Go™").
- Keep this device out of the reach of children. Do not use Flash&Go™ on children or allow children to use it.

## 7. Reducing the risk of injury

As with any electrical device, certain precautions must be taken in order to ensure your safety when using Flash&Go™.

**No modification of this equipment is allowed!**

**Keep Flash&Go™ away from water!**

Flash&Go™ is an electrical device. As such it should always be kept away from water.

Do not place or store Flash&Go™ where it can fall or be pushed into a tub, sink or any other vessel containing water. Do not place in, or drop into water or any other liquid.

- This may cause severe electrocution.
- Do not use Flash&Go™ while bathing.
- Do not use Flash&Go™ if it becomes damp or wet.
- Do not reach for Flash&Go™ if it has fallen into water.
- Unplug Flash&Go™ immediately if it has fallen into water.
- Keep Flash&Go™ away from water! Unplug Flash&Go™ immediately if it has fallen into water.

### **Never open Flash&Go™!**

Do not attempt to open or repair your Flash&Go™ device. Opening Flash&Go™ may expose you to dangerous electrical components and to pulsed light energy, either of which may cause bodily damage and/or eye injury.

**Do not attempt to open or repair your Flash&Go™ device. Only authorized Flash&Go™ repair centers are permitted to perform repairs.**

- Trying to open Flash&Go™ may also damage the device and will void your warranty.
- Please contact Flash&Go™ Customer Service if you have a broken or damaged device in need of repair.
- Use Flash&Go™ only for its intended use and as described in its manual.
- Flash&Go™ should never be left unattended when plugged into an outlet.
- Do not operate Flash&Go™ if it has a damaged cord or plug and keep the power cord away from heated surfaces.
- Do not use Flash&Go™ if you see or smell smoke when it is in use.
- Do not use Flash&Go™ if it is not working properly or if it appears damaged.
- Do not use Flash&Go™ if the fan vent in its treatment surface is cracked, coming off or missing.
- Do not use Flash&Go™ if the Skin Color Sensor is cracked, or broken.
- Do not use Flash&Go™ if the outer shell is cracked or coming apart.
- Always unplug Flash&Go™ from the electrical outlet immediately after use.
- Unplug Flash&Go™ before cleaning.
- Do not use Flash&Go™ with any attachments or accessories not recommended by Home Skinovations Ltd.
- Flash&Go™ must be operated only with its designated adaptor (for details on the right adaptor to use please see chapter 14 - "Specifications")

## 8. Risks of using flash&Go™

Side effects and complications associated with the use of Flash&Go™ according to its instructions are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore it is important that you understand and accept the risks and complications that can occur with pulsed light hair removal systems designed for home use.

Adverse Event	Likelihood of the Adverse Event	Adverse Effect	Likelihood of the Adverse Effect
Stacking or Overlapping of multiple pulses on the same skin spot	Minor	Minor Skin Discomfort	Minor
		Skin Redness	Minor
		Increased Sensitivity of the Skin	Minor
		Skin Wounds and Burns	Rare
		Scarring	Rare
		Pigment Changes	Rare
		Excessive Redness and Swelling	Rare
		Infection	Negligible
Bruising	Negligible		

#### ■ Minor Skin Discomfort

Although home pulsed light hair removal is generally very well-tolerated, most users do feel some mild discomfort during use, usually described as being a mild stinging sensation on the treated skin areas. The stinging sensation usually lasts during the application itself or for a few minutes thereafter. Anything beyond this minor discomfort is abnormal and means that either you should not continue to use Flash&Go™ because you are unable to tolerate the hair removal application, or that the energy level setting is too high.

#### ■ Skin Redness

Your skin may become red right after using Flash&Go™ or within 24 hours of using Flash&Go™. Redness generally clears up within 24 hours. See your doctor if redness does not go away within 2 to 3 days.

#### ■ Increased Sensitivity of the Skin

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

#### ■ Skin Wounds and Burns

Very rarely, burns or wounds to the skin can occur following the application. The burn or wound can require a few weeks to heal and, extremely rarely, may leave a noticeable scar.

#### ■ Scarring

Although very rare, scarring may occur. Usually when scarring occurs, it is in the form of a flat and white lesion on the skin (hypotrophic). However, it can be large and red (hypertrophic) or large and extend beyond the margins of the injury itself (keloid). Subsequent aesthetic treatments may be required to improve the appearance of the scar.

#### ■ Pigment Changes

Flash&Go™ targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless, there is risk of temporary hyperpigmentation (increased pigment or brown discoloration) or hypopigmentation (whitening) to the surrounding skin. This risk of changes in skin pigmentation is higher for people with darker skin tones.

Usually discoloration or changes to skin pigment are temporary and hyperpigmentation or hypopigmentation rarely occur.

#### ■ Excessive Redness and Swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but one should avoid exposure to sun.

#### ■ Infection

Infection of the skin is exceedingly rare but is still a possible risk following a skin burn or wound caused by Flash&Go™.

#### ■ Bruising

Very rarely, use of Flash&Go™ may cause blue-purple bruising that can last 5 to 10 days. As the bruising fades, there may be a rust-brown discoloration of the skin (hyperpigmentation).

## 9. Long Term Hair Removal the Flash&Go™ Way

A typical full hair growth cycle may take 18-24 months. During this time multiple Flash&Go™ sessions may be required in order to achieve permanent hair removal.

The efficiency of hair removal varies from person to person according to body area, hair color, and how Silk'n Flash&Go™ is used.

Typical Silk'n Flash&Go™ hair removal plan during a full hair growth cycle:

Treatments 1-4 – plan two weeks apart

Treatments 5-7 – plan four weeks apart

Treatments 8 + – treat as needed, until desired results are achieved.

### Expected Treatment Outcome for Hair Removal

Fitzpatrick skin type	Hair Color	Light Fluence [J/cm <sup>2</sup> ]	Anatomic area of the body <sup>2</sup>	Average Number of Treatments <sup>1</sup>	Effectiveness (% of hair removal/regrowth) <sup>1</sup>
I-IV	Light brown to brown	3-4.5	Lower or upper legs	Session 10-12, approx. four weeks apart.	60% Hair reduction
			Arms		50% Hair reduction
			Bikini Line		60% Hair reduction
			Underarms		55% Hair reduction
I-IV	Brown to black	3-4.5	Lower or upper legs	Session 8-10, approx. four weeks apart.	70% Hair reduction
			Arms		60% Hair reduction
			Bikini Line		70% Hair reduction
			Underarms		65% Hair reduction

- Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.
- One cannot expect permanent epilation in a single or even double treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

## 10. How To Use Flash&Go™

### 10.1. Flash&Go™ Setup

- Remove the Flash&Go™ hair removal device, and other components from box.
  - Plug the power cord into the Flash&Go™ Device socket.
  - Plug the other end of the power cord into an electrical outlet.
  - All six LED lights will blink on and off, one after the other, in a circle.
- Your Flash&Go™ is now ready to start.

### 10.2. Treating with Flash&Go™ for the first time

- The skin should be **shaved, clean, dry** and **free** of any powders, antiperspirants or deodorants.
- Press the central button briefly. The Energy level 1 led will turn on.  
About 1 second later, the Ready Indicator Light will turn on. The device is now set to deliver the lowest energy level pulses. The device is ready for you to trigger the first pulse.
- If you wish to change the energy level, press the CENTRAL BUTTON quickly and repeatedly to increase the energy level from 1 to 5 and from 5 directly back to 1, cyclically, until the desired energy level is set. The Energy Level indicator Lights will indicate the energy level setting. (See "Energy Level" frame).

#### Energy Level

Energy level determines the intensity of the Flash&Go™ light pulse delivered to your skin, from the lowest level (1) to the highest level (5). Corresponding indicator lights on the control panel represent the increases in energy. As energy level increases, so do hair removal results as well as the risk of possible side effects and complications.

Always start your first Flash&Go™ hair removal session at the lowest energy setting. If you experience little or no discomfort during and after the hair removal session using the lowest energy level, raise the energy level by one indicator light at the next session, and so on for each subsequent session.

**Do not raise the energy level if you experience abnormal discomfort during or after treating with Flash&Go™ (See section 8 - "Risks of Using Flash&Go™").**  
Do not raise the energy level during a hair removal session even if you experience no discomfort. Discomfort may also appear sometime after the session.

Whenever Flash&Go™ is turned on; its energy level will automatically be reset to the lowest energy level. Only one Energy Level indicator Light will be on.

To set the energy level, press the central Button repeatedly to increase the energy level. After reaching level 5, an additional press will bring the level back down to Level 1, and so on, cyclically. The number of Energy Level indicator Lights will coincide with the change in energy level.

- Apply the TREATMENT SURFACE to the skin, making sure the skin is spread evenly and smoothly and the TREATMENT SURFACE is in full contact with the skin. Make sure the Ready Indicator Light starts blinking slowly.
- Press the Pulse Trigger. The device will first determine the color of your skin.  
If the color of the skin is light enough for safe application, the device will flash a pulse of light onto your skin. You will see a bright flash of light and simultaneously hear a subtle pop sound, which is a normal noise for the device. You may feel a mild sensation of warmth and tingling.  
Flash&Go™ will immediately recharge for the next pulse. After a few seconds (depending on the set energy level) the Ready Indicator Light will turn on again and start blinking.  
If you see no light pulse, and the Ready Indicator Light blinks rapidly, this is now the SKIN COLOR WARNING INDICATOR LIGHT alerting that your skin tone is too dark for safe application. Try using the device on a different body part or contact Silk'n support.
- Move the TREATMENT SURFACE to another area of skin.  
Administer one pulse at a time at your chosen energy level. Use the pressure marks the device just made in your skin to guide you for proper positioning of the next pulse, avoiding both gaps and overlaps between pulses. Flash&Go™ can be applied in two ways: "Pulsing" and "Gliding".

#### PULSING

When Pulsing, administer one pulse at a time at your chosen energy level. Use the pressure marks the device just made in your skin to guide you for proper positioning of the next pulse, avoiding both gaps and overlaps between pulses.

#### Covering skin areas

Treatments should be administered in rows, starting at one end of each row and progressing sequentially towards the other end. This technique allows better control of skin coverage, and helps you avoid treating the same area more than once or overlapping skin areas. Ideal body areas for Flash&Go™ use include the underarms, bikini line, arms and legs and below cheekbones line.

- Repeat the process starting again with Step 5.

## GLIDING

Alternatively, if you choose the Gliding method, hold the PULSE TRIGGER down continuously, and slide the device smoothly along the treated area. Please note that this method is best applied when working with low energy pulses. This way, the device takes only 1-1.5 seconds to recharge between pulses and provides a nonstop sequence of flashes.

In both methods of treatment you are cautioned not to treat the same area of skin more than once per hair removal session! Treating the same area of skin more than once per session increases the likelihood of adverse effects.

Gliding should be applied while holding the PULSE TRIGGER down continuously, and using a light sliding motion, smoothly moving the device along the treated area. The treatment is then administered quickly and more efficiently. Even though Gliding is usually done using low energy pulses for fast recharge, it is still essential to control skin coverage and avoid double treatments.

When applied to the skin, the Flash&Go™ TREATMENT SURFACE is designed to create temporary pressure marks on the treated area. These visible marks can be used for proper positioning of the device for treatment sessions.

If your skin blisters or burns, STOP USE IMMEDIATELY!

### 10.3. What to Expect when treating with Flash&Go™?

For many people, using Flash&Go™ may be their first experience with a light-based device designed for home use. Flash&Go™ is simple to use, and hair removal sessions go by quickly. During a Flash&Go™ session it is normal to experience and feel:

- **A Flash of Light** – The bright light of Flash&Go™ will not harm the eyes when applied to non-facial sites, and special eye-protection is not needed when using Flash&Go™.
- **A Fan Noise** – The cooling fan in Flash&Go™ sounds similar to a hairdryer. This is normal.
- **A Pop Sound with Each Pulse** – When a pulse of HPL™ light is activated, it is normal to hear a subtle pop sound simultaneously with the flash of light.
- **Moderate Pressure of the Treatment area** – This is necessary and helpful for placement of adjacent pulses of light, and is part of the unique safety features of Flash&Go™.
- **A Sensation of Warmth and Tingling** – During each pulse of light, it is normal to feel a mild sensation of warmth and tingling from the light energy. Remember it is important to always use low energy settings for initial hair removal sessions. You may feel some warmth for up to an hour after your Flash&Go™ session.
- **Some Mild Red or Pink Color** – During and just after your Flash&Go™ session it is not uncommon to see some very mild, pink-like color of the skin. This is usually most noticeable around the hairs themselves. However if you see full redness of the skin, blistering or burns, stop use of Flash&Go™ immediately.

### 10.4. After treating with Flash&Go™

- When a Flash&Go™ session has been completed, turn Flash&Go™ with an extended press of the CENTRAL BUTTON. Be sure to remember the last energy level setting you used, as it will not be restored when turning Flash&Go™ on again.
- Unplug the power cord from the electrical outlet.
- After each hair removal session, it is recommended that you clean your Flash&Go™ device, especially the Treatment surface (See section 11.1: "Cleaning Flash&Go™").  
After cleaning, it is recommended you store your Flash&Go™ device in its original box, and keep it away from water.

#### Skin care following hair removal session

Do not expose the treated areas of skin to the sun. Be sure to carefully protect the treated skin with sunscreen, throughout the hair removal period and for at least 2 weeks following the last Flash&Go™ session.

### Side effects and complications

Some patients may experience pigmentation changes resulting from treating with Flash&Go™. These effects, if they occur, are generally mild and transient.

If you experience any complication (See section 8 – "Risks of using Flash&Go™") please contact your physician immediately.

## 11. Flash&Go™ Maintenance

### 11.1. Cleaning Flash&Go™

After each hair removal session, it is recommended that you clean your Flash&Go™ device, and especially the treatment surface.

- Unplug Flash&Go™ before cleaning.
- Use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe the Flash&Go™ outer shell, and especially the treatment surface.

Never immerse Flash&Go™ or any of its parts in water!

This device is to be orderly disposed of when reaching End of Life stage. In this way the device's valuable materials are recycled and the environmental burden avoided.

Deliver your old equipment to the nearest electronic waste collection point for safe and efficient disposal.

## 12. Troubleshooting

### 12.1. “My Flash&Go™ does not start.”

- Make sure the power cord is properly connected to the Flash&Go™ device.
- Make sure the power cord is plugged into an electrical outlet on the wall.

### 12.2. “A light pulse is not emitted when I press the PULSE TRIGGER”

- Make sure the Ready Indicator is blinking.
- Make sure you have good contact with the skin and that the treatment surface is evenly and firmly pressed to the skin. For your safety, the Pulse trigger will activate a pulse only if the treatment surface is firmly pressed against the skin.
- Check the Indicator Light. If it blinks slowly, the Skin Color sensor has stopped the light pulse; the skin color window was not properly covered by the skin or your skin may be too dark for you to use Flash&Go™. Try to apply Flash&Go™ again on a different part of your body.
- If the Ready light is off, together with any number of blinking Energy Level light indicators, you have a system error. If this problem persists, contact your local Flash&Go™ Customer Service Center.
- If these problems persist, contact your local Flash&Go™ Customer Service Center.
- Check if the System Warning Indicator Light is on, and all other lights are blinking. If so, unplug the device from its power source. The device has reached its End of Life stage. Please refer to section 11 – “Flash&Go™ Maintenance” for further instructions on how to safely dispose of the device.

Do not attempt to open or repair your Flash&Go™ device. Only authorized Flash&Go™ repair centers are permitted to perform repairs.

Opening Flash&Go™ may expose you to dangerous electrical components and to pulsed light energy, either of which may cause bodily damage and/or eye injury.

Trying to open Flash&Go™ may also damage the device and will void your warranty.

Please contact Flash&Go™ Customer Service if you have a broken or damaged device in need of repair.

## 13. Customer Service

For more information about Silk'n™ products please enter your local Silk'n™ website, [www.silkn.com](http://www.silkn.com). If your device is broken, damaged, in need of repair, or for any other user assistance, please contact your local Silk'n™ Customer Service:

For USA and Canada customers: 1-877-DO-SILKN / 1-877-367-4556, [contact@silkn.com](mailto:contact@silkn.com)

## 14. Specifications

Spot size	0.9cm x 3cm [2.7cm <sup>2</sup> ]	
Technology	New Home Pulsed Light™	
Pulse Duration	500 – 800 µsec	
Wavelength	475-1200nm	
Power Output	3-5 J/cm <sup>2</sup>	
Charging time / Power source	Continuous operation	
Time needed to treat lower legs	30 minutes	
Operation and Safety	Skin Color Sensor seamlessly ensures use only on appropriate skin types.	
Classified standard	IEC 60601-2-57:2011	
Risk Group	Exempt	
Package size	Height 9 inch, Width 9 inch, Depth 5.2 inch   80x180x100mm	
System weight	200 gr.	
Temperature	Operating	10°C to 35°C
	Storage	-40 to +70°C
Humidity	Operating	30 to 75% RH
	Storage	10 to 90% RH
Atmospheric pressure	Operating	700 to 1060hPa
	Storage	500 to 1060hPa
Adaptor Rated	Model: KSA24A   200150HU 100-240V ; 50-60Hz 0.5A	KSAPV036   200300DS 100-240V ; 50-60Hz 0.9A

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
## 15. Labeling

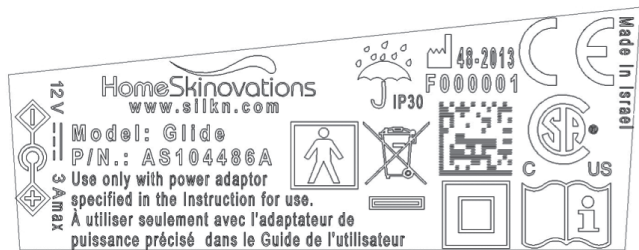Class II equipment

WEEE - Waste Electrical and Electronic Equipment

CE Mark

Degree of protection against ingress water: ordinary

Follow operating instructions

Degree of protection against electric shock: type BF applied part

CSA Mark for USA and Canada

This device is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

## 15. IEC 60601-1-2 Declaration

### 1. Summary of Test Results

Test	Standard	Class/ Severity level	Test result
<b>Documentation</b> (IEC 60601-1-2 sections 4 and 5)			
General requirements General requirements for EMC	section 4.1.1.	--	Complies
External labels	section 5.1	--	Complies
Conformity of Users' Manual	section 5.2.1	--	Complies
Accuracy of Technical Description	section 5.2.2	--	Complies
<b>Emission</b> (IEC 60601-1-2 section 6.1)			
Conducted emission Freq. range: 150 kHz - 30 MHz	sec. 6.1.1 & CISPR 11	Group 1 Class B 230 VAC mains; 120 VAC mains	Complies
Radiated emission Freq. range: 30 - 1000 MHz	sec. 6.1.1 & CISPR 11	Group 1 Class B	Complies
Harmonic current emission test	sec. 6.1.3.1 & IEC 61000-3-2	230AC mains	Test NA
Voltage changes, Voltage fluctuations and Flicker test	sec. 6.1.3.2 & IEC 61000-3-3	230AC mains	Complies
<b>Immunity</b> (IEC 60601-1-2 section 6.2)			
Immunity from Electrostatic discharge (ESD)	sec. 6.2.2 & IEC 61000-4-2	6 kV contact discharges & 8 kV air discharges	Complies
Immunity from radiated electromagnetic fields	sec. 6.2.3 & IEC 61000-4-3	3.0 V/m 80 MHz $\pm$ 2.5 GHz, 80% AM, 1kHz	Complies
Immunity from Electrical Fast transient (EFT)	sec. 6.2.4 & IEC 61000-4-4	$\pm$ 2.0 kV Tr/Th - 5/50 ns, 5 kHz on 230AC mains	Complies
Immunity from Surge	sec. 6.2.5 & IEC 61000-4-5	$\pm$ 1.0 kV DM/ 2.0kV CM Tr/Th - 1.2/50 (8/20) $\mu$ s on 230 VAC mains	Complies
Immunity from conducted disturbances induced by radio-frequency fields	sec. 6.2.6 & IEC 61000-4-6	3.0 V <sub>RMS</sub> 0.15+ 80 MHz, 80% AM 1kHz	Complies
Immunity from Voltage interruptions	sec. 6.2.7 & IEC 61000-4-11	100% - 10 ms; 60% - 100 ms; 30% - 500 ms; 100% - 5sec on 230/ 120VAC mains	Complies

Approved: Eng. Yuri Rozenberg  
Position: Head of EMC Branch

Written by: Galit Gorodetsky  
Position: Technical Writer

Electronics and Telematics  
Laboratory  
April 7, 2013

Name: Alexander Karlinas  
Position: Testing Engineer

Name: Michael Feldman  
Position: Testing Engineer