Copyright 2013 © Home Skinovations Ltd. All rights reserved.  
Print date: August 2015

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability.

Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd.

Data is subject to change without notification. Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

## Table of Contents

<b>1. Before You Start</b> .....	<b>2</b>
1.1. Silk'n Silhouette™ device description .....	2
1.2. Package contents .....	3
<b>2. Silk'n Silhouette™ intended use</b> .....	<b>3</b>
<b>3. Safety with Silk'n Silhouette™</b> .....	<b>3</b>
<b>4. Contrarindications</b> .....	<b>4</b>
<b>5. When not to use Silk'n Silhouette™</b> .....	<b>5</b>
<b>6. Get started!</b> .....	<b>5</b>
<b>7. Treatment regimen</b> .....	<b>7</b>
<b>8. Maintenance of Silk'n Silhouette™</b> .....	<b>10</b>
8.1. Cleaning Silk'n Silhouette™ .....	10
<b>9. Troubleshooting</b> .....	<b>10</b>
<b>10. Customer Service</b> .....	<b>11</b>
<b>11. Specifications</b> .....	<b>11</b>
<b>12. Labeling</b> .....	<b>133</b>

# 1. Before You Start

Before using Silk'n Silhouette™ for the first time, please read this Instructions for Use booklet in its entirety, and watch the Silk'n Silhouette™ training video.



Pay particular attention to sections on device use procedures, device operation, and after-use procedures. We recommend you re-familiarize yourself with this Instructions for Use booklet before each use of Silk'n Silhouette™.

Silk'n Silhouette™ is a powerful electrical device. As such, it should be used with special attention to safety. Please read all warnings and safety precautions before use, and strictly follow them when using Silk'n Silhouette™.

## 1.1 Silk'n Silhouette™ device description

The Silk'n Silhouette™ has a TREATMENT SURFACE with four ELECTRODES and a THERMAL SENSOR, a TREATMENT MODE SELECTION BUTTON, and an ADJUSTABLE STRAP. (To enlarge pop out inner rubber section located in center of strap).

- ① Adjustable strap
- ② Treatment mode selection button
- ③ Treatment mode indicator lights
- ④ Treatment surface with four electrodes
- ⑤ Thermal sensor
- ⑥ Treatment led lights
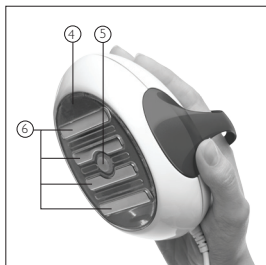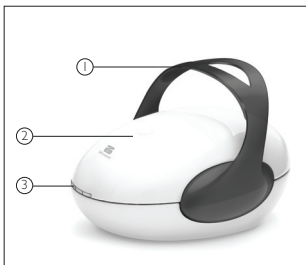
## 1.2. Package contents

Upon opening the Silk'n Silhouette™ package, you will find the following parts:

- Silk'n Silhouette™ Body Contouring Device
- A bottle of Slider Gel
- An AD/DC wall mount adaptor
- This Instruction for Use booklet

## 2. Silk'n Silhouette™ Intended Use

Silk'n Silhouette™ is intended for cellulite and circumferential reduction.

## 3. Safety with Silk'n Silhouette™

**Silk'n Silhouette™ protects your skin:**

Silk'n Silhouette™ comes with a THERMAL SENSOR designed to measure skin temperature and stabilize the emitted heat energy.

To prevent overheating, a unique built-in heat stabilizer inside the Silk'n Silhouette™ measures the treated skin temperature during session. As a rise in the treated skin temperature is detected (42°C-43°C), the stabilizer will automatically stop emitting heat, until skin temperature lowers enough for the stabilizer to resume heating.

**Note: Make sure the device is tightly pressed to your skin before starting treatment!**

## 4. Contraindications

### Important Safety Information – Read Before Use!

Silk'n Silhouette™ is not designed for everyone. Please read and consider the information in the following section before use. For further information and personalized advice you may also visit [www.silkn.com](http://www.silkn.com) or your local Silk'n™ domain.

- Do not use Silk'n Silhouette™ if you have a pacemaker or internal defibrillator, or any other active electrical implant anywhere in the body.
- Do not use Silk'n Silhouette™ if you have a permanent implant in the treated area such as metal plates and screws, silicone implants or an injected chemical substance.
- Do not use Silk'n Silhouette™ if you have a current history or condition of cancer, or pre-malignant moles.
- Do not use Silk'n Silhouette™ if you have received radiation therapy or chemotherapy treatments within the past 3 months.
- Do not use Silk'n Silhouette™ if you are pregnant or nursing.
- Do not use Silk'n Silhouette™ if you have a history of diseases stimulated by heat, such as recurrent Herpes Simplex in the treatment area. If so you may treat only following a prophylactic regime.
- Do not use Silk'n Silhouette™ if you suffer from poorly controlled endocrine disorders, such as diabetes and thyroid conditions.
- Do not use Silk'n Silhouette™ if you suffer from epilepsy
- Do not use Silk'n Silhouette™ if you have a disease related to photosensitivity, such as Porphyrria, Polymorphic Light Eruption, Solar Urticaria, Lupus, etc.
- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated.
- If your skin was treated with a physician-prescribed drug within the past 6 months, please consult with your physician.
- Do not use Silk'n Silhouette™ if you have any active condition in the treatment area, such as sores, psoriasis, eczema, and rash.
- Do not use Silk'n Silhouette™ if you had any surgery in the treated area within 3 months prior to treatment.

## 5. Do not use Silk'n Silhouette™:

- With attachments or accessories are not recommended by Home Skinovations Ltd.
- If there is a damaged cord or plug
- If you see or smell smoke when it is in use.
- If it is not working properly or if it appears damaged.
- If its THERMAL SENSOR is cracked, or broken.
- If the outer shell is cracked or is coming apart.

Please contact Silk'n Customer Service.

## 6. Get Started!

1. Remove Silk'n Silhouette™ DEVICE and other components from box.
2. Plug the AC/DC ADAPTOR cord into Silk'n Silhouette™ socket.
3. Plug the AC/DC ADAPTOR into an electrical outlet. All three indicator lights will turn on and off repeatedly one after the other:
4. The device is now in Standby mode and is ready to start.

### **Treat with Silk'n Silhouette™ safely and properly!**

5. Spread a thin layer of Slider Gel on the designated treatment area.  
Alternatively, spread a moderate amount of Slider Gel on the device's TREATMENT SURFACE. The Slider Gel is designed to better conduct the energy into your skin, for optimal results.
6. Press the TREATMENT MODE SELECTION BUTTON. TREATMENT MODE I INDICATOR LIGHT will turn on.
7. Attach the TREATMENT SURFACE to the skin. Please note that the skin must be clean and dry before applying slider gel. Skin should be free of other moisturizers and lotions that can harm the device.

## Treatment Modes

Silk'n Silhouette™ offers 3 types of treatment represented by 3 TREATMENT MODE INDICATOR LIGHTS.

The treatment mode determines the nature of Silk'n Silhouette™ energy delivered to your skin. Modes differ but complement each other.

Setting the Treatment Modes

- To set a TREATMENT MODE, press the TREATMENT MODE SELECTION BUTTON repeatedly to move between modes:
  - LOW: Energy Level 1
  - MEDIUM: Energy Level 2
  - HIGH: Energy Level 3
- When the desired treatment is set, the TREATMENT MODE INDICATOR LIGHTS will indicate the mode chosen.

**LOW**

Energy Level 1

**MEDIUM**

Energy Level 2

**HIGH**



Energy Level 3

8. Attach the device to the designated skin area. You will feel a mild sensation of warmth.
  9. Move the device in gentle massage- like linear motions on the designated area (HIPS, BUTTOCKS, THIGHS, CALF, BACK, THORAX, BELLY AND ARMS). Progress slowly and carefully.
  10. If you wish to change treatment modes in mid-treatment, press the BUTTON again. This will change the treatment mode on the spot and allow you to treat the same, or different areas with different modes according to your needs.
- Please note: the TREATMENT SURFACE cannot be in contact with the skin when changing mode.

## 7. Treatment Regimen

Treating with Silk'n Silhouette™ is easy! Just follow the illustrations below, and you are on your way to a whole new look!

It is important to preheat treatment area before moving to higher levels. Always begin treatments on low as indicated in illustration, for best results be sure to increase energy level to medium or high for remaining session. As energy level increase so do results.

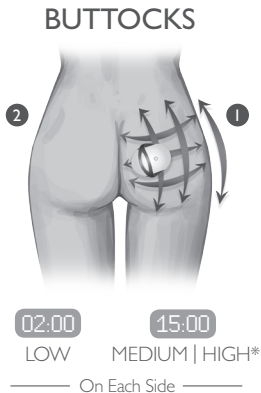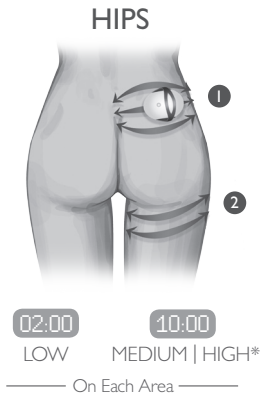\* Depending on user tolerance

## THIGH

02:00

LOW

15:00

MEDIUM | HIGH\*

## CALF

02:00


LOW

10:00

MEDIUM | HIGH\*

On Each Side

## BACK

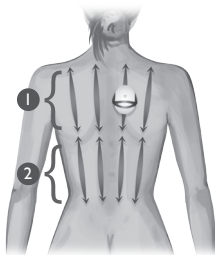02:00


LOW

12:00

MEDIUM | HIGH\*

On Each Area

## THORAX

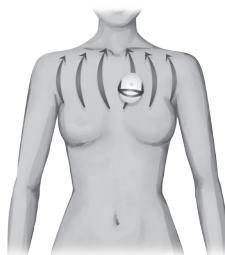02:00

LOW


10:00

MEDIUM | HIGH\*


\* Depending on user tolerance

## BELLY
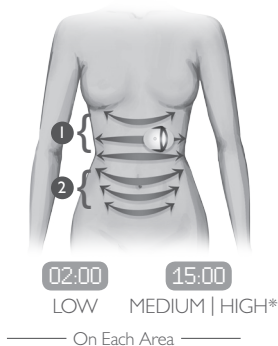
## ARM

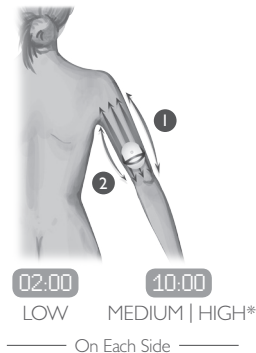\* Depending on user tolerance

### 7.1. What to expect when treating with Silk'n Silhouette™?

For many people, using Silk'n Silhouette™ may be their first experience with a heat-based device designed for home use. Silk'n Silhouette™ is simple to use, and sessions go by quickly. During a Silk'n Silhouette™ session it is normal to experience and feel:

- **A Sensation of Warmth** – During each session it is normal to feel a mild sensation of warmth from the heat energy.
- **Some Mild Red or Pink Color** – During and just after your Silk'n Silhouette™ session it is not uncommon to see some very mild, pink-like coloration of the skin.

However if you see full redness of the skin, or blistering, stop using Silk'n Silhouette™ immediately.

#### For best results keep the following schedule

Start with sessions twice a week for 10 weeks to give your body an intensive treatment. By then you will be able to see a change in your body contours. Maintain results by treating once every 2 weeks.

Please note: optimal results depend on accurately following the treatment regimen!

## 7.2. After treating with Silk'n Silhouette™

- When Silk'n Silhouette™ session has been completed turn Silk'n Silhouette™ off by unplugging the AC/DC ADAPTOR cord from the electrical outlet.
- After each session it is recommended that you clean your Silk'n Silhouette™ device, especially the TREATMENT SURFACE (See: 8.1 "Cleaning Silk'n Silhouette™").
- After cleaning, it is recommended to store your Silk'n Silhouette™ device in its original box, and keep it away from water!

# 8. Maintenance of Silk'n Silhouette™

## 8.1. Cleaning Silk'n Silhouette™

After each treatment session, it is necessary to clean your Silk'n Silhouette™ device, and especially the TREATMENT SURFACE. Unplug Silk'n Silhouette™ before cleaning

- Unplug Silk'n Silhouette™ before cleaning.
- Use absorbent paper (such as toilet paper or hand drying napkins) to gently wipe the gel residue off the TREATMENT SURFACE.
- Never immerse Silk'n Silhouette™ or any of its parts in water!

# 9. Troubleshooting

## “My Silk'n Silhouette™ does not start.”

- Make sure the POWER ADAPTOR is properly connected to the Silhouette™ device.
- Make sure the POWER ADAPTOR is plugged into an electrical outlet on the wall.

## “Heat energy is not emitted and I don't feel any heat emitted from the device even though I pressed the TREATMENT MODE SELECTION BUTTON”

- Heat will not be emitted if the TREATMENT SURFACE is not touching to your skin. This is a safety feature. Make sure the device is tightly pressed to your skin before starting treatment.

Do not attempt to open or repair your Silk'n Silhouette™ device. Only authorized Silk'n™ personnel are permitted to perform repairs.

- Trying to open Silk'n Silhouette™ may also damage the device and will void your warranty.
- Please contact Silk'n™ Customer Service if you have a broken or damaged device in need of repair.

## 10. Customer Service

For more information about Silk'n™ products please enter your local Silk'n™ website, [www.silkn.com](http://www.silkn.com).  
If your device is broken, damaged, in need of repair, or for any other user assistance, please contact your local Silk'n™ Customer Service:

For USA and Canada customers: 1-877-DO-SILKN / 1-877-367-4556, [contact@silkn.com](mailto:contact@silkn.com)

For customers in other countries: [info@silkn.com](mailto:info@silkn.com)

## 11. Specifications

Technology		HT™
Wavelength Red light		620nm±20nm
Wavelength IR		850nm±20nm
Radio Frequency		1 MHz ±30%; 24W max
Operation and Safety		Heat stabilizer preventing skin overheating
Package Size		22cm X 22cm X 14cm
System weight		200gr
Temperature	Operating	10°C to 35°C
	Storage	-40 to +70°C
Relative Humidity	Operating	30 to 75%rH
	Storage	10 to 90%rH
Atmospheric pressure	Operating	700 to 1060hPa
	Storage	500 to 1060hPa
Adaptor		Model: KSAPVO361200300D5
■ Rated		100-240V; 50-60Hz; 0.9A
■ Output		12Vdc; 3A

Copyright 2013 © Home Skinovations Ltd. All rights reserved.  
Print date: August 2015

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability.

Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd.

Data is subject to change without notification. Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

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