mē™

User guide
Mode d'emploi

User manual - original English instructions

These instructions, warnings and safety information within this user manual are intended to help ensure proper use and optimal results. It is important that you read this entire user manual before using your mē for the first time.

Table of contents

Introduction	3
Intended use of the mē	3
Getting to know your mē.....	4
Warnings	5
Using your mē for the first time	7
Skin sensitivity test.....	8
How to use your mē	8
What to expect during treatment.....	11
Treatment schedule	11
How to clean and store your mē.....	11
Disposal/recycling	12
Traveling with mē.....	12
Troubleshooting	13
Warranty.....	14
Manufacturer information.....	15

Congratulations!

mē is the only in-home compact hair removal device using patented elōs technology that is clinically proven to be safe and effective for all skin tones for permanent results on the body and face.

Introduction

Developed by the leading global aesthetic device company with the same technology used by dermatologists around the world, mē treatments disable the hair follicle so that unwanted hair does not grow back. Gentle, compact and easy to use, get ready to experience silky, smooth skin without the hassles of shaving and waxing!

Intended use of the mē

mē is an over-the-counter device intended for permanent hair reduction on the body for men and women and on women's facial hair below the cheek bone line.


Skin Tones

mē is intended for use on ALL skin tones (the skin tone you were born with).

Hair Colors

mē is clinically proven to be effective on the widest range of hair colors.

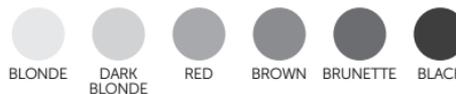Results may vary based on individual's skin tone and hair color being treated and may require additional treatments beyond the full protocol. For more information, please contact Customer Care at www.mepower.com

Getting to know your mē mē package includes:

- mē device
- Power supply
- User manual

Device components

1. Control panel

- ① Power button
- ② elōs energy setting button with 3 energy level indicators: LOW- MEDIUM- HIGH
- ③ Power supply port

2. Cooling vents

Cool air is continuously circulated to keep the device cool while in use.

3. Treatment window



The light energy flashes out through the treatment window.

4. Radio frequency (RF) contact bars

Both RF bars must be in complete contact with the skin for the light energy to flash.

5. Power supply port

Connects the device to an electrical source during use.

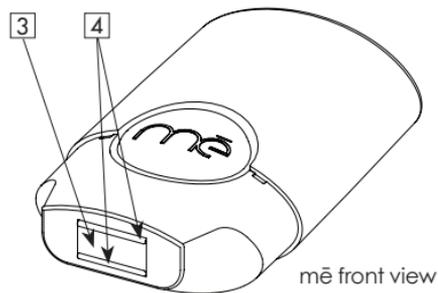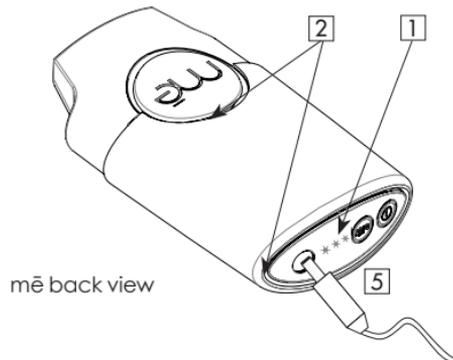Warnings

Who should not use (contraindications)

DO NOT use if you have epilepsy or light-induced seizures or chronic migraines. The energy flashes might cause severe headaches and migraines.

DO NOT use if you have an active implant, such as a pacemaker, incontinence device, insulin pump, etc. Use of the mē has not been tested in individuals with active implants, and, therefore, the effects are unknown.

KEEP OUT OF REACH OF CHILDREN. The mē should not be used by anyone, or on anyone, under the age of 18. The use of the mē on or by children has not been studied and could result in serious damage.

The mē is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.

Warnings

DO NOT use if you are pregnant, trying to become pregnant (conceiving) or breast-feeding. The mē has not been tested on pregnant or breast-feeding women. Hormonal changes could increase skin sensitivity and the risk of skin damage.

DO NOT treat with the mē immediately after or before an

active tanning session - in natural light, on a tanning bed or when using sunless tanning products. If you are not certain whether you are tanning or sunburned it is recommended to do a skin sensitivity test (page 8) and wait 24 hours.

DO NOT use the mē on areas that have been sunburned as this can make the skin more susceptible to damage.

DO NOT use directly on parts of your skin where you have tattoos, dark brown or black spots (such as large freckles, birthmarks or moles). These areas may absorb too much light energy which can damage the skin.

DO NOT use on the ears, neck, nipples, genitalia or around the anus. These areas may have more sensitive skin and using the mē there may cause serious damage.

DO NOT use on or around the eyes, eyebrows or eyelashes. Doing so could cause permanent eye damage.

DO NOT use on a man's face, jaw or neck. Using the mē in these areas may result in skin damage.

DO NOT use the mē on any area where you may want the hair to grow back. The hair removal results may be permanent.

DO NOT use the mē if you are sensitive to light (photosensitive).

- Some diseases can be related to photosensitivity

including porphyria, polymorphic light eruption, solar urticaria and lupus.

- Some medications can be related to photosensitivity including: NSAIDS; steroids; antibiotics; phenothiazine; thiazide, diuretics, sulfonyluraes, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), Retin-A, Accutane® and/or topical retinoids.
- Undergoing radiation or chemo-therapy treatments.

These conditions can make the skin more sensitive to light and more susceptible to damage. Consult your doctor.

DO NOT use if you have an active or known skin condition in the treatment area including acne, sores, psoriasis, eczema, sunburns, melasma, herpes simplex, open wounds or a rash. These conditions can make the skin more sensitive to light and more susceptible to damage.

DO NOT use if you have a history of skin cancer or areas of potential skin malignancies. Use of the mē has not been tested in areas with skin cancer, and therefore, the effects are unknown.

DO NOT hold the treatment window still over one spot on the skin and emit multiple light flashes on the skin. Doing so may generate too much heat, which can damage the skin.

DO NOT cover the cooling vents during use, as they help

to keep the device cool during operation. Covering the cooling vents could result in the device overheating.

DO NOT get the mē damp or wet.

DO NOT use the mē while bathing or put it next to a bath or sink.

DO NOT operate the mē with a damaged plug or cord.

For means of isolation – disconnect the device from the mains.

Use only the power supply provided with the mē.

DO NOT open the outer covers of the mē or touch any inner parts.

DO NOT attempt to repair your mē device.


ALWAYS keep the cord away from heated surfaces. Failure to follow these precautions can cause electrocution, electric shock or result in damage.

ALWAYS attach the AC plug from the power supply into the device first, then into the wall socket.


DO NOT use the mē in an environment where a spark could cause an explosion, such as near a gas source.

Using your mē for the first time


Setting up your mē

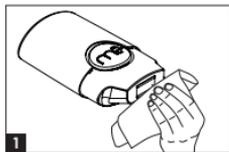
1. Clean the treatment window with a non-alcohol based wet wipe or soft tissue. This is done in order to remove any debris or dust from the treatment window.

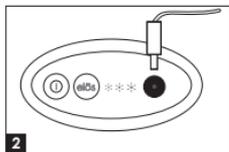
2. Connect the power cord into the power supply port found on the control panel of the mē.

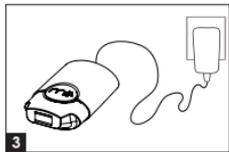
3. Connect the power cord to an AC outlet.

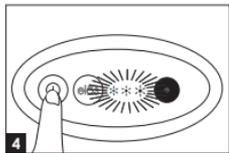4. Press the power button to turn the device on. All 3 energy level indicators will light up blue and start blinking.

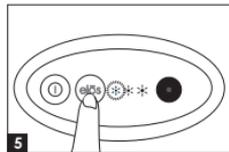5. Press the elōs button once to select the LOW energy level. When the energy indicator lights up steady blue, the mē is ready for use.

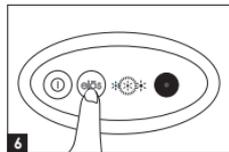6. To increase the setting to MEDIUM or HIGH, press the elōs button until the indicator displays the desired level.

Note: When moving from HIGH to LOW, the device will reset. The energy indicator lights will cycle momentarily before the LOW energy indicator light appears a steady blue. Proceed treating on LOW only after the the LOW indicator light is steady.

Skin sensitivity test

If you are using the mē for the first time or treating a new body area, choose a small area close to the area you intend to treat. The Skin Sensitivity Test should be performed on LOW energy level:

- Prepare your skin for the test and follow the treatment instructions.
- Wait 30 minutes and check your skin for any reaction. If the skin appears normal, then you are ready to perform a full treatment.
- If the treated spot appears red or irritated, please wait 24 hours and repeat the test. Otherwise consult your doctor.

How to use your mē

If this is your first time using the device, be sure to conduct the Skin Sensitivity Test.

Step 1: Preparation

To prepare for your mē treatment:


- First remove hair from the area you want to treat.
- Clean your skin and dry it with a towel.
- This will ensure that the elōs technology is focused on the hair follicles beneath the surface for optimal efficacy.

If you have light hair, it is recommended that you wax or epilate immediately before treating with the mē in order to get the best results. Make sure to clean your skin of any wax or hair residue before starting treatment. Do not treat if you feel your skin is irritated from waxing or epilation.

Step 2: Activate

Press the power button to turn on your mē device and select your energy level (LOW, MEDIUM or HIGH):

- **LOW:** Press the elōs button once to select the lowest energy level. The left energy level indicator light will turn blue (image A) and you will hear the fan inside the device.
- **MEDIUM or HIGH:** Use the elōs button to select higher energy levels for greater efficacy (image B, C). Always start treating on LOW and if you are comfortable, move to MEDIUM or HIGH energy level for greater efficacy.

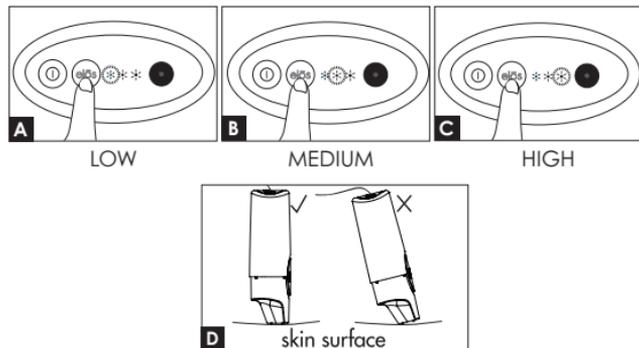Step 3: Treating your body

Gently place the device on your skin, so that both RF contact bars are in complete contact with the skin (image D). The device will not work unless both RF contact bars are touching your skin.

- The mē will emit flashes of bright light - this is the elōs energy at work.
- Slowly glide the device over the treatment area, without lifting the device move to the next spot only after the treatment flash has been emitted.
- Continue gliding the device over the area so that you pass over the same area 2-3 times.

Step 4: Treating your face

Prepare your treatment area (women's faces below the cheek line):


- You can shave (dark hair), wax or tweeze (light hair) before using the mē.
- Make sure to clean your skin of any wax or hair residue before starting treatment.
- Do not treat if you feel your skin is irritated from waxing or tweezing.
- Gently place the device on your skin, so that both RF bars are in complete contact with the skin. The device will not work unless both RF contact bars are touching your skin.
- Allow one flash of elōs energy on your first target area.
- Lift the device and place it on the the next adjacent area and wait for each flash of elōs light to be emitted before moving to the next spot.
- Repeat this stamping motion until you have completed 2-3 passes of the entire area.

Treatment tips

When using the mē on difficult to treat areas (e.g. underarms) apply the device in two directions: left/right and up/down for full coverage (image A).

- If you find the light from the elōs flash too bright, you can wear sunglasses or goggles for your comfort.
- DO NOT hold the treatment window still over one spot and flash multiple light flash on the skin. Doing so may generate too much heat, which can damage the skin.
- The mē will automatically shut off if it is not in use for more than 3 minutes. Press the power button to turn the device back on.

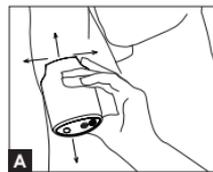What to expect during treatment

During treatment you may feel a warming sensation. You should always start your treatments at the LOW energy level and if you are comfortable, increase to MEDIUM or HIGH for even greater efficacy. If your skin looks slightly red after treatment, this should disappear within 24 hours.

Treatment schedule

Hair removal is a gradual process. Actual results will vary from person to person and depend on the body site being treated. Hair re-growth may fluctuate over time, depending on the person, body area, hair growth cycle and whether the instructions were followed correctly.

For best results, follow the full protocol of treating the target area once a week for 7 treatments. Use your temporary hair removal solution in between treatments. After your initial protocol, use your mē for touch-ups or as needed.

Note: Facial areas can be treated more frequently up to 2-3 times per week for 6 weeks if desired. Additional treatments may be needed for optimal results.


How to clean and store your mē

Turn the mē off by pressing the power button on the control panel. Unplug the power supply from AC outlet and power supply port (image A).

Clean the treatment window and outside surface of the mē using ONLY a non-alcoholic based wet wipe or soft tissue (image B).

Store the mē in a cool and dry place AWAY from the reach of children.

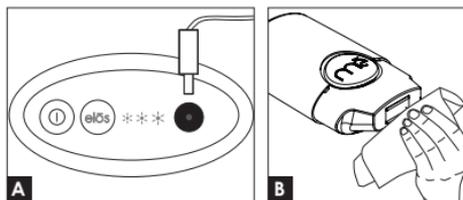Disposal/Recycling

Please follow your local regulations for disposal of electric systems when recycling or disposing of your mē. If you have questions, please contact us at www.mepower.com.

Waste Electrical and Electronic Equipment (WEEE).

Traveling with mē

The mē is designed with universal voltage. That means that you can use the device all over the world. Interchangeable plugs are sold separately and vary by country. These are designed to work with the mē power supply. Make sure that the mē does not get wet while travelling, and do not use it in a wet environment.

Troubleshooting

Problem	Possible cause	Solution
No indicator lights are on.	<ul style="list-style-type: none"> Power supply is not properly connected to the device or to an electrical source. Device wasn't turned on. 	<ul style="list-style-type: none"> Ensure the DC plug of the power supply is securely plugged into the device and the other end is connected to an AC outlet. If the power supply has been properly connected, press the power button located on control panel of the device.
All 3 energy indicator lights are blinking blue.	<ul style="list-style-type: none"> An energy level was not selected after turning on the device. 	<ul style="list-style-type: none"> Press the eIōs button once to select the LOW energy level. To select MEDIUM or HIGH energy levels, use the eIōs button until the indicator light displays the energy level you desire.
No energy light flashes from treatment window.	<ul style="list-style-type: none"> Device wasn't turned on. eIōs energy level was not selected by pressing the eIōs button. RF skin contact bars not in contact with the skin in treatment area. Device was left idle for more than 3 minutes and turned off. 	<ul style="list-style-type: none"> In order for the mē to emit flashes, make sure all 3 conditions (device on, energy level selected, RF skin contact bars in contact with skin) are met. If the device was idle, press the power button to reactivate the device for treatment again.
All 3 energy indicator lights are steady blue (not blinking).	<ul style="list-style-type: none"> Device flashes consumed. Cooling vents covered Room temperature too high. 	<ul style="list-style-type: none"> Turn off and disconnect the device from the power supply and contact our Customer Care Specialists for assistance. Let device cool. Turn the device back on to re-attempt your treatment. Move to a cooler room.
3 energy indicators are blinking in a continuous cycle.	<ul style="list-style-type: none"> When the device was set to HIGH energy level the eIōs button was pressed to lower the energy level to LOW. 	<ul style="list-style-type: none"> Wait for the indicator lights to stop cycling and the LOW energy level indicator is a steady blue. Then treat the desired area.

TO THE EXTENT PERMITTED BY LAW, ILUMINAGE SHALL NOT BE RESPONSIBLE FOR ANY DIRECT OR INDIRECT DAMAGES RESULTING FROM THE USE OR MISUSE OF THIS PRODUCT INCLUDING ANY GENERAL, INCIDENTAL, SPECIAL, OR CONSEQUENTIAL DAMAGES - TO THE MAXIMUM EXTENT PERMITTED BY THE LAW. THE WARRANTY PROVIDED IN THIS DOCUMENT IS THE SOLE AND EXCLUSIVE WARRANTY PROVIDED BY ILUMINAGE TO THE PRODUCT AND IT REPLACES ANY REPRESENTATION AND/OR WARRANTY INCLUDING FITNESS FOR ANY PARTICULAR USE AND MERCHANTABILITY.

This warranty does not cover products determined by Iluminage to be damaged by any of the following: accident, misuse, abuse or alteration, servicing by unauthorised personnel, use with unauthorised accessories or use other than as instructed, connecting to an incorrect current or voltage and/or any other conditions beyond Iluminage's control.

Manufacturer information

Labels And Symbols

The mē device shows information pertaining to the device's part number, serial number, manufacturer information and the electrical and regulatory standards the mē complies with.

This device is not suitable for use in the presence of a flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

Definition of Symbols

WEEE - Waste Electrical and Electronic Equipment

Follow operating instructions

Class II equipment

Do Not use in wet environment

Degree of protection against electric type BF applied part

CE Mark

Illuminage Ltd.
Yokneam Illit 2069200, IL

mē specifications*

Light spectrum: 550 - 1200 nm
Max optical: 4J/cm²
RF wave frequency - 6.78 MHz
Max radio frequency power - 5 Watt
Pulse light repetition rate: up to 1.1 Hz (1 pulse every 0.9 second)
Electrical rating - 100-240~0.4A, 50-60 Hz
Dimensions - 6 x 10 x 16 cm (W x L x H)
Weight - 0.3 kg

*Patented technology US6702808 and US5405368

Applicable Standards

Safety Standards

EN 60335-1:2012 / AC:2014;
EN 60335-2-23:2003+A1:2008+A11:2010 / AC:2012;
EN 60335-2-27:2013; EN62233:2008
IEC 62301 Edition 2.0 2011-01; EN 50564:2011

EMC Standards

EN 61000-6-1:2007
EN 61000-6-3:2007
IEC 60601-1-2:2007

Acceptable Environmental Conditions

Storage/Transport

Temperature - 10°C - 55°C (14°F - 131°F)
Humidity - 90% at 55°C (131°F)
Pressure - 700-1060 hPa

Operating

Temperature - 10°C-30°C (50°F-86°F)
Relative Humidity - Up to 80% at 37°C (99°F)
Altitude - Up to 2000m (6,500ft)

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© 2015 All rights reserved. mē is a trademark of Illuminage Beauty Ltd.
The mē is manufactured by Illuminage Ltd.