

In-Use Screen

After a Massage program is selected, In-use screen displays all settings to that program.



- 1 Time left
- 2 4D intensity level
- 3 Active airbags
- 4 Active rollers
- 5 Bluetooth connected
- 6 Indicators
- 7 Back-roller position & direction

Nature Sounds

From the SOUNDS button on the main menu you can select your choice of soothing sound therapy. Adjust volume to one of 5 levels also found under SOUNDS.

- White Noise
- Ocean Surf
- Crickets
- Rain
- Nature
- Airplane
- Car Ride
- Fireplace
- Thunderstorm
- Windy Night

Share a photo and tag us:

@Infinitymassagechairs

@infinity_chairs

@Infinitymassagechairs

Don't forget to use **#lovemyinfinitychair**

INFINITY®

72 Stard Road, Seabrook, NH 03874
603-910-5000
www.infinitymassagechairs.com

INFINITY®

Quick Start Guide

Genesis Max 4D Massage Chair

INFINITY®

Multi-Functional Screen

Use the directional buttons below to navigate these menu screens. Then press OK to select.

Power – Powers up the remote or returns the chair to the exit position and powers down remote


Home – Shortcut to the home page of all menu screens

OK – After navigating to the command desired on the display, push OK to select

Zero Gravity – Press for Zero Gravity 1, press again for Zero Gravity 2, press again to restore recline

Foot Roller – Press to turn off Foot Roller, press again to resume Foot Roller

Footrest – Press and hold up arrow to bring foot rest up, press and hold down arrow to bring foot rest down

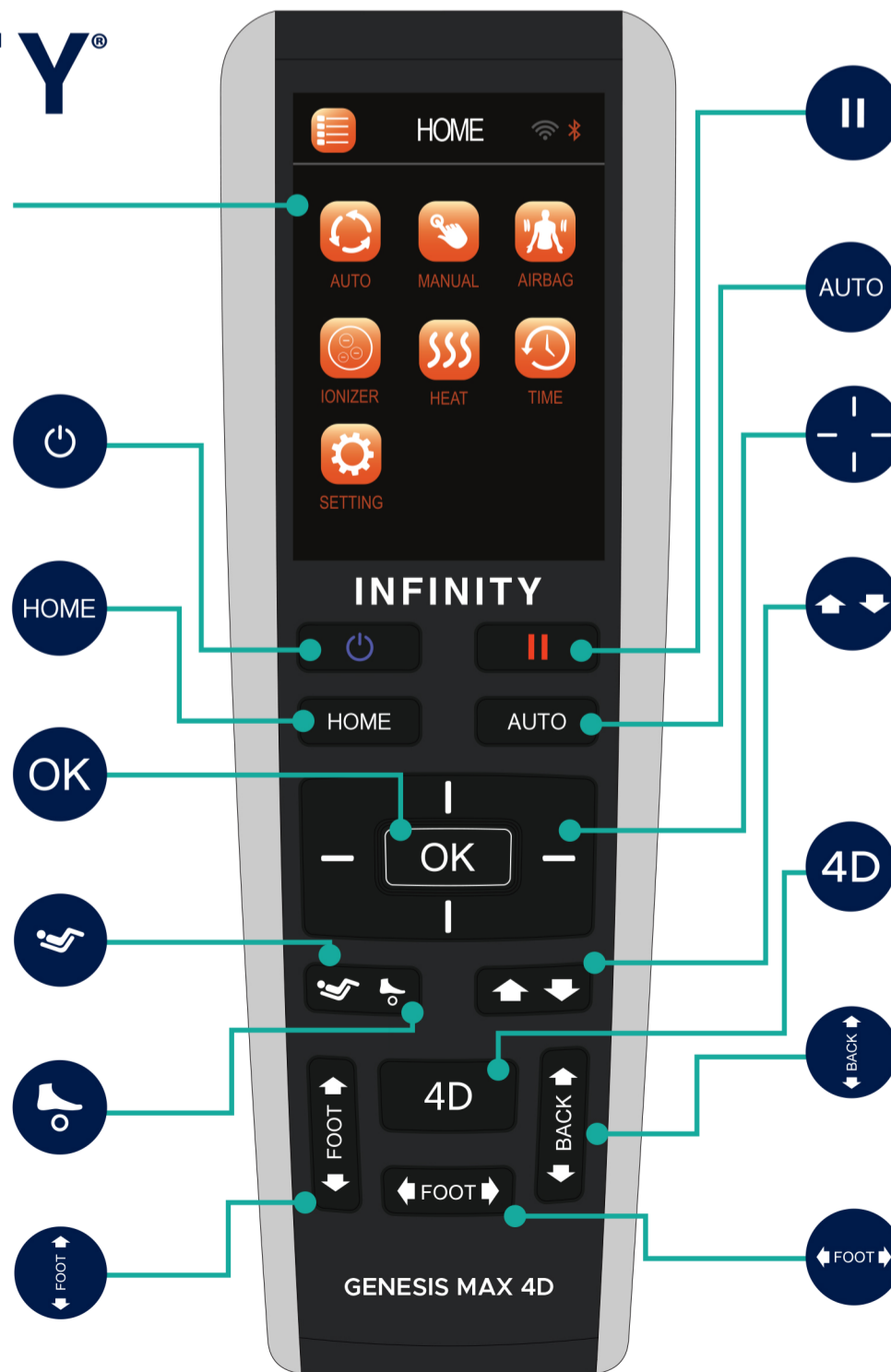

Pause – Press pause to immediately stop all massage chair activity, press again to resume

Auto – Shortcut to the Auto Programs menu screen

Directional Button – Press up, down, left or right to navigate the display screen

Back Adjustment – Back mechanism position adjustment; press and hold up to raise, press and hold down to lower. Use during initial body scan fine tuning or in manual point zone mode

4D – Press to increase or decrease 4D intensity of back mechanism. In-use screen displays intensity dial

Seat Back Position – Press and hold down arrow to recline the seat back, press and hold the up arrow to upright the seat back

Leg Length – Press and hold right arrow to extend the leg length adjustment, press and hold the left arrow to retract the leg length adjustment