

The brand name 'sonicsmooth' is written in a purple, lowercase, sans-serif font. The letter 'o' is replaced by a stylized gear or sunburst icon. A trademark symbol (TM) is located at the top right of the text.

by MICHAEL TODD BEAUTY

World's First

TRUE SONIC  
DERMAPLANING SYSTEM

INSTRUCTIONS  
FOR USE





# CONGRATULATIONS

Dear Customer,

Thank you and congratulations on choosing Michael Todd Beauty's SONICSMOOTH, one of the safest and most effective dermaplaning systems available for home use. We are confident that you will adore your instantly smoother, brighter and "fuzz free" skin for years to come.

SONICSMOOTH is the world's first sonic dermaplaning system specifically designed to expertly trim facial vellus hair (peach fuzz) while helping to gently remove dead skin cells and smooth away rough, uneven textured skin.

The SONICSMOOTH sonic dermaplaning system uses advanced sonic wave technology paired with Michael Todd Beauty's Safety Edge to expertly provide safe, instant and painless professional-grade dermaplaning results in the comfort of your own home.

3 Speed settings allow you dial up or dial down the power/intensity setting on the SONICSMOOTH so that you are able to customize your dermaplaning session to your specific speed preferences.

SONICSMOOTH is cordless and rechargeable with a 3 hour max. running battery life. The device handle is coated in a soft rubber for comfort, to prevent slipping and for ease of use. SONICSMOOTH has an ergonomically designed handle that comfortably sits securely in the palm of your hand.

SONICSMOOTH will quickly become your favorite weekly skin smoothing treatment. You are going to love the silky soft feeling of "fuzz free" skin and will immediately notice how much brighter your complexion looks. Your serums and creams will be able to work better after a SONICSMOOTH session because the dead skin cells and build-up will be gently sloughed away by state-of-the-art sonic technology. In addition, your makeup is going to sit flawlessly on your skin after removing the peach fuzz and layers of dead skin. ENJOY!

READ ALL INSTRUCTIONS BEFORE USE  
AND SAVE FOR FUTURE REFERENCE

## REGULATORY INFORMATION

U.S. This Class B digital product complies with Part 15 of the U.S. Federal Communications Commission (FCC) rules. Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Canada. This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference and (2) this device must accept any interference including interference that may cause undesired operation of the device.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes : (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement."

Mfg: Michael Todd Beauty  
Model: **sonicsmooth**

Device complies with  
CAN ICES-3 (B)/NMB-3(B)

For home or office use. Not intended for use in machinery, medical or industrial applications.

Radio and TV Interference Regulations.  
This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency (RF) energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications devices (such as televisions, baby monitors, cordless phones, etc). However, there is no guarantee that RF interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antennae of the other radio-communications device
- Move the equipment farther away from the radio or tv
- Connect the product into an outlet on a circuit different from that to which the other device is connected
- Consult the manufacturer's customer service or an experienced radio/tv technician for more suggestions.

# IMPORTANT SAFETY INFORMATION

BASIC SAFETY PRECAUTIONS. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER. When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF, YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

## DANGERS. To reduce the risk of electrocution:

1. **DO NOT** place or store the Charging Cradle where it can fall or be pulled into a tub or sink.
2. **DO NOT** place the charging cradle or plug into water or other liquid. The charging cradle should not be immersed or used in shower or bath.
3. **DO NOT** reach for the Charging Cradle or the plug that has fallen into water. If unplug while in water, might get electrocuted.
4. **DO NOT** use sonicsmooth in a location where the power cord can come in contact with water while it is plugged in.
5. Keep sonicsmooth away from open flame or flammable aerosol products. sonicsmooth is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

**WARNINGS. To reduce the risk of burns, electrocution, fire or injury to persons:**

1. Close supervision is necessary when this product is used by, on, or near children or physically challenged persons.
2. Use this product only for its intended use as described in these Instructions for Use. **DO NOT** use attachments other than those recommended and provided by MICHAEL TODD BEAUTY.
3. Keep the power cord away from heated surfaces.
4. Never use **sonicsmooth** while sleeping or drowsy.
5. **DO NOT** use **sonicsmooth** outdoors or where oxygen is being administered.
6. Never operate **sonicsmooth** if it is damaged.
7. **DO NOT** use the Charging Cradle for recharging if the power cord or plug is damaged, the Charging Cradle is not working properly, or if the Charging Cradle has been dropped on the ground or into water.
8. The Handle and Charging Cradle have no serviceable parts. Never attempt to open the Handle for any reason. If the Charging Cradle cord is damaged, it must be unplugged and replaced. Contact customer service at 772.800.7654 or email us at [customerservice@michaeltoddbeauty.com](mailto:customerservice@michaeltoddbeauty.com) for instructions on how to obtain a replacement.

9. Use **sonicsmooth** with the supplied power adapter only.

#### IMPORTANT CAUTIONS.

1. **sonicsmooth** is not a medical device and is not intended to treat medical conditions. If you are under a doctor's care for any skin condition, if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using **sonicsmooth**.
2. If you experience any pain or discomfort from using **sonicsmooth**, discontinue use and consult a physician.
3. To avoid risk of injury, we **DO NOT** recommend using **sonicsmooth** over the eyelid, directly on the nose or on the lips.

# CONTENTS

Your sonicsmooth comes with the following:

- A. Sonicsmooth Dermaplaning Tool
- B. Pre Treatment Cleanser 25 ml
- C. Post Treatment Gel 25 ml
- D. Medical Grade Steel Single-Use Safety Edges 8 wks supply
- E. US Wall Plug
- F. Charging-Storage Cradle





# GETTING TO KNOW YOUR SONICSMOOTH



# OPERATING INSTRUCTIONS

1. We recommend that you first prep your skin with the enclosed degreasing cleanser in the lavender tube. Sonic Dermaplaning is most effective on skin that is clean, dry with all traces of residual surface oil removed. If your skin feels a bit tight and dry after using the degreasing cleanser, do not worry. You will be applying a hydrating soothing serum after your session.
2. Before using your sonicsmooth for the first time, charge the handle for 12 hours. See charging instructions.
3. We do not recommend using the sonicsmooth directly over the eye lid, on the nose or lips.
4. Attaching/Removing Safety Edges:
  - sonicsmooth must be in OFF mode to change Safety Edges
  - To remove a Safety Edge from the Handle, with the device turned OFF, grip the outer edges of the Safety Edge with your thumb and pointer finger and slide upward away from the device handle
  - To attach a Safety Edge, with the Safety Edge still in the Cartridge Case, slide the device Handle through the opening in the Cartridge Case along the Safety Edge tracks until the Safety Edge "clicks" into place.
  - Tip: Use your pointer finger on the hand not holding the device Handle to hold the Safety Edge in place when attaching the device.
  - When the Safety Edge is secured into place, lift the device and attached Safety Edge from Cartridge




5. Starting: Press the ON/OFF Button to turn on
6. Selecting Speeds: To customize your dermaplaning session, sonicsmooth provides a choice of 3 power/intensity settings (low, medium and high). We recommend that you start on a lower speed setting and gradually increase power/intensity as you get used to sonic dermaplaning.

After pressing the on/off button one time the unit will power on in low speed. When you press the button again, the unit will change to medium speed. Pressing the power button a third time will change the speed to high and pressing the button one more time will shut off the device.

Through smart technology, if you hold the power button down for 2 seconds in your favorite speed setting, the unit will power off and remember your favorite setting the next time you power on the device. For example, when the unit is powered on in medium speed setting, hold the button down for two seconds to power off the device. When you press the button one time to power on the device, the handle will power on in the medium speed setting

7. Sonic Dermaplaning. After prepping your complexion with the degreasing cleanser in the lavender tube, rinsing clean and patting skin dry, you are ready to begin your session. There is no wrong way to dermaplane but if you would like to mimic the pattern your aesthetician uses to perform a dermaplane session, follow the tips below.

Caution: Pull hair back away from face and take precaution when using around hairline and eyebrows.

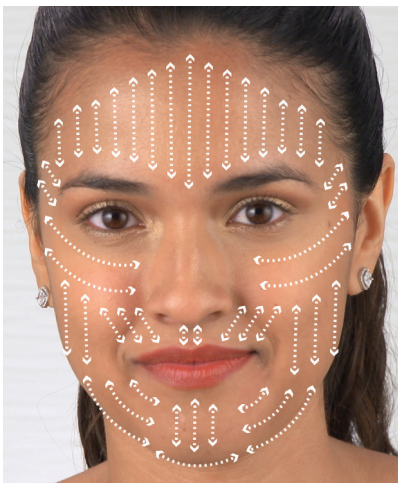

## Sonic Dermaplaning Tips

This is a graphic that shows the aesthetician recommended directions to dermaplane in different areas of the face for the best possible result.

8. Gliding movements. Hold skin taut with the second and third finger of the hand not holding the **sonicsmooth**. Starting near the two fingers holding skin taut, gently glide the **sonicsmooth** along the skin moving away from fingers.

Treatment areas should be between 1"-2" long and the width of your Safety Edge. You may gently glide over skin in a treatment area 1-3 times positioning the **sonicsmooth** at a 45° angle.

9. When you are finished, power off the **sonicsmooth** and return it to the charging cradle. Discard your Safety Edge.
  
10. After a sonic dermaplaning your skin may be slightly pink or red in color. This is normal after the complexion. Apply the **Sonicsmooth Soothe After Treatment Gel** to the skin. Any sensation of post treatment redness will fade in the next couple of hours.

## Expert Tips from an Aesthetician

1. Best Time to use **sonicsmooth**. Sonic Dermaplaning is one of the safest and most effective at home dual anti-aging / hair removal treatments. It is important to space treatments apart every 7-10 days.
2. When you perform sonic dermaplaning with **sonicsmooth** the complexion is excited resulting in a refreshed plumping effect and an increase in skin luminosity. As you dermaplane you are removing unwanted facial hair and also the top layer of dead dull skin as well as pore clogging debris. It's important to protect your skin from sun exposure, harsh chemicals and exfoliants around your dermaplaning sessions.
3. Ideally, performing **sonicsmooth** treatments at night before bed will allow your skin to refresh itself while you sleep. Keep in mind that after dermaplaning skin is more readily able to absorb serums and moisturizers because the dead layer of skin and vellus hair that blocks absorption has been safely removed. You may notice a dramatic improvement in the look and feel of your skin when your serums and creams are able to perform properly.
4. It is not recommended that you apply makeup directly after a **sonicsmooth** session. Allow freshly revealed skin to receive soothing nutrients from the **sonicsmooth Soothe After Treatment Gel** as well as revitalizing oxygen. Allowing the pores to "breathe" while you sleep will yield a refreshed glowing complexion in the morning.

5. Did you know that your makeup will look better when applied to skin that has been dermaplaned? It's true, liquid, cream and powder makeup will have a smooth fresh palette elevating the appearance of makeup when applied the morning after or hours after a **sonicsmooth** session.
6. We do not recommend using the **sonicsmooth** in the eye socket area, on the nose or directly on the lips. These areas are not flat in surface and are not suitable for the flat Safety Edge on the **sonicsmooth**.
7. Holding your **sonicsmooth** at a 45° angle is the safest and most effective way to sonic dermaplane. It is not recommended that you hold the **sonicsmooth** perpendicular to skin.
8. Unlike shaving other parts of your body, the best possible result from **sonicsmooth** is achieved by using light feathery strokes in a downward motion which is the same direction as hair growth. When we shave other parts of our bodies we typically shave against the direction of hair growth but when dealing with soft vellus facial hair, dermaplaning in the same direction as hair growth yields a better result.
9. When working in areas better accessed in a slightly sideways angle, it is suitable to angle the **sonicsmooth** to fit the contour of the face. An example of this is when working on and around the cheek bone area. It may be easier to successfully dermaplane this area by directing strokes from the hairline towards the center of the face. Just remember to hold skin taut with the opposite hand.

# CHARGING

1. Your **sonicsmooth** is powered by a rechargeable battery. Recharge the battery only as needed in the Charging Cradle. It is not necessary or recommended that the **sonicsmooth** be recharged after each use. Depending on the power setting used, the fully charged **sonicsmooth** should provide 3 hours of continuous use before needing to be recharged.
2. **First Use.** Before you use your **sonicsmooth** for the first time, you need to charge it fully for 12 hours. For subsequent recharging allow 8-10 hours.
3. **Charging Cradle.** Your **sonicsmooth** comes with a universal Charging Cradle to store and recharge the Handle. The universal Charging Cradle is compatible with international and domestic outlets (100-240 Volts AC; 50/60 Hz) with the respective country plug adapter. If the shape of the plug does not fit the power outlet, use an attachment plug adapter of the proper configuration for the power outlet.
4. To charge the **sonicsmooth**, plug the power adapter into a wall outlet and plug the DV jack into the back of the Charging Cradle. Place the handle securely into the Charging Cradle. The Charge Indicator will light up in red indicating the unit is charging. When the unit is finished charging the Indicator will light up in green. Unplug the adapter from the wall and Charging Cradle when not in use. **DO NOT** unplug the adapter from the wall outlet when it is still connected to the unit as this will cause discharging of the battery.

5. Infrequent Use. The Handle needs to receive a charge occasionally. If you rarely use the Handle, recharge it in the Charging Cradle every 3 months to maintain its charging capability.
6. Charge Indicator. When the charge is good the Power Light will be green. When the charge is low the Power Light will be red. After 10 minutes of use with a low charge the Power Light will Flash Red and the device will power off.

If the sonicsmooth is not sitting in the Charging Cradle Properly it will fall out and will not receive a charge.

If your sonicsmooth stops working, charge it in the Charging Cradle for a minimum of 5 hours. If the issue is not resolved, please contact Customer Service.





## CARE AND MAINTENANCE

1. Handle. To periodically clean the Handle, use an alcohol wipe. Using water on the handle may damage the unit.
2. Safety Edges. The Safety Edges are designed for single use. Dermaplaning with a fresh Safety Edge everytime is a good hygienic practice as well as eliminates issues associated with using a dull blade.
3. Charging Cradle. **DO NOT** immerse the Charging Cradle in water. Unplug the Charging Cradle before cleaning it. After unplugging the Charging Cradle, wipe with a soft, damp cloth.
4. **DO NOT** clean sonicsmooth or any of its components using a dishwashing machine.
5. **DO NOT** use abrasive cleaners or chemicals to clean your sonicsmooth.
6. Storage. When not in use for prolonged periods, store your sonicsmooth in a dry place at a temperature between 32°F and 140°F (0°C to 60°C).
7. sonicsmooth does not require any other regular maintenance. If the unit is not working properly, discontinue use and contact Customer Service at [customerservice@michaeltoddbeauty.com](mailto:customerservice@michaeltoddbeauty.com) or call (772) 800-7654.

## LIMITED WARRANTY

Please use this as your reference for any customer service or warranty needs you may have in the future.

You must register your product within 90 days of purchase in order to obtain the warranty. Your warranty is not effective unless you register your product within the time provided. In order to register your product go on-line to [www.warrantymichaeltoddbeaauty.com](http://www.warrantymichaeltoddbeaauty.com) and complete the information required for registration.

LIMITED 2 YEAR WARRANTY. This limited warranty gives you specific legal rights and you may have other rights which vary from jurisdiction to jurisdiction. Any claim for breach of this limited warranty must be commenced within 2 years following the original purchase.

For a period of two years from the date of purchase, Michael Todd Beauty LP ("MICHAEL TODD BEAUTY") warrants the SONICSMOOTH (the "Product") against Covered Defects. If there is a Covered Defect during the warranty period MICHAEL TODD BEAUTY agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty.

The limited warranty begins on the date of purchase and expires 2 years thereafter. In no event is the warranty period tolled or extended - including any period during which the Product (or any parts thereof) is in repair or if the Product is replaced by MICHAEL TODD BEAUTY during the warranty period.

This warranty is limited and applies only to a purchase of the Product from MICHAEL TODD BEAUTY or an authorized reseller. The warranty is for the original end user of the Product only and cannot be transferred under any circumstances.

Except as otherwise provided by law, to obtain service under this warranty, return the Product with a Covered Defect to us with your purchase receipt and, if you are returning the Product after one (1) year from the date of purchase, \$30 for postage and handling. For shipping instructions, contact our Customer Service department at [customerservice@michaeltoddbeauty.com](mailto:customerservice@michaeltoddbeauty.com) or call (772) 800-7654.

WHAT IS COVERED? "Covered Defects" means damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instruction manual included with the Product. Covered Defects does NOT include using the Product with a voltage converter; damage, defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; improper storage; improper or incorrectly performed maintenance or repair of the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer's instructions; normal wear and tear; alteration or overuse of the Product; as well as any professional or commercial use of the Product.

WARRANTY LIMITATIONS. THIS LIMITED WARRANTY IS EXCLUSIVE OF ALL OTHER WARRANTIES, WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED. THERE ARE NO IMPLIED WARRANTIES CREATED BY THE MANUFACTURE, SALE OR USE OF THE PRODUCT.

IN NO EVENT SHALL MICHAEL TODD BEAUTY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, STRICT PRODUCT LIABILITY, TORT OR ANY OTHER LEGAL THEORY.

THESE WARRANTY TERMS, CONDITIONS AND LIMITATIONS ALLOCATE THE RISKS OF PRODUCT FAILURE BETWEEN MICHAEL TODD BEAUTY AND YOU, THE END USER. THIS ALLOCATION IS RECOGNIZED BY BOTH PARTIES AND IS REFLECTED IN THE PRICE OF THE PRODUCT. BY USING SONICSMOOTH YOU ACKNOWLEDGE THAT YOU HAVE READ THESE TERMS, CONDITIONS AND LIMITATIONS, UNDERSTAND THEM AND ARE BOUND BY THEM.

SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IN SUCH STATES, THE LIMITS CONTAINED HEREIN AND THE DISCLAIMER OR EXCEPTIONS FROM WARRANTY MAY NOT APPLY. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEREIN.

PRODUCT CHANGES. MICHAEL TODD BEAUTY reserves the right to make changes and improvements to its products without incurring any obligation to incorporate such improvements in products previously sold or to notify you of these changes or improvements.

DISPUTES. This limited warranty is governed by the laws of the State of Florida, without regard to principles of conflicts of laws, and applicable United States federal law. The United Nations Convention on Contracts for the International Sale of Goods shall not apply. YOU AGREE TO SUBMIT ANY CLAIM, DISPUTE, OR CONTROVERSY (REGARDLESS OF THE FORM OR TYPE OF ACTION) AGAINST US TO ARBITRATION ADMINISTERED BY THE NATIONAL ARBITRATION FORUM UNDER ITS CODE OF PROCEDURE THEN IN EFFECT (available at <http://www.arb-forum.com> or via telephone at 1-800-474-2371). Except to the extent otherwise required by law, you agree that the arbitration proceedings will take place exclusively where the principal office of MICHAEL TODD BEAUTY is located at the time of institution of the proceedings. You may be entitled to file suit under state laws without submitting your claim to arbitration first. Information on arbitration may be obtained from, and a claim may be filed with, the National Arbitration Forum at P.O. Box 50191, Minneapolis, MN 55405-0191 USA.

Distributed by Michael Todd Beauty LP

Port St Lucie FL 34953

For more information, updates, reviews frequently asked questions and more visit us on the web at [www.michaeltoddbeauty.com](http://www.michaeltoddbeauty.com)

© 2018 MICHAEL TODD BEAUTY LP. All rights reserved.

