



nuovaluce
BEAUTY

Nuovaluce
ANTI AGING DEVICE
WITH MICROCURRENT & LIGHT THERAPY
TECHNOLOGY



ULTIMATE DEVICE TO REDUCE
WRINKLES AND TIGHTENING SKIN

Table of Content

	PAGE
INTRODUCTION	1
What is the Nuovaluce Anti-Aging Device?	1
Nuovaluce Anti-Aging Microcurrent & Light Therapy Device Includes	1
Nuovaluce Anti-Aging Microcurrent & Light Therapy Device Benefits	2
HOW TO USE	2
A. First-Time Usage	2
B. Need Recharging	3
C. Changing the Treatment Head	3
MICROCURRENT TREATMENT & HOW TO USE IT	4
What is the Microcurrent Facial Massage Treatment?	4
So How Does Microcurrent Facial Massage Treatment Work?	4
Intensity Level Settings	6
THE RED PHOTON LIGHT THERAPY TREATMENT & HOW TO USE IT	7
What is the Red Photon Light Therapy Treatment (Not applicable for eye area)?	7
So How Does Microcurrent Facial Massage Treatment Work?	7
FAQ	9
ANYTHING ELSE	16

INTRODUCTION

What is the Nuovaluce Anti-Aging Device?

Nuovaluce Anti-Aging Device combines the power of Microcurrent treatment and red photon light therapy to reduce the signs of aging.

It puts the power of light and Microcurrent technology into your hands in an ergonomically designed hand held device that stimulates the facial muscles and promotes the production of collagen and elastin to treat periorbital wrinkles.

Thus, giving you skin that is firmer, smoother, and more radiant looking.

As such, it is a perfect pairing of science and technology.

Nuovaluce Anti-Aging Microcurrent & Light Therapy Device Includes:

- Nuovaluce Anti-Aging Device
 - Microcurrent treatment stainless steel head - Good for firming and toning your skin.
 - **Red light** treatment stainless steel head - Good for healing skin tissue and reducing wrinkles.
- Product Charging Cable



Nuovaluce Anti-Aging Microcurrent & Light Therapy Device Benefits:

- Light - You can fit it in the palm of your hand no problem!
- Easy to use - It's as simple as choosing a head and pressing a button to get the results that you want on your skin.
- Easy to carry around - You can take on the go wherever, whenever you want to.
- Multifunctional - The device comes with two interchangeable heads: one for Microcurrent tightening and one for **LED Red Light Therapy**.
 - Microcurrent therapy stimulates the facial muscles and promotes the production of collagen and elastin to treat periorbital wrinkles.
 - **Red LED light therapy** repairs the tissue on your skin in order to reduce periorbital wrinkles and imperfections on your skin.
 - Altogether, the Microcurrent and **Red Light Emitting Diode** technology stimulate collagen production, thereby giving you that fuller, more youthful face you always dreamed of having. **NO ADDITIONAL SURGERY PROCEDURES, CREAMS OR SERUMS REQUIRED!**

HOW TO USE

A. First-Time Usage

After receiving the device, please make it fully charged by charging it at least 8 hours each and every single day (when needed). When all the **blue lights** on the device stop blinking, it means that the device is fully charged.

NOTE: You can charge it easily by using your mobile charger, laptop, desktop, pad etc.

B. Need Recharging

Situation 1: When the device is in use, all **blue LED lights** flash 5 times and the device beeps 5 times before shutting down.

Situation 2: When charging, **Blue LED lights** will indicate the charging status.

Situation 3: Fully charged battery. No lights appear!

C. Changing the Treatment Head

Step 1: Press the release button  on the back of the Nuovaluce Anti-Aging Beauty device to release the treatment head.

Step 2: Select the Nuovaluce Anti-Aging Beauty treatment head of your choice and attach it to the body of the device.

NOTE: The device does NOT vibrate when in use. **DO NOT** change the treatment head when the power is on.

Please check the pictures below for clarification:



Attached
Microcurrent Head

Attached
LED Light Head

MICROCURRENT TREATMENT & HOW TO USE IT

What is the Microcurrent Facial Device Treatment?

The Nuovaluce microcurrent treatment head is used to define and contour larger facial areas including the jowl, jawline, cheeks and forehead and eye area.

It goes 3 layers deep into your skin and reactive your natural collagen production. This way, your skin tones and tightens all over your face.



So How Does Microcurrent Facial Device Treatment Work?

Step 1: Thoroughly cleanse skin, removing all makeup prior to using the Nuovaluce Anti-Aging Beauty device.

In order to ensure that the skin is receptive for maximum benefits, gels, serums, moisturizers or creams should not be applied until the treatment is completed.

If you have dry skin and facing problem to glide the microcurrent head, we recommend that you use aloe vera gel or any moisturizer cream.

Step 2: Attach the facial treatment head to the body of the Nuovaluce Anti-Aging Beauty device.

Step 3: Press the power button  once to turn on the device.

Step 4: Adjust the unit according to the *intensity level (starting at 0 and ending at 4) by pressing the up (+) or down (-) button located on the front of the device.

The default intensity level is set at level 0 (no microcurrents). We advise that initially start at the lowest setting at 1-3 and increase the setting as you get comfortable using the device.

**Please see page 7 for more details*

Step 5: While lifting and holding the skin in place for 5 seconds, gently press both spheres of the light head onto your face, or glide the device in a gentle upward motion along your face, according to the image provided.

The Nuovaluce Anti-Aging Beauty device will shut down automatically after 10 minutes of continuous use.

Step 6: Apply your routine daily skin-care regimen when treatment is completed.

NOTE: You can charge it easily by using your mobile charger, laptop, desktop, pad etc.



Gently pressing the microcurrent spheres onto your face

For treatment around eyes like under eye bags, wrinkles around eye and to tone your eyelids, hold the device in place for 5 seconds.

Intensity Level Settings

There are five intensity levels for the Nuovaluce Anti-Aging Beauty microcurrent facial treatment. facial treatment (levels 0 -4). While 0 (the default level setting), contains no currents, 4 contains the most.

It is recommended to begin with the first level and gradually increase levels over time, depending on your individual skin therapy needs. Be aware that the default level setting is level 0 (no current).

1. Press to power on and intensity level 0 begins (1 blue light appears).
2. Press (H) once: Switch to intensity level 1 (2 blue lights appear).
3. Press (H) twice: Switch to intensity level 2 (3 blue lights appear).
4. Press (H) three times: Switch to intensity level 3 (4 blue lights appear).
5. Press (H) four times: Switch to intensity level 4 (5 blue lights appear).
6. Press (L) to lower the intensity level.

The use of intensity level will be as follows :

- First two weeks intensity level: 0-1
- Next two weeks level: 1-2
- Next two weeks level: 2-3

In other words, every two weeks, the intensity level will gradually increase.

NOTE: The Nuovaluce Anti-Aging Beauty device will automatically shut down after 10 minutes of continuous use. To shorten the usage time, press the power button  for one second. The device will beep three times before shutting down.

THE RED PHOTON LIGHT THERAPY TREATMENT & HOW TO USE IT

What is the Red Photon Light Therapy Treatment (Not applicable for eye area)?

The Nuovaluce Anti-Aging Beauty treatment's **red photon light therapy** is designed to target fine lines and wrinkles and discoloration on your skin while healing skin tissue in order to improve skin tone and promote collagen production.



Our Red Photon Light Therapy Treatment has a wavelength of 620 nanometers, which is strong enough for you to see results in as little as a couple of sessions without damaging your skin.

So How Does Red Light Therapy Treatment Work?

Step 1: Thoroughly cleanse skin and remove all makeup prior to applying the Nuovaluce Anti-Aging Beauty treatment.

Step 2: Attach the **red photon light** head to the body of the device.


Step 3: Turn on the Nuovaluce Anti-Aging Beauty device to begin the treatment.

Step 4: Place the Nuovaluce Anti-Aging Beauty lights in direct contact with the skin and hold for 5-8 seconds before moving to the next area, using this technique to treat each area of the face. See images below!

Please note that the device will automatically shut down after 10 minutes of use.

Step 5: Apply your routine daily skincare products after treatment is complete.

NOTE: For best results, use Nuovaluce Anti-Aging Beauty treatment for 10 minutes each therapy 5 minutes, 3 to 4 times each week.

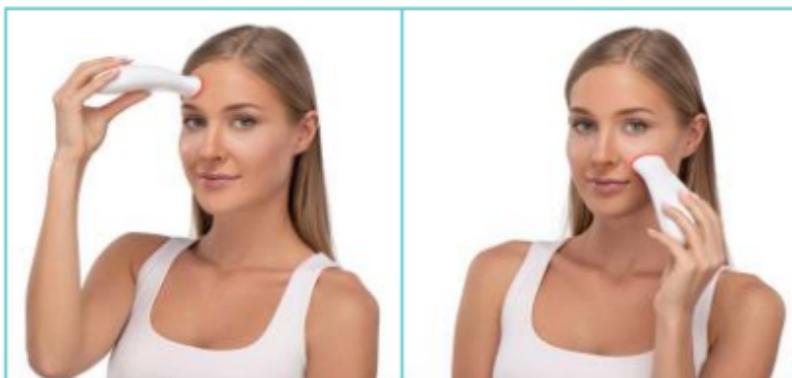

Gently pressing the Red Light Therapy onto your face

FAQ

Q: How does the Nuovaluce Anti Aging Device tone my skin? How much current will I feel?

A: There are five levels with the Nuovaluce Anti Aging Device (0-4) and with each level the intensity of the microcurrent increases, starting at 0, the lowest.

Depending on which level you choose, you will experience a low-level electrical current that stimulates muscle tone at your dermis level.

It mimics the microcurrent of your body, stimulating collagen, thereby making your skin look more beautiful and youthful as you sit back and relax to the feel of a light tingling sensation.

Q: Am I able to use the Nuovaluce Anti Aging Device with my existing skin care products?

A: Yes! Absolutely! The Nuovaluce Anti Aging Device works harmoniously with many different skin care lines. So you can use any existing skin care products with the device.

Q: How fast does the Device work in and when will I see results?

A: On average, you should be able to see results within 2 weeks. However, the speed of results defer from person to person. Most of our customers are seen satisfactory and noticeable results after 60 days of use.

Q: Am I able to use the Device if I plan on getting fillers, botox, or plastic surgery?

A: Yes, The Nuovaluce Anti-Aging Device is FDA-cleared and as such is safe to use on your skin. If you plan on getting fillers, botox, or plastic surgery, you should be able to use the Device without any problems.

Q: I am facing problem in gliding the device on face. What should I do now?

A: If you have dry skin, we recommend to use aloe vera gel or any brand moisturizing cream before using our device. It will help for smooth gliding on skin.

Q: What kind of results will I start to see when I use the Nuovaluce Anti-Aging Device?

A: When using the microcurrent toning treatment, the device will rebuild protein such as collagen, making the face tone and tighten. Areas that have seen the most effect have been bags under the eyes, eyelids, upper lips, jawline, and the entire face in general.

When using the light therapy treatment, you will see the appearance of fine lines

and wrinkles go away. It will also remove discoloration, dark spots, and reduce the size of your pores.

Q: How often can I use my Nuovaluce Anti Aging Device?

A: Your Nuovaluce Anti Aging Device should be used 3-4 times a week for the first 90 days to the setting of your choice as long as it is comfortable for you.

After the initial 90 day, you may scale down to using it 2-3 times a week to maintain your results.

During one session, we recommend that both the microcurrent toning treatment and the light therapy treatment be done for 10 minutes each.

Q: I have seen people using the device differently then what you guys are suggesting. Who should I follow?

A: While we do understand that others might learn better from observing users, we, at Nuovaluce Beauty, suggest that you first

read the directions and watch the videos to see how different people use the device in ways that have worked for them.

This will allow you to figure out which way is best for you, giving you the best results for your beautiful skin.

Q: When shouldn't I use the Nuovaluce Anti Aging Device?

A: If you...

- Have a pacemaker
- Have seizures
- Are pregnant
- Under the age of 18

Then you shouldn't use the device.

Q: In what order should I use the Microcurrent Toning and the Red Light Therapy ?

A: We recommend that you use the Nuovaluce Microcurrent Toning attachment first, then the Nuovaluce Red Light Therapy attachment head, and finish off with the regular serums and moisturizers you use.

Q: Do I have to wear goggles when using the Nuovaluce Red Light Therapy ?

A: No, you do NOT need to wear goggles because the coverage of the red LEDs is not broad enough to reach the eye area. It is designed to focus only on the coverage of the problem areas.

Q: Should I use the Nuovaluce Microcurrent Toning Device or the Red Light Therapy more often?

A: Neither! You can use both or one more than the other, whichever works best for you. We recommend using both treatments every time you use the device.

Q: How do I register my warranty?

A: To register your warranty, follow these steps below:

1. Visit <https://nuovalucebeauty.com/pages/warranty> Click on the “warranty” tab located at the top of the Nuovaluce website. Through this tab, you find the warranty registration form.
2. Fill out all of the contact information and product details.
3. Hit the “submit button” located at the bottom of the form.

That’s it!

Q: What Payment Methods do you accept?

A: You can make payments on our website using PayPal, Google Pay, Afterpay or by Credit Card.

Q: Do you offer a payment plan?

A: Yes, we offer a payment plan! You are able to pay for your Nuovaluce device over four equal installments using Afterpay. Simply select, “Afterpay,” at payment method page.

Q: How does your Device compare with Others in the Market?

A: Our device offers both microcurrent and red light therapy to maximize results. Thus giving you tighter and more beautiful looking skin compared to other leading products in the market at a fraction of the price.

Q: Can I use the device if I have broken capillaries?

A: Yes, absolutely! Nuovaluce Anti Aging Device is not only safe to use on skin with

broken capillaries, but helps reduce their

appearance on the skin. Additionally, it helps to reduce skin discoloration and dark spots, as well as shrink pores.

Q: How does this Device Compare to Professional Microcurrent Treatments?

A: The Nuovaluce Anti Aging Device works just as well as professional microcurrent treatments you would find at a dermatologist or esthetician office, only at a fraction of the cost and available for use in the privacy of your own home!

Q: What do I do if I am not feeling the current in the microcurrent attachment head?

A: Some users feel it more than others. As such, you may adjust the microcurrent level to your comfort.

Q: Can I use the device more often to get better results?

A: No, it is not recommended to use your Nuovaluce Anti Aging Device more often than 3-4 times per week. Besides, by using our product for just 10 minutes, we promise that you will see the results that you want it on time.

Q: Are there any benefits to using the red light therapy on its own?

A: Yes, the Nuovaluce Anti Aging Device red light therapy attachment on its own reduces the appearance of fine lines on skin with its tissue healing abilities. In addition, it removes discoloration, dark spots, and capillaries.

Q: Are there any benefits to using the microcurrent on its own?

A: Yes, microcurrent therapy on its own has many benefits such as toning, tightening,

and lifting the skin and reduces puffiness under the eyes. Also, the microcurrent helps rebuild skin proteins, including collagen.

Q: I have rosacea. Is this device ok to use still?

A: Yes, it is okay to use the Nuovaluce Anti Aging Device if you have rosacea ! Our product is safe for use on all skin types.

Q: Do you run a credit check for the AfterPay payment plan?

A: No, there is no credit check for the payment plan offered. It is simple and very fast. Simply select the, "Afterpay Payment Option" at checkout.

Q: Where does your product ship from?

A: Our product ships from San Diego, California. We offer " Free 2 Day Shipping ".

Q: Do you ship internationally?

A: Yes, absolutely! We ship Internationally!

ANYTHING ELSE?

In that case, feel free to contact us by going to <https://nuovalucebeauty.com> click the “contact us” button and fill out the form.

We’ll be more than happy to answer any additional questions or comments you have about any of our products.