

## Shiatsu Foot Massager with heat



**Item no: SFM630  
Instruction Manual**

### IMPORTANT SAFEGUARDS

WHEN USING ELECTRICAL PRODUCTS, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING

#### DANGER: TO REDUCE THE RISK OF ELECTRIC SHOCK:

- ALWAYS unplug this appliance from the electrical outlet immediately after using or before cleaning.
- DO NOT reach for an appliance that has fallen into water. Unplug it immediately.

#### WARNING: TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, FIRE OR INJURY TO PERSONS:

- An appliance should NEVER be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
- When children are present, close supervision is necessary.
- DO NOT use outdoors.
- Use this appliance only for its intended use as described in this manual.
- DO NOT use any attachments not provided with the unit.
- DO NOT carry this appliance by the power cord or use the cord as a handle.

2

- NEVER operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped into water.
- Keep cord away from heated surfaces.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- To disconnect, turn off the machine by pressing the Power button, then unplug from outlet.
- Use heat function carefully. Heated surfaces may cause serious burns. DO NOT use on sensitive skin areas or if you have poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- Do not stand on the appliance. Use only while seated.

3

### SAVE THESE INSTRUCTIONS

CAUTION: PLEASE READ ALL INSTRUCTIONS CAREFULLY BEFORE OPERATING.

- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Consult your doctor prior to using this product, if
  - You are pregnant
  - You have a pacemaker
  - You have any concerns regarding your health
- NOT recommended for use by diabetics.
- NEVER leave the appliance unattended when plugged in, especially if children are present.
- NEVER cover the appliance when it is in operation.
- DO NOT use this product for more than 15 minutes at a time.
- Extensive use could lead the product to overheat and shortened lifespan of product. Should this occur, turn the unit off and allow the unit to cool down before operating.
- NEVER use this product directly on swollen, inflamed area, rashes or wounds.
- DO NOT use this product as a substitute for medical use.
- DO NOT use this product before bed. The massager has a stimulating effect and can delay sleep.
- NEVER use this product while in bed.
- This product should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- This unit should NOT be used by children without adult supervision.
- For home use only. Not for commercial use.

NOTE: This unit is equipped with a protective thermostat that will shut off the motor if it overheats. If this happens, unplug the unit and let it cool down for 30 minutes before plug in again and continue use.

4

### Instructions for Use

1. Plug into a 120-volt outlet.
2. It is highly recommended to wear sock while using the unit. Do not use the unit while wearing shoes.
3. Place massager on the floor, in front of your chair or seat.
4. To activate massage action, press Power button once (Fig.1).  
Tip: By changing the pressure of your feet on the massager, it will change the feel of the massage. If the unit tickles you, then you may need to apply more pressure or you may need to become accustomed to the sensation. If the massage feels too strong, you may need to lighten the pressure of your feet on the unit.
5. To activate soothing heat with the massage, press the POWER button one more time (Fig.1). Allow a few moments for the heating function to warm up.
6. When finished, press the POWER button one more time to turn the unit off and unplug.



Fig. 1

Note: As a safety feature, the soothing heat function can be used only with the massage feature.

5

### Maintenance

#### To Store

- Place massager in its box or in a safe dry cool place.
- Avoid contact with sharp edges or pointed objects.
- To avoid breakage, DO NOT wrap the power cord around the unit.
- DO NOT hand the unit by the cord.

#### To Clean

- Unplug the unit and allow it to cool down before cleaning.
- Clean only with a soft damp cloth.
- NEVER allow water or any other liquids to come into contact with the unit.
- DO NOT immerse in any liquid to clean.
- NEVER use abrasive cleaners, brushes, or strong chemicals that might or might not be flammable and/or damaging to the product to clean.
- DO NOT attempt to repair the massager. These are no user serviceable parts. For customer service, please contact service@aurorahhealthbeauty.com.


### Technology Parameter

Name: Shiatsu Foot Massager  
Item NO: SFM630  
Voltage: DC12V, 2A  
POWER Consumption: 24W  
Made in China



6

## Shiatsu Foot Massager with heat

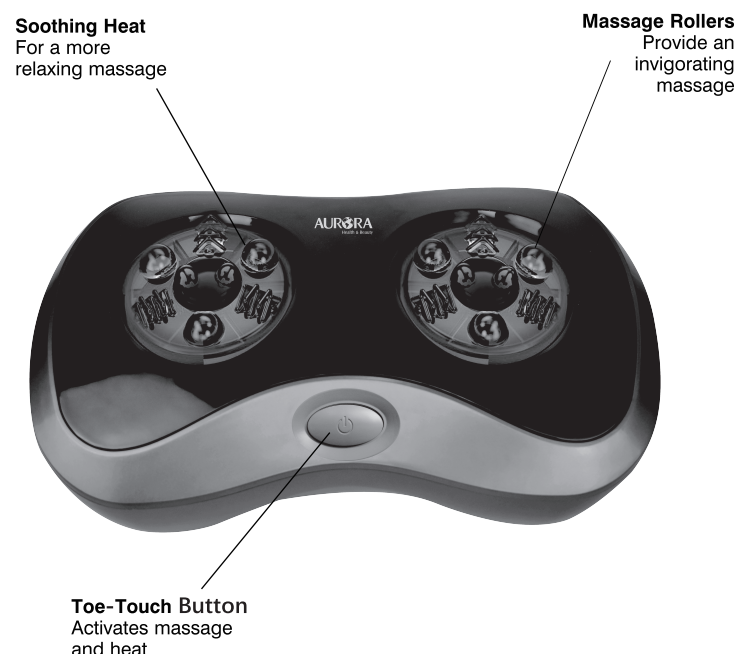

7