


MASSAGER SEAT

Item no: MSS610

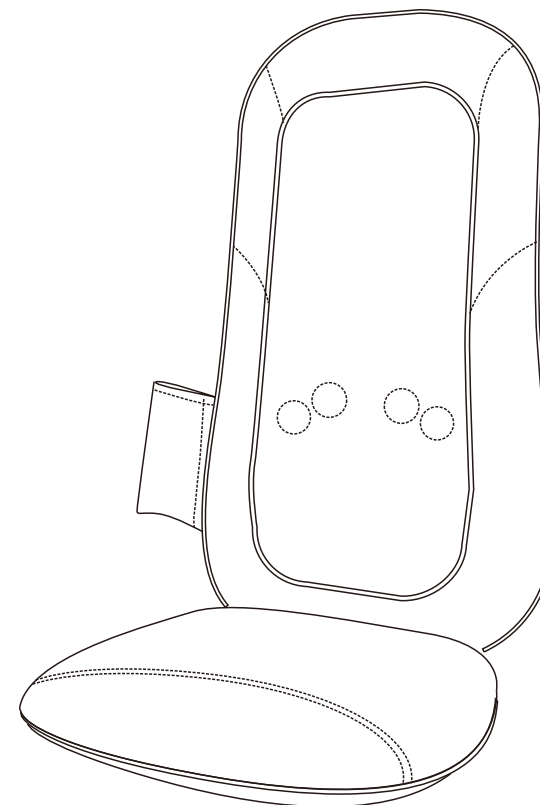

IMPORTANT SAFETY INSTRUCTIONS

When using the product, all instruction and basic precautions should always be followed, including the following:

DANGER – To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use or before cleaning.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not in use, before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this product is used by, or near children, invalids, or disabled persons.
- 4) Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) If the supply cord is damaged, it must be replaced by the manufacture, its service agent or similarly qualified persons in order to avoid a hazard.
- 6) Do not carry this product by its supply cord or use the cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never drop or insert any object into any openings.
- 9) Do not use outdoors.
- 10) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 11) To disconnect, switch the ON/OFF power off and remove the plug from the outlet.
- 12) The product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or the lack of experience and knowledge - unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
- 13) Children should be supervised to ensure that they do not misuse the product.
- 14) Use heated surface carefully. May cause serious burns. Do not use over insensitive skin areas or if you have poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

This manual should be kept on hand
for easy reference.

It is your guide to a safe and
efficient operation.

SAVE THESE INSTRUCTIONS

This product is intended for household use only

HOW TO OPERATE

Plug the power into a power source.

When the power is on, the massager nodes will restore into the lowest position.

Press the ON/OFF button and the massager will work in AUTO mode with shiatsu massage functions. The auto mode setting is programmed to 15-minutes interval.

SUITABLE CROWD




Officer/ Student

Using the computer for a long time leads to muscle stiffness.



Elderly people

Muscle tension and backache.

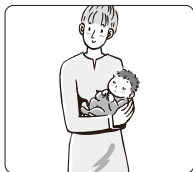

Housewife

Muscle soreness after busy house work.

HOW TO OPERATE



ON/OFF

Press the ON/OFF button to start or stop the device.



AUTO

Press the AUTO button and the massage nodes will work up and down. This cycle will last 15 minutes.



UPPER BACK

Press UP button to adjust the positioning to the upper back.



LOWER BACK

Press DOWN button to adjust the positioning to the lower back.



HEAT

Press the HEAT button to turn the heating function on or off. The massager node will light up in red with the heat is on.




VIBRATION

Press the VIBRATION button to turn the vibration function on or off.



PART LISTS

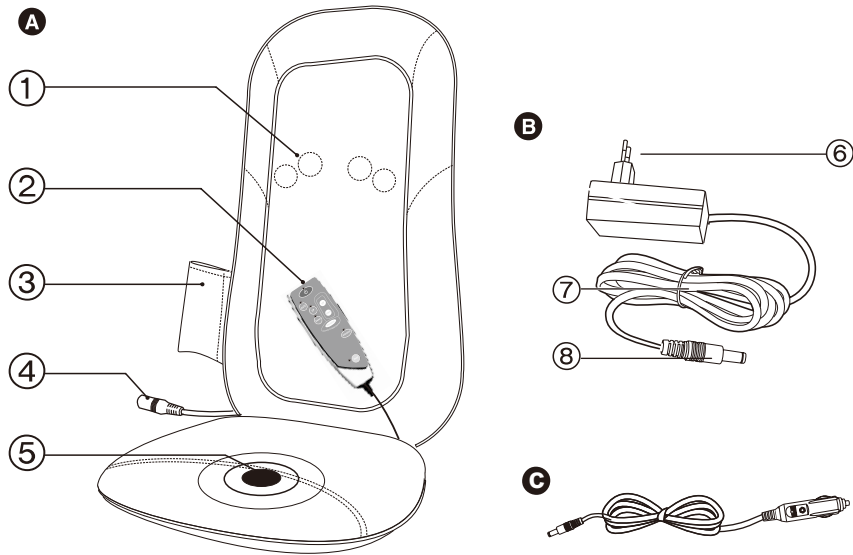

- A** Main unit
- ① Massager Nodes
 - ② Remote Control
 - ③ Remote Control Pocket
 - ④ Dc Connector
 - ⑤ Seat Vibration

- B** Power adaptor
- ⑥ Plug
 - ⑦ Cord
 - ⑧ Appliance Plug
- C** Vehicle Power Adaptor

MAINTENANCE INSTRUCTIONS

1. Cleaning

Please ensure the ON/OFF button is off and power is not plugged into the outlet prior to cleaning. Use dry cloth to clean the product.



ATTENTION

- Do not use benzene, thinner, or any other aggressive detergents to clean this product in order to avoid the colors fading!

2. Maintenance

Ensure the ON/OFF button is off and power is not plugged into the outlet after each use. Keep the product clean and store in a safe, dry and cool area.

TECHNOLOGY PARAMETER

Name: Massager Seat
 Item No: MSS610
 Voltage: DC 12V, 2A
 Power Consumption: 24W
 Made in China

