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HEAD MASSAGER

Manual

Tips: Please read this manual completely before use, and keep this manual for future reference


Product Introduction

Multi-dimensional and multi-point massage according to the characteristics of human head acupoints. By massaging head dot-matrix distribution, different vibration frequencies, multi-speed mode to have different massage experience massage. It can effectively relieve the symptoms of insomnia, dreaminess, dizziness, mental stress, anxiety and tension caused by study pressure and work fatigue.

Product parameter

Product name:	Multifunctional Bionic Head Massager
Product model:	IKR-HM-210/IKR-HM-214/IKR-HM-310/IKR-HM-314
Product size:	162 X 162 X 152mm
Packing size:	175 X 175 X 167mm
Material:	ABS + Silicone
Input:	5V
Battery capacity:	1100mAh


Product structure and operating instructions

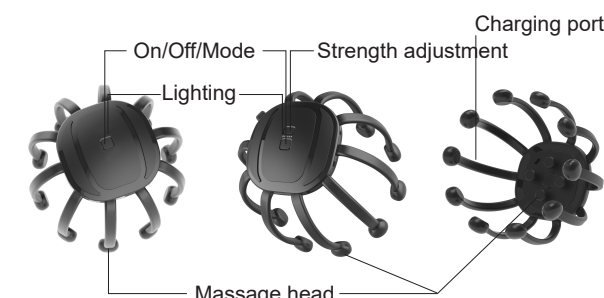- Press the power button for 3 seconds to turn on.
- Double-click the button to turn on/off the music function (voice and music start by default when open).
- Single click the button to switch modes:
① Dynamic mode; ② Happy mode; ③ Health mode; ④ Decompression mode; ⑤ Soothing mode, each mode has different massage methods and vibration intensity to meet different needs.
- IKR-HM-310/314 series short press the strength key to increase the strength in turn, 5 strengths can be adjusted.

- It will automatically shut down after 20 minutes of work by default, or press and hold the power button for 3 seconds to shut down
- Bluetooth function:

According to the voice prompt, open the phone settings, turn on the Bluetooth, search for "BT-IKR" and click to connect, as shown in the figure.

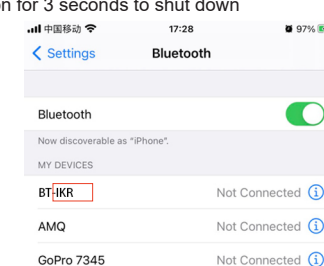Tips

- The indicator light turns red and flashes for 5 seconds to automatically shut down when the battery is low. Please charge in time.
- The red light flashes slowly when charging, and the blue light is on when it is fully charged.
- Cannot be used while charging.

Packing

- Product*1
- USB charging cable*1
- Product Manual*1
- Product warranty card*1

Precautions

- Do not use this product in places with high heat, flammability, electromagnetic radiation, and humidity;
- Do not use this product under thunder or lightning conditions;
- Do not disassemble, repair and modify this product privately, as it may cause accidents or malfunctions;
- This product is non-waterproof design, it is strictly forbidden to enter the product, if it enters water, please do not use it;
- Do not use for infants, head trauma, or people who have had recent head surgery;
- Patients with deep stroke, hypertension and hypotension should be used with caution, please consult a professional physician. If you feel unwell during use, please stop using it immediately;
- Use with caution, slow movement, physical and mental discomfort, or mental retardation, please use it under the supervision of an escort. If you feel unwell during use, please stop using it immediately;
- Please avoid using this massager when drinking alcohol or being unwell;
- This product is a personal health care product and cannot replace medicines and medical treatment periods;
- This product is a household electronic health care product, not suitable for diagnosis, treatment and other commercial purposes.

Warning

Do not use for infants, head trauma, or people who have had recent head surgery.

Troubleshooting

If there is a failure during the applicable operation, please perform a check. If the failure cannot be eliminated or the product cannot operate satisfactorily.

Basic failure	Possibly caused by	Basic elimination method
Can not turn on	Weak battery	If yes, please charge
Vibration intensity is too weak	Weak battery	If yes, please charge
Head vibration too weak	Whether the massage head touches head	If yes, the massage head need to closer to head

Maintain

Please pay attention to moisture-proof, rodent-proof, and pollution-proof. Please be sure to turn it off when cleaning the product. The product should be placed in a safe place after use, do not put heavy objects on the product. When there is a malfunction, you should ask a professional to repair it.

Storage and transportation

- Temperature: $-20^{\circ}\text{C} \sim +40^{\circ}\text{C}$.
- Atmospheric pressure: $860\text{hpa} \sim 1060\text{hpa}$
- Relative humidity: no more than 80%, no corrosive gas and well ventilated indoor
- The room should be well ventilated to avoid strong sunlight and corrosive gas intrusion
- Not close to the ground, walls and roof
- Transportation: Simple earthquake-proof facilities for product roommates, suitable for air, railway, road and ship transportation. Rain, snow, splashing, upside-down and collision should be avoided.
- If the storage period exceeds 3 months, you should open the box for inspection, recharge the machine, and use it again.

Warranty

- The warranty period of this product is one year;
- During the warranty period, the failure occurred under normal use in accordance with the instruction manual (judged by the company's official work night wish) free repair;
- During the warranty period, if one of the following situations occurs, it will be charged with repair:
* Cannot provide the essential warranty and valid proof of purchase;
* Failure and damage caused by misuse and improper maintenance;
* Think of damage caused by transportation, handling, or falling;
* Improper use causes damage to the equipment caused by water or other solutions;
- If the product does not meet the warranty, the company will charge the consumer for repairs;
- The warranty card is no longer issued separately, please keep it properly.

Warranty Card

Product name	Mode	
User information	Name	Contact
	Address	
Distributor	Malfunction	