



# INSTRUCTIONS MANUAL

TYL-TB9033

## SONIC ELECTRIC Toothbrush & UV Sanitizer

For a dentist-clean feeling every brush



[www.tylrhome.com](http://www.tylrhome.com)



## Table of Contents

<b>IMPORTANT SAFETY INSTRUCTIONS .....</b>	<b>2</b>
<b>KNOW YOUR SONIC TOOTHBRUSH .....</b>	<b>5</b>
Brushing Modes .....	6
Intelligent Timer .....	6
<b>CHARGING YOUR SONIC TOOTHBRUSH .....</b>	<b>7</b>
<b>USING YOUR SONIC TOOTHBRUSH .....</b>	<b>8</b>
<b>CLEANING YOUR SONIC TOOTHBRUSH .....</b>	<b>9</b>
<b>KNOW YOUR SONIC TOOTHBRUSH UV SANITIZER .....</b>	<b>11</b>
<b>CHARGING YOUR UV SANITIZER .....</b>	<b>11</b>
<b>USING YOUR UV SANITIZER .....</b>	<b>13</b>
<b>CLEANING YOUR UV SANITIZER .....</b>	<b>14</b>
<b>WARRANTY TERMS AND CONDITIONS .....</b>	<b>15</b>

## **IMPORTANT SAFETY INSTRUCTIONS**

When using electrical appliances, basic safety precautions should be followed, including the following:

### **DANGER!**

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc.
- Do not immerse the charge in water or any other liquid. After cleaning, make sure the charger is completely dry before connecting to wall socket.
- Do not plug or unplug the charger with wet hands.
- Keep the sanitizer away from water. Do not place or store the sanitizer over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the sanitizer in water or any other liquid. After cleaning, make sure the sanitizer is completely dry before you connect it to the mains.

### **WARNING**

- The main cord cannot be replaced. If the main cord is damaged, discard the charger. Always have the charger replaced within one of the original types to avoid a potential hazard.
- Do not use the charger outdoor or near a heated surface.
- If the appliance is damaged in any way (brush head,

toothbrush handle or charger), discontinue use.

- This appliance contains no user-serviceable parts. Contact customer service for repair.
- This appliance can be used by children and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Discontinue use of the sanitizer if the UV light bulb remains on when the top cover is open. UV light can be harmful to the human eye and skin. Always keep the sanitizer out of the reach of children.

## **CAUTION**

- Do not clean the brush head, handle, charger, and sanitizer in the dishwasher.
- Consult your dentist before using the toothbrush if you have had oral or gum surgery in the previous 2 months.
- Consult your dentist if excessive bleeding occurs after using this toothbrush, or if bleeding continues to occur after one week of use, or if you experience discomfort or pain when you use the toothbrush

- This toothbrush complies with safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the manufacturer of the implanted device prior to use.
- If you have any medical concerns, consult your doctor before using the toothbrush.
- This product has been designed for cleaning teeth, gums, and tongue. Do not use if for any other purpose.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace brush head every 3 months, or sooner if signs of wear appear. Do not use other brush heads other than the ones recommended by the manufacturer,
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpaste), thoroughly clean the brush head after each use.
- Check and make sure that the product has no signs of damage before use.
- Do not look directly at the UV light when it is turned on. Always close the top cover when sanitizer is in use.

**SAVE THESE INSTRUCTIONS. FOR HOUSEHOLD USE ONLY.  
NOT FOR COMMERCIAL USE.**

## KNOW YOUR SONIC TOOTHBRUSH



## **Brushing Modes**

**Clean:** For everyday cleaning

**White:** For whitening and removal of surface stains

**Polish:** For polishing and teeth brightening

**Gum Care:** Combines cleaning with additional action to gently massage and stimulate gums

**Sensitive:** Extra-gentle mode for sensitive teeth and gums

## **Intelligent Timer**

Your sonic toothbrush has a built-in timer for 2 minutes, with a brief pause at 30 second intervals. The 2-minute timer helps you brush for the dentist-recommended 2 minutes while the pauses every 30 seconds guide you when to move to the next section of your mouth. After the 4<sup>th</sup> interval, the toothbrush will automatically turn off.



## CHARGING YOUR SONIC TOOTHBRUSH

**Charge for 12 hours before first use.**

1. Place the charging base on a level non-skid surface, then plug into a power supply.
  2. Place the handle vertically on the charging base (see Fig.1). The charging indicator lights up in red while the brushing mode indicators alternately flash in blue. The charging indicator turns to blue when the toothbrush is fully charged.
- A fully charged handle should last around 180 minutes or 45 days (based on 2 brushings a day at 2 minutes each brushing).

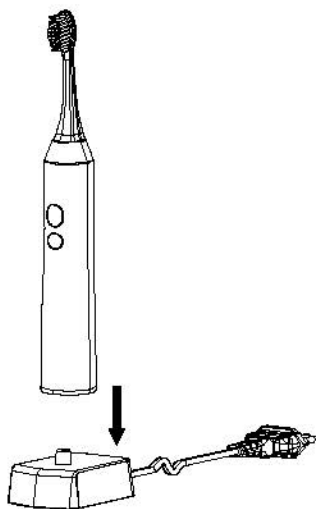


Fig.1

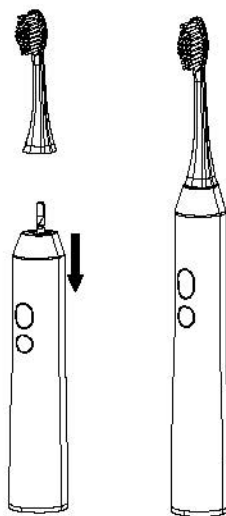


Fig.2

## USING YOUR SONIC TOOTHBRUSH

1. Align the brush head so the bristles face the same direction as the front of the handle. Firmly push the brush head down on the metal shaft of the handle until it locks in place. (See Fig. 2)

**NOTE:** *It is normal to see a slight gap between the brush head and the handle.*

2. Wet bristles with water, then apply toothpaste. To avoid splattering, place the toothbrush inside your mouth before turning the toothbrush on. (See Fig. 3)

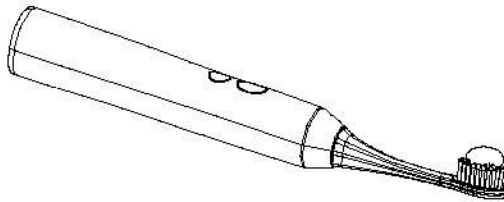


Fig.3

3. Hold the toothbrush horizontally. Place the toothbrush bristles against the teeth at a 45° angle so that bristles are directed towards the margin between the teeth and gums. (See Fig.4)

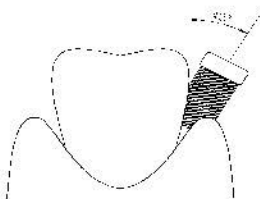


Fig.4

4. Press the Power button to turn on the toothbrush. Press the Mode button to select preferred mode.

Move the toothbrush horizontally in a back-and-forth motion, using short strokes – less than ¼". Be sure to keep the bristle points in contact against the gum margin. Brush each area for 30 seconds until the brief pause, then move brushing to the next area of your mouth. After the 4<sup>th</sup> interval, the toothbrush will automatically turn off.

**NOTE:** *When the handle is low in battery, the charging indicator will light up in red. Place the handle on the charging base to charge.*


## **CLEANING YOUR SONIC TOOTHBRUSH**

1. Rinse the brush head handle under running water after every use. Wipe dry with a clean cloth or towel.
2. At least once a week, detach the brush head from the handle. Rinse the metal shaft of the handle under running water, making sure to remove residual toothpaste.
3. Unplug the charging base from the power supply. Use a clean damp cloth to wipe the surface of the charging base.

**NOTES:**

- *Do not use abrasive cleaners or any sharp objects to clean the toothbrush.*
- *Replace the toothbrush head every 3 months or sooner if signs of wear appear.*




## KNOW YOUR SONIC TOOTHBRUSH UV SANITIZER

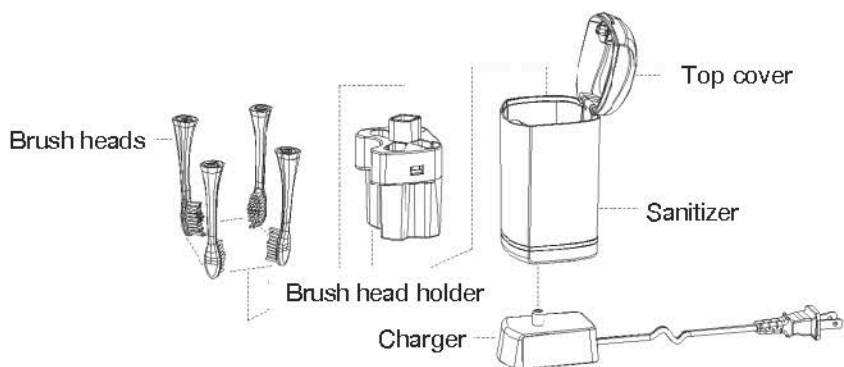

## CHARGING YOUR UV SANITIZER


**Charge for 12 hours before first use.**

1. Place the charging base on a level non-skid surface, then place the sanitizer on the charger. Make sure the charging hole at the bottom of the sanitizer aligns with the charging pole on the charger.
2. Connect the charger to a power source. The charging indicator lights starts blinking in red while charging, then turns to blue when charging is complete.

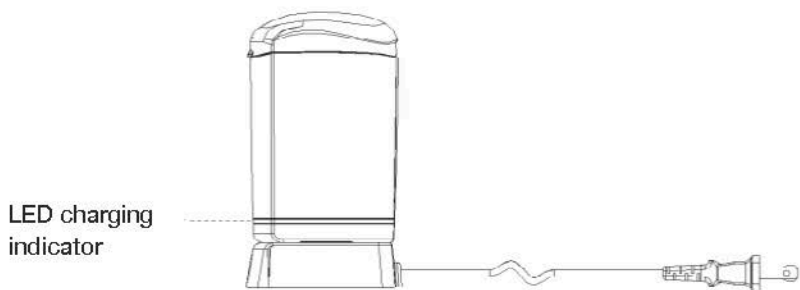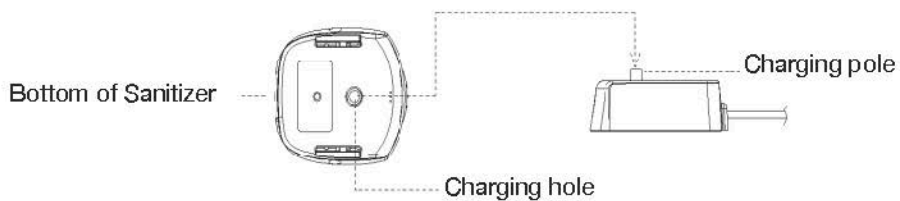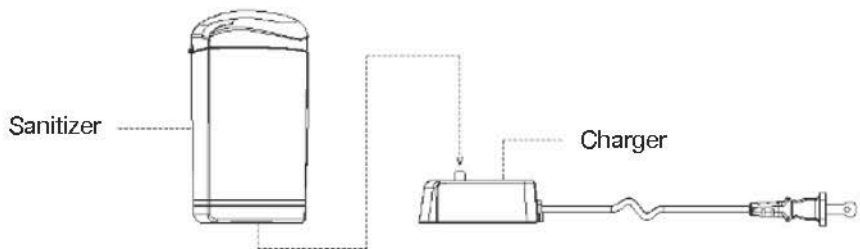


## USING YOUR UV SANITIZER

1. Place the sanitizer on a flat, level surface. Open the top cover and place the toothbrush head(s) head side down into the toothbrush head holder.

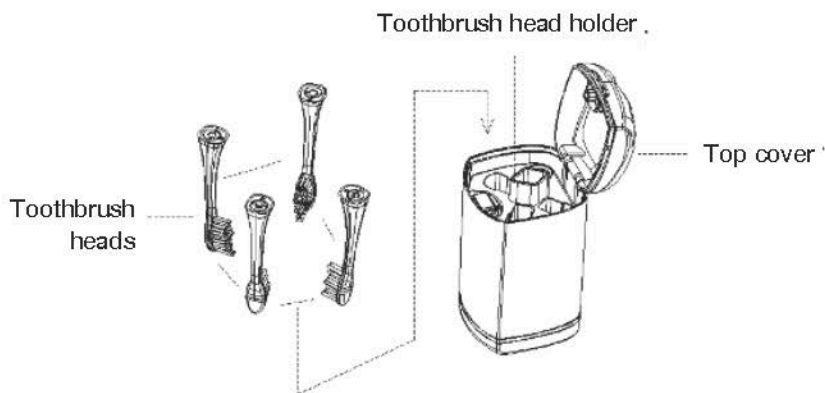


2. Close the top cover, then press the Power button. The LED charging indicator flashes in blue, signifying that the UV lamp is on and is sanitizing. The sanitizer will automatically shut off after 8 minutes, or when the top cover is opened.

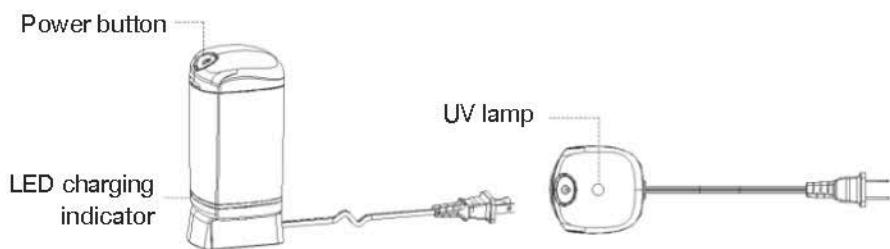

**NOTE:** When the sanitizer is low in battery, the LED charging indicator will flash in red. Place the sanitizer on the charging base to charge.

## CLEANING YOUR UV SANITIZER

1. Unplug and disconnect sanitizer from power supply before cleaning.
2. Open the top cover and take out the toothbrush heads out of the sanitizer.
3. Detach the upper cover from the lower cover and rinse both parts with water. Wipe dry with a soft cloth or towel.
4. Wipe the surface of the sanitizer with a soft cloth or towel. Do not rinse or immerse in water.

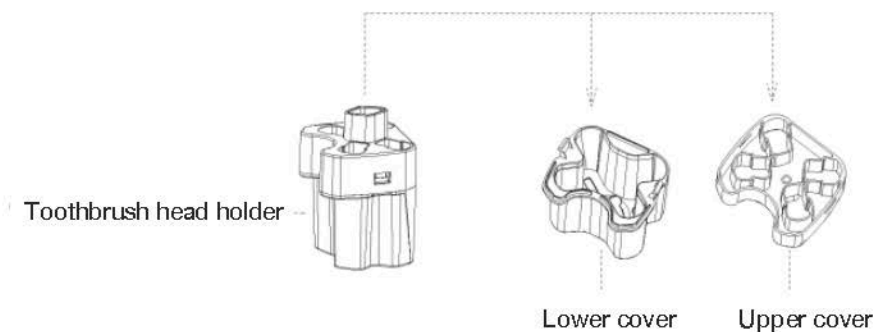



## **WARRANTY TERMS AND CONDITIONS**

### **ONE-YEAR LIMITED WARRANTY ON PARTS AND LABOR**

TYLR provides a warranty to the original purchaser of this product against defects in materials and workmanship for one year. If a product covered by this warranty is determined to be defective within the warranty period, TYLR will repair or replace the defective part with a new or remanufactured part free of charge subject to terms and conditions stated herein.

To obtain warranty service, proof of purchase in the form of a sales invoice or copy thereof is required to show that a product is within the warranty period. If any one of the documents cannot be produced, corresponding fees for labor and replacement parts will be charged.

#### **Parts and Labor**

There will be no charge for parts and labor during the warranty period. Replacement parts or units may be new or recertified and are upon TYLR's option and sole discretion. Replacement parts and units are warranted for the remaining portion of the original warranty. Our decision on all questions relating to complaints as a result of defects, either workmanship or materials, shall be conclusive and you shall agree to abide by such decision. Any unit or defective part that have been replaced shall become TYLR's property. After the lapse of this warranty, a charge will be made on all labor and replacement of parts.

## **Types of Services**

Defective Products must be sent to any **TYLR** service center to obtain warranty service. **TYLR** is neither responsible for transportation costs to the service center nor **TYLR** will cover return shipping to the customer. Product returns to **TYLR** service centers must utilize either the original carton box and shipping material or packaging that affords an equal degree of protection.

For home service, the transportation fee of **TYLR** service personnel will be charged to the customer, depending on the location. Our service personnel reserve all rights to refuse to attend to any equipment installed in an inaccessible location or any other hazardous situation.

## **LIMITATIONS AND EXCLUSIONS**

**TYLR** one-year limited warranty only covers defects in materials and workmanship; however, this warranty does not cover the following:

1. Damage resulting from accidents, transport, misuse, abuse, alteration, tampering, or failure of the purchaser to follow normal operating procedures outlined in the user's manual.
2. Damage caused by negligence and poor maintenance.
3. Defects or damage due to spillage of food/liquids or improper usage of electrical supply or electrical circuit, major fluctuations in voltage, exposure to rain or moisture, wrong fuel or oil used.
4. Damage, losses, and defects as a result of fire, flood, or other Acts of God.
5. Normal wear-and-tear, corrosion, rusting, or stains.

6. Defects & damage arising from improper testing, operation, usage of the wrong component, demonstration, maintenance, installation, adjustment, or any alteration or modification of any kind.
7. Scratches & damage to the outer surface areas and externally exposed parts that are due to normal customer use.
8. General maintenance & routine servicing.
9. Claims for damaged/missing parts after 7 days from the date of purchase.
10. If any part or parts of the unit are replaced with a part or parts not supplied or approved by us, or if the unit has been dismantled or repaired by any person other than a **TYLR** authorized technician.
11. Any equipment/product which has its serial number removed or made illegible/tampered with.
12. If the warranty card is altered, defaced, or erased in any manner whatsoever.
13. If the unit is used for purposes other than household use, such as commercial use.

The six (6) months warranty period covers the rechargeable batteries and the one (1) year warranty period covers the mainboard and control PCB, motor and power supply. However, this warranty does not cover parts (i.e. accessories, chargers base, bristles, UV lamp, fuse, gears, interior and exterior casing) that are subject to wear and tear, rust, stains, or corrosion.

This one-year limited warranty does not cover Products sold **“AS IS”**, **“FACTORY RECTIFIED”**, or by a non-authorized reseller.

THERE ARE NO EXPRESS WARRANTIES OTHER THAN THOSE LISTED OR DESCRIBED ABOVE. ANY IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE PERIOD OF TIME SET FORTH ABOVE. **TYLR'S** TOTAL LIABILITY FOR ANY AND ALL LOSSES AND DAMAGES RESULTING FROM ANY CAUSE WHATSOEVER INCLUDING **TYLR'S** NEGLIGENCE, ALLEGED DAMAGE, OR DEFECTIVE GOODS, WHETHER SUCH DEFECTS ARE DISCOVERABLE OR LATENT, SHALL IN NO EVENT EXCEED THE PURCHASE PRICE OF THE PRODUCT. **TYLR** SHALL NOT BE RESPONSIBLE FOR LOSS OF USE, COMMERCIAL LOSS, LOST REVENUE OR LOST PROFITS, OR OTHER INCIDENTAL OR CONSEQUENTIAL DAMAGES. THIS WARRANTY IS SUBJECT TO CHANGE WITHOUT NOTICE.

**IMPORTANT:** Please present sales invoice as proof of purchase whenever you require our service during the warranty period.

*[www.tylrhome.com](http://www.tylrhome.com)*